# Bluetooth Baby Scale

## Instruction Manual

English / Nederlands / Deutsch / Français

**moa**

color your life

[www.moacolors.com](http://www.moacolors.com)

## BEFORE USE


Before use, please read this manual carefully.

### BATTERY

Battery included before shipped, take away battery's insulating film before using.


### APP DOWNLOAD

Scan the QR code to download App or search for FITBRICK in the Appstore.
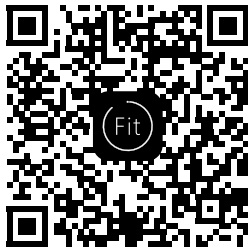
## DISPLAY

1. Display Screen 2. Hold button 3. Operating 4. Tare off 5. Switch On/Off
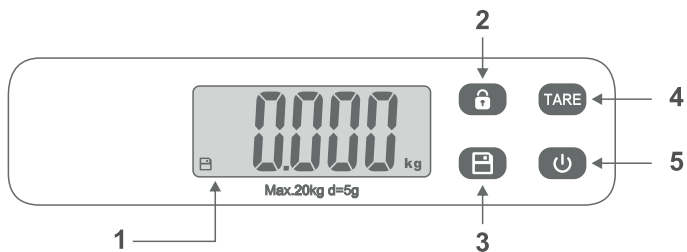
## INSTRUCTIONS (1/2)

### SETTING DATA IN THE APP

1. Place your scale on a hard and flat surface to ensure the accuracy.
2. Turn on the App, add right item and choose "Baby Scale".
3. Identification device, click "Turn on Bluetooth" button. (Skip this step if Bluetooth already turned on)
4. Press "On/Off" button to turn on the scale about 2 second When App on your phone identify the scale, please keep your baby scale turn on then click "confirm add " ,add scale successfully.
5. Data setting,please make sure set the real data(age and gender).

### WEIGHTING AND DATA

Weight baby after finished setting data (Weight data will show on the scale and the App simultaneously).

1. Weight directly:

Press "On/Off" button, place the baby on the scale carefully when "0.000kg" appears on the display.

2. Weight with container:

Press "On/Off" button, when "0.000kg" appears on the display, put the container on the scale (shows container's weighing data), then press "Tare weight" button, when "0.000kg" appears on the display again, place the baby into the container( shows baby's weighing data).

3. Hold Function:

**AUTO HOLD:** When the baby is placed on the platform, he/she may move around. When he/she is still enough, the AUTO HOLD function will establish the weight on display. The sign "HOLD" flashes and then weight is established.

**MECHANICAL HOLD:** If the baby is moving around too vigorously, you may select the mechanical hold function by pressing "HOLD" button. The sign "HOLD" flashes and weight is established on display. The sign "HOLD" remains on display until you press "HOLD" button. Then it shows "0.000kg" which means the scale is ready for next weighing.

## INSTRUCTIONS (2/2)

### 4. Data storage in scale function

When the data displayed on the screen, press "OPERATING" button, the weighing data will flash for some times, when stopped flashing, it means data storage successfully in the scale.

## PROMPT SYMBOL

### 1. Overload indication

When the scale is overloaded (over 20kg), the symbol will appear on the display, please put down the scale in case of damage to the scale.

### 2. Low Battery Indication

When the batteries need replacing, the symbol will appear on the display, please change to new batteries.

### 3. Unit Switch

The scale can only change unit on the App. For example, change unit from kg to lb on the App, by pressing the Synchronize button, please wait a few seconds, the scale unit will change to lb.

## TIPS

1. When use the APP, please open the Bluetooth and GPS
2. When first time use this scale, please download the App first, then put baby's data (Gender/Age). For other baby, please add new user name and fill the data, just choose the user name when next time use.
3. If the scale would not be used in a long time period, please take the battery out, and make sure there is nothing on the scale.
4. User need to add baby's height and head size in the App
  1. Open the App and click the user's photo in upper left corner.
  2. Enter the unit selection page.
  3. Choose the right unit (Weight: kg, lb; Height: cm, in).

## DISPOSAL

Do not dispose of the device in normal domestic waste.  
Dispose of the device over a registered waste disposal firm or through your communal waste disposal facility.  
Observe the currently valid regulations. In case of doubt, consult your waste disposal facility.

## WARRANTY

Dear Valued Customer,

Thank you very much for purchasing a MOA product.

We wish you to inform you that this product is covered by a warranty which complies with all legal provisions concerning existing warranty and consumer rights in the country where the product was purchased.

Should you find any defect or malfunction of your MOA product, please contact the appropriate Customer Care Center.

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