Meltazo®

Incredible Possibilities

MELTAZO Pro 6 IPL Hair Removal
Model: Ramavi Pro 6

PRODUCT FEATURES

■ **MELTAZO PRO 6 IPL HAIR REMOVAL DEVICE** is a new version of a permanent painless hair removal device with the latest semiconductor refrigeration technology that has 3-in-1 functions; HR (Depilation), SC (Skin Rejuvenation), and RA (Acne Treatment). It is ice-cooling hair removal that uses an advanced laser system to remove body hair painlessly and can give a full body hair treatment. It has up to 999,999 light pulses, and two head lamps, and is intended for HR, SC, and RA functions. It effectively uses light pulses to work on the hair roots and achieve permanent hair removal and offers you a similar experience to salon hair removal. It has 9 adjustable energy levels and offers the most comfortable, comprehensive, and effective treatment experience for both men and women. It has both manual and automatic flash modes designed to meet the needs of hair removal on different body parts. The manual mode is mainly for removing small areas of hair such as; armpits, fingers, bikini area, lips, and other small sensitive areas. The automatic mode is used for relatively large hair removal areas such as; arms, legs, back, and stomach. It's safe, easy to use at home, and ideal for a wide range of skin types.

- **The unique refrigeration technology makes the whole process of hair removal feel comfortable and removes the burning sensation**

MELTAZO PRO 6 IPL HAIR REMOVAL DEVICE has high strength, high-temperature resistance, high light transmittance, high thermal conductivity, and other valuable physical properties that can be combined with semiconductor refrigeration technology which can achieve a safe and comfortable ice-cooling painless hair removal experience. In the depilation process, it can reach 10°C when contacting the skin, effectively cooling the epidermis, reducing the skin burning damage, and the whole process is comfortable and painless.

- **The intensity pulsed light technology is perfectly applied to the home hair remover field**

MELTAZO PRO 6 HAIR REMOVAL DEVICE uses intense pulsed light (IPL) technology which is the most effective method known on the market to inhibit hair regrow over years, through close cooperation with skin experts, we have improved the optical technology originally used in professional beauty salons, so that you can use this technology safely, conveniently and effectively at home. The light of hair remover is soft, the irradiation is convenient, and the effective light intensity makes people feel very comfortable. After repeated use, excess hair will eventually disappear. Let you enjoy the pleasure of getting rid of hair, and glow with charm every day.

- **Long-life integral lamp holder**

This product has a long-lasting integrated lamp head that does not need to be replaced and is suitable for multiple hair removal areas such as; face, lip hair, underarms, bikini line, arms, legs, etc., saving you expensive lamp head replacement costs, and no need of purchasing accessories frequently.

SAFETY PRECAUTIONS

To better protect your skin, make sure to use the device close to the skin:

- Do not apply any gel or care products (aloe vera gel, etc.) before use.
- Incorrect operation or misusing the product can cause product damage or personal injury.
- ⓘ Dangerous-related content that can cause personal injury due to incorrect operation.
- Be sure to assemble and use the product according to the guidelines indicated in the "User Manual".
- ⓘ Warning-related content that can cause danger due to incorrect operation.
- 📖 Please read the "User Manual" carefully and keep it for future reference.
- ⓘ Attention-related content that can cause loss of personal safety or property due to incorrect operation.

ⓘ DANGEROUS

- ⓘ **Warning:** Do not use this product near bathtubs, shower basins, washbasins, or other utensils that hold water. Please keep the product and the power adapter dry.

Do not use it in the bath or shower.

- ⓘ If the product has been damaged, do not touch any internal parts to avoid electric shock.
- ⓘ Please do not insert any object into this product and do not replace any part of this product unintentionally.
- ⓘ Do not use this product around your eyes. Wearing protective glasses is a common operation to protect the user. However, if you operate around your eyes, the protective glasses may not be able to completely block the energy, and the light will be absorbed by the iris melanin in your eyes at close range.
- ⓘ We recommend that you wear protective glasses for every operation. Although the light does not directly affect the eyes, the high-frequency flash light may affect your vision.
- ⓘ If you need a back hair removal or ice compress operation, you may need an assistant or partner, so it may be risky for you to use it alone.
- ⓘ If your hair is longer than 1mm, please shave first. We recommend that you use a shaving knife/trimmer. Excessively long hair will cause the skin surface to absorb a lot of heat and cause skin redness and swelling, and burned hair will produce a peculiar smell that will affect the user's product experience.
- ⓘ If the skin of the irradiated part is not dry and clean, it will affect the expected outcome, and may also cause the product to not flash properly.
- ⓘ Please pay attention to preventing children under 15 years old from touching this product, and also please keep away children under 15 years old when using this product.
- ⓘ Please use different light intensity for different parts of the body as required because the skin color of different parts of the body is not the same.
- ⓘ If the power cable is damaged, it must be replaced by professionals to avoid danger.

REMARK

- ❗ Please do not apply antiperspirant immediately after irradiating the armpit.
- ❗ Please check and clean the dirt of the irradiation port regularly to avoid scald and malfunction.
- ❗ Please do not place the irradiation port near items that are easy to degauss, in order not to cause failure and damage.
- ❗ Please do not use it when the vent is blocked, in order not to cause smoke and fire.
- ❗ Please do not point the irradiation port toward a black object to avoid machine malfunction.
- ❗ Please do not irradiate clothes, hair, and other items to avoid blackening, discoloration, and breakage.
- ❗ Please do not use it for pet depilation to avoid injury.
- ❗ Please do not use this product around the eyes to avoid eye damage.
- ❗ Please wipe the hair and skin care products left on the hair remover device after use to avoid malfunction and skin damage.
- ❗ Please do not use it around children. The hair remover device will emit strong light and cause eye injury.
- ❗ Please set the irradiation intensity according to your skin condition. Please make sure to use first gear for the first time.
- ❗ It is a normal phenomenon if the skin becomes slightly hot and red after using it. Please turn on the ice/cold compress function to relieve skin discomfort.
- ❗ The product can't be left unattended after the boot, after use, please shut down in time.
- ❗ Cleanse the skin to ensure that the skin surface is free of hair residue, completely dry, and free of oily substances and blemishes.
- ❗ If the ambient temperature difference reaches more than 20°C (from a low-temperature environment to a high-temperature environment, or vice versa), please turn off the device and let it stand for about 2 hours before using it again.
- ❗ Please do not use it near flammable and explosive materials (including sprayers and gas stations), which may cause fire, electric shock, and other hazards.
- ❗ The ambient temperature of this epilator is 0-40°C. Please check the room temperature before use.
- ❗ Please do not touch the irradiation port after use, because the irradiation port produces a high temperature and can easily cause burns.
- ❗ Please avoid direct ultraviolet radiation on the irradiated part within 24 hours after irradiation.
- ❗ Please do not place the product in an environment with high temperature and humidity, as it may cause deterioration and malfunction of the product.
- ❗ If the built-in fan does not operate after turning on the power supply (can't hear the buzzing sound), please do not use it, so as not to cause smoke or fire.

CONTRAINDICATIONS

- Please do not apply antiperspirant immediately after irradiating the armpit.
- Severe heart disease, moderate to severe hypertension, malignant tumor, folliculitis, infectious skin disease, burned skin, and anticoagulant disorders are prohibited to use it.
- Patients with abnormal collagen metabolism, including a history of scarring and poor wound healing are prohibited to use it.
- People with vascular diseases such as; varicose veins and vasodilation at the irradiation site are not allowed to use it.
- Patients with photosensitive dermatitis are not allowed to use it.
- Please do not use it if the irradiation site on the skin is infected, herpes simplex, skin lesions, hematoma, etc.
- Patients with diseases that suppress the immune response (including HIV infection, lupus erythematosus, porphyria, etc.) are not allowed to use it.
- It is not allowed to be used for any patient with hemorrhagic disease. People who use isotretinoin, fruit acid (AHA), salicylic acid (BHA), or azelaic acid on the skin cannot use this product.
- Please use caution if you have used chemical hair removal within 30 days, skin whitening products, or cosmetics containing vitamin A in the short term.
- Do not use essential oils at the same time. Make sure the depilated area on the skin is clean, dry, and free of creams, lubricants, and other cosmetics.
- If the skin surface is not clean and dry, it may affect the usage experience of the product, and even may damage the product.
- This product is forbidden to use while taking aspirin, immunosuppressive drugs, and isotretinoin.
- Those who have recently used light-sensitive drugs (such as; tetracycline, sulfa, etc.) are not recommended to use this product.
- Light damage (sunburn), mucous membranes, and other sensitive parts are prohibited from using this product.
- Do not use this product on unhealed wounds.
- It is not recommended for patients with hirsutism and other diseases-induced skin hirsutism.
- If you take painkillers, it will reduce the skin's sensitivity to heat. It is recommended to use it carefully before the effect is relieved.
- Please do not use it if your skin or hair color is not suitable.

DO NOT USE FOR THE FOLLOWING BODY PARTS OR THE BELOW CONDITIONS

1. Oral cavity, lips, vagina, mucous membrane, etc.
2. Plastic surgery and the surgical site.
3. Adam's apple.
4. Eyeball, eyelids, and around the eyes.
5. The head.
6. Acne area.
7. Spot and mole sites.
8. Eyebrows, hair, forehead.
9. The ear.
10. Nipples and areola.
11. Inside your belly button.
12. Large veins.
13. Tattoos and body painting.

BE CAREFUL

Do not use this product for naturally dark skin.

Please do not insolate before and after depilation.

Do not use in a dark environment.

Please always wear protective glasses when using this device.

Before hair removal, please shave first (do not use chemical hair removal method).

For the first time, make sure you start in first gear.


Do not use in a direct sunlight environment in order not to cause device failure.


Do not slam and shake the product.

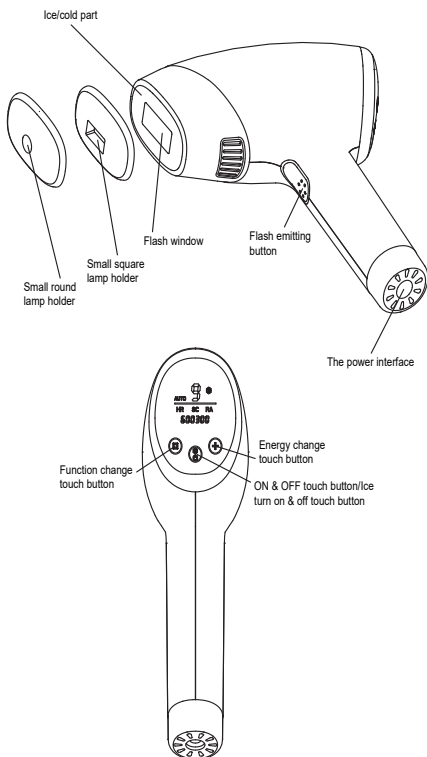PRODUCT STRUCTURE

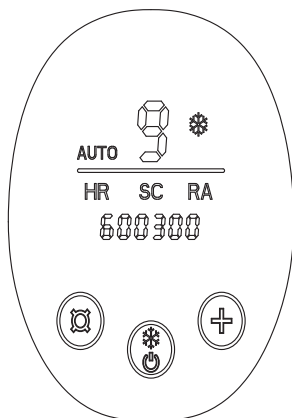- ※ It is strictly forbidden to block the air inlet and outlet with debris, otherwise, it will cause serious damage to the product.
- ※ Do not slam and shake the product.
- ※ You can use the flash window for all 3 functions; HR, SC, and RA.

CONTROL PANEL

LCD Screen	Function	Parameter
HR	Hair remover	999,999
SC	Skin care	
RA	Remove acne	
	Ice	
AUTO	Auto light mode	
1-5	Light wave intensity gear display	1-5
Six figures	Number of remaining	
	Function change button, available in "HR", switch between "SC" and "RA"	
	Long press this button for 2 seconds to turn on or off the device, or short press to turn on or off the ice function.	
	The energy change can be adjustable from 1 to 9 gears.	

PRODUCT PARAMETERS

Type of light source	Intelligent pulse light	Power supply	100-240V,50-60Hz
Light area	4 square centimeters	Control panel	Sense of convex button
Main unit size	162*64*207	Package size	362*252*72
Total weight	1.3kg	Work environment	Temperature 20 - 60°C, Relative temperature 20 - 90%, No condensation environment
Main unit weight	0.38kg		

ACCESSORIES

Main unit: is automatically affixed with 3,6c m ² Cap which can be used for larger areas such as; legs, arms, back, stomach, and other larger areas.	1 Piece
Main unit + square lamp holder: is 1,8 cm ² Cap which can be used for medium areas such as; bikini, armpit, and other medium areas.	1 Piece
Main unit + round lamp holder: is 1,4 cm ² Cap which can be used for small areas such as; fingers, chin, face, upper lips, and other small areas.	1 Piece

PREPARATION BEFORE USE

- Clean skin entirely to ensure that skin is hair-free, completely dry, and free of the residue of oily substances (perfume, lotion, sunscreen, etc.).
- Make sure that there is no dirt left on the irradiation port, no damage or abnormality. If there is dirt, please make sure to clean it, so as not to cause product failure, deterioration, and skin problems.
- Please do not apply any gel or care products (aloe vera gel, etc.) before hair removal. This product has the function of ice/cold compress and hair removal at the same time. It can take care of the skin by making it more comfortable and painless simultaneously.

THE IRRADIATION INTENSITY SETTING

This product has been specially developed to be used at home for hair removal, so its light energy is softer and offer 5 different light energy. You could adjust the light energy to suit your skin and body hair color and keep yourself at a comfortable level. The recommended light energy scale below will help you determine the appropriate light energy.

1. Please refer to the table below to determine which light intensity is suitable for your skin and hair color.
Tip: You could put the skin tone card on your skin and compare it with your skin tone.
2. If you feel any pain or discomfort, lower the energy gear. People with darker skin are more likely to develop skin allergies because the skin absorbs more energy. Therefore, it is recommended that people with darker skin tones use a lower energy scale. When the pulse of light shines on the skin, it may feel warm or hot, but it should not be painful.

Suitable skin color, hair color, and gears contrast table						
SKIN COLOR HAIR COLOR	White	Yellow	Light brown	Medium brown	Dark brown	Black dark
White/gray	1-3	1-3	1-3	1-3	×	×
Red	1-3	1-3	1-3	1-3	×	×
Light gold	3-5	3-5	3-5	1-3	×	×
Light brown	5-9	5-9	5-9	1-3	×	×
Medium brown	5-9	5-9	5-9	1-3	×	×
Dark brown	5-9	5-9	5-9	1-3	×	×
Black	5-9	5-9	5-9	1-3	×	×

Note: X is not applicable.

Note: Please refer to this table to select the appropriate gear.

Note: Due to various reasons, your skin may react differently due to different exposure times/situations.

Note: no matter what skin color, it is recommended to start from the first gear and gradually increase the gear.

SKIN SENSITIVITY TEST


If you are using this product for the first time, please select a small area of your skin first and conduct skin sensitivity according to the following requirements.

Testing:

- 1: Refer to the above-mentioned information to start the device.
- 2: Choose gear 1 or 2.
- 3: Follow the instructions and precautions to operate the selected skin area.
- 4: After 30 minutes of the skin test, check skin reaction.
- 5: If the tested area of your skin is normal, then you can continue to use it. If abnormal, please stop using it.

INSTRUCTIONS

1. Before use, please shave, clean and wipe the hair in the depilated area. Make sure that the hair follicles under the skin absorb the maximum energy for optimal hair removal.
2. Connect the power adapter, then connect the power, and wear the protective eyeglasses equipped with the product. The background light of the LCD screen of the hair removal device is ON, and the screen state is "OFF".
3. Long press the start button for 2 seconds to turn on or off the device, and short press to turn on or off the ice/cold function, then the hair removal device will start, and the current working state will be confirmed by referring to the previous "screen presentation".
4. Press the "Function change" touch button to enter the depilation mode you want when the LCD screen shows HR, SC, and RA. Then you can use any of these lamp holders according to your skin area for depilated.
5. The machine is equipped with 1-9 gear, short press the "Gear selection" touch key to adjust the intensity of the light wave. Please always start using the first gear, and increase gradually the gear according to your own skin experience. The higher the gear, the better the effect.

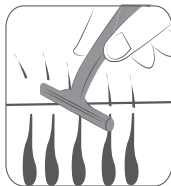6. Manual light setting (suitable for depilation of small areas).

The default mode of the hair removal machine is to manually emit light. Press the flash emission button once and emit light once.

Automatic light setting (suitable for large area depilation).

After the hair removal device is started up, long press the flash emitting button, and the device will switch to automatic light output, and the "AUTO" on the LCD screen will keep flashing. You only need to make the lamp holder vertically close to the skin, and the device will automatically emit light.

7. Apply the light output clinging to the area where hair removal is required and use it vertically.

- When the device is in manual flash mode, the content on the screen will start flashing, then press the light emitting button to flash, and then the hair removal is completed in the area covered by the light output. Then, you move the lamp head to the next hair removal area and cling to the skin, when the content on the screen starts flashing, press the flash emitting button to emit the flash, then change over to the next hair removal area, and so on. It is recommended to irradiate each part 2-3 times.

- When the device is in automatic flash mode, long press the light emitting button, and the flash will be emitted, the area covered by the light output has been depilated, then move the light output to the next depilation area and cling to the skin. The device will emit flashing light and then change over to the next hair removal area and so on. It is recommended to irradiate each part 2-3 times.

Tip: With the automatic flash mode, the light output is completely clinging to the skin. While you are sliding the light output on the skin surface, it can flash continuously for 8 seconds.

8. If you want to use the ice-cooling/cold compress function, just press the "ON/OFF" button, "☒" will appear on the LCD screen, and once again press the "ON/OFF" button, the "❄" on the LDC screen will disappear, then turn off the ice-cooling/cold compress function.

9. After finishing the depilated process, long press the "ON/OFF" button for 2 seconds to shut down or the device will shut down automatically if it has no action for 5 minutes.

10. Unplug the adapter from the device and remove the adapter from the power socket.

11. After 10 minutes, when the surface of the lamp cap cools, wipe the lamp cap light window with a non-alcoholic soft towel or cloth.

12. After cleaning, put the device and adaptor back in the packaging box. Place the device in a ventilated, dry place, and place it out of reach of children.

13. Care of the depilated area on the skin. After 20 minutes of depilation, moisturizing skin care products can be selected to moisturize and repair the depilation part. If you feel uncomfortable, you can apply the ice-cooling/cold compress for 30-60 minutes until the skin feels no burning pain.

FUNCTION SWITCHING

1. Follow the procedure of the device shown on the LCD screen.
2. Ice-cooling/cold compress function: press the "ON/OFF" button, and the snowflake icon appears on the LCD screen, the ice-cooling/cold compress function is turned on.
3. Hair removal & skin care & acne removal function: simply press the "Function change" button, and the LCD screen will switch between "HR", "SC" and "RA". Then you can use one of the lamp holders according to the area to be depilated on the skin.

3 FUNCTIONAL SUGGESTIONS FOR HAIR REMOVAL / SKIN REJUVENATION / ACNE REMOVAL


Only one functional photon can be irradiated in each part of the body in one day. If other functional photons are needed in the same part of the body, the depilation process shall be carried out the next day. This machine is equipped with an extra round and square small lamp holder, which can be used for medium and small parts of the body for specific hair removal. When you want to use one of the lamp holders, just put the lamp head on the flash exit of the main machine, and it will automatically adsorb on the main unit, and you can use it directly.

SKIN CARE

- If the skin becomes red after hair removal by light irradiation, please use the skin care mode for nursing. It is recommended to avoid touching the water on the depilated area within 6 hours after hair removal.
- Do not use perfume, facial masks, and cosmetics containing abrasive pearls on the skin after the hair removal process.
- Avoid UV exposure for 24 hours after hair removal. Apply sunscreen if necessary.
- Consult your doctor if you feel unwell.
- Sunbathing is likely to tan the skin or cause sunburn on the top of the skin. IPL irradiation at this time will cause redness, inflammation, and other conditions. Therefore, please be careful not to expose yourself to the sun for a long-term during the period that you still need to undergo the hair removal process. The recovery time from sunburn should not be less than 1 week.
- It is strictly forbidden to expose the depilated area to the sun for a long-term within 48 hours of the depilation operation. Even if the depilation area needs to be exposed after 48 hours, always pay attention if the skin feels a bad condition. Please use SPF30+ sunscreen for sun protection after 72 hours.
- The use of artificial light is not allowed to tan the skin (Change the skin color to bronze color through physical or chemical formula.) during hair removal.
- If you need to undergo plastic surgery after hair removal (only the hair removal area), please consult your doctor first.

Note: Even if you follow all the instructions, you may still experience skin reactions. In this case, please stop using the device and consult your doctor.

RECOMMENDED USE PERIOD

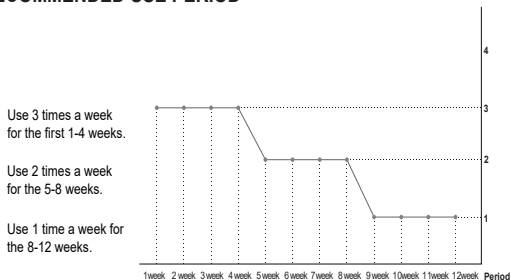- ▲ Please use it multiple times according to the growth of your hair. For intense hair, you can increase the gear and frequency accordingly. The specific situation varies from person to person.

CLEANING AND MAINTENANCE

To ensure the effectiveness and durability of the product, please clean the product before and after each hair removal process and even if necessary during the hair removal period. If you don't clean it properly, this product will not get the right effect.

- ◆ Inspect the product, especially the transparent window and body parts.
- ◆ Make sure they are all clean and tidy.
- ◆ Do not wash the product and its accessories under the faucet or in the dishwasher. Do not use steel wire balls, or abrasive cleaners (such as; toothpaste or corrosive liquid gasoline, oil, acetone, etc.) to clean the product.

Note: Please stop using this product if you can no longer clean the transparent window.

- ◆ Turn off the product, unplug the power plug and allow it to cool.
- ◆ When using this product, the flash window will get hot, make sure it cools down before cleaning.
- ◆ If necessary, clean the outside of the product with a clean soft towel or cloth.
- ◆ A wet soft towel can be used to clean the flash window.

Note: If the wet soft towel is not effective, a few drops of high concentration alcohol (no more than 96%) can be used to clean the above-mentioned areas.

ENVIRONMENTAL PROTECTION

When disposing of this product, please do not dispose of it with normal household waste. You should give it to the designated recycling center. This will be good for our environment.

ELECTROMAGNETIC FIELD (EMF)

This product complies with all applicable standards and regulations relating to exposure to electromagnetic fields.

HAIR REMOVAL ESSENTIALS

Before using the machine, you should shave the visible hair first and then use a toner or wet towel to wet the hair removal area to avoid discomfort caused by dry skin.

■ The gist of each part's depilation

Please understand the following points when removing hair with the IPL hair removal device:

• Underarm hair removal

- ※ Before the hair removal process, please wipe the anti-perspirant, perfume, cream, and other cosmetics applied to the armpit.
- ※ For areas with loose skin, stretch the skin with your fingers and then remove the hair. It is normal for the skin to have a slight burning sensation at the radiation port.

• Bikini line hair removal

- ※ For areas with an accumulation of pigment and dark colors, please start using the low gear first, and keep using gear 1 until the skin adapts.

• Hands and feet hair removal

- ※ Joints and uneven areas are more prone to redness and pain than other parts, but this is normal.

▲ Do not use it for the following parts:

There are burn marks, dark colored parts, eyeballs, eyelids, around the eyes, chilblain, moles and other parts, eyebrows, forehead, sunburns, and other parts.

- ※ After hair removal, avoid the place with ultraviolet rays for 24 hours to prevent the sun from burning the skin.
- ※ Be careful not to look directly at the flash because this will damage your eyes.

PRECAUTIONS

■ Possible Side Effects

- Slight redness of the skin or a slight itching and warmth are normal, harmless, and will soon go away.
- After hair removal, the reaction of the depilated area for a short time is similar to that of exposure to the sun. If it does not disappear after 48 hours, you are advised to consult a doctor.
- Shaving and hair removal process together at the same time will cause dry and itchy skin, this phenomenon is a normal reaction, the discomfort will gradually reduce, and you can apply an ice-cooling compress to that part. After 24 hours, you can apply a mild, non-irritating moisturizer.

■ RARE SERIOUS ADVERSE REACTIONS

1. **Burns, severe redness, and swelling pain:** These reactions are rare, usually caused by the use of light intensity that is too high for skin color. If these reactions do not disappear within 3 days, it is recommended that you consult a doctor. Please wait until the skin is healed before performing the second irradiation, and make sure to use a lower light intensity.
2. **Skin discoloration:** This phenomenon occurs rarely. The skin discoloration usually appears darker or lighter than the surrounding area. This is caused by the intensity of the light used being too high for the skin tone. In rare cases (Usually caused by the use of light intensity that is too high for skin color), the effect of pigment changes may be longer than normally observed. If the skin discoloration does not disappear within two weeks, it is recommended that you consult a doctor. Please irradiate the discolored area after the discoloration disappears and the skin returns to normal skin tone.
3. **Skin infections and inflammations:** This is rare and can be caused by using this product on skin with wounds created while shaving the hair.
4. **Extreme pain:** If you use this product on unshaved skin, the light intensity you use is too high for your skin color; flashing on the same area many times; having open wounds, inflammation, infection, tattoos, or using this product on skin with burns and other areas may cause this reaction during or after hair removal.

PROBLEMS AND SOLUTIONS

This section outlines the problems, possible reasons, and solutions that you encounter with the product. If the following information is not available to resolve the problem, contact the customer service center.

PROBLEMS	POSSIBLE REASONS	SOLUTIONS
Unable to boot	The power supply is not properly connected	Check if the power is plugged in
No light	Perhaps the main unit head and body temperature are too high, or the main unit head is close to the skin	Wait for the device to cool down naturally before using it again. Then, place the main unit against the skin & wait for the light indicator to turn on
The indicator light is not ON	Bad device	Check if the power supply is working properly, repair it immediately or return the product
There is a peculiar smell when flashing, and the hair removal area turns red	Filters are dirty, and contaminated, and hair remains on the body surface	Wipe the filter lens, shave the hair, and reduce light intensity
Unbearable pain during hair removal	The hair on the area that needs to be removed is not shaved	Before using this product, please shave the hair on the area to be irradiated. If shaving causes skin irritation, please cut the hair as short as possible and use the first gear and gradually increase to the next appropriate gear
	The area being depilated is not suitable for using this product	Do not use this product on the inner labia, vagina, anus, nipple, and areola. Men should not use this product on the scrotum
	Excessive light intensity	Reduce the light intensity to a comfortable level
	The light window is damaged	Repair immediately or replace the product
The product becomes hot and humid during use	Perhaps because the device worked for a long time	Turn off the device to rest, and dissipate the heat completely before use

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MELTAZO is registered trademark with premium products.