
Lifeproducts

Live a Better Life

Shiatsu Nek- en Rug Massagestoel

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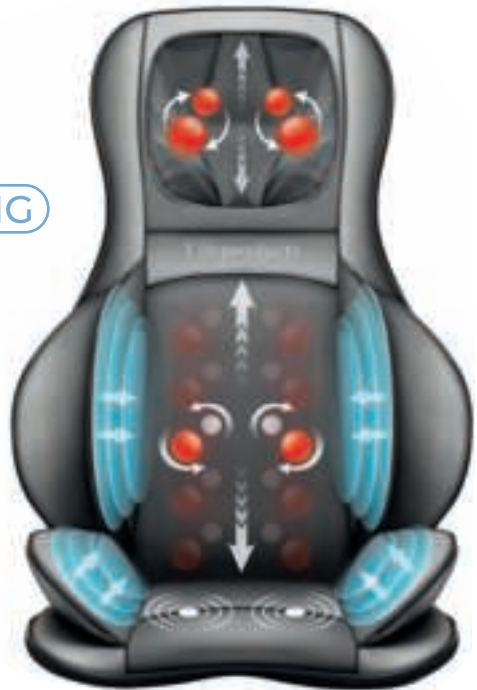FC CE RoHS
MADE IN CHINA

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1 Veiligheidsinformatie

2 Handige Informatie

3 Bediening

4 Diversen

5 Garantie

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Lifeproducts

Live a Better Life

Air Compression Shiatsu Neck and Back Massager

MANUAL

FC CE RoHS
MADE IN CHINA

Contents

1 Safety Instructions

2 Useful Information

3 Operating

4 Miscellaneous

5 Warranty

IMPORTANT INFORMATION! RETAIN FOR FUTURE USE!

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

Explanation of symbols

This instruction manual belongs to this device. It contains important information about starting up and operation. Read the instruction manual thoroughly. Non-observance of these instructions can result in serious injury or damage to the device.

WARNING

These warning notes must be observed to prevent any injury to the user.

CAUTION

These notes must be observed to prevent any damage to the device.

NOTE

These notes give you useful additional information on the installation or operation.

Protection category II

LOT

Lot number

Manufacturer

1 Safety Instructions

Safety Instructions

about the power supply

- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- Only connect the mains plug to the wall socket when the device is switched off.
- Only use the mains adapter that is provided with the massage cushion.
- Keep the power supply, the cable, the control unit and the massage cushion away from heat, hot surfaces, moisture and liquids. Never get hold of the power supply with wet hands.
- Do not take hold of a device if it falls into water. Pull out the power plug and the power supply unit from the socket immediately.
- Always unplug the equipment immediately after use.
- To disconnect the unit from the power supply, always pull the plug out of the socket. Do not pull on the cable!
- Do not carry, pull or twist the massage cushion on the cable, the power supply unit or control unit.
- You must not continue to use the unit if the lead or the mains adaptor are damaged. For safety reasons, these parts may only be replaced by authorized service personnel. In order to prevent hazards, always send the unit in to the service centre.
- Ensure that the cable will not cause tripping. The cable must not be kinked, trapped or twisted.

specific groups of people

- The surface of the appliance is hot. People who are insensitive to heat must be careful when using the appliance.
- Children must be supervised to ensure that they not play with the device.
- Do not use this unit to supplement or replace medical treatment. Chronic complaints and symptoms may be made worse.
- You should refrain from using the AIR COMPRESSION SHIATSU NECK AND BACK MASSAGER, or should con-sult your doctor before using it, if
 - you are pregnant,
 - you have a cardiac pacemaker, artificial joints or electronic implants,
 - you suffer from one or more of the following illnesses and complaints: Circulation disorders, varicose veins, open wounds, bruises, broken skin, inflammation of the veins.
- If you experience any pain or discomfort during the massage, stop immediately and consult your doctor.
- Never leave the unit unattended if it is connected to the power supply.
- Take care when using the heater. Do not use if you have circulatory problems or on areas not sensitive to heat. Risk of burning!
- Consult your doctor if you have any questions concerning the therapeutic use of massage seats.
- If you have unexplained pain, are receiving treatment from a doctor and/or use medical equipment, consult your doctor before using the massage cushion.

1 Safety Instructions

before using the unit

- Always carefully inspect the power supply, the cable, the control unit and the massage cushion for damage before use. Do not attempt to use a defective unit.
- Do not use the unit if there are visible signs of damage on the unit or the cable components, the unit does not work properly or the seat or the control unit have been dropped or become damp. In order to prevent hazards, always send the unit in to the service centre.

operating the unit

- The massage cushion must only be used for its intended purpose as described in the instruction manual.
- Using the massage cushion for any other purpose invalidates the warranty.
- Never leave the unit unattended if it is connected to the power supply.
- The maximum operating time for one session is 15 minutes.
- Do not allow the massage cushion to come into contact with pointed or sharp objects.
- Never put the unit down next to an electric stove or other sources of heat.
- Never use the massage cushion if it is folded up or crumpled.
- Do not stand on the unit.
- Do not touch the power supply unit when you are standing in water and always take hold of the plug with dry hands.

**The device is not to be used for commercial use or medical purposes.
If you have health concerns, consult your doctor before using the massage cushion.**

Only use the massage cushion indoors.

Do not use the massage cushion in damp rooms (e.g. when bathing or showering).

maintenance and cleaning

- The unit is maintenance free. If a fault does occur, just check whether the mains power adapter is properly connected.
- You may only clean the device yourself. If the device malfunctions please do not repair it yourself as this leads to the automatic expiry of all guarantee rights. Ask your specialist dealer and have repairs carried out by the service.

1 Safety Instructions

- Do not use the unit if there are visible signs of damage on the unit or the cable components, the unit does not work properly or the massage cushion or the control unit have been dropped or become damp. In order to prevent hazards, always send the unit in to the service centre.
- Do not immerse the device in water or other fluids.
- If, nonetheless, liquid should penetrate the device, unplug the mains cable immediately.
- Cleaning and user maintenance shall not be made by children without supervision.

Do not wash!

Do not dry clean!

2 Useful Information

Thank you very much


Thank you for purchasing Comfier Air Compression Shiatsu Neck and Back Massager .
With normal care and proper treatment, it will provide years of reliable service.
Please read all the instructions carefully before using this product.
Please retain this manual for future reference.

2.1 Features

- Electrically height-adjustable Shiatsu neck massage
- Rolling massage with adjustable width control
- Customize Spot Massage
- With optional heat function on back massage
- Massage areas can be selected (upper, lower, full back)
- Seat vibration massage (3 intensity levels for options)
- Air massage on back and seat (3 intensity levels for options)
- Optics PU leather and breathable mesh

2.2 Contents

- Air Compression Shiatsu Neck and Back Massager
- 12V home adpater

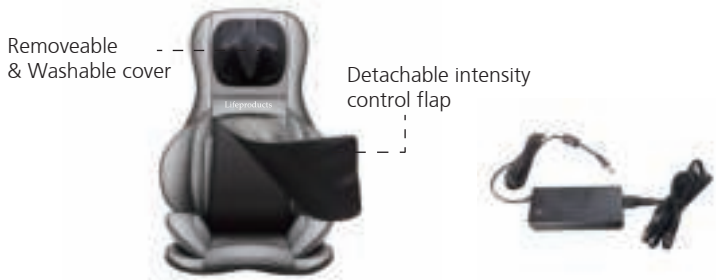
WARNING

Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!

3 Operating

1. Attach the massager to a chair with elastic straps or place it on other support.
2. Connect the adapter cable to the corresponding cable on the cushion.

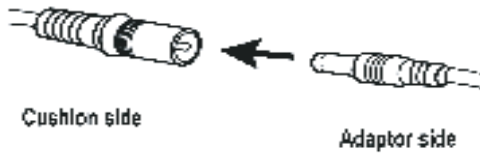3. Plug the home adapter into an electric outlet.
4. Turn on the device using the controller (instructions on page 18).
5. When finished place the controller in the pouch at the side of the cushion.
6. Detachable intensity control flap for back lets you choose a softer or more intense massage.
7. Detachable cover for neck lets you choose a softer or more intense massage.


Note :

If the Shiatsu is too strong,

- ① Install the flap onto the cushion to get a softer massage.
- ② Place a blanket or towel between you and the massager to achieve a far softer massage intensity if you still feel too strong shiatsu after installing the flap onto the cushion.

If the shiatsu can't reach your neck,shoulder,lower back area

- ③ If you want to massage your neck, shoulder,lower back, but the shiatsu nodes can not reach the area in need. You may adjust the height of your back when sitting, by placing a seat cushion or a pillow on the sitting part.

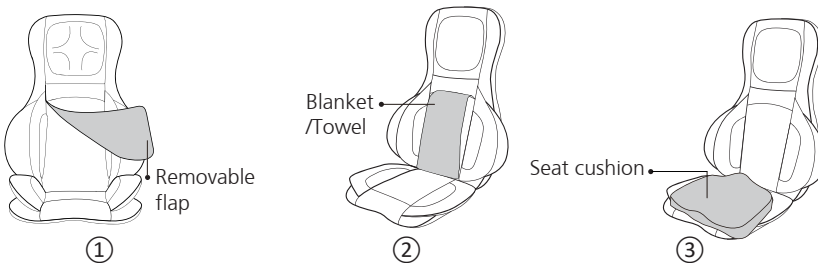- The device will automatically switch off after 15 minutes timer run out
- Remove the adapter from the outlet when not in use.
- Do not setup or use the massager in a bathroom or similary wet / damp areas.

!Attention

- You should allow your body some rest periods. To avoid straining your muscles too much, we recommend to not exceed a continuous massage of 15 Minutes per time and not exceed 3 times use one day.

3.1 Controller Instructions

Vibration Button:
Press to turn on the seat Vibration, you may cycle through High, Medium, Low or off settings

Back Shiatsu Button:
press the direction button to switch the rotation

Spot Button:
Press to keep shiatsu function in one area. Press **arrow** buttons to adjust the shiatsu nodes position up or down.

Air Intensity Button:
Press to set both air massage intensity at Low, Medium or High level

Air Back Button:
Press to turn on or off air massage on the back

Shiatsu Area Button:
Select Full, Upper or Lower back to designate the area covered by the shiatsu massage.

Auto Button:
Quick demonstration of the different functions. The L.E.D. lights will blink to indicate the active functions. Once complete, the device will automatically turn off.

Neck Shiatsu Button:
Press to turn on Shiatsu neck massage, Press again to shift rotation direction, Press once again to turn off
Press **arrow** buttons to adjust the shiatsu nodes position up or down

Rolling Button:
press to activate Rolling function, the massage nodes travel up & down across the muscle gently. Continue pressing the button to adjust the width of the massage nodes to custom fit the body.

Air Seat Button:
Press to turn on or off air massage on the seat

Back Heat Button:
Press to turn on or off Heat
The heat can't be activated if back shiatsu function is not turned on

Power Button:
Press to turn on or off the massager

4 Miscellaneous

4.1 Cleaning and maintenance

1. Do not wash in washing machine or immerse the product or the adapter in water.
2. Always switch off before cleaning and remove the adapter from the socket with clean, dry hands.
3. The massager should be wiped with a slightly moist cloth and dried thoroughly. Do not wipe the product with alcohol or harsh chemicals such as thinner or bleach.
4. The Cover for neck is detachable and washable.
5. Do not overstretch the elastic back straps.
6. Keep away from heat, fire or exposure to direct sunlight.
7. Do not apply excessive pressure or impact onto massager/massage heads. Application of a direct force in excess of 110 lbs. may damage your massager and invalidate your warranty. Refrain from over usage of the device.
8. Store the massager in a cool, dry place out of reach of children.
9. Avoid contact with sharp edges or pointed objects which might cut or puncture the cover of massager.

4.2 Disposa

This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

4.3 Technical Data

Name and model	: Air Compression Shiatsu Neck and Back Massager CF-2307A
Voltage supply	: 100-240 V~ 50/60 Hz
Power supply unit	: 12V = 4000 mA
Power consumption	: 48 W
Timer	: 15 minutes
Operating conditions	: only use in dry rooms
Storage conditions	: in a cool, dry place
Dimensions	: approx. 81 x 56 cm (Back) 53 x 38 cm (Seat)
Weight	: approx. 9.5 kg (without power supply)
Adapter cable length	: approx. 1.8 m

In accordance with our policy of continual product improvement, we reserve the right to make technical and optical changes without notice.

5 Warranty

If you have any issue about the product, please feel free to contact us by sending email to Lifeproducts@hotmail.com. We will strive to provide the best service possible within 24 hours.

30 days unconditionally Return

Lifeproducts product can be returned to receive a full refund for any reason within 30days. Please contact our customers service (Lifeproducts@hotmail.com), our staff will contact you within 24 hours.

90 days refund/replace

Comfier product can be returned / replaced within 90 days if the product break down in the period of proper use.

24 months Warranty

If the product break down within 24 months in the period of proper use, customers can enjoy the relevant warranty to get them replaced.

Attention!

No warranty will be given to any force majeure and man-made causes to a defective product, such as improper care, personal tear down and willful damage,etc.

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