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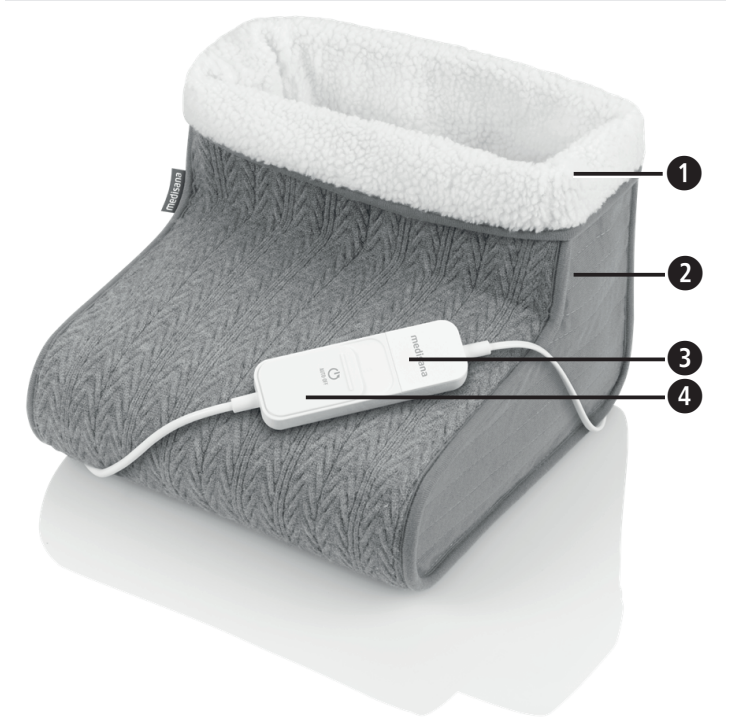
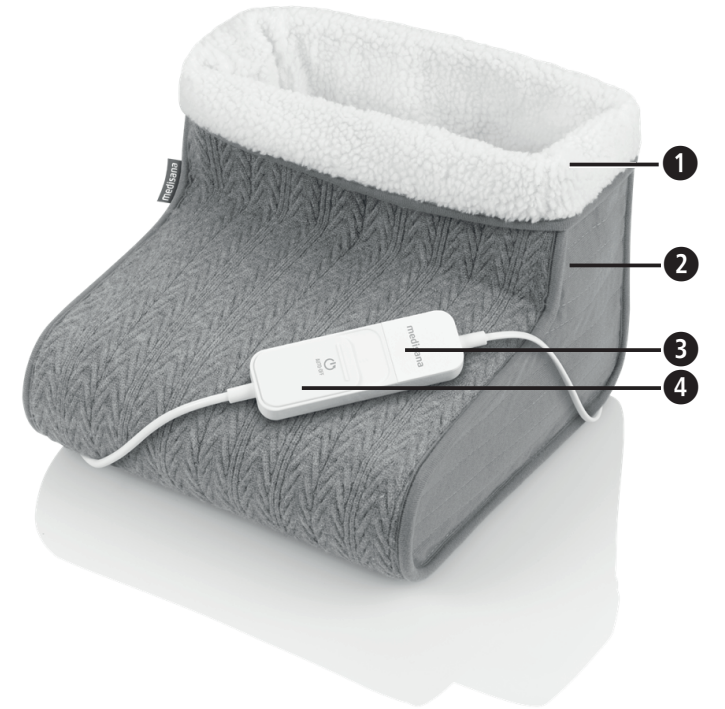
Explanation of symbols

 **IMPORTANT**  
Non-observance of these instructions can result in serious injury or damage to the device.

 **WARNING**  
These warning notes must be observed to prevent any injury to the user.

 **CAUTION**  
These notes must be observed to prevent any damage to the device.

 **NOTE**  
These notes give you useful additional information on the installation or operation.

 Only use the massage cushion indoors.

 **Protection category II**

 **LOT** **Lot number**

 **Manufacturer**

**IPX1** Drip-proof (does not apply for control unit)  
IPX 0 for control unit

 Do not puncture the foot warmer!

 Do not use the foot warmer when it is folded!

 Not to be used by very young children (0-3 years)!

 Do not wash the foot warmer!

 Do not chlorinate!

 Do not dry the foot warmer in a tumble drier! The cover can be dried in a drier!

 Do not iron the foot warmer!  
The cover can be ironed at a low temperature!

 Do not dry clean!

 The cover (interior) can be washed in a normal wash at max. 30 °C!

**GB Safety instructions**  
 **Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.**

- 
- The device must only be used for its intended purpose as described in the instruction manual. Using the appliance for any other purpose invalidates the warranty.
  - The device is not to be used for commercial use or medical purposes.
  - Prolonged use of the foot warmer at a high setting could lead to skin burns.

- about the power supply**
- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
  - Only insert the plug into a socket when the appliance is switched off.
  - Only use the control unit supplied for the appliance.
  - Keep the plug, cable and appliance well away from heat, hot surfaces, moisture and liquids.
  - Do not touch the plug if you are standing in water or if your hands are wet or damp.
  - Do not take hold of a device if it falls into water. Pull out the power plug from the socket immediately.
  - Only plug the device into a socket where the plug is freely accessible.
  - Always switch off the appliance after use with the sliding switch on the control unit and unplug the power plug.
  - To disconnect the unit from the power supply, always pull the plug out of the socket. Do not pull on the cable!
  - Do not carry, pull or twist the device by the cable.
  - You must not continue to use the unit if the cable or plug are damaged. For safety reasons, these parts may only be replaced by authorized service personnel. In order to prevent hazards, always send the unit in to the service centre.
  - Ensure that no one can trip over the cable. The cable must not be kinked, trapped or twisted.

- specific groups of people**
- Do not use the device on children or persons who are disabled, asleep or insensitive to heat (persons who cannot react on overheating).
  - Children below the age of 3 years are not allowed to use this device, because they are not able to react on overheating.
  - The device must not be used by young children over the age of 3 years, unless the control unit has been adjusted accordingly by its parents or legal guardians or unless the child has been well instructed how to use the control unit properly.
  - This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
  - Children shall not play with the appliance.
  - Do not use this unit to supplement or replace medical treatment. Chronic complaints and symptoms may be made worse.
  - Before using the device, consult your doctor regarding any unexplained pain, if you are undergoing medical treatment or using medical equipment.

- before using the unit**
- Check the power plug, cable and massager carefully for damage before each use. Do not attempt to use a defective unit.
  - Do not use the appliance if it, or the power cord show signs of damage, if it is not in perfect working order, or if it has been dropped or become damp. In order to prevent hazards, always send the unit in to the service centre.
  - Do not keep or place any electrical device on a surface where it is in danger of falling into a bathtub or sink.
  - Only use the appliance in the proper place, as described in this instruction manual.
  - Do not use the appliance while you are asleep or when you are lying in bed.
  - Prolonged use of the foot warmer at a high setting could lead to skin burns.
  - Outdoor shoes must be removed before use.
  - The device must not be used for warming animals.

- Never store or use the appliance next to an electric stove or other sources of heat.
  - Never leave the unit unattended if it is connected to the power supply.
  - Do not allow the device to come into contact with pointed or sharp objects.
  - Do not cover up the appliance when it is switched on. Never use it under blankets or cushions. There is a risk of fire, electric shock or injury.
  - Do not use an external time switch or a separate remote control system to operate the appliance.
  - Do not use the device whilst using aerosol sprays or simultaneously using oxygen.
- maintenance and cleaning**
- The unit is maintenance free. If a fault does occur, just check whether the mains power adapter is properly connected.
  - You may only clean the device yourself. If the device malfunctions please do not repair it yourself as this leads to the automatic expiry of all warranty rights. Consult your specialist dealer and have repairs carried out by the service centre.
  - Do not immerse the device in water or other fluids.
  - If, nonetheless, liquid should penetrate the device, unplug the mains cable immediately.
  - Cleaning and user maintenance shall not be carried out by children without supervision.
  - Avoid bending the pad too much as objects placed on it during storage may cause damage.

**Intended use**  
The device is used to warm your feet for relaxation and well being. It is not intended for medical therapeutic use. It is not suitable for clinics, massage practices, wellness studios etc. The device is not a medical device and is not intended for the treatment of pain. In case of doubt or uncertainty, consult your doctor before use.

 **WARNING**  
**Please ensure that polythene packing is kept away from the reach of children! Risk of suffocation!**

- Device and controls**
- ❶ Cover (interior)
  - ❷ Foot warmer
  - ❸ level indicator light
  - ❹ Sliding switch (switches device on/off, level 0 - 1 - 2 - 3)

- Scope of supply**  
Please check first of all that the device is complete and is not damaged in any way. In case of doubt, do not use the it and contact your dealer or your service centre. The following parts are included:
- 1 **medisana** foot warmer **FW 150** including control unit
  - 1 Instruction manual
- If you notice any transport damage during unpacking, please contact your dealer without delay.

**Operation**  
Always make sure that the foot warmer is not creased or folded while it is in use. You must always use the foot warmer with the cover.

- Insert the power plug into the mains socket. Make sure that it remains easily accessible.
- Switch on the heat function of the unit by moving the slide switch ❹ to position “1”. If necessary, select one of the higher levels 2 or 3.
- Slide both switches to position “0” if you no longer wish to use the unit. Disconnect the power supply unit from the socket.

**After 90 minutes the unit will automatically switch off.**

- Cleaning and maintenance**
- Before cleaning the unit, make sure it is switched off and disconnected from the mains. Clean the unit once it has cooled down using a soft, damp cloth.
  - Under no circumstances should you use aggressive detergents, strong brushes, solvents or alcohol, as they may damage the surface of the unit.
  - When cleaning, never immerse the unit in water and make sure that no water gets into the unit.
  - Wait until the device is completely dry before using it again.
  - Undo any kinks in the cable.

- You should ideally store the device in its original packaging and keep it in a clean, dry place.

- Insert cover (interior)**
- Before cleaning the foot warmer, the insert cover must be taken out.
  - The cover can be cleaned according to the care instructions on the labels.

 **Disposal**  
This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a  municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

**Technical specifications**  
Name and model : **medisana** Foot warmer **FW 150**  
Current supply : 220 - 240 V~ 50/60 Hz

Consumption : 100 Watt  
Dimensions LxWxH : approx. 30 x 30 x 24 cm  
Weight : approx. 0.7 kg  
Storage conditions : clean and dry  
Item number : 60259  
EAN number : 40 15588 60259 7


**In accordance with our policy of continual product improvement, we reserve the right to make technical and visual changes without notice.**

The current version of this instruction manual can be found under [www.medisana.com](http://www.medisana.com)

**Warranty and repair terms**  
Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.  
The following warranty terms apply:

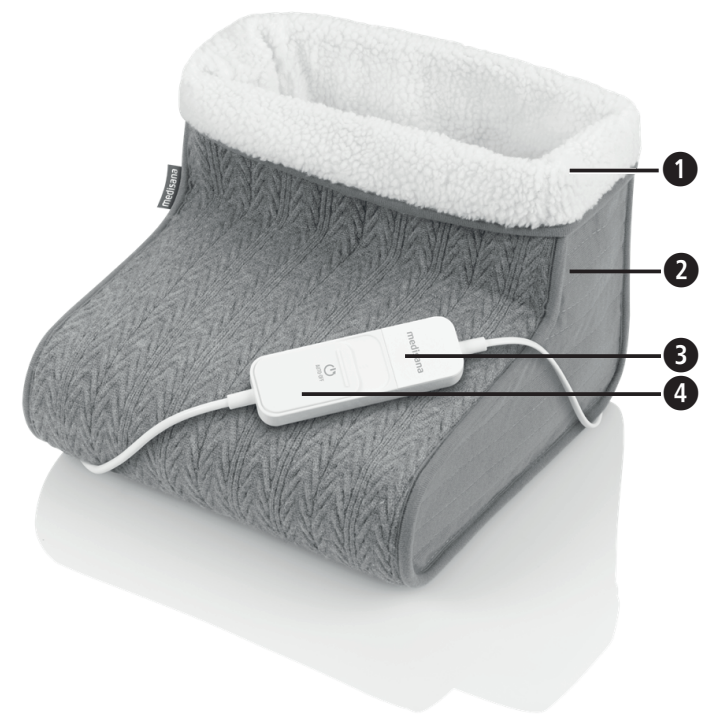1. The warranty period for **medisana** products is three years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
  - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
  - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
  - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
  - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

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**medisana GmbH**, Jagenbergstr. 19, 41468 NEUSS, DEUTSCHLAND