Lifeproducts
Live a Better Life

Massage Gun Elite

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Index

Dutch	2-8
English	9-16
French	17-24
German	25-32
Polish	33-40
Spanish	41-48
Swedish	49-56

1.	Introduction	10
2.	Warnings + safety instructions	10
3.	Received from shipment	12
4.	Contents packaging + product features	12
5.	Instructions for use	13
6.	Massage heads explanation	14
7.	Problems and solutions	16
8.	Maintenance and cleaning	16
9.	Dispose and recycle	16
10.	Service + Warranty	16

1. Introduction

Thank you for choosing one of our Lifeproducts products! **Please follow the safety instructions at all times for safe use.** Should there be any questions about the product, please contact us via email: klantenservice@lifeproductsstore.nl

The Lifeproducts Massage Gun Elite is suitable for the entire body. The massage gun has 8 different massage heads, 9 different speeds and 2 different settings. It charges in 3 to 4 hours and then can be used for 5-6 hours.

2. Warnings + safety instructions

- ⚠ Make sure that the packaging material does not get into the hands of children. They may suffocate!
- ⚠ Read the instructions carefully before use and keep them in a safe place!

General

- Be aware of the risks and consequences of improper use of the massage gun. Misuse of the massage gun can result in personal injury and damage to the unit or products connected to it.
- Use the massage gun only for the purposes indicated in this manual. If the massage gun does not function properly, is dropped or damaged, DO NOT use it and contact us by e-mail.
- The massage gun is designed and intended for private use. Do not use it for commercial purposes.
- DO NOT leave the massage gun unattended when it is in use.
- Inspect the massage gun and battery carefully before each use.
- The massage gun is suitable for indoor use only.

Persons and usage

- This device is NOT suitable for persons under 18 years of age or for persons with physical or mental disabilities. Use of the appliance is permitted only when a supervisor is present, who is responsible for the safety of the user.
- Use and store the massage gun out of the reach of children under the age of 18. Supervise children so that they do not play.
- Do not use the device in the same place for more than 180 seconds and not more than 30 minutes at a time. Stop using immediately if abrasions or other discomforts occur.
- Keep the massage gun away from your hair.

- Keep your fingers away from the unit's telescopic pliers and the cooling holes on the back of the motor to avoid entrapment.
- Keep the cooling holes on the back of the motor free of obstructions. Do not insert objects into the cooling holes of the motor.
- Use high speed and pressure only on the soft tissue of the body to avoid personal injury.
- Do not use the device if you are tired, have taken medication, or have consumed alcohol.

Medical Notes

- This device is NOT a substitute for treatment from a physician or a physical therapist.
 - Consult a physician before using the massage gun if you have unexplained pain, are undergoing medical treatment, or are using medical equipment.
 - Do NOT use the device: If you are pregnant; If you have a pacemaker or other electronic implants, artificial joints, high blood pressure, heart disease, cerebral atherosclerosis, cerebral thrombosis, broken bones, malignant tumors, circulation problems, varicose veins, open wounds, bruises, broken skin, inflammation of the veins, jaundice, diabetes, nervous disorders, acute inflammation or other health problems. If you have unexplained pain; On the abdomen, bones (such as joints and spine), head, on or around the throat, on or around the genital area, on or around the kidneys or other sensitive or weak areas;
 - On burnt or sensitive skin, swollen or inflamed areas, areas of poor circulation, broken skin, or any other area that is sensitive or weak.
 - On unconscious or sleeping people; In the same place for a long time, this may cause symptoms or pain; In bed.
 - Stop using immediately and consult your doctor if you experience pain or discomfort while using the massage gun, or find the massage unpleasant.
- DO NOT fall asleep while using the Massage Gun.

Electrical

- DO NOT disassemble, repair or restore the massage gun. This may cause damage to the unit. For your own safety, these parts should only be replaced by authorized service personnel. Always send the massage gun to a service center for repairs to avoid damage and personal injury.
- Before connecting the charger and the massage gun, make sure that the mains voltage indicated on the type plate matches that of your mains supply.
- DO NOT use the massage gun near a bath, a shower, a swimming pool, or near a (wash) tub filled with water. No liquids, heat sources, dust or dirt should be allowed to penetrate the motor's cooling vents.

- Do not use the device near highly flammable substances, The massager may generate a spark which could lead to fire or explosion.
- Do NOT allow the device to come into contact with water. DO NOT immerse the device in water. If water does get into or on the unit, turn the unit off immediately and contact our customer service department.
- Do not use the massager on wet skin. DO NOT touch the appliance or the power cord with wet hands.
- Keep the device away from moisture, heat, rooms where the temperature is relatively high and heat sources such as ovens.
- Do not cover the massage gun when it is on. Never place pillows or blankets on or against the unit. This can lead to personal injury or damage to the unit.
- Charge the device only with the Lifeproducts charger provided.

Other

- Choose the right massage head according to the massage needs of different- the parts of the body.
- Do not expose the device to direct sunlight.
- Do not use the massage gun in combination with massage oil, lotion, ointments, etc.

3. Received from shipment

Check that the contents of the package match the packing slip when you receive the shipment. Inform your supplier if any parts are missing. If the equipment appears damaged, file a complaint with the carrier immediately and provide the supplier with a detailed description of the damage. Keep the damaged packages to substantiate your complaint.

4. Contents packaging

Lifeproducts Massage Gun Elite - MA-08-NL/EU - EAN 8719327429662

Accessoires

- 8x Attachments
- 1x EU-adapter
- 1x Carrying bag
- 1x Instruction manual

Dimensions packaging	290x251x79mm
Weight	0,81 KG
Color	Black
Material	ABS
Battery	Lithium-Ion 2600mAh
Battery life	5-6 hours
Charging time	3-4 hours
Power supply	Frequency power: 100V 240V AC 50/60Hz 30W (charger) 16,8V, 0,8A
Speeds	9 speeds: 20/30/40/50/60/70Hz from 1200 - 3000 rpm
Different modes	2
Overheat protection	Ja

5. Instruction for use


- ⚠ • DO NOT use the massage device on wet skin.
- ⚠ • DO NOT use the device on the head, abdomen, joints on bone parts or on sensitive or weak parts of the body.
- ⚠ • Stop using immediately if you experience pain or discomfort.

Use

1. Press the power button on the back of the unit for 3 seconds. The device is in position 1.
2. To increase or decrease the speed, briefly press the same button you use to turn the unit on/off.
3. To turn the device off, press and hold the power button for 3 seconds. The device will also turn off after pressing the power button at the 9th position.
4. The battery level is displayed on the right side on the back of the massage gun. Charge the massage gun for 4 hours for a full battery.


The Massage Gun Elite has a special **medi-function**. This is shown on the left side of the massage gun. The medi-position has 2 speeds. You can switch to the medi-position by pressing the on/off button 2x quickly.

What is the medi-mode? The medi-mode is a special program where the gun automatically speeds up and slows down. So you experience a wonderful therapeutic massage. You have 2 positions, position 2 is faster and therefore more intense. The fitness mode has 9 speeds from 1400 to 3200 RPM! In this mode the speeds remain the same.

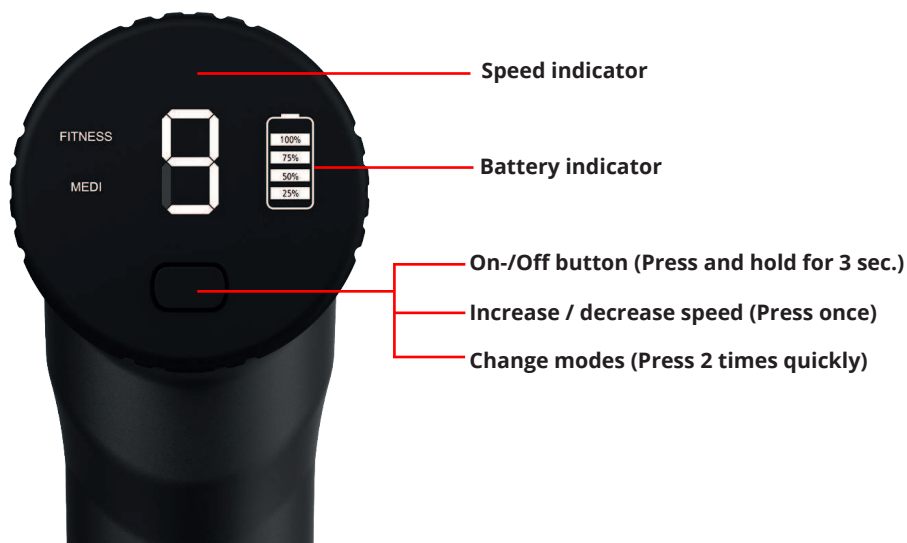
Change massage head

1. Pull the current massage head from the massage gun

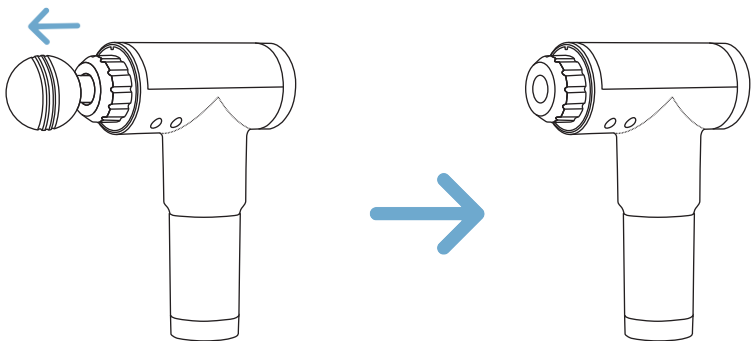2. Click the desired massage head into the massage gun

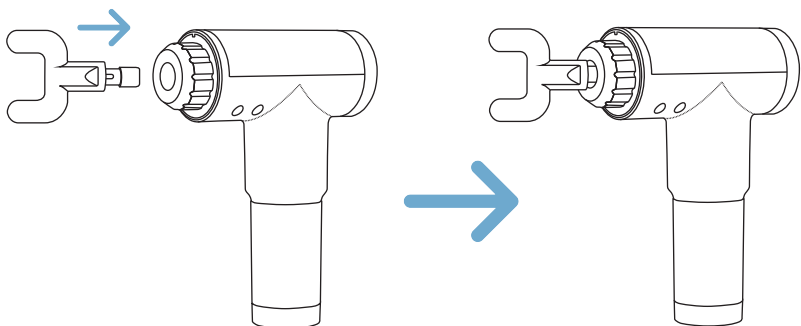6. Massage head explanation


Round head

Suitable for all muscle groups


Flat round head

Suitable for the whole body

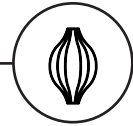U-shaped head

Suitable for the back

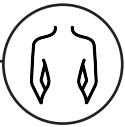Spherical head

Suitable for meridian acupressure points and soles of the feet

Air cushion head

Suitable for the whole body.
Reduces muscle tremors.

Flat head with massage points

Suitable for the whole body

in addition to the above massage heads, there are 2 more special metal heads:

Ball shaped metal head

Suitable for people with sensitive skin. Can be used with massage oil. This stimulates the massage allowing you to enjoy ultimate relaxation!

Flat round metal head

Suitable for people with sensitive skin. Can be used with massage oil. This stimulates the massage allowing you to enjoy ultimate relaxation!

7. Problems and solutions

1. The device stops working after 20 minutes of operation

The overheating protection has been activated. Restart the device after 5 minutes.

2. The device turns off during use

The massage heads are pressed too hard on the skin, so the device temporarily stops working. Restart the device by pressing the power button for 3 seconds.

8. Maintenance and cleaning

- Wipe the surface of the device with a slightly damp cloth and dry the device with a clean and dry cloth.
- Store the massage gun in a dry and clean place.
- Do not store the device in extremely high or low temperatures.

9. Dispose and recycle

Think of the environment and contribute to a cleaner living environment! At the end of its life, do not throw away the massage gun with normal household waste, but hand it in at a government designated collection point to be recycled.

10. Service & warranty

If you need service or information about your product, do not hesitate to send an e-mail to klantenservice@lifeproductsstore.nl.

We provide a 1 year warranty on our products. The customer must be able to show proof of the purchase and the warranty must be valid within a period of 1 year from the date of purchase.

Please note that the warranty does not apply to: Damage caused by improper use or repair; for parts subject to wear; for defects of which the customer was aware at the time of purchase; for defects caused by the customer's own fault, for damage caused by third parties.

Disclaimer

Specifications may change without reason.

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Lifeproducts B.V.