
**Rechargeable Sonic Power Toothbrush**  
**Ultrasoon Elektrische Tandborstel**  
**Operating Instructions**  
**Handleiding**  
Model: DT20

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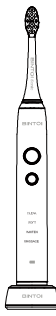**IMPORTANT SAFEGUARDS**

The Power toothbrush is safe and effective when used according to directions. Use all electrical appliances, however, basic safety precautions should be followed, especially when children are present.

**DANGERS-WATER**

To reduce the risk of burns, electrocution, fire, or injury to persons:  
 • Do not place or store the Charger base where it can fall or be pulled into a bathtub, toilet, or sink.  
 • Do not place or store the Charger into water or other liquid.  
 • Do not use the product that has been taken into water or other liquid, including immediately.  
 • Do not use any hair-follicle.  
 • Unplug charger base before cleaning and assure that it is dry before plugging into an electrical outlet.

**DANGERS-VOLTAGE**


To reduce the risk of electric shock, electrocution, fire, or injury to persons:  
 • Do not plug the charger base into an outlet with a voltage other than specified on the label of the Charger.  
 • This product is designed to operate within a range of 100-240 volts.  
 • Voltage conversion is not guaranteed compatibility.  
 • Do not wire the plug into an outlet.  
 • The plug does not fit easily into the outlet, do not force it.

**DANGERS-WARNINGS**

To reduce the risk of burns, electrocution, fire, or injury to persons:  
 • Do not use the product if it is not working properly, if it has been dropped or damaged or dropped into water, immediately discontinue use of the product especially damage in any way.  
 • Do not use the product if it is not working properly, if it has been dropped or damaged or dropped into water, immediately discontinue use of the product especially damage in any way.  
 • Do not use the product if it is not working properly, if it has been dropped or damaged or dropped into water, immediately discontinue use of the product especially damage in any way.

**CONFIGURATION**

Professional brush head  
 Color ring  
 Power button  
 Mode button  
 Charging indicator  
 Charger base  
 Power cord  
 Power plug

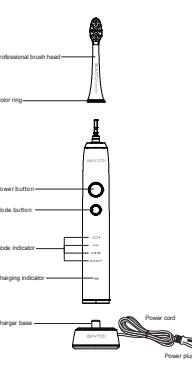**SAVE THESE INSTRUCTIONS**

**GETTING STARTED**

**Charging the toothbrush**

Be sure to turn off the toothbrush before charging. Plug the power plug into the outlet. Then insert the toothbrush into the charging. The charger indicator light will be on, indicating the toothbrush is charging. After the battery is fully charged, the charger indicator light will turn off, indicating the toothbrush is fully charged.

**HOW TO USE**

Slide the brush head onto the handle. Push brush head down until it locks into place on handle. There is a small gap between color ring and the handle (Fig 1).  
 Use the toothbrush and apply a small amount of toothpaste.  
 Apply toothpaste to the brush head (Fig 2).  
 Push the mode button to select mode.  
 If the toothbrush has not been used for a long time, it may need to be recharged. There are 4 modes: a. CLEAN: Standard brushing with polishing. b. SOFT: Gentle brushing with polishing. c. WHITE: Whitening brushing with polishing. d. MASSAGE: Pulse with gentle gum stimulation.

**TRUBLESHOOTING**

Problem	Possible cause	Solution
The toothbrush does not work.	You have not purchased the correct model for your country.	Change the main unit for the correct model.
The toothbrush does not work properly.	The toothbrush is pressed against your teeth.	Press the brush head gently against your teeth.
The toothbrush does not work properly.	The battery has reached its capacity.	Contact service center for support.
The toothbrush does not work properly.	Changing the time of use.	Contact service center for support.

**WARRANTY**


Your toothbrush is guaranteed for 1 year against manufacturer's defects. The warranty is not valid if the product is damaged due to misuse, accident, or other causes. The warranty is not valid if the product is damaged due to misuse, accident, or other causes. The warranty is not valid if the product is damaged due to misuse, accident, or other causes.

**HOW TO USE**

Apply light pressure and gently move the brush head across the teeth in a semi-circle motion.  
 • Apply light pressure and gently move the brush head across the teeth in a semi-circle motion.  
 • Apply light pressure and gently move the brush head across the teeth in a semi-circle motion.

**MAINTENANCE**

To clean the plug from the electrical outlet before cleaning the product:  
 • To clean the plug from the electrical outlet before cleaning the product:  
 • To clean the plug from the electrical outlet before cleaning the product.

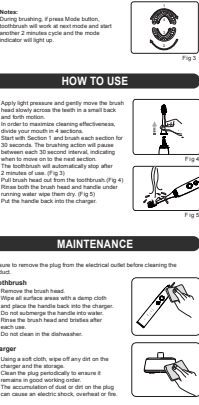
**DISPOSAL**

The general advice for this product is to recycle. Please do not dispose of the product in the trash. Please do not dispose of the product in the trash. Please do not dispose of the product in the trash.

**BINTOI**

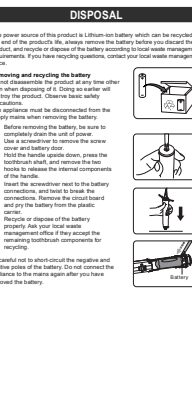Remove and recycle the battery.  
 • Remove and recycle the battery.  
 • Remove and recycle the battery.

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