
melodii

UV-C & Ultrasonic Denture Cleaner Clean Like Never Before!

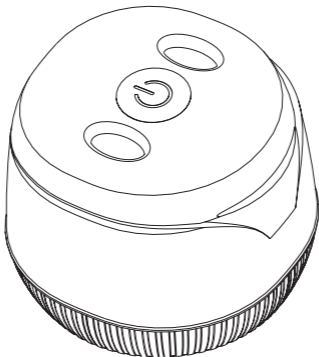
English French German Italian
Dutch Spanish Swedish

EN	3
UV-C & Ultrasonic Denture Cleaner	
FR	8
UV-C & Nettoyeur à Ultrasons pour Prothèses Dentaires	
DE	13
UV-C & Ultraschall-Prothesenreiniger	
IT	18
UV-C & Pulitore Ultrasonico per Dentiere	
NL	23
UV-C & Ultrasone Kunstgebitreiniger	
ES	28
UV-C & Limpiador Ultrasónico para Dentaduras	
SE	33
UV-C & Ultraljudsprotesrengörare	

melodii

UV-C & Ultrasonic Denture Cleaner

Clean Like Never Before!

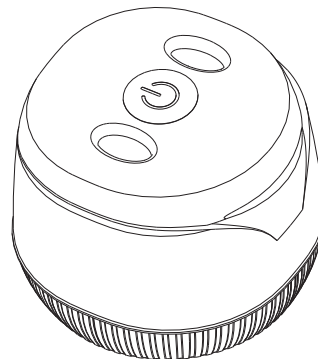Introduction

Dive into the world of ultimate cleanliness. Here's your easy guide to make the most of your cleaner.


Cool stuff your cleaner does

1. **Ultrasonic Magic:** Gets rid of dirt using a powerful 45,000Hz frequency.
2. **Saves Energy:** Only uses 24W DC power. It's greener than many battery-run ones out there!
3. **Easy-Peasy:** One button. Two cleaning modes:
 - Quick Clean: 5 minutes with ultrasonic and UV-C.
 - Deep Clean: 10 minutes - first 5 for cleaning, next 5 for special UV light magic.
4. **Germs? Bye-Bye!** The UV-C light can zap up to 99.99% of bacteria.

Get to Know Your Cleaner

- * **Sleek Look:** It has a roomy 200ml water tank.
- * **Not Just For Dentures:** Clean your jewelry, makeup brushes, toothbrushes, aligners, and more!
- * **Travel Buddy:** Handy size for travel and a fashionable gift pick.

Product details

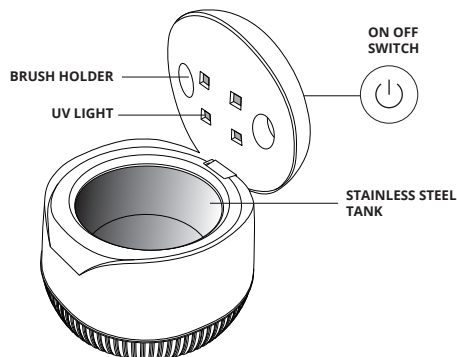Specifications

POWER SUPPLY	DC 100-240V 50/60Hz 0.5A max.
RATED VOLTAGE	12V/2a
WATT HOUR	≤24W
ULTRASONIC FREQUENCY	45khz
TANK SIZE	Ø80*39mm
ITEM SIZE	120*120*100mm
MODE 1	A 5-minute combination of ultrasonic cleaning & UV light action.
MODE 2	A 10-minute cycle, split between 5 minutes of ultrasonic and UV cleaning, followed by 5 minutes of UV pulsing.
TANK CAPACITY	200ml

How to use

1


Get ready: Pour in water and a cleaning solution that's safe to use.

2

Place Items: Pop your items into the water tank.

3

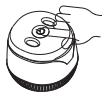Pick a Cleaning Mode:

- Push once for a **Quick Clean**.
- Push twice for a **Deep Clean**.
- Push three times if you want to turn it off.

Pro Tip: A Quick Clean is great for everyday use. For stubborn dirt, go with the Deep Clean mode 2-3 times.

Safety First

- * Start with water always. Your cleaner likes it that way!
- * Use it 3 times in a row? Let it chill for 10 minutes.
- * Fill water between the 'Min' and 'Max' signs.
- * Be careful! Don't let water touch the electronic bits.
- * No strong chemicals or bleach, please.
- * Pull the plug after you're done. Look over the cord and plug now and then.
- * **Kids Alert:** This isn't a toy. Keep it out of their reach.

Things to remember

- * **Electric Bits:** Every now and then, check the power cord and plug.
- * **Water Tips:** Keep it in a dry spot and don't dunk it underwater.
- * **Cleaning Fluid:** Stick to suggested solutions. Bleach or harsh stuff is a no-no.

Take good care of your cleaner by following these tips. Got questions or issues? This guide's got your back. Or, chat with our friendly customer support team.

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