HOME HUID SCHOONHEID INSTRUMENT

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HOME SKIN BEAUTY INSTRUMENT

INSTRUCTIONS FOR USE

Directory

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Directory

- Safety Precautions.....01
- Product structure.....02
- Product Profile.....03
- Preparation before use.....04
- Use Operation.....05
- User Care.....06
- Cleaning and Maintenance.....07
- Warrant.....08

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SAFETY PRECAUTIONS

- ⚠️ Do not apply any or care products (like vera gel, etc.) before use. Please be sure to assemble and use the product according to the regulations contained in the instruction manual. Please refer to the instructions for use carefully and keep it for inquiry if you have any questions during use, please contact our customer service in the specified language.
- ⚠️ Please do not use this product near a bath, shower, bath, or near others. Please do not use the instrument near children, pregnant women, people with a pacemaker or other electrical devices.
- ⚠️ It is not recommended to use the instrument on the face of the eyes, the neck, the arms, the hands, or the feet.

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
⚠️ BE CAREFUL


- ⚠️ Do not use this instrument for natural black skin.
- ⚠️ Do not use this product near a bath, shower, bath, or near others. Please do not use the instrument near children, pregnant women, people with a pacemaker or other electrical devices.
- ⚠️ It is forbidden to use in dark light environment.

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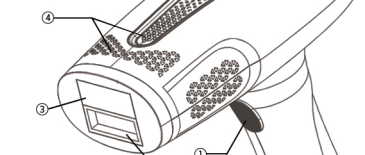
- ① flash knob
- ② Flash Venster
- ③ Pannelen
- ④ DC socket interface

STRUCTURE (ICE FREEZING PL HAIR REMOVAL) (BHRL-04)

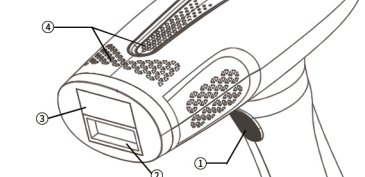
- ① flash button
- ② flash window
- ③ Ice panel
- ④ Air inlet
- ⑤ Display screen
- ⑥ DC socket interface

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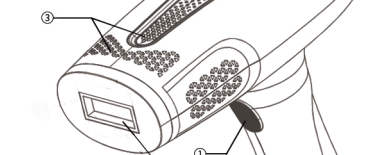STRUCTURE (NON-ICE FREEZING PL HAIR REMOVAL) (BHRL-04B)

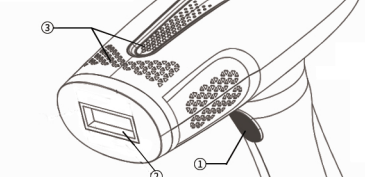- ① flash button
- ② flash window
- ③ air inlet
- ④ Display screen
- ⑤ DC socket interface

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PRODUCT INTRODUCTION

ICE freezing point removal instrument, first use a shaver to clean the hair on the surface of the body, skin, and surrounding tissues due to temperature rise, thereby inhibiting hair growth and pain, and ensuring the experience of using the instrument. It is not recommended to use the instrument on the face of the eyes, the neck, the arms, the hands, or the feet.

Model	Specification parameters	
BHRL-04	Weight: 0.32kg	Atmospheric pressure: 86kPa-106.0kPa
	Voltage: AC220V 50Hz	Environment temperature: 5-40°C
	Power: 30W	Light intensity: 5-stage adjustable
BHRL-04B	Light intensity: 5-stage adjustable	Temperature: 1-100°C
	Charging time: approximately 2 hours	Charging mode: 100%
	Relative humidity: $5\% \sim 95\%$	Implementation standard: GB4706.1

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Procedures for using ice compress care mode

- 01 Please wear protective eye caps on the hair that is being shaved after the screen is closed.
- 02 Clean the skin and dry it for 3 seconds, turn on the power by pressing the On/Off key, and the screen will be closed after 3 seconds.
- 03 Clean the skin and dry it for 3 seconds, turn on the power by pressing the On/Off key, and the screen will be closed after 3 seconds.
- 04 Clean the skin and dry it for 3 seconds, turn on the power by pressing the On/Off key, and the screen will be closed after 3 seconds.
- 05 It is not recommended to use the instrument on the face of the eyes, the neck, the arms, the hands, or the feet.

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NOTICE

- 01 Please do not use any or care products (like vera gel, etc.) before use. Please be sure to assemble and use the product according to the regulations contained in the instruction manual. Please refer to the instructions for use carefully and keep it for inquiry if you have any questions during use, please contact our customer service in the specified language.
- 02 Please wear protective eye caps on the hair that is being shaved after the screen is closed.
- 03 Clean the skin and dry it for 3 seconds, turn on the power by pressing the On/Off key, and the screen will be closed after 3 seconds.
- 04 Clean the skin and dry it for 3 seconds, turn on the power by pressing the On/Off key, and the screen will be closed after 3 seconds.
- 05 It is not recommended to use the instrument on the face of the eyes, the neck, the arms, the hands, or the feet.

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Uncomfortable groups and conditions

- Do not use the instrument on the face of the eyes, the neck, the arms, the hands, or the feet.
- Do not use the instrument on the face of the eyes, the neck, the arms, the hands, or the feet.

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WARRANTY

The warranty period of this product varies according to the country in which it is sold, and will vary from country to country. It is not recommended to use the instrument on the face of the eyes, the neck, the arms, the hands, or the feet. It is not recommended to use the instrument on the face of the eyes, the neck, the arms, the hands, or the feet.