# auronic

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## Introduction

- Thank you for choosing our Auronic product! Please always follow the instructions for safe usage. If you have any questions about the product or experience any problems, please contact us by email: [service@Auronic.nl](mailto:service@Auronic.nl)
- The Auronic Massage Pillow is suited for different body parts: your legs, upper back, neck, lower back, shoulders, feet, calves, and arms. The Massage Pillow will turn off automatically after 15 minutes of usage to ensure user safety. Choose between three speed modes and a heat mode. Easy to carry with you and including a carry-on bag.

## Warnings and Safety Instructions

- ⚠ Please ensure that the packaging is kept out of the reach of children! Risk of suffocation!
- ⚠ Make sure to keep this user manual. Read the manual carefully before using your Auronic Massage Pillow

### General

- If you have not used the Auronic Massage Pillow before, it is important that you read and follow the instructions as stated in the user manual. Be aware of the risks and consequences resulting from misuse of the Massage Pillow. Misuse of the Massage Pillow can result in personal injury, or destruction of the instrument/appliance/device and of the installations.
- Use the Massage Pillow only for purposes described in this user manual.
- If the Massage Pillow malfunctions, is dropped or (appears to be) damaged in any manner, do NOT use the product and contact Customer Service.
- The Massage Pillow is designed and intended for private, non-commercial use only.
- Do NOT leave the Massage Pillow unattended when in use. When you are not using the Massage Pillow, turn the device off and disconnect the device from the power socket
- Do NOT disconnect the Massage Pillow from its power supply by pulling the cable.
- Do NOT carry, pull or twist the device by its cable.
- Make sure that you do NOT trip over the cable. The cable should not be kinked, trapped or twisted.
- Do NOT use the Massage Pillow when you are driving.

### Persons and Use

- This device should NOT be used by persons (including children) with physical, sensory or mental disabilities, or persons (including children) with lack of experience with and knowledge of the device, unless they are supervised or have been given instruction concerning the use by a person responsible for their safety.
- Do NOT use the Massage Pillow: If you Are pregnant; If you have a cold with high temperatures, varicose veins, thromboses, jaundice, diabetes, nerve disorders (e.g. sciatica) or acute inflammation; If you have unexplained abdominal pain or pain in the calves; If you are sensitive to heat; On bones (e.g. joints, spine) the head, the throat area, the genital area, around the kidneys, or other sensitive areas; On burned or sensitive skin, swollen or inflamed areas, areas with poor circulation, broken or injured skin, bruises, or areas with skin rashes; On unconscious or sleeping people; On the same body part for a long time, this can cause complaints and burns; In bed.
- We recommend to only use the Massage Pillow if you are 16 or over. Keep the Massage Pillow away from children when using and storing.
- Do NOT use the Massage pillow for more than 15 minutes at a time and not more than 30 minutes a day.
- The Massage Pillow will turn off automatically after 15 minutes of continuous use. Let the device cool down before you use it again.
- If you feel any pain or discomfort from use of the Massage Pillow, or find the massage unpleasant, stop the treatment and consult your doctor.

### Medical remarks

- Persons with cardiac pacemakers, artificial joints or electronic implants should consult their doctor before using the Massage Pillow and should never massage the region of and around the heart.
- This device is NOT a substitute for treatments administered by doctors or physiotherapists.
- Before using the Massage Pillow, consult your doctor regarding any unexplained pain, if you are undergoing medical treatment or using medical equipment.

### Electronics

- Do NOT reconstruct, repair or disassemble the Massage Pillow. This can damage the device. For safety reasons, these parts may only be replaced by authorized service personnel. In order to prevent hazards, always send the Massage Pillow to the service centre for repairs.
- Make sure your outlet supplies the correct power voltage to avoid damage to the Massage Pillow.
- Always check the power cord and plug before using the Massage Pillow. Do NOT operate the Massage Pillow with a damaged cord or plug. A damaged cord or plug must be replaced by a qualified person in order to avoid damages and personal injuries.
- Only use the adapter provided by Auronic.
- Do NOT use the Massage Pillow in or near a bath, shower, swimming pool, or over a basin filled with water.
- Make sure to prevent any contact with water. Do NOT submerge the Massage Pillow in water. If water does get inside the casing, turn the Massage Pillow off immediately and contact our customer service.
- Do NOT operate the device with damp or wet hands and do NOT touch the power cord with damp or wet hands.
- If the power fails during the massage, turn off the Massage Pillow immediately and pull the plug from the socket.
- The Massage Pillow heats up during use. Do NOT use the device if you are sensitive to heat.
- Keep the Massage Pillow away from moisture, heat and in or close to areas where the surrounding temperature is relatively hot (e.g. stoves or other sources of heat).
- Do NOT cover the Massage Pillow when it is switched on. Never put blankets or cushions on top of or against the Massage Pillow. This poses a risk of fire, electric shock or injury.
- Do NOT insert any objects into the openings of the device. Make sure that the moving parts can always move freely and no objects are jammed.
- Make sure to keep your fingers and hair away from the moving parts, to avoid injuries.
- Only use the adapter, massage heads and accessories provided by Auronic, to avoid overheating, explosions, electric shocks, fire, other damages and personal injuries.
- We recommend NOT using this device directly before retiring - the massage that it delivers has a stimulating effect and may delay sleep

### Receiving Your Shipment

- When receiving your shipment, make sure that the contents are consistent with the packing list. Notify your distributor of any missing items. If the equipment appears to be damaged, file a claim immediately with the carrier and notify your distributor at once, giving a detailed description of any damage. Save the damaged packing container to substantiate your claim.

### Contents of Package

- Auronic Massage Pillow - Model AU3032 - EAN 8720195254812

#### Included Accessories


- 1x User Manual
- 1x Power Adapter
- Carry-on bag

### Product Features


Suitable for


Neck & Shoulders

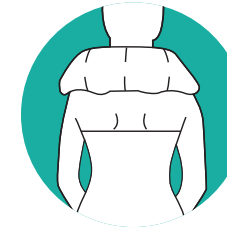
Neck

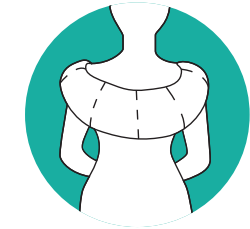Back

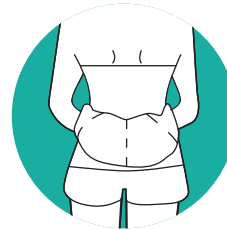Waist

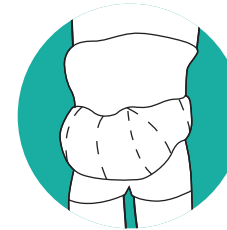Abdomen

Thighs

Legs

Feet

Calves

### Technical Data

Size	41 x 17 x 50 cm	Speed modes	3
Weight	1,2 kg	Automatic shut-down	After 15 min.
Colour	Black/Grey	Overheat protection	Yes
Material	PU leather	Massage Type	Shiatsu massage
Input - Frequency - Output	12V - 50 Hz - 5A	Number of Massage Nodes	8
Power	24W	Heat function	Yes
Suitable for (body part)	Legs, upper back, neck, lower back, shoulders, feet, calves, arms	Size Carry-on bag	40,00 x 18,00 x 12,00 cm
		Material Carry-on bag	Leather

## Operating Instructions

	<b>On/Off</b> Turn the massager on or off. The massager turns off automatically after 15 minutes of continuous use
	<b>Speed</b> Increase/decrease the rotation speed of the massagenodes
	<b>Direction</b> Choose whether you want the massage nodes to turn clockwise or counterclockwise
	<b>Heat mode</b> Turn the heat mode on or off

## Maintenance and Cleaning

- Unplug the Massage Pillow prior to cleaning. Use a soft damp cloth and non-abrasive cleaners. Do NOT submerge the Massager in water or any other liquids and keep it away from solvents and harsh detergents.

## Storage

- Store the Massage Pillow in a clean and dry place. Do not store the device in extreme high or extreme low temperatures.

## Disposal - Recycle

- Contribute to a cleaner living environment! This product is subject to the European Directive (EU) 2012/19/EU. Do not dispose of the Massage Pillow with the household waste. For proper recycling, please contact your local authority or your household waste disposal service for further details on your nearest designated collection point.

## Service and Warranty

- If you need service or information regarding your product, please contact the Auronic Customer Support at [service@auronic.nl](mailto:service@auronic.nl)
- Auronic provides a 2 year warranty on its products. To obtain service during the warranty period, the product needs to be returned with a proof of purchase. Product defects have to be reported within 2 years from the date of purchase. This warranty does not cover: damages caused by misuse or incorrect repair; parts subject to wear; defects that the customer was aware of at the time of purchase; damages or defects caused by customer neglect; damages or defects caused by third parties.

## Disclaimer

- Changes reserved; specifications are subject to change without stating grounds for doing so.

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© **LifeGoods B.V.**