Yunmai Fascia Massage Gun Extra Mini

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The product user guide contains product features, how to use, and the operating procedure. Read the user manual carefully to get the best experience and avoid unnecessary damage. Keep this manual for future reference. If you have any questions or comments about the device, please contact the customer line.

✉ www.alza.co.uk/kontakt

☎ +44 (0)203 514 4411

Importer Alza.cz a.s., Jankovcova 1522/53, Holešovice, 170 00 Praha 7, www.alza.cz

Safety Instructions

Please read this manual carefully before using the product and keep it properly.

Warning

Attention:

- It is forbidden to use this product in the head, face, elbow, knee and other joints, spine, and other bone parts with less muscle coverage.
- It is forbidden to use this product in the lateral, anterior and posterior cervical vertebrae of the neck, so as to avoid damaging the carotid sinus, resulting in fainting and cardiac arrest.
- It is forbidden to use this product in the armpit, the medial side of upper arm, the medial side of clavicle, the medial side of elbow joint, the popliteal fossa at the back of knee, abdomen, kidney area, genitalia and other vulnerable parts with dense nerves and blood vessels.
- Patients with kidney stones, ureteral stones, bladder stones, urinary tract lithiasis, urinary tract infection and other urinary diseases are prohibited to use this product below the abdomen.
- It is forbidden to use this product in the body area with bruise, pain, surface wound, local fever, swelling discomfort, acute inflammation, muscle spasm, muscle laceration, acute sprain, contusion, and strain.
- In case of fever, nausea, vomiting, muscle pain and weakness after strenuous exercise, it is forbidden to use this product and you should seek medical advice immediately.
- Patients with the following diseases are prohibited to use this product: cardiovascular and cerebrovascular diseases, blood diseases, coagulation dysfunction, thrombocytopenic purpura, neck plaque, varicose veins, venous thrombosis, cancer, tumour, mental diseases.
- Patients with fracture or internal nail are prohibited to use this product.
- It is forbidden to use this product in precise operation or driving.
- Please do not use this product continuously for a long time in order to prevent muscle damage caused by excessive massage.

Yunmai Fascia Massage Gun can help the body relieve pressure and fatigue, relieve muscle pain, and discomfort, and improve your sports performance. However, it cannot replace the function of professional medical equipment. It has no effect on weight loss, cannot cure any disease, and cannot solve medical problems. There is no compensation for medical accidents caused by your own disease or improper operations.

Warnings

Attention:

- For indoor use ONLY!
- During charging the unit must be placed in well ventilated area.
- The battery is to be disposed of safely.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory, or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Children should be supervised to ensure that they do not play with the appliance.
- This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased.
- They can take this product for environmentally safe recycling.

Attention

Dos and Don'ts

- The product has been charged in advance with about 60% to 80% of the power before leaving the factory and can be used without being charged after opening.
- When the machine is running, do not put fingers, hair or other parts of the body close to the interface between the massage head and the machine so as to avoid injury.
- Please do not block or insert foreign matters into the installation groove of massage head and the charging interface, otherwise it may cause scald, electric shock, or fire due to short circuit.
- Please do not soak the product in water or other liquids.
- Please do not use this product in humid environment.
- Please use and charge it at 0~ 40°C.

- Keep away from water, fire and corrosive environment, so as not to affect battery performance.
- Please do not repair, disassemble, or modify the product by yourself.
- Please do not use non-original charger for charging, and only use the very original charger provided by our company.
- Please do not use this product when it is being charged.
- After finishing charging, please take away the charger from the product in time.
- If you have not used this product for a long time, please charge it in advance before you reuse it someday.

Unsuitable people

This product is not allowed to be used by people over 80 years old or underage in order to avoid accidents. If any of the following conditions are met, please contact the doctor before use, and only use with permission.

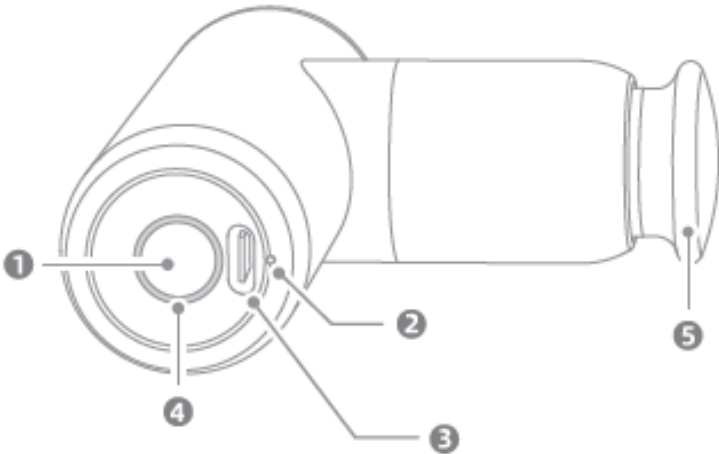
- Pregnant and puerperal women.
- People over 60 years old.
- People with rhabdomyolysis syndrome.
- People with limb loss or/and high paraplegia.
- People with pacemaker and stent.
- People who took some surgery within 12 months.
- People with orthopaedic diseases such as abnormal spine, joints and bones, osteoporosis, cervical spondylosis, lumbar disc herniation and others.
- People with abnormal limb perception and limited mobility caused by disease or/and age.
- People with complications caused by other diseases in the body, such as neuropathy or retinal injury.
- Patients with kidney stones, ureteral stones, bladder stones, urinary tract lithiasis, urinary tract infection and other urinary diseases.

In case of abnormality and failure:


- Avoid beating, shaking, or hitting the product. Strong impact may cause damage to parts of the machine and lead to abnormal work. In case of the above situation, please consult the after-sales service centre in time.
- Long term continuous use of this product may cause the surface temperature of the machine to rise and damage the machine, please do not use it continuously for a long time.
- In case of fever, smoking, excessive buzzer of the machine, please stop using it immediately and contact the after-sales service centre for solutions.
- If high temperature, deformation, expansion, and battery leakage occur in the process of using the machine, please stop using it immediately and contact the after-sales service centre.

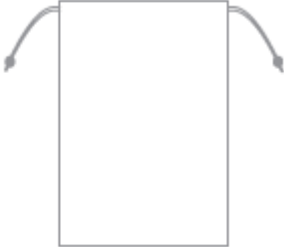Product Instruction

Product Diagram

Product Accessories List

			
Storage Bag	Manual	Charging Cable	
			
Round head	Oblate Head	Mushroom-shaped head	U-shaped head

Product Function Introduction

Host

1. power key (gear switch and power on/off)
2. power indicator light
3. charging interface
4. gear indicator light (breathing light)
5. massage had (silicone head can be installed)

Functions of Fascia Massage Gun

- To relieve the delayed muscle soreness after exercise and relax the tight muscles.
- To relieve the back and shoulder muscle pain and stiffness caused by sedentary and bad posture.
- Daily massage to relax the body.
- Silicone massage heads
- Round Bobo head: Soft and elastic, suitable for most of the body muscle relaxation.
- Mushroom-shaped head: deeper and stronger, suitable for impacting deep tissues and local sore muscles, such as palm, sole and shank.
- U-shaped head: suitable for muscles on both sides of spine and trapezius of shoulder to create beautiful back and right-angle shoulder.
- Oblate head: beating muscles as the palm, suitable for most of the body muscle relaxation.

Tips:

The illustrations of products, accessories and user interface in the manual are schematic diagrams for reference only. Due to the updating and development of the product, there may be some differences between the actual product and the schematic diagram. Please refer to the actual product.

Intelligent protection:

This machine has intelligent protection function. In order to prevent muscle injury caused by prolonged massage, the gear indicator light goes out after 10 minutes of continuous work, and the power indicator light remains in the state. Click the power key to enter the first gear within 5 seconds to time it again. If there is not any operation, the machine will shut down automatically.

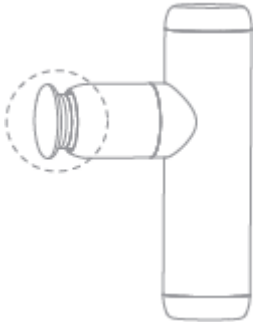
Locked-Rotor protection:

This machine has locked-rotor protection function. If locked-rotor state lasts for more than 3 seconds, it will enter standby mode. Click power key to reactivate in 10 seconds. If there is not any operation, it will shut down automatically.

Installation Process

1. Installation

The Fascia Massage Gun has its own massage head (flat head), which can be used independently or combined according to the demand. The installation effect is as follows

Installation display

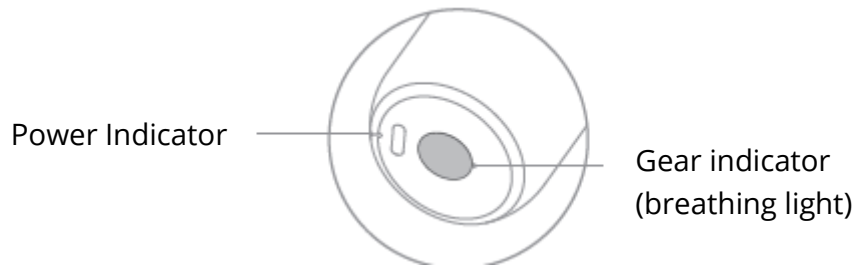
Warning

Please install and remove the massage head in the state of power off and ensure that the head is installed firmly before starting the machine. when using it, the massage head cannot be adjusted or loaded so as to avoid pinch.

Usage

Indicator light

Gear indicator (gear shift)


After power on, click the power key to change gears in turn.

Gear Level	Gear indicator (breathing light)	Rated Speed (RPM)	Function
First gear	Blue light breathing	1500	Arousing muscles
Second gear	Green light breathing	2400	Massage to relax
Third gear	Red light breathing	3200	Deep massage
Fourth gear	Yellow light breathing	350	Gentle massage
Standby gear (shutdown after 2s)		0	

Power indicator (power indicator in use)

- The green light is on in a stable state, which means full power.
- The green light flashes slowly, which means medium power.
- The green light flashes quickly, which means low power. Please turn it off and charge it in time.

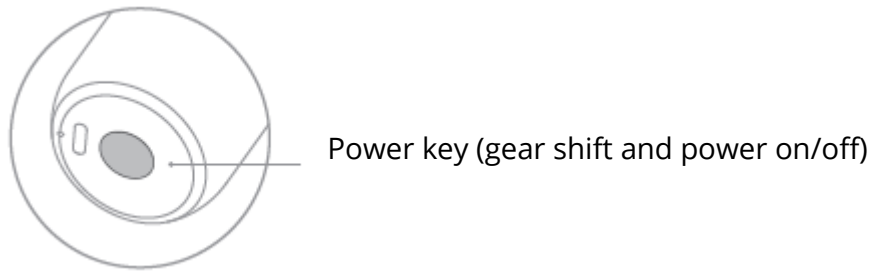Power key


- Long press the power key to turn it on or off.
- After power on, single click the power key to change gears. It could cycle in the order of one, two, three and four gears in turn.

Attention

It is recommended to run it for 5 seconds after starting the machine before using it. It is forbidden to use it when it is being charged!

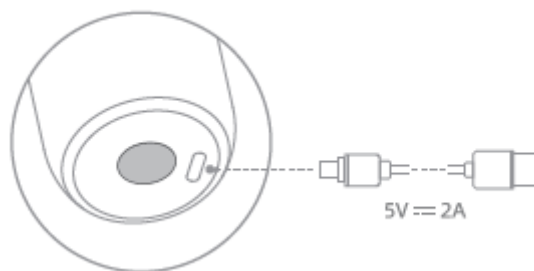Charging

The charging port is located at the bottom of the handle of the Fascia Massage Gun. Plug the charging cable into the charging port and connect the power supply to start charging.

Power indicator (charging indicator)

- It is being charged if the red light is on.
- The massage gun is full charged when the green light is on.

Attention

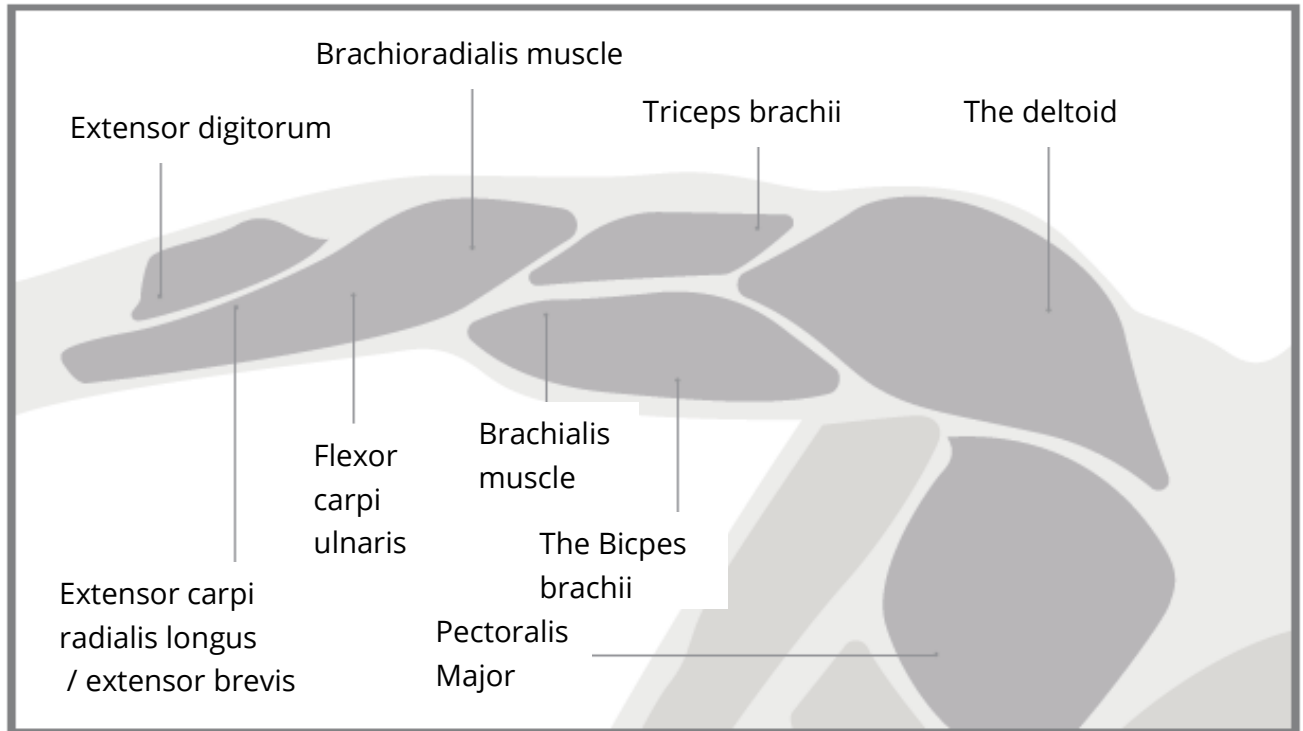
Do not use it when it is being charged!

How to use

- Choose the appropriate massage head for installation, open the Fascia
- Massage Gun, then you can change to the appropriate gear according to your needs. And it is recommended to try it from the low gear level.
- Gently press the massage head on the muscle part to be relaxed, and slowly slide back and forth. If there is an extra sore part in the local area, apply appropriate pressure at the tenderness point, and stay for 15s 30s.
- The Fascia Massage Gun itself has impact force, so it is not appropriate to use excessive force when applying pressure so as to avoid damaging muscles and bones.
- It is recommended to use it in the part with more muscle. Do not stay in the same part for more than 60s.
- The recommended duration is 10 minutes. If you still feel sore, extend the massage time appropriately.
- It is recommended that the single use time should not exceed 15 minutes and the continuous use time should not exceed 30 minutes.
- Please operate the product in strict accordance with the massage position and operation methods written in this manual.

Massage Position


Tip: do not use it in neck, armpit, inside of upper arm and inside of clavicle.

Arm


Position	Duration (s)	Frequency
Brachioradialis muscle	10 – 15 s	3 times
Flexor carpi radialis	10 – 15 s	3 times
Palmaris longus	12 – 15 s	3 times
Flexor carpi ulnaris	10 – 15 s	3 times
Extensor digitorum	05 – 10 s	Twice
Extensor carpi longus radialis	05 – 10 s	Twice
Extensor carpi radialis brevis	05 – 10 s	Twice
The biceps brachii	12 – 20 s	3 times
Brachialis muscle	12 – 20 s	3 times
The deltoid	12 – 20 s	3 times
Pectoralis major	12 – 20 s	3 times

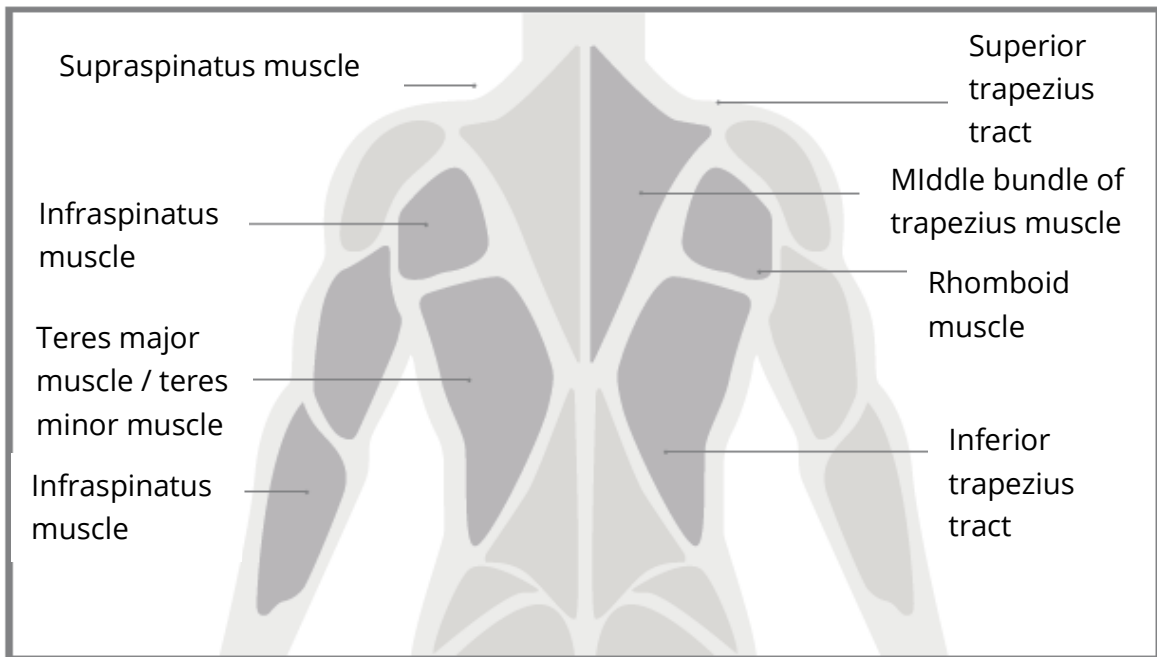Tip: it is forbidden to use in the neck, spine, and other skeletal parts with less muscle coverage.

Position	Duration (s)	Frequency
Triceps brachii	15 – 20 s	3 times
Rhomboid muscle	20 – 30 s	3 times
Superior trapezius tract	15 – 20 s	3 times
Middle bundle of trapezius muscle	15 – 20 s	3 times
Inferior trapezius tract	15 – 20 s	3 times
Supraspinatus muscle	10 – 15 s	3 times
Elevator scapulae	10 – 12 s	Twice
Teres major muscle	10 – 12 s	Twice
Teres minor muscle	10 – 12 s	Twice
Infraspinatus muscle	10 – 12 s	Twice

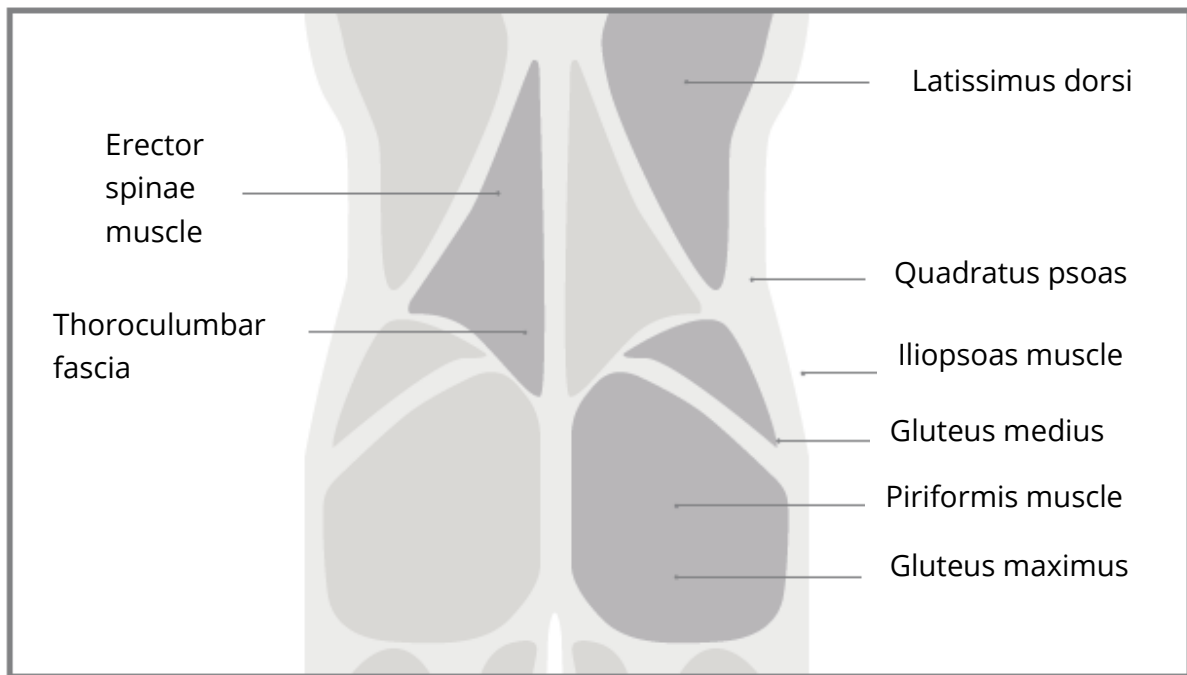
Tip: massage the muscle from the middle to both sides, from top to bottom; Massage along the muscle, and pay attention to avoid protruding bones.

Waits

Position	Duration (s)	Frequency
Erector spinae (lumbar)	30 – 60 s	4 times
Latissimus dorsi	30 – 60 s	4 times
Thoracolumbar fascia	40 – 60 s	3 times
Psoas major	30 – 60 s	3 times
Internal and external oblique muscle of abdomen	10 - 15 s	Twice
Erector spinalis (thoracic)	15 – 20 s	Twice
Quadratus psoas	20 – 30 s	Twice
Iliopsoas muscle	20 – 30 s	3 times
Gluteus maximus	30 – 60 s	4 times
Gluteus medius	30 – 60 s	4 times
Piriformis muscle	60 – 90 s	Once

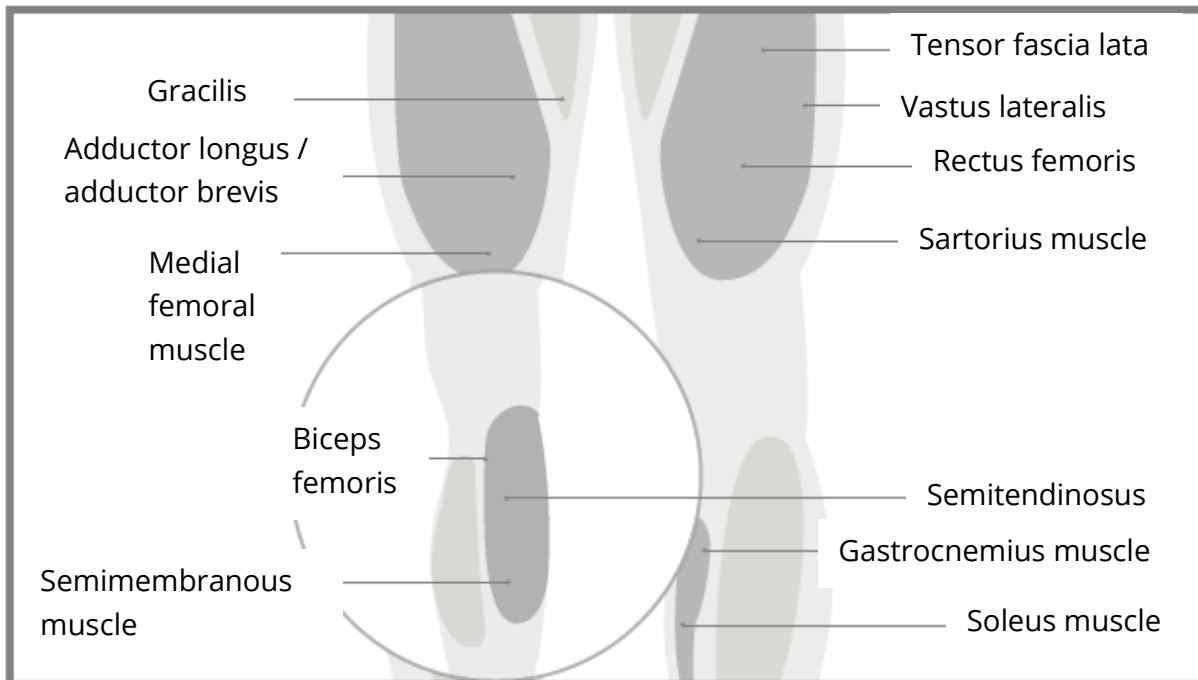
Tip: do not use it in the popliteal fossa at the back of the knee and pay attention to avoid protruding bone joints.

Legs

Position	Duration (s)	Frequency
Semi membranous muscle	20 – 30 s	3 times
Gracilis	10 – 15 s	Twice
Adductor brevis	10 – 15 s	Twice
Adductor longus	10 – 15 s	Twice
Rectus femoris	30 – 45 s	Twice
Vastus latarelis	30 – 45 s	Twice
Medial femoral muscle	30 – 45 s	Twice
Sartorius muscle	30 – 45 s	Once
Tensor fascia lata	15 – 20 s	Twice

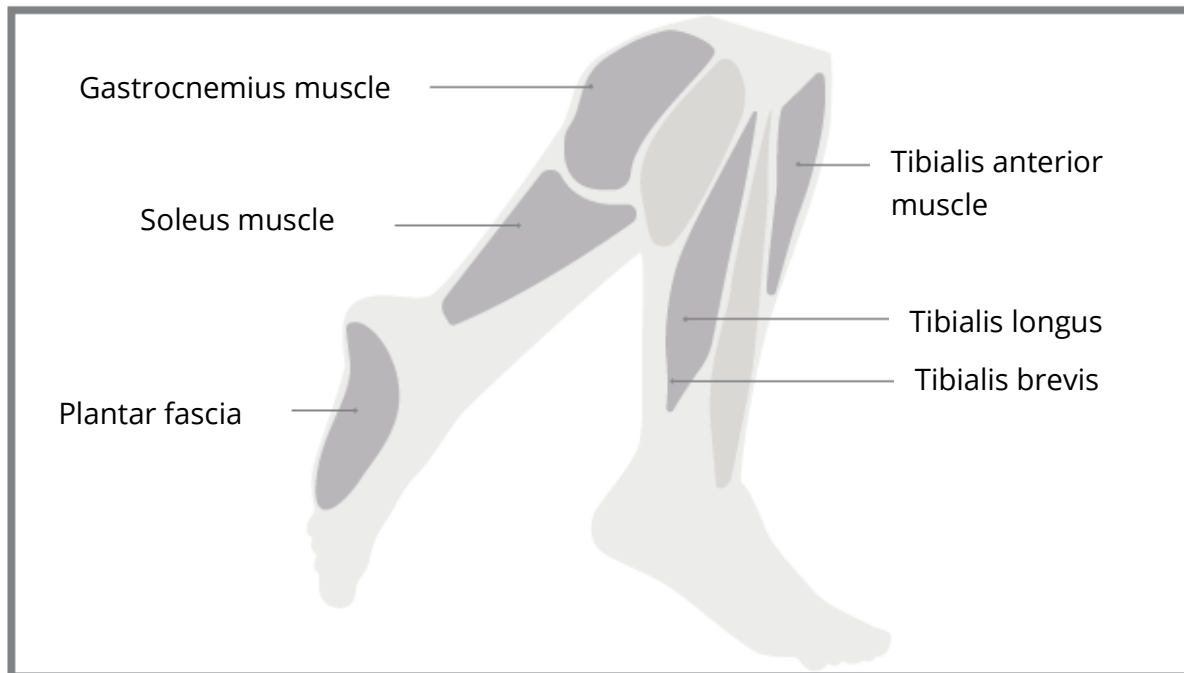Tip: pay attention to avoid protruding bone joints.

Position	Duration (s)	Frequency
Gastrocnemius muscle	20 – 30 s	3 times
Soleus muscle	20 – 30 s	3 times
Tibialis anterior muscle	15 – 20 s	Twice
Peroneus longus	15 – 20 s	Twice
Peroneal brevis	15 – 20 s	Twice
Plantar fascia	20 – 30s	3 times

Basic Parameters

Product Name: Yunmai Fascia Massage Gun extra mini

Product model: MVFG-M281

Net weight: 240 g

Battery capacity: 1300 mAh

Rated input: 5V-2A

Strength level: four gears

Power supply mode: rechargeable lithium battery

Product size: 88 x 35 x 126.6 mm

Ages for use: for adults over 18 years old; and pregnant women should not use.

Maintenance

1. Cleaning

Please wipe the machine body and massage heads with clean and slightly wet soft cloth and wet towel regularly; then dry them with dry cloth and paper towel, and finally put them in the recycling bag or packing box.

Please do not use water to wash or expose the massage heads to the sun.

When cleaning, make sure that the power switch is turned off and the charging plug is removed.

Please do not use corrosive solution to wipe the machine body, otherwise it may cause product failure or product cracking and discoloration.

2. Storage

Please store the product in a cool and dry place.

Please do not store the machine body in high temperature environment above 40°C and low temperature environment below -5°C. Keep it away from humid, water, fire, and corrosive environment, so as not to affect the battery performance.

When it needs to be stored or carried out, the Fascia Massage Gun can be stored in the attached portable storage bag.

Please do not touch, squeeze, drop, throw, or/and scratch the product.

The product is suitable for transportation by car, train, boat, plane, etc. During transportation, the product should be fixed steadily and put up gently. We should prevent the product from extrusion, violent vibration, collision, rolling and falling.

If you have not used this product for a long time, please charge it in advance before reusing it.

If high temperature, deformation, expansion, and battery leakage occur in the process of using the machine, please stop using it immediately and contact the after-sales service centre.

Warranty Conditions

A new product purchased in the Alza.cz sales network is guaranteed for 2 years. If you need repair or other services during the warranty period, contact the product seller directly, you must provide the original proof of purchase with the date of purchase.

The following are considered to be a conflict with the warranty conditions, for which the claimed claim may not be recognized:

- Using the product for any purpose other than that for which the product is intended or failing to follow the instructions for maintenance, operation and service of the product.
- Damage to the product by a natural disaster, the intervention of an unauthorized person or mechanically through the fault of the buyer (e.g. during transport, cleaning by inappropriate means, etc.).
- Natural wear and aging of consumables or components during use (such as batteries, etc.).
- Exposure to adverse external influences, such as sunlight and other radiation or electromagnetic fields, fluid intrusion, object intrusion, mains overvoltage, electrostatic discharge voltage (including lightning), faulty supply or input voltage and inappropriate polarity of this voltage, chemical processes such as used power supplies, etc.
- If anyone has made modifications, modifications, alterations to the design or adaptation to change or extend the functions of the product compared to the purchased design or use of non-original components.

EU Declaration of Conformity

Identification data of the manufacturer's / importer's authorized representative:

Importer: Alza.cz a.s.

Registered office: Jankovcova 1522/53, Holešovice, 170 00 Prague 7

IČO: 27082440

Subject of the declaration:

Title: Massage Gun

Model / Type: Yunmai Extra Mini

The above product has been tested in accordance with the standard(s) used to demonstrate compliance with the essential requirements laid down in the Directive(s):

Directive No. (EU) 2014/53/EU

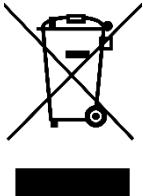
Directive No. (EU) 2014/35/EU

Directive No. (EU) 2011/65/EU as amended 2015/863/EU

WEEE

This product must not be disposed of as normal household waste in accordance with the EU Directive on Waste Electrical and Electronic Equipment (WEEE - 2012/19 / EU). Instead, it shall be returned to the place of purchase or handed over to a public collection point for the recyclable waste. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in fines in accordance with national regulations.

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