

Xiaomi scales with the function of measuring body fat


Precautions

1. Do not place weight on wet and slippery surfaces, and should not stand on the scale with wet feet or stand on a wet surface weights to avoid slipping.
2. When weighing both feet must be located in the center of scale to avoid tipping of device.
3. Operational environment must be dry, it is forbidden to put weight in water.
4. Use a damp cloth and neutral detergent to clean the surface of the weights, contact with organic solvents and other chemicals is forbidden.
5. When battery charge is over, immediately remove a battery to avoid leakage of electrolyte and damage of scales.
6. Persons with pacemakers and other medical equipment is not recommended to use this device. In body composition analysis of professional athletes, pregnant women or others, who are a few special provisions, result can be error.


Main specifications

Dimensions: 300mm x 300mm x 20 mm;
Net Weight: 1.6 kg;
Limit of measurement: 5 kg - 150 kg;
Units of measurement: Jin / kg / pound;
Step of measurement: 50 g;
Case Material: ABS plastic;
Wireless: Bluetooth 4.0;
Operating temperature 0 ° C ~ 40 ° C;
Battery: 4 AAA batteries;
Hardware Requirements: OS Android 4.4 / iOS 8.0;
and above, supports Bluetooth 4.0.

Installing of battery

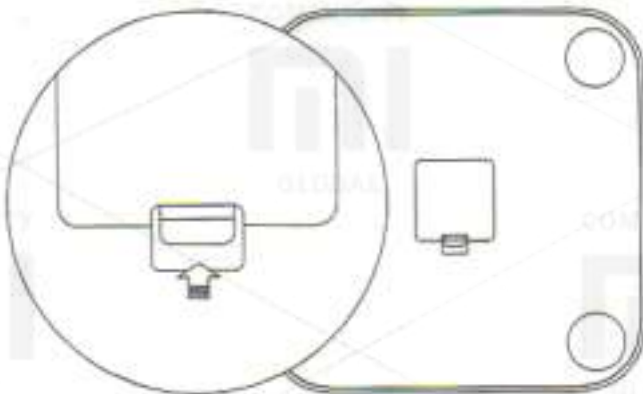

Open a battery compartment cover, under direction from the picture, insert four batteries are included.

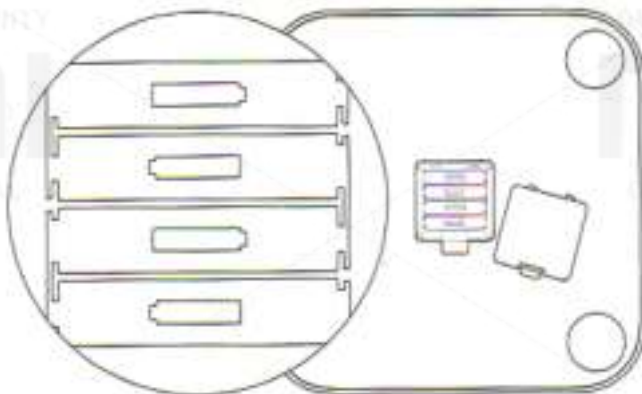

Installing the App

Enter a search «Mi Fit» at Xiaomi app store or App Store, or scan the QR-code below to download and install the latest version.



Adding the scales

Open MiFit Xiaomi app and add scales, follow the prompts to enter personal data. After adding scales view a history of weighing and information on the state of the body in a mobile application.

Usage

Place the scales on a hard, flat surface. Barefoot stand on the surface of the scales as indicated in the picture, provide good contact of skin with two metal electrodes.

Once weight indicator stabilized, starts measuring a composition of body. When a progress indicator below the display weight flashes, which means that a measurement is completed, then you can enter at MiFit application to view the information of body composition.



If the icon «OL» flashes while weighing, it means, that the measurements exceeded the maximum limit of 150 kg weight.



List of faults

Failure in the functions, mentioned in the manual. You can not turn on the device, it is impossible to measure weight, fat mass.

No image on the display / present errors / broken pixels.

Can not sync with these devices.

Crack the shell, incurred due to the structure or materials used.

Body composition, that the device shows, used for informational purposes only and is not a medical opinion.

Setting units of weight

In factory settings of Xiaomi scales installed «jing», please open MiFit, menu settings wt replace equipment unit weights. Scales supports three measurement units, «Jin», «kg» and «pound».

Literal tips

If the icon «L» flashes while weighing, it means low battery, please replace the battery.



Warranty

After-service of Xiaomi floor scales with function of measuring body fat is obliged to: perform warranty after-sales service under the Law of the PRC «Consumer Protection» and the Law of the PRC «Quality of the product» warranty service involves: During the warranty period you have a right to repair, exchange or return of goods, according to the regulations. Repair, exchange and return of goods carried if the sales receipt.

1. Within 7 days of purchase, in case of failure indicated in the «List of faults», which was approved by Xiaomi service centers, you are entitled to a free refund, exchange or repair the goods.

2. Within 8-15 days after the purchase, in case of a malfunction, indicated «List of faults», which was approved by Xiaomi service centers, you are entitled to free repair or exchange of goods.

3. Within 12 months from the date of purchase, in case of failure, indicated in the «List of faults», which was approved by Xiaomi service centers, you are entitled to free repair.

The following conditions are not covered by the warranty:

1. The quality faults, resulting from unauthorized repair, misuse, falling into a liquid, accident, modification, misuse of components, not designed for this device or install the system software; break, fix labels or serial number of the device.
2. Expiration of warranty.
3. Damage, caused by force majeure.
4. Do not meet the faults, listed in the «List of faults»

The title and content of harmful substances in the product

Name of part	Harmful substances					
	Pb	Hg	Cd	Cr(VI)	PPB	PPDF
Housing	O	O	O	O	O	O
Circuit board	X	O	O	O	O	O
Support of snakes	O	O	O	O	O	O

This table has been compiled in accordance with SJ / T 11364.

O: means, that the content of harmful substances in homogeneous material of this component is within the requirements of the provisions of GB / T 26572.

X: means, that the content of harmful substances in at least one of the homogeneous materials of this component, exceeds the threshold requirements set GB / T 26572.

This product complies with RoHS environmental regulations; Currently there are no technology, that can replace or reduce the led content in electrotechnical ceramics and copper alloys.

The term environmentally sound use of the term, refers to a safe period within which under normal operating conditions contained in this product are harmful substances or elements pose no danger.

Before using the product, please read this manual and keep it properly.

Xiaomi scales with the function of measuring body fat

Model: XMTZC02HM;

CMIIT ID: 2016DP6264;

Applicable standards: Q / HM02-201;

Producer: Ltd. «Anhui Huami Informeysn Technologies»

Address: Hefei Economic Zone and scientific development, ul.

Vantsyansi 800 Dunman state database, build. A4, Suite 1201;

Visit web page after-service Xiaomi: www.mi.com/service
or phone number service support Xiaomi
and learn more about after-sales service.