



Vibrating Massage Ball



User manual



Thank you
for purchase!

Designed in The Netherlands

Made in China

VOTIQ B.V.

Charging method:

1. Connect the power adapter charging plug to the product charging jack.
2. Plug the power adapter into a 100-240V-, 50/60Hz power socket.
3. When charging starts, the red LED flashes during charging.


*If the massage ball is continuously charged for 3 hours, the three battery indicator lights are still not fully illuminated (full state), please unplug the power adapter and contact the maintenance personnel to avoid danger.

*The vibration function cannot be used while charging.

4. After the charging is completed, the LED light continues to illuminate to indicate that the battery is fully charged.

*It takes about 1.5 hours to send it from empty to full charge.

5. When charging is complete, disconnect the AC adapter and charging plug.

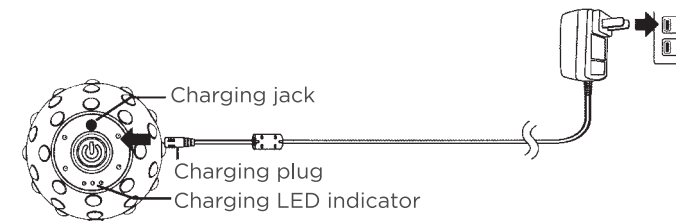

1

Specifications:

Power Specifications

Input AC 100-240V-, 50/60 Hz
DC9V two 1000 mA

Executive Standard: GB4706.1-2005/ GB4706.10- 2008

The company reserves the right to modify the product and the manual.
Notice that the appearance of the product is subject to the actual product.



6

Specifications:

| | |
|-----------------------|---|
| Product name | Vibrating massage ball |
| Product number | J F.005 |
| Battery specification | DC 7.4V 1200mAh |
| Charging time | About 1.5 hours (full charge) |
| Quiescent Current | Less than 80uA |
| Full discharge time | 50 minutes (5 cycles) for 10 minutes, 10 minutes for 10 minutes |
| Working time | 50 minutes |
| Vibration revolution | 1st gear motor speed: 2000RPM soil 200RPM |
| Automatic shut-down | 2 speed motor speed: 2500RPM 200RP 3-speed motor speed: 2100-3600RPM 4-speed motor speed: 3600RPM soil 400RPM |
| Ambient temperature | About 10 minutes |
| Bearing weight | 10°C-4D°C |
| Product Size | 80kg |
| Product weight | Φ95 xL 90mm0.3 kg (body weight) |
| Material | Silica gel |
| Product accessories | Power adapter, manual |
| Manufacturing origin | China |

*The actual time and number of uses, it depends on the use environment.

5

Charging method:

The vibrating massage ball offers 1 automatic mode and 3 massage strengths. You can make a selection by pressing the function button.

Boot

Press **⏻** button once, the massager is turned on; the automatic vibration massage starts, the vibration speed is one gear, and an intensity indicator light is displayed;

Press the **⏻** button for the second time, the vibration speed is increased to the second gear, and two intensity indicators are displayed.


Pressing the **⏻** button for the third time, the vibration enters the automatic mode, the vibration speed is automatically weak to strong, and then from strong to weak, and the cycle is repeated, displaying three intensity indicators;

Pressing the **⏻** button for the fourth time, the vibration speed becomes the highest level, showing four intensity indicators;

Shut down

Press the **⏻** button for the fifth time, the massager is turned off, the movement stops all massage vibration function, enters the shutdown state, and the intensity indicator does not light.

Press and hold for 2 seconds to turn off: Press and hold the **⏻** power button for 2 seconds to turn off the power when the product is turned on. The indicator light is off.

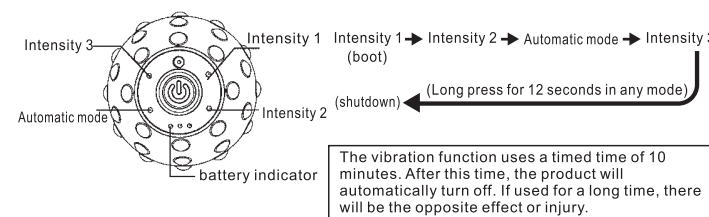

2

Charging method:

How to maintain

Before cleaning the product, make sure the power switch is turned off and unplug the power cord.

Wipe the cleaning product with a clean, soft cotton cloth. Do not use corrosive cleaning agents such as benzene or thinner to clean the product. Do not wipe the switch control surface with a damp cloth to avoid machine malfunction. Wipe it with a dry cloth.

Do not dip the product into water or other liquids. If water is spilled on the product, it may cause electric shock.

Storage method

To keep the product clean and clean, keep the product away from high humidity.

You should ideally store the equipment in its original packaging and keep it in a clean, dry place.

Do not place in direct sunlight or high temperatures.

If not used for a long time. Cover with a dust cover and remove the power plug from the power outlet.

Correct handling:



This symbol indicates that this product cannot be disposed of with household waste in any EU country. In order to avoid any possible damage to the environment and human health, it should be recycled to increase the reusability of the remaining available substances. Recycle used equipment, follow the recycling process, or contact the original retailer. They can recycle products for environmental safety.

3

Common problems and how to deal with them:

| Phenomenon | Possible reason |
|---|---|
| This product does not vibrate even if the power switch is pressed | <ul style="list-style-type: none"> • The body and power adapter are still connected. • The battery is exhausted. Please charge. (The vibration function cannot be used during charging). |
| Vibration intensity weakens | <ul style="list-style-type: none"> • Insufficient power to use. Please recharge. |
| Even if you just rushed through the electricity, This product can only be used for a few minutes. | <ul style="list-style-type: none"> • Extremely low or high ambient temperatures can cause the product to malfunction. Please use it within the recommended ambient temperature range of 10 ° C - 40 ° C. • Battery life is exhausted. Please consult the dealer repair center. (Do not replace the battery yourself. Failure to do so may result in damage to this product. |
| The product suddenly stopped during the process. | <ul style="list-style-type: none"> • The product stops automatically at 10 minutes. Please reboot. |
| There is abnormal noise during the use of the product. | <ul style="list-style-type: none"> • This product structure will emit an action sound when vibrating. Not an exception. |

If it is done, there is still an abnormality in the above processing. Please contact your local dealer for maintenance.

4