auronic

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Introduction

- Thank you for choosing a Auronic product! Please always follow the instructions for safe usage. If you have any questions about the product or experience any problems, please contact us by e-mail: service@auronic.nl.

Intended use

- This product is only intended for the purposes described in the manual.
- The supplier is released from liability for any injury or damage caused by using the product outside of its intended use.
- Always observe the following precautions when using 'the water flosser', hereinafter referred to as 'the product'.

Warnings and Safety Instructions

- ⚠ Please ensure that the packaging is kept out of the reach of children! Risk of suffocation!
- ⚠ Make sure to keep this user manual. Read the manual carefully!

General

- If this is your first time using the product, it is important to read the manual and instructions carefully.
- If the product is not working properly, has been dropped or damaged, do not use the product and contact our customer service.
- Before each use, inspect the parts and accessories for damage or wear (cracks/breaks) that could impair function and if so, replace them with original parts and accessories recommended by the manufacturer.
- Do not use accessories not offered by the supplier. They may pose a danger to the user and damage the product.
- Use the product only for the purposes indicated in this manual.

- Do not disassemble, repair or repair the product. This may damage the product. For your own safety, these parts should only be replaced by authorized service personnel. Always send the product to the service center for repairs to avoid damage and personal injury.
- Be aware of the risks and consequences associated with incorrect use of the product. Improper use of the product can lead to personal injury and damage to the product.
- The product is designed and intended for private use. Do not use the product for commercial purposes.
- Keep the product away from heat, areas where the temperature is relatively high and heat sources such as ovens.
- Do not leave the product unattended when in use. When not using the product, turn off the product and unplug the power cord.
- Do not cover the product, even when it is turned off.
- Fill the product with clean water. Do not put any other liquids in the product.

Persons

- This product may be used by children aged 8 years and above if they are properly supervised and informed about the safe use of the product and understand the potential hazards involved. Cleaning and maintenance of the product should not be carried out by children unless they are older than 8 years and supervised. Keep the product and cord out of the reach of children under 8 years old.
- The product can be used by people with physical or mental disabilities, provided they are properly supervised and informed about the safe use of the product and understand the potential hazards involved.
- Cleaning and maintenance of the product should not be carried out by children unless they are older than 8 years and supervised.
- Do not let children play with the product. Use and store the product out of the reach of children.

Use

- Do not cover the product when it is on. Never place pillows or blankets on or against the product. This could result in a fire hazard, electric shock, personal injury or damage to the product.
- Do not clean the product while it is still plugged in.
- Do not use the product for longer than 2 minutes in a row.
- Do not use the product to spray water into your ears, nose, or other sensitive areas. The water pressure may cause personal injury in these areas.
- Do not use the product near oxygen or aerosol sprays.
- Do not use the product if you wear oral jewelry.

Medical notes

- If you experience pain or discomfort while using the product, stop immediately and consult your doctor.
- This product is not a substitute for treatment from a physician or physiotherapist. Chronic complaints can worsen.
- Do not use the product if you have a wound or inflammation on your tongue or in your mouth.
- Consult your doctor before use if you have dental conditions or have had oral surgery.
- This product does not serve as a substitute for treatments from a (dental) doctor.

Electronics

- If the power fails during use, turn off the product immediately and unplug the power cord.
- If the surface of the product is cracked, please turn off the product to avoid electric shock.
- Periodically inspect both power cord and plug for signs of damage or wear, and if any, do not operate the product and have the power cord and/or plug replaced or repaired by authorized service personnel.
- To avoid overheating, explosion, electric shock, fire, other damage and personal injury, use only the adapter and/or accessories supplied by the supplier.

- To avoid electric shock, do not use the product if you suspect that the power cord or plug is damaged. A damaged power cord or plug must be repaired or replaced by a professional (an authorized distributor).
- Unplug the product and turn off the product when not in use and before cleaning the product.
- Do not unplug the product by pulling on the power cord. This may cause fire or electric shock.
- Do not let the cable hang over the edge of a table or counter.
- Do not drop the product or subject it to large shocks.
- Do not place the cable on hot surfaces or near oil.
- Do not tuck the power cord under rugs or blankets. Make sure that the cable is never in a walkway.
- Do not move the product by pulling or twisting the cord.
- Before first use, remove all packaging material and any stickers.
- Before connecting the product, make sure that the mains voltage indicated on the type rating plate corresponds to that of your mains supply.
- Make sure that you do not trip over the cable and that the cable is not bent, twisted, or caught.
- If black smoke comes out of the product, unplug the power cord immediately. Wait for the smoke to clear before removing the contents of the product and contact our customer service.
- Do not insert objects into the openings of the product. Make sure that the moving parts are not blocked and that nothing is pinched.

Receiving your shipment

- When receiving your shipment, make sure that the contents are consistent with the packing list. Notify your distributor of any missing items. If the equipment appears to be damaged, file a claim immediately with the carrier and notify your distributor at once, giving a detailed description of any damage. Save the damaged packing container to substantiate your claim.

Contents of package


- Auronic Water flosser - Model AU3041 - EAN 8720195255192

Included accessories

- 1x User Manual
- 4x Attachment
- 1x USB-charging cable
- 1x USB adapter (euro plug)
- 1x Travel Bag

Product features

Product overview

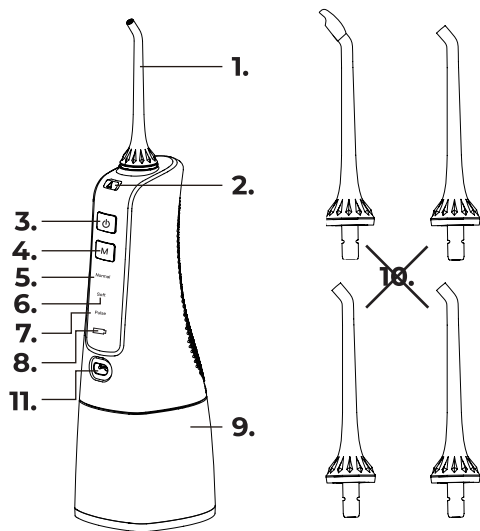- | | |
|---|---|
| 1. Attachment | 7. Indicator light massage mode (Pulse) |
| 2. Release button | 8. Indication light battery |
| 3. On/off button | 9. Water tank |
| 4. Mode button | 10. Attachments |
| 5. Indicator light normal mode (Normal) | 11. Charging point |
| 6. Indicator light soft mode (Soft) | |

Technical data

Dimension	7.6 x 6.7 x 20.4 cm
Weight	305 g
Color	Black
Material	ABS Plastic
Cord length	80 cm
Adapter input	100 - 240V~
Adapter output	5V 1A
Frequency	50/60Hz
Current	5W
IP number	IPX7
Battery	3.7V Lithium 1400mAh
Pressure	30 - 100 psi
Pulses	1400 - 1800 times/minute
Water tank capacity	300 ml
Charging time	4 hours
Battery life	20 days (if fully charged)
Modes	3: normal, soft, massage

Other features

- Removable water tank;
- 360° rotatable attachment;
- Cordless use;
- Memory function to remember your preference mode;
- 4 attachments including 1 orthodontic.

Operating instructions

Before use


- Fully charge the product before use. This takes about 4 hours. Do not use the product while it is charging.

Attach/remove attachment

- Attach: place the attachment at the top of the product's opening and turn until you hear a click.
- Remove: press the release button and remove the attachment from the opening by gently pulling it upwards.


Using

1. Make sure the product is off. Open the lid of the water reservoir and pour in water. Note: Make sure not to exceed the 300 ml line.
2. Choose the desired mode by pressing the mode button: standard, soft or massage. You can choose the desired mode both when the product is on and when the product is off.
3. Holding the product at an angle downward, place the end of the attachment in your mouth and aim at the teeth. Position the product so that the water sprays at a 90° angle to your gums. You can rotate the attachment to find the right position. Close your mouth slightly to prevent the water from squirting out.
4. Turn the product on by pressing the power button.
5. Start on the outside of your upper teeth (1). Then continue with the inside (2). Then continue with the inside of your lower teeth (3). Finally, start working on the outside of your lower teeth (4).

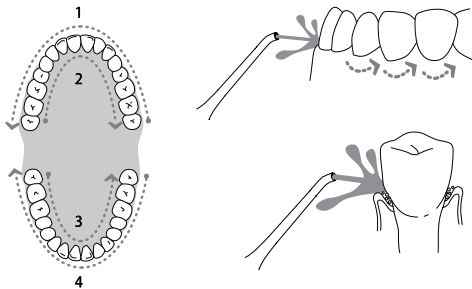
- The product turns off automatically after 2 minutes of continuous use.
- If the product is in massage mode, the battery life is displayed after 3 seconds.
- The battery is low if all three indicator lights flash during use. In this case, charge the product.
- The product has a memory function. When you turn on the product, it is automatically on the last used mode.

After use

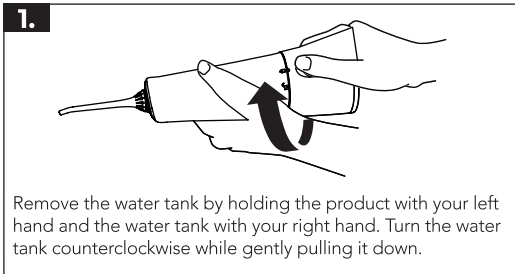
Remove the water tank by holding the product with your left hand and the water tank with your right hand. Turn the water tank counterclockwise while gently pulling it down.

2.

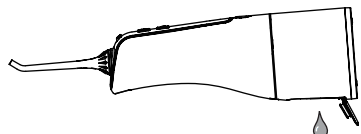Empty the water tank by opening the water tank lid and pouring it out.

3.

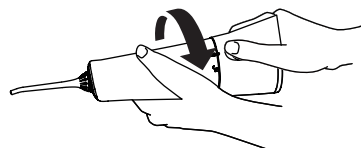Place the water tank back on the product by rotating it clockwise back onto the product.

Maintenance and cleaning

- If you have not used the product for a long time, let it fully charge before using it again. Charge the product in a timely manner to extend the battery life.
- Clean the product with a damp, clean cloth and dry it thoroughly. Do not use abrasive or harsh cleaning agents.

Storage

- Store the product in a clean and dry place. Do not store the device in extreme high or extreme low temperatures.

Disposal - recycle

- Think about the environment and contribute to a cleaner living environment! This product is subject to the European Directive (EU) 2012/19/EU. At the end of its service life, do not dispose of the product with normal household waste, but deliver it to a collection point designated by the authorities for recycling.

Service and warranty

- If you need service or information regarding your product, please contact our Customer Support at service@auronic.nl.
- We provide a 2-year warranty on its products. To obtain service during the warranty period, the product needs to be returned with a proof of purchase. Product defects have to be reported within 2 years from the date of purchase. This warranty does not cover damages caused by misuse or incorrect repair; parts subject to wear; defects that the customer was aware of at the time of purchase; damages or defects caused by customer neglect; damages or defects caused by third parties.

Disclaimer

- Changes reserved; specifications are subject to change without stating grounds for doing so.

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© **LifeGoods B.V.**