MINI

HYPERSPHERE

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READ ALL INSTRUCTIONS BEFORE USING THE HYPERSPHERE MINI

DANGER

To reduce the risk of electric shock: Do not reach for an appliance that has fallen into water. Unplug immediately. Do not use while bathing or in a shower. Do not place or store appliance where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.

WARNING

To reduce the risk of burns, fire, electric shock, or injury to persons: Do not charge unattended or leave charging overnight. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons. Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons. Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer. Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair. Do not carry this appliance by supply cord or use cord as a handle. Keep the cord away from heated surfaces. Do not use outdoors. Do not operate where aerosol (spray) products are being used or where oxygen is being administered. Do not immerse unit in water. Keep liquids away from ventilation ports, buttons and charging port. Do not remove screws or attempt to disassemble. Do not operate continuously for more than one hour. After one hour use, allow unit to rest for 30 minutes before re-using. Unplug the unit after charging or prior to use. Do not use directly on a hard surface. Use only on a carpeted floor or a fitness mat. Persons weighing above 250 lbs./113 kgs may damage or break the unit during normal intended use. This is not a toy. For adult use only. Do not use if injured. Consult your doctor before using this product. Use product only according to Instructional Use materials available at hyperice.com.

CAUTION

The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not disassemble, heat above (manufacturer's maximum temperature limit), or incinerate. Replace battery with (battery manufacturer's name or end product manufacturer's name and part number) only. Use of another battery may present a risk of fire or explosion. Dispose of used battery promptly. Keep away from children. Do not disassemble and do not dispose of in fire. Please do not

use the Hypersphere Mini, or any vibration high-intensity exercise device without first obtaining approval from your doctor if any of the following apply: Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pace-makers, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements or IUD's, metal pins or plates or any concerns about your physical health. Frail individuals and children should be accompanied by an adult when using any vibration device. These contra indications do not mean that you are not able to use a vibration or exercise device but we do advise you to consult a doctor first. Ongoing research is being done to the effects of vibration exercise for specific medical disorders. This will likely result in the shortening of the list of contra indications as shown above. Practical experience has shown that integration of vibration exercise into a treatment plan is advisory in a number of cases. This must be done on the advice of and in accompaniment with a doctor, specialist or physiotherapist.

OPERATING INSTRUCTIONS

READ ALL INSTRUCTIONS BEFORE USING THE HYPERSPHERE MINI

CHARGING

Fully charge for two hours before first use. To charge, connect the USB cable from the charging port to any standard USB charger (not included). LED lights illuminate to show battery level and indicate active charging. LED lights will remain illuminated when fully charged while adapter is plugged in. Unit may be recharged at any time and at any battery level. It is not recommended to fully deplete the battery to the red LED level. Recharge the unit monthly to best prolong battery life. The average run-time is 2 hours at level three, and longer at lower levels.

OPERATING

To turn on, press down on the power button for three seconds until LED lights illuminate. Press again once for level one, twice for level two and three times for level three. The blue LED lights will indicate the corresponding power level. Do not use directly on a hard surface. Use only on a carpeted floor or a fitness mat. To turn off, press down on the power button for three seconds. The unit is fully turned off when all the LED lights are off. The unit will turn off after 10 minutes of use. To continue use, turn unit on again and set to desired level.


1. Power Button
2. Battery level indicator LED (green, red)
3. Frequency level indicator LED (blue)
4. USB Charging Port
5. USB Charging Cable (not shown)


SPECIFICATIONS

- Rated Charging Voltage: 5V DC- 1A
- Rechargeable Lithium Ion battery
2400mAh
- Frequency: Low 42Hz, Medium 52Hz,
High 62Hz
- Motor RPM: Low 2500, Medium 3100,
High 3700
- Weight: .75 lbs / .34 kg
- Designed in USA. Assembled in China

MAINTENANCE, CLEANING & STORAGE

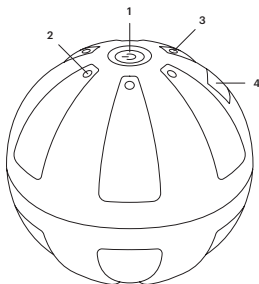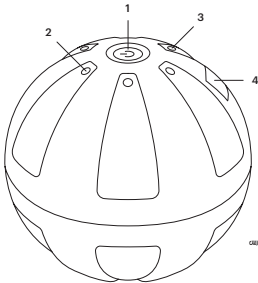- Use a towel or soft cloth to wipe any moisture from unit.
- Keep unit away from heat sources, liquids, and children.

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