

Frequently asked questions(FAQ)

1.Can a minor (person under age 19) use this product?

A: This products is only recommended for use by an adult

2. Does the light of the product have impact on the eyes?

A: Because of the unique surface design of this product, it illuminates light sources evenly on the face and is structurally designed not to harm near and around the eyes. In the case of sensitivity, please use product with eyes closed.

3. Is there a limit to how often this product can be used?

A: This device is safe to be used everyday

5. Does this product have any side effect to the skin?

A: The periods of LED optic therapy or LED photom treatment is suitable for non-invasive painless

Repairing

Problem	Cause of trouble	Action for trouble
Not turning on the light	Disconnection of plug	Connect to the plug
	Short of electricity in auxiliary batery	Connect to the plug after charging it
	Cannot push the button	Push the switch button
	Out of order	Contact customer e center



7 Color LED Masker User Manual

Please read the manual carefully to ensure safe and correct use.

Caution

Only use mask wireless


1. Use on clean skin, do not use harsh products such as acids, spot treatments and retinols under the mask.
2. Electric circuit insider the device should not be pressed or folded and it should be placed flat to protect the shield after every use.
3. Store shield back in box to ensure the shield does not damage.

Principles of the product

This product uses 3 types of light sources including red, blue, and orange colors; lowers the acne development by illuminating the light sources evenly with different wavelengths using the principle of bio-ray therapy and light material with nanotechnology; and accelerates the microcirculation of the cells by stimulating collagen generation, therefore it is helpful to skincare effect.

Name of each component

- 1.USB charging line
- 2.Frame of nose ridge
- 3.Switch

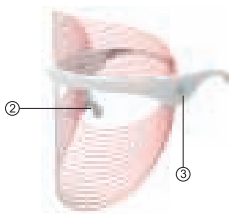

Function of the product

The processes of LED optic therapy or LED photon treatment belong to non-invasive painless skincare, and their effects are significant without harmful light. Hence, it is used widely by dermatologists. Wavelength of red light varies from 620 to 750nm and it facilitates blood circulation, accelerates metabolism, generates collagen, helpful to the wrinkle care and the skin improvement.

Wavelength of blue light varies form 476 to 495nm and it removes propionibacterium acnes, is helpful to the acne care b improving the environment in the hair follicle.

Wavelength of orange light varies from 590 to 620nm and it improves the function of oxygen exchange in the cells, accelerates the blood circulation together with fading he pigments therefore, is helpful to improves to soften skin.

Product specifications

Model	LJ-106
Product name	Three Color LED Mask
Product size	210*165*125mm
Product weight	73g
Material	ABS
Wavelength	Red 620-750nm Blue 476-495nm Orange 590-620nm
Output voltage	5V 1A
Accessories	LED mask, Manual and cable line

How To Use

1. Cleanse skin
2. Apply soothing and hydrating products
3. Push button for 3 seconds to TURN ON the shield
4. Push the button until desired LED light colour appears
5. Place the shield on face for 10 to 20 minutes
6. Hold button down until shield POWERS OFF

LED Light Therapy can be used daily for optimal results

Maintenance and storage

Maintenance:

1. Polish the main body with a dry cloth after use, and do not wash with water directly
2. Always unplug after use and do not disassemble the product.

Storage:

1. Store this product in a cool and dry place avoiding humidity and direct sun light.
2. Do not submerge the product in water or liquid.
3. Keep out of reach of children.