Handleiding

Manual

**Massagegun XG10 – met 10
massagekoppen**

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table of contents

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1. Massage Gun XG10

1.1 Specifications

- Charging time: 3-4 hours
- 10 massage settings
- 10 different massage heads
- Usage time: 4-6 hours
- Amplitude of 8mm
- With an LED display
- 60W motor
- Lithium battery
- Weight: 850g
- Complies with CE marking and EU declaration of conformity

1.2.1 General Safety Warnings

- ⚠ This massage gun is designed for massaging muscles in adults. Do not use the device for other purposes and never on pets.
- ⚠ Do not use the massage gun if there is visible damage to the casing, battery, or charger.
- ⚠ Keep the massage gun away from liquids and heat sources. Never submerge the device in water.
- ⚠ Do not disassemble the massage gun or remove any screws. This may lead to injury or damage.
- ⚠ Allow children or individuals with limited cognitive abilities to use the device only under direct adult supervision.

1.2.2 Safe Use on the Body

- ⚠ Do not use the massage gun without a doctor's permission if you have: pregnancy, heart problems or pacemakers, diabetes with complications, epilepsy, migraines, hernias, recent surgeries, joint prostheses, metal implants, or other health issues.
- ⚠ Avoid use on open wounds, bruises, inflammations, fractures, eczema, rashes, or other skin conditions.
- ⚠ Move the massage head slowly over the muscle group. Do not press too hard and never keep the device in one spot for too long.

- ⚠ Do not use the massage gun on the head, throat, spine, or directly on bony protrusions.
- ⚠ The massage gun automatically shuts off after 20 minutes. Use it for a maximum of one hour continuously and allow the device to cool down for at least 30 minutes afterwards.
- ⚠ Bruising or slight muscle tenderness may occur. This is normal, but stop immediately if you experience pain, dizziness, nausea, or tingling.

1.2.3 Electronic and Charging Warnings

- ⚠ Charge the massage gun only with the supplied charger and cable.
- ⚠ Charge the device at room temperature, in a dry environment, and never unattended.
- ⚠ Never use the massage gun while it is charging.
- ⚠ Disconnect the charger as soon as the battery is fully charged.
- ⚠ If the device becomes unusually warm, produces strange noises or smells, or emits smoke: turn it off immediately and stop using it.

1.2.4 Maintenance and Storage

- ⚠ Clean the casing with a dry, soft cloth. Do not use moisture, solvents, or aggressive cleaning agents.
- ⚠ Keep the ventilation openings free from dust and dirt.

- ⚠ Store the massage gun in the supplied case, in a dry and cool place. Do not expose the device to direct sunlight, moisture, or extreme temperatures.

1.2.5 Environment and Responsible Disposal

- ⚠ The massage gun contains a lithium battery and must not be disposed of with regular waste. Return the device to a recognized e-waste collection point.
- ⚠ Separate cardboard and plastic packaging for recycling.

2. Installation

2.1 Product diagram

2.2 Associated accessories

Round head
Larger muscle groups – Back and thighs

Fork head
Spine & neck

Flat head
Sensitive and bony areas

Pneumatic head
For sensitive areas or beginners

Thick flat head
Wide muscle groups

Thumb-shaped head
for deep, precise pressure like hands and feet

Bullet head
Deeper tissue massages

Shoulder head
For the shoulder area

Conical head
For joints and pressure points

Spiral head
For acupressure and deeper muscles

USB-C charging cable

Attachment rubbers

Carrying case

2.3 Use of the massage gun

Attention! Check if the massage gun is fully charged before using it.

2.3.1 Turning on the massage gun

Turn on the massage gun using the slide switch at the bottom.

2.3.2 Adjusting speed

Click on the fingerprint on the LED display to adjust the speed.

2.3.3 Turning off the motor

By holding down the fingerprint, the motor of the massage gun will turn off.

2.3.4 Turning off the massage gun

Completely turn off the massage gun using the slide switch at the bottom.

2.4 Changing the massage head

Attention! Never change the massage head while the massage gun is on.

2.4.1 Confirming

Press the pin of the massage head as far as possible into the massage gun.


2.4.2 Removing

Pull on the massage head to remove it.

2.4.3 Replacing rubbers

With frequent changes of massage heads, the rubber rings on the pins can wear out. Replace them in time with the included spare rubbers to maintain proper function.

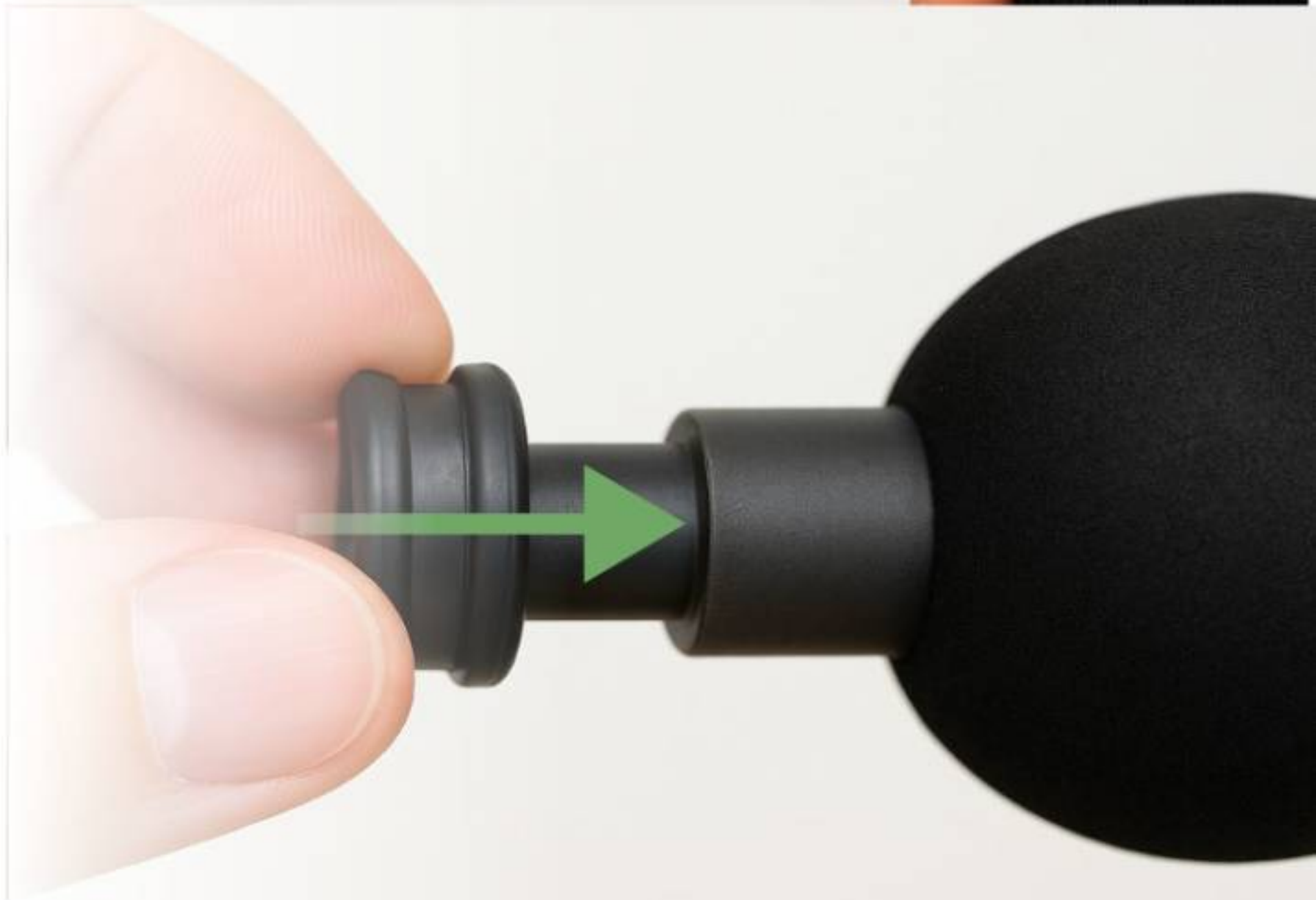3. Sustainable recycling

Your contribution to a better environment

We value sustainability highly and advise our customers to recycle products responsibly at the end of their life cycle. Do not throw products in the general waste, but hand them in at a recognized recycling organization or take them to the recycling center in your municipality.

Together we can ensure a more sustainable future!

4. Warranty provision

4.1 Quality guarantee

2BEHOME guarantees its customers a high standard of Dutch quality, and we take pride in that. Therefore, all our products come with a standard warranty period of 2 years against defects in materials or workmanship.

4.2 Warranty expiration

Unfortunately, we cannot provide a warranty in the following cases:

- a. After the expiration of the above warranty period.
- b. The defect is due to improper use, neglect, or due to falls and bumps.

4.3 Claiming warranty

Do you want to claim warranty or are you unsure whether your situation falls within the warranty? Please contact us so we can discuss the warranty option together!

4.4 Right of return

Are you not completely satisfied or do you have another reason to return the product? You can return the product within 30 days unless otherwise stated in the terms of specific platforms. The terms of the relevant platform determine whether a product can be returned and within what timeframe.

Make sure the product is undamaged and returned in its original packaging. After receipt and acceptance of the return, you will receive a full refund, minus any return costs (if applicable).

4.5 About 2BEHOME

2BEHOME has been the specialist in everything related to home since 2020. At 2BEHOME, you are at the right place for the best quality aroma diffusers, toilet brushes, dryer balls, and more!

2BEHOME is complemented by the associated brands 2BEFUN & 2BEHEALTHY, which stand out in products including the best massage gun, yoga mats, PS5 cases, and Spikeball sets!

Each product is extensively tested for quality, durability, and sustainable production before being included in the range. This ensures you are always guaranteed the best products in the Netherlands and Belgium! What product are you looking for? Below you will find an explanation of all products!

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2BEHOME B.V.