beper

IT	TERMOMETRO MULTIFUNZIONE AD INFRAROSSI	pag.2
EN	MULTIFUNCTIONAL INFRARED THERMOMETER	pag.9
FR	THERMOMÈTRE INFRAROUGE MULTIFONCTIONNEL	pag.15
DE	INFRAROT-MULTIFUNKTIONSTHERMOMETER	pag.22
ES	TERMÓMETRO INFRARROJO MULTIFUNCIÓN	pag.29
GR	ΠΟΛΥΛΕΙΤΟΥΡΓΙΚΟ ΘΕΡΜΟΜΕΤΡΟ ΥΠΕΡΥΘΡΩΝ	pag.36
RO	TERMOMETRU MULTIFUNȚIONAL CU INFRAROȘU	pag.43
CZ	MULTIFUNKČNÍ INFRAČERVENÝ TEPLOMĚR	pag.49
NL	MULTIFUNCTIONELE INFRAROODTHERMOMETER	pag.56
LV	DAUDZFUNKCIONĀLS INFRASARKANAIS TERMOMETRS	pag.63

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General safety advice

READ THE OPERATING INSTRUCTIONS CAREFULLY BEFORE FIRST USE OF THE APPLIANCE.

Read the following instructions carefully, since they give you useful safety information about installation, use and maintenance and help you avoid mishaps and possible accidents.

Remove the packing and make sure that the appliance is intact, with special attention to the supply cable. The elements of the packing (plastic bags, polystyrene, etc) do not have to be accessible to the children, in order to preserve them from danger sources.

It is recommended not to throw these packing material in the domestic waste, but to deliver them to the appropriate station of collection or destruction refusals, asking eventually information to the manager of the Sanitation Department of your own municipality. This appliance is to be destined to its special purpose only.

Any other use is to be considered as improper and thus dangerous. The manufacturer cannot be responsible for any damages deriving from improper, wrong or incautious use. To avoid any heating risk, it is recommended to unroll completely the supply cord and to unplug the appliance, when not used.

To protect against electric shock, do not immerse cord, plug or any parts of the appliance into water or other fluids.

Do not hold the appliance with wet hands or feet.

Do not hold the appliance with naked feet.

Do not expose the appliance to atmospheric agents such as rain, wind, snow, hail.

Do not use the appliance near explosives, highly flammable materials, gases, burning flames, heaters.

Do not allow children, disabled people or every people without experience or technical ability to use the appliance if they are not properly watched over. Close supervision is necessary when any appliance used by or near children.

Make sure that children do not play with this appliance.

If you decide not to use the appliance any longer, make it unserviceable by simply cutting the supply cable after disconnecting the plug from the socket.

Also make any dangerous parts unserviceable, especially for children who could play with the appliance.

Introduction

This infrared thermometer is specifically designed for safe measurement of forehead temperature. This thermometer can measure body temperature by detecting the intensity of infrared light emitted by the forehead; it then converts the detected heat in temperature measurement that is shown on the LCD display. This thermometer is designed for intermittent measurement of human body temperature through the skin surface of the forehead, and can be used on people of all age.

Warnings

Carefully read this manual before using the item.

There are no genre nor age limitations for the use of this device.

This device is exclusively designed for home use. The use of this device does not replace medical consult. Keep the device and its batteries out of the reach of children. Do not let children play or use the device without supervision; some parts are very small and can be easily swallowed. Do not immerse the device in water or other liquids. Do not try to modify and/or tamper with the device in any way.

Do not expose the device to extreme temperatures (below $-25^{\circ}\text{C}/-13^{\circ}\text{F}$ or above $55^{\circ}\text{C}/131^{\circ}\text{F}$), or to high humidity ($>95\%$ RH). Remove batterie when the device is not used for long periods of time.

Description of product Fig. A

1. Infrared probe
2. START button
3. button 
4. button 
5. Battery cover

Description of LCD display Fig. B

Mounting and/or replacing of batteries

This thermometer works with 2 alkaline batteries, type AAA (not included). During first use, remove battery cover by sliding it downwards, and insert the batteries following the polarity indicated in the battery case; then place the battery cover back into position. Proceed in the same way for replacing batteries, when the symbol appears on the display. 

Date and time setting

When thermometer is off, keep  pressed for about 3 seconds. Select time format (12H/24H) pressing START button; once selected format appears on the display, press  time starts flashing on the display.

To set the hours, press START button until the correct hour appears; press  to confirm selection.

To set the minutes, press  again; then press START button until the correct minutes appear. Press  to confirm selection.

Once hour and minutes have been confirmed, press  again: year will start flashing on the display. Press START button until the correct year appears; then press  to set the month, and press START button until the correct month appears; press  once again to set the day, and press START button until the correct day appears. Press  one last time to confirm selection and to exit date and time setting.

Measurement scale setting

This thermometer can measure temperature both in Celsius (°C) and Fahrenheit (°F) degrees.

To select measurement scale, when thermometer is off press START button for 3 seconds until the measurement scale appears on the display. Press START button again until the desired scale appears on the display; at this point, press  to exit setting.

Hints for correct assessment of temperature

To ensure most accurate reading, it is necessary to consider different factors that can affect the measurement, and to follow some basic hints.

It's important to know each individual's normal temperature when healthy and not ill; this is the only way to correctly diagnose fever.

It is advisable to stay inside for at least 30 minutes before assessing body temperature.

Note: this thermometer and the user need to be in the same room with the same room temperature for at least 10 minutes before reading

It's better not to eat nor drink and to avoid hot baths, showers, etc. before taking the temperature.

Remove all hair, headbands, hats etc. from forehead at least 10 minutes before taking the temperature.

Lotions, oils, make up and other cosmetic preparations might affect the reading of the temperature, which can be lower than actual one. After washing forehead, wait at least 10 minutes before taking the temperature.

Do not take temperature on scar tissues, open sores or abrasions. Do not use thermometer on sweating or wet forehead.

Do not take temperature during or immediately after nursing a child.

Do not use this thermometer outdoors.

Do not use this thermometer to assess temperature in proximity of surfaces/places with high temperatures (e.g. fireplaces or stoves).

The probe window is the most delicate part of the thermometer; do not touch nor damage it in any way.

Keep the thermometer in the same room where reading takes place.

If the thermometer is stored in a significantly different environment than testing location, make sure to place it in the testing location for at least 30 minutes.

Do not use this device in oxygen-rich environment, and in environments where flammable anesthetic mixtures with air, oxygen or nitrous oxide can be found.

Instructions for use

Measurement of body temperature (forehead)

Press ; the thermometer emits a beep sound, and the display activates, showing all the available functions and backlight colors. After a few seconds, the thermometer makes another beep sound, confirming that it is ready for reading.

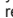
Aim the probe of the thermometer at the center of the forehead, with a distance of 2-3 cm, and press START button.

Note: do not remove the thermometer from the forehead before hearing the final beep sound.

The result of the measurement is shown on the display; backlight will be in 3 different colors, based on temperature assessed:

green  when the thermometer is ready for reading, and when temperature assessed is below 37.3 °C (99.1°F)

yellow  when temperature assessed is below 37,8 °C (100.0°F)

red  when temperature assessed is equal to or above 37,8 °C (100.0 °F)

Press  to turn thermometer off.

Waiting at least 3 seconds between each reading, to ensure complete accuracy.

To turn the thermometer off after taking the temperature, press ; display will show the OFF writing.

Thermometer will automatically shut off if not used for 1 minute, to preserve battery life.

Measurement of surface temperature

To select surface mode, when the thermometer is on press /. Aim the probe window at the center of the object or surface, with a distance of 1-2 cm, and press START button.

Wait for the second beep sound before removing the thermometer.

Press  to turn thermometer off.

Memory mode

Memory mode can be used both in forehead and surface mode.

When thermometer is on and ready to be used, or after taking the temperature, long press  for 3 seconds: the symbol M appear on the display (to the right). At this point, the thermometer will show last 10 readings; to see them, keep pressing  (each reading is preceded by its correspondent number from 1 to 10, where reading number 1 represents the latest and reading number 10 represents the oldest).

Memory mode can be used for both body temperature reading and surface/objects/liquids temperature reading.

Cleaning and maintenance

The probe window must be kept clean and dry, to ensure accuracy of reading; to clean this extremely delicate part, use a soft cloth slightly dampened with a 75% isopropyl alcohol disinfectant solution. Do not use abrasive detergents. Do not use chemicals other than isopropyl alcohol to clean the probe window.

After cleaning, wait at least 10 minutes to allow the thermometer to dry completely before using it. Use a soft, dry cloth to clean outer body and display.

⚠ Warning: do not immerse thermometer in water.

Store the thermometer in a clean, dry environment, away from dust, dirt and other contaminants, and away from direct sunlight.

Strong electromagnetic fields may compromise the correct operation of the thermometer. It's better to store the thermometer in its original box when not in use.

Specifications

Measuring range	Forehead mode: 34.0°C – 43.0°C (93.2 °F – 109.4 °F) Object mode: 0°C – 100°C (32°F – 212 °F)
Measuring site	Forehead (in forehead mode)
Reference body site	Oral (this thermometer converts the forehead temperature to display its "oral equivalent")
Operation mode	Forehead mode (adjust mode)
Laboratory accuracy	Forehead mode: ± 0.2°C (0.4°F) during 35.0°C ~ 42.0°C (95.9°F ~ 107.6 °F) at 15°C ~ 35°C (59.0°F ~ 95.0°F) operating temperature range ± 0.3°C (0.5°F) for other measuring and operating temperature range Object mode: ± 4% or ± 2°C (4°F) whichever is greater
Display resolution	0.1°C or 0.1°F
Measure time	approx 1 second
Operating temperature range	10°C ~ 40°C (50°F ~ 104°F), 15% ~ 85% RH, non-condensing Atmospheric pressure: 700 hPa ~ 1060 hPa
Transport and storage temperature range	-25°C ~ 55°C (-13°F ~ 131°F), 15% ~ 95% RH, non-condensing Atmospheric pressure: 700 hPa ~ 1060 hPa
Clinical accuracy	Clinical bias: -0.09°C (-0.16°F) Clinical repeatability: 0.13°C (0.23°F) Limits of agreement: 0.87°C (1.57°F)
Shock	Withstands drops from 3 feet (0,90 m approx)

Dimension	153*41*44 mm
Weight	approx 84 g (with batteries)
Batteries	DC3V (2ñ batterieAAA)
Battery life	approx 1 year / 6000 readings
Ingress protection rating	IP22

Troubleshooting

Error	Problem	Solution
Er1	Measurement taken before thermometer is ready	Take temperature when the symbol  or appears on the display
Er2	The room temperature is not within the range between 10°C and 40°C (50°F ~ 104°F)	Place the thermometer in a room for at least 30 minutes with a temperature between 10°C and 40°C (50°F ~ 104°F)
Er3	The thermometer is not correctly placed.	Refer to "Instructions for use" paragraph for correct use and proceed with a new reading
Er4	The thermometer detects a rapid change in the room temperature	Allow the thermometer to rest in a room for at least 30 minutes at room temperature between 10°C and 40°C (50°F ~ 104°F) before taking the temperature
Er5	The thermometer is not functioning properly	Remove batteries, wait 1 minute and place batteries back in the battery case. If the error message still appears, contact distributor
Hi	In forehead mode: measured temperature is higher than 43.0°C (109.4°F) In object/surface mode: measured temperature is higher than 100°C (212°F)	Refer to "Instructions for use" paragraph for correct use and proceed with a new reading
Lo	In forehead mode: measured temperature is lower than 34.0°C (93.2°F) In object/surface mode: measured temperature is lower than 0°C (32°F)	Refer to "Instructions for use" paragraph for correct use and proceed with a new reading
	Low battery	Replace the batteries

The European directive 2011/65/EU on Waste Electrical and Electronic Equipment (WEEE), requires that old household electrical appliances must not be disposed of in the normal unsorted municipal waste stream. Old appliances must be collected separately in order to optimize the recovery and recycling of the materials they contain, and reduce the impact on human health and the environment. The crossed out "wheeled bin" symbol on the product reminds you of your obligation, that when you dispose of the appliance, it must be separately collected.

GUARANTEE CERTIFICATE

This appliance has been checked in the factory. From the date of original purchase a 24 month guarantee applies to material and production defects. The purchase receipt and the guarantee certificate must be submitted together in case of claim to guarantee.

The guarantee is valid only with a guaranty certificate and proof of purchase (fiscal receipt) indicating the date of purchase and the model of the appliance

For any technical assistance, please contact directly the seller or our head office in order to preserve the appliance's efficiency and NOT TO void the guarantee. Any intervention on this appliance by non-authorized persons will automatically void the guarantee.

GUARANTEE CONDITIONS

If the appliance shows defects as a result of faulty material and/or production during the period of guarantee, we guarantee repair free of charge on condition that:

- The appliance has been used properly and for the purpose for which it has been intended.
- The appliance has not been tampered with, otherwise it is not maintainable.
- The purchase receipt shall be presented.
- The appliance showing a fair wear and tear shall not be covered by this guarantee.

Therefore, any part that could be accidentally broken or having visible signs of use in consumable products (such as lamps, batteries, heating elements...) the aesthetic parts are excluded from the guarantee, and any defect whatsoever resulting from non-respect of the rules for use, negligence in usage and/or maintenance of the appliance, carelessness, wrong or improper installation, damage during transport and any other damage not attributable the supplier.

For every defect that could not be repaired within the guarantee period, the appliance will be replaced free of charge.

In any case, if the part to be replaced for defect, breakage or malfunction is an accessory and/or a detachable part of the product, Beper reserves the right to replace only the very part in question and not the entire product

For technical assistance and/or repairs after the guarantee period you can directly contact the below address:

Contact your distributor in your country or after sales department beper.

E-mail assistenza@beper.com which will forward your enquiries to your distributor.

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Beper Group S.p.A. - Via Saliere, 30 - 37050