



# Portable Oral Irrigator User Manual

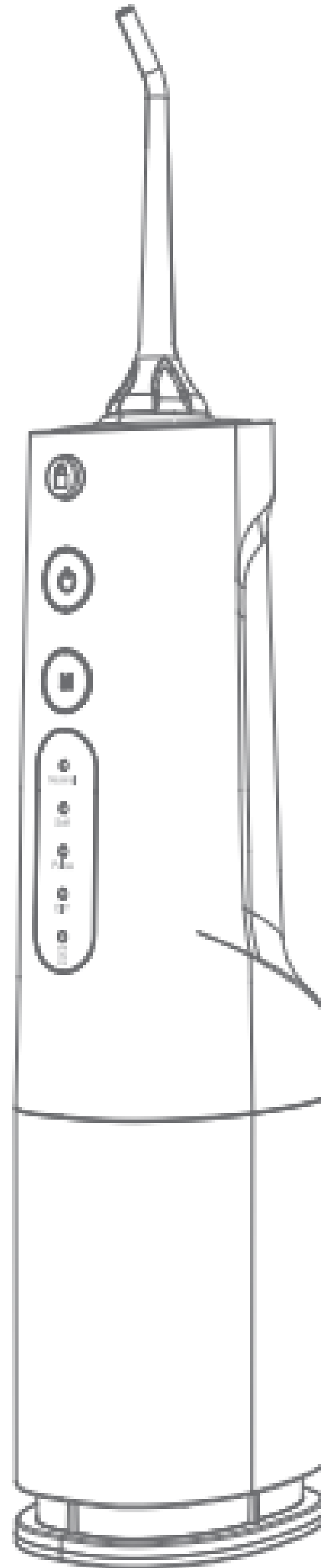

## Special Tips

- ◆ **Check Battery Volume :** When the product is in standby, long press the mode button to show the battery volume . If above 75% battery volume , 4 mode lights will be bright ; 50%-75% battery volume, 3 mode lights will be bright ; 25%-50%, battery volume, 2 mode lights will be bright ; below 25% battery volume, Only DIY mode light will be bright.
- ◆ **DIY Mode Setting:** When switching to the fourth "DIY" working mode, the system remembers the last use of the fourth mode status to start working. In the fourth mode working status, press the ON/OFF Button again, the working status will raise from the lowest water pressure to the highest water pressure and stay at the highest water pressure mode temporarily , it will circulate according to this mode. Release the ON/OFF Button to indicate the selection of the working status (new status is locked), the system maintains and work under the newest working status.

## Product Accessories

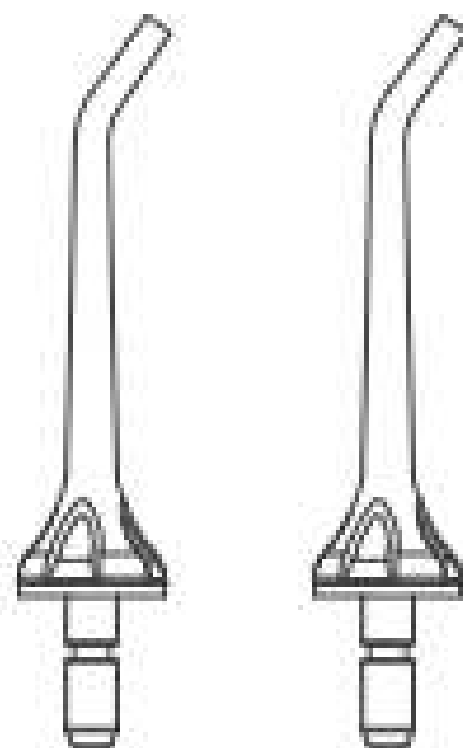

2 Jet tips  
1 USB Cable




## Notice for usages

- ◆ In the first time usage or not used for a long time, you shall wait for a moment before the Oral Irrigator's jet tip can discharge water after pressing the ON/OFF Switch. Under the circumstances, follow these steps: Under such circumstances, follow the steps as below:
  - ◇ Pull the inner water tank at full expanding to the bottom end first with hearing the click sound, it means the inner water tank has been pulled in right place (water will leak from the inner water tank if not do so), then open the water inlet cover to fill water into the water tank.
  - ◇ After adding adequate amount of water into the jack of the jet tip seat, install the jet tip again, then switch ON/OFF Button.
  - ◇ Kind notices:
    1. set 'soft' mode for using in the first time.
    2. keep the jet tip pointing to your tooth gaps closely during using. Farther you keep, stronger the waterflosser pressure will be.

1. Do not disassemble and repair products' internal components.
2. This model's jet tip can be replaced and universal for same products from our company.
3. Pls check the USB Cable if damaged or not before and after usage, keep the USB Cable away from heated surfaces.
4. Do not put or insert any strange objects into any opening jacks or hoses
5. Do not touch the product's power supply or charge the product in a non-dry environment.
6. Do not place or store the product where it may fall or be pulled into a bathtub, shower or sink.
7. The product can only be filled with water or other non-corrosive liquid recommended by professional doctors.


## Warning

1. If you have had oral or dental surgery in the past 2 months, please consult your doctor before using this product and follow the doctor's advice and guidance.
2. When using this product for oral cleaning in the first time, you may experience minor gum bleeding even under the case of healthy gums. This is because your gums are stimulated by the jet pressure for the first time, and usually stop bleeding after 1--2 weeks of adaptation. If your gum bleeding lasts for more than 2 weeks, it may mean abnormal conditions such as periodontal disease in your gums. In this case, please stop using the product and consult your dentist.
3. Do not use this product if the product's USB cable or plug is damaged, or it cannot work, falls, is damaged, or falls into any liquid.
4. Do not connect this product to a voltage system other than this device or charger.
5. This product is not suitable for persons (including children) who have poor physical, sensor or mental abilities or lack of experience, without the knowledge of the guardian.
6. Do not inject water into the bottom of tongues, ears, noses or other sensitive body parts under working status. The pressure generated by this product may cause serious damages.

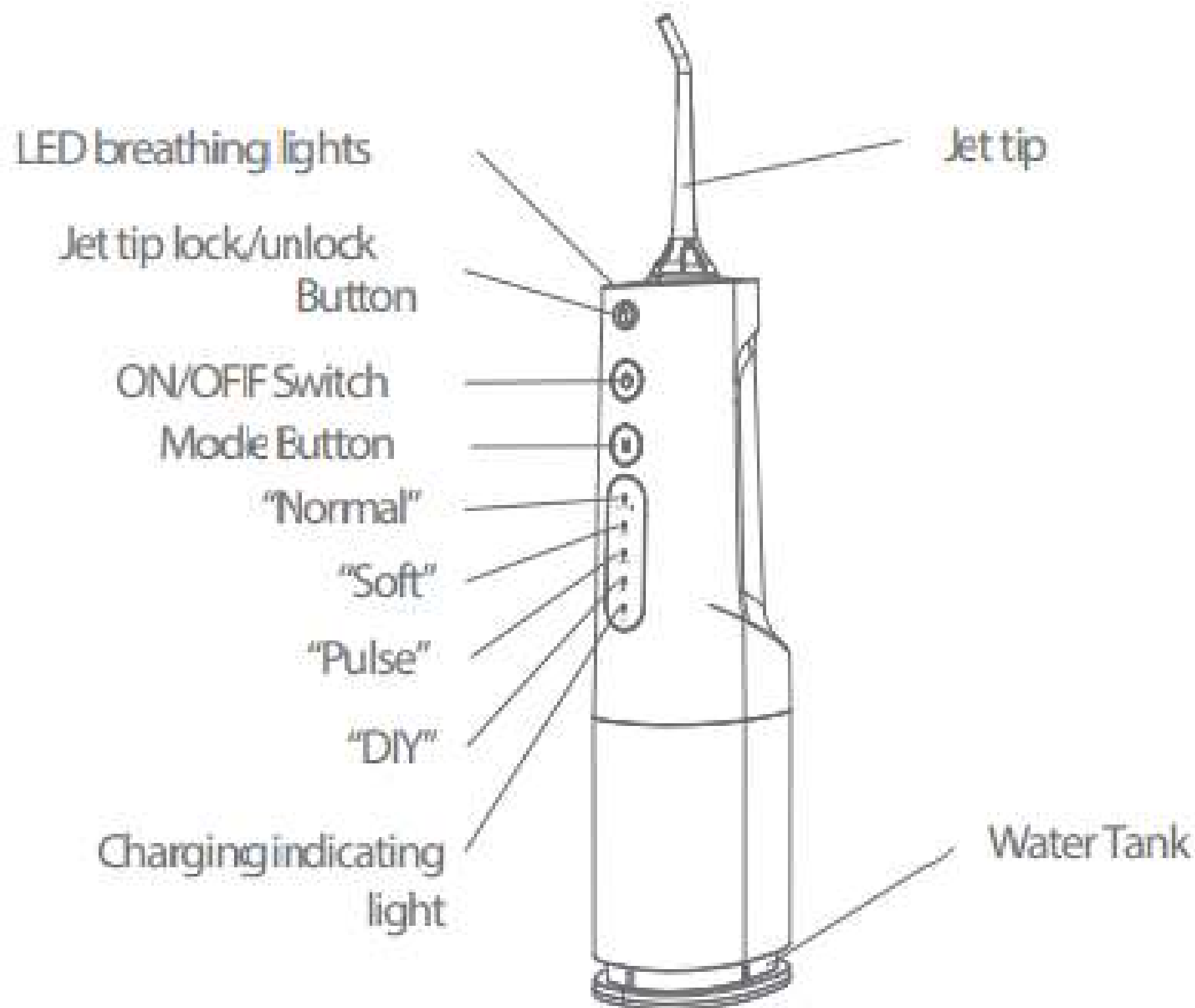

## Product Features

1. 360ml detachable expandable double water tank, small body, portable, the water tank is easy to clean.
2. Combining different using circumstance requirements We add colorful breathing lights, soft lighting effects and three-color indicator status for bringing you different using experience.
3. IPX7 full body waterproof so it can be used in the shower and supports full body with water washing.

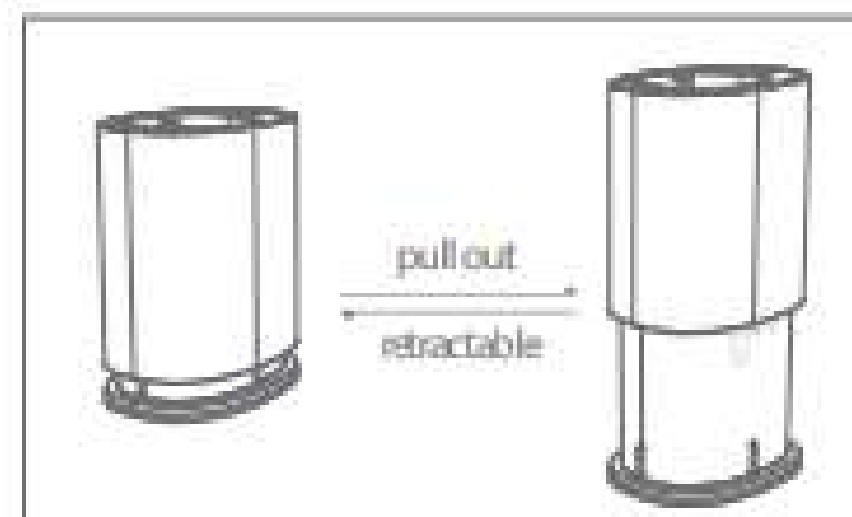

1. Start working: Be sure of fully charging the device before using for the first time. First press the On/Off Button, then press mode button to select working mode.
2. Filling the Water Tank: Open the water inlet of the water tank, fill the water tank with water, close the inlet window, and press tightly  
(Notice: The water tank is detachable and can be disassembled for filling water. Do not start working during such operation)
3. Install / Replace the jet tip: Insert the jet tip into the jet seat in the middle of the Oral Irrigator's top head. To replace or remove the jet tip, press the jet tip lock/unlock Button.
4. Select the working mode: You can press the mode button on the handle to select the appropriate cleaning mode for the first time usage. You can also choose the cleaning mode according to professional dentists' guidance. (Normal mode: standard cleaning mode, it can effectively remove the food residue in the gap of the teeth, deeply clean teeth, keep mouth fresh; Soft mode: gentle massage to gums; Pulse mode: pulse massage to gums; promote blood circulation; DIY mode: users can customize the cleaning intensity according to personal preference, and set various cleaning strengths within the range.)

### Charging Precautions

1. Charge the product for 8 hours before the first time usage. Please use 5V 1A charger to charge the product.
2. Remove silicon cover from the charging seat, insert the USB pin into the charging socket of the oral irrigator (make sure there is no residual liquid in the socket before operating)
3. Plug another side USB connector into the charger and then plug the charger into a power socket.
4. When charging, the charging indicator lights is red and slowly flashing. When fully charged, the charging indicator will remain green till the charger is removed.

## Recommended Skills

- ✓ Keep unit upright with the jet tip towards the teeth  
Keep the product perpendicular to the mouth for best effects when using
- ✓ Slightly open mouth so that the inlet water can flow out.  
Make waterflow spray to teeth and gums at a certain angle. Move slowly the water flosser along the teeth

## Cleaning Guide

Never immerse the device into any liquid.

### Water Tank

- It can be washed with soapy water or placed on the top dish rack of the dishwasher. Unplug the power cord before cleaning.
- Always drain and dry the water in the tank before exposing it to a freezing temperature environment.
- Wash off the dirt with a brush or similar item, and then wipe it off with a cloth.
- Dry the water tank when not in use for a long time.

### Nozzle:

- Press the jet tip lock/unlock Button to remove the jet tip. Clear it with water and wipe with a cloth.
- It is recommended to replace to the new jet tip every 3 months. Otherwise, the jet tip will become too large to effectively clean the teeth.

### Main Body:

Wipe the dirt with a cloth dampened with a mild detergent, then wipe with a dry cloth. If it is very dirty, it can be washed with water.

- Do not use alcohol, benzene, thinner, etc. to wipe the product. Do not wash with hot water above 50 °C.
- Wipe off the water after cleaning.

### Suction Pipe

- Wipe with a cloth after cleaning.
- Do not bend the suction pipe. Otherwise, the suction pipe may work in failure such as stopping water absorption during working status.
- Clean the filter net of the suction pipe in time to prevent clogging.

### Water Tank Cover

- Wipe off the dirt, otherwise it will cause water leakage.

## Troubleshooting Guide

Problems	Reason	Method of handling
Nozzle cracking	Severe nozzle damage	Replace the new jet tip in the same model
Water tank leaks water	The seal is damaged	Replace the seal and return to the supplier for repair
	Water or other liquid is left in the water tank	Make sure to empty the water tank after using every time.
	The water tank cover is not closed	Close the water tank cover
Leakage between jet tip and its sockets	The jet tip is not properly installed	Remove the jet tip and reinstall
No liquid flow or insufficient pressure	Equipment is inverted	The nozzle must be higher than the water tank
	Insufficient water or liquid in the tank	Filling the tank with liquid
	Wrong working mode	Select the "Normal" mode



1. This warranty card is valid only in the People's Republic of China (except Hong Kong, Macao and Taiwan).
2. Within one year since the purchase date, free repairing warranty due to manufacturing process, components default or under the condition that our authorized maintenance station personnel confirm that it is damaged in the normal usage.
3. This free warranty does not include brittle parts (such as jet tips).
4. Please note the following situations, even within the warranty period, you can not get free warranty service:
  - 1) In case of failure to present this warranty card and purchase receipt.
  - 2) In case the device is repaired and damaged by personally that has not been authorized by our company's authorized service station.
  - 3) Other damage caused by force majeure factors (such as natural disasters, voltage abnormalities, etc.)
  - 4) In the case of failures caused from usage of outside the family (such as industrial or commercial use) or exposing the device into extreme hot or cold environments.
5. Please fill the following warranty form carefully and send it to the company's designated after-sales service point for repairing along with the device to be repaired.

### Maintenance and repairing

Customer Name	Contact Number
Product Model	Product Color
Purchase Date	Purchase Platform/Channel
Contact Address	
Maintenance Date	