Table des matières

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Table of contents

1. Before you start	134
1.1. The Professional depilation solution for use at Home	134
1.2. Description of the Homelight appliance	135
1.3. Package contents	135
1.4. The replaceable flash lamp	135
2. Intended use of your Homelight appliance	136
3. Using the Homelight epilator safely	136
4. Contra-indications	137
5. When should I not use the Homelight appliance?	138
6. Precautions – How to use your Homelight appliance safely	139
7. Reduce the risk of injury	141
8. Possible side effects with the Homelight epilator	143
9. Long-lasting hair removal with the Homelight epilator	145
9.1. How does the light eliminated the hair?	146
9.2. What effect does the growth cycle of the hair have on photo-epilation?	146
9.3. Optimising the results by planning your Homelight depilation sessions	146
10. Taking the first steps with your Homelight epilator	148
10.1. Installing your Homelight appliance	148
10.2. The first session with your Homelight epilator	148
10.3. Epilation of the face with the Homelight epilator	150
10.4. What to expect from the Homelight epilator	154
10.5. After a treatment with your Homelight epilator	154
11. Maintenance of your Homelight epilator	155
11.1. Cleaning your Homelight epilator	155
11.2. Changing the replaceable flash lamp	155
12. Trouble shooting	156
12.1. My Homelight epilator won't start	156
12.2. There are no pulses when I press the release mechanism	156
13. Customer service	157
14. FAQ	157
15. Specifications	160
16. Packaging - Labelling	160
17. Quick start-up	161

The BaByliss Homelight epilator is a powerful electrical appliance and so for reasons of safety it should be used with great care.

Before using the Homelight epilator for the first time, please carefully read this user manual completely.

Keep the device away from your eyes and do not look at the area being treated when the appliance flashes!

Do not use the Homelight appliance on tanned skin or after recent exposure to the sun! Using it under these conditions may cause serious burns or severely damage your skin.

For epilation of the face, please use only the precision head especially designed for the face and sensitive areas.

Do not use the Homelight on naturally dark skin!

Before using you Homelight appliance, please read thoroughly and follow precisely all safety warnings and directions given by BaByliss.

Note: The Homelight appliance is not as effective on white or light-coloured hair.

For best results with the Homelight depilator, plan your sessions well in advance.

Below is an example of a depilation schedule over a full growth cycle:

- The first 3 or 4 sessions with the Homelight depilator should be 2 weeks apart.
- Sessions 5 to 7 should be 4 weeks apart.
- After that, you need only use the Homelight epilator when hair regrows until epilation is lasting. Leave at least 2 weeks between sessions.

You can plan your personal sessions to suit your rhythm and still achieve satisfactory results. However, we recommend you follow the schedule suggested here because the results will be excellent.

Note: Do not treat the same patch more than once in two weeks as this will not improve results and it will risk causing undesirable effects.

You can plan maintenance depilation sessions with your Homelight depilator once you have achieved your objectives.

Because of hormonal factors and physiological changes, dormant follicles can become active. Maintenance sessions may thus sometimes be necessary.

Before using the Homelight epilator for the first time, please read this user manual fully. Pay particular attention to the paragraphs on the use and the operation of the unit. We recommend that you reread this user manual before every use of your Homelight epilator.

The BaByliss Homelight epilator is a powerful electrical appliance and so for reasons of safety it should be used with care.

Please read all notices and safety precautions before using and make sure you apply them strictly when using the Homelight epilator.

1. Before starting

1.1. The professional depilation solution for use at home

For more than 25 years, BaByliss has placed itself at the leading edge of technology to better respond to the increasingly demanding and pointed desires of women in depilation. Thanks to the pulsed light depilation technique, BaByliss's Homelight (Home Pulsed Light) epilator lets you gently and painlessly remove unwanted hair for silky smooth skin that lasts.

This exclusive depilation procedure produces a light that, pulsed then absorbed by the melanin (the hair's pigment), is carried to the bulb and transformed into heat. This thermal effect painlessly destroys the hair and roots and prevents regrowth. This exclusive, patented technology offers a long-lasting, scientifically-proven solution that provides fast, long-lasting results.

The BaByliss Homelight epilator guarantees a painless treatment for all areas of your body (do not use on your face) and has a unique safety mechanism based on a skin colour sensor.

The BaByliss Homelight epilator offers long-lasting depilation worthy of a professional in the privacy and comfort of a treatment at home.

1.2. Description of the Homelight appliance

The Homelight epilator is a lasting photo-epilation device that can be used at home.

Your Homelight unit includes an epilator which has a control panel and a replaceable flash lamp. The release mechanism is on the appliance handle and the replaceable flash lamp is located in the end of the appliance.

1.3. Package contents

- Homelight epilator and its flash lamp with a capacity of 5000 impulses
- Safety glasses
- Mains adapter
- Precision head with a capacity of 1000 pulses for epilation of the face and sensitive areas
- Instruction manual
- Accompanying DVD.

1.4. The replaceable flash lamp

The flash lamp of the Homelight epilator has a capacity of 5000 pulses (irrespective of light intensity level). It is estimated that on a person of average size the flash lamp's 5000 pulses are sufficient for a complete and durable treatment (11 sessions) of the lower legs, the bikini line and the underarms.

2. Intended use of your Homelight appliance

The Homelight epilator is designed to eliminate unwanted hair.

The Homelight epilator is an epilatory appliance that serves to eliminate undesirable body hair (in particular, the underarms, the bikini line, the legs and arms) and facial hair (cheeks, lips, chin) thanks to its precision flash lamp, designed especially for the face and sensitive areas.

The Homelight epilator can be used by men and women aged 18 and older.

3. Using the Homelight epilator safely

Safety is a priority with the Homelight epilator

- **Home Pulsed Light (HPL) technology provides maximum safety and minimum risk**
Depilation using the Home Pulsed Light technology lets you achieve long-lasting results with less power intensity compared to other photo-epilation methods. The epilator's low energy use reduces the chances of injury or complications and contributes to your overall safety.
- **The Homelight epilator has a protective mechanism.**
The BaByliss Homelight depilation appliance should not be used by people whose skin is naturally dark or tanned. The appliance has a safety mechanism, a sensor that checks the pigmentation of the skin being treated and detects skin types with which it is compatible. This exclusive safety procedure will not allow you to treat areas where the skin is too dark or too tanned. Moreover, the optical lens through which the energy pulse is emitted and which has a surface of 4 cm² is fitted into the flash lamp. It is this that ensures optimal protection by preventing direct contact between the skin and the optical lens.
- **The Homelight epilator has a mechanism to protect your eyes**
The Homelight epilator is designed with a built-in mechanism that will protect your eyes. This was designed so that the pulse cannot be emitted if the epilator is aimed at open space. The safety mechanism only allows the appliance to work when the application surface is touching the skin. Do not look at the area being treated when the unit flashes.

It is nevertheless strongly advised that you wear the safety glasses included during treatment sessions with the Homelight epilator.

4. Contra-indications

Important information on safety – Read before using

Using the Homelight epilator cannot be used by everyone. Before using, read the following guidelines. For more information or personal advice, please visit the website www.babyliss.com or the website for the distributor in your region.

Finally, if you are unsure about using the Homelight epilator, please consult your dermatologist or your GP.

- DO NOT USE the Homelight epilator on naturally dark skin.

The Homelight epilator will remove unwanted hair by selectively attacking the pigments in the hair shaft.

There is a wide range of pigments in the surrounding skin tissue. The amount of pigment in a person's skin determines not only the colour, but also the level of risk to which the user is exposed when using a form of depilation that uses light, regardless of which one.

The treatment of dark skin with the Homelight epilator may lead to side effects, like burns, blisters and changes in skin colour (hyper or hypopigmentation).

The Homelight epilator was designed with a sensor that will detect the colour of the skin to be treated at the beginning of each session and from time to time during use. This sensor checks that the energy is only administered to the appropriate skin types.

DO NOT USE the Homelight appliance on naturally dark skin! It will not work!

- DO NOT USE the Homelight appliance on tanned skin or after recent exposure to the sun!

Using the Homelight appliance on skin recently exposed to the sun may have side effects, like burns, blisters and changes in complexion (hyper or hypopigmentation).

DO NOT USE the Homelight appliance on tanned skin or after recent exposure to the sun! Such exposure may cause serious burns or severe damage to your skin. Avoid exposing yourself to sunlight for 4 weeks before treatment with the Homelight epilator, and for 2 weeks following treatment.

- **Keep away from your eyes!**

Pointing the BaByliss Homelight appliance at your eyes when activated can cause serious eye damage.

Keep away from your eyes!
Use the safety glasses included when using the Homelight epilator.

- **NOTE:** The Homelight appliance is not effective on fair hair.

The Homelight appliance is not effective on white, grey or blond hair. If your hair is one of these colours, the Homelight appliance may not work.

5. When should I not use the Homelight appliance?

Using the Homelight epilator is temporarily contra-indicated under certain conditions. **DO NOT USE** the Homelight device if:

- You are pregnant or breast feeding.
- You have had lengthy exposure to the sun or UV equipment over the previous 28 days.
- You have tattoos or permanent make-up on the area to be treated.
- There are abnormalities on the skin in the area to be treated: brown or black marks, large freckles, birthmarks, warts or nevus pigmentosus.
- You suffer from eczema, psoriasis, lesions, open wounds or evolutive infection in the area to be treated. Wait until the affected area has healed before using the Homelight appliance.
- You show antecedents to the genesis of keloids, you have a known sensitivity to light (photosensitivity) or one caused by medication, including non-steroidal anti-inflammatories (e.g.: aspirin, ibuprofen, paracetamol), tetracycline, phenothiazine, thiazide, diuretics, sulphonylureas, sulfamides, DTIC, fluorouracil, vinblastine, griseofulvin, Alpha Hydroxy Acids (AHAs), Beta-Hydroxy Acids (BHAs), Retin-A®, Accutane® or topical retinoids.

- You have recently been treated or are being treated with Alpha-Hydroxy Acids (AHA), Beta-Hydroxy Acids (BHA), Retin-A®, topical retinoids or azelaic acid.
- You suffer from an abnormal skin condition caused by diabetes, for example, or another systemic or metabolic disease.
- You have been undergoing treatment with Accutane® (isotretinoin) for the previous 6 months.
- You have been following a steroid regime for the previous 3 months.
- You show signs of pre-herpes growth in the area to be treated, unless you have consulted a physician and been given preventive treatment before using the Homelight epilator.
- You suffer from epilepsy.
- You have an active implant, like a pacemaker; an implant for urinary incontinence, an insulin pump, etc.
- You have a pathology linked to photosensitivity, like porphyria, polymorphic eruption to light, solar urticaria, lupus, etc.
- You have an antecedent to skin cancer or areas subject to dermal malignancies.
- You have been undergoing radio or chemotherapy for the last 3 months.
- You suffer from any other condition your physician deems a risk for the treatment.
- **If you have any doubts about the risks to you of the use of the Homelight appliance, consult your GP or your dermatologist.**

6. Precautions – How to use your Homelight appliance safely

Test the appliance on a small patch of skin on your leg or arm (a few flashes) at the lowest power intensity to ensure that you have no adverse reaction to the treatment in the following days.

■ Choose the intensity level with CAUTION!

The intensity level determines the power of the pulsed light that is projected onto your skin during the treatment, ranging from the lowest (-) to the highest (+). The power intensity indicator lights on the control panel indicate the level of intensity selected. By increasing the power, you will obtain better results, but you also increase the risk of side effects (see 'Possible Side Effects' below).

The first time you use your appliance, always use the lowest setting (one indicator light at '-').

Only increase the power in your next session if you have not felt no, or only slight, discomfort during your first session or afterwards while your appliance is set to the lowest setting. In the next session, increase the power, one notch at a time, to a higher indicator light, and continue in this way in the following sessions.

For detailed instructions on the selection of the power, see the box 'Power Intensity Level' in the section 'The First Session with the Homelight epilator'.

■ **Avoid side effects!**

Do not treat the same patch of skin more than once per session!

Avoid using a pulse on the same patch!

If burns or blisters appear, STOP IMMEDIATELY!

■ **Avoid complications after using the Homelight epilator!**

Avoid exposure to sunlight.

"Exposure to sunlight" means continuous exposure to direct sunlight without protection for about fifteen minutes or continuous exposure to indirect sunlight without protection for an hour.

Protect your treated skin with sun screen that has a protection factor of 30 or more throughout the treatment period and for at least two weeks after the last session of depilation.

Always shave the area to be treated and make sure you have clean skin before using the Homelight appliance.

Make sure birthmarks and tattoos are covered during the Homelight epilator treatment and also cover brown or black marks, as well as freckles, nevus, warts or growths before using the Homelight epilator.

- Never look directly into the light emitted by the replaceable flash lamp. Use the safety glasses included whenever you use the Homelight epilator.
- Never use the Homelight epilator on the tip of the breasts or the genital organs (male or female).
- Never use the Homelight epilator on any part of the body where you might later want hair.
- Never use the Homelight epilator for anything but depilation.
- Never point the application surface of the Homelight epilator to emit a light pulse into open space. Always ensure you point the treatment surface at the skin and that it is in full contact with your skin before releasing a pulse.
- Never use flammable liquids, like alcohol (including perfume, deodorant or any other preparation containing alcohol) or acetone, to clean you skin before using the Homelight epilator.
- Using your Homelight appliance may cause a temporary change in pigmentation (see below 'Possible Side Effects').
- Keep the appliance out of the reach of children. Do not use the Homelight epilator on children nor allow children to use it.
- This appliance is not intended for use by people (adults and children) whose physical, sensory or mental capacities are reduced, or by persons who lack experience or knowledge, unless they are supervised or have been given instructions beforehand on the use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not use the appliance like a toy.

-  CAUTION!: The plastic bags used to wrap this appliance or the packaging may be dangerous. Keep these bags out of the reach of babies and children. Do not use them in cribs, children's beds, pushchairs or playpens. The thin film will stick to the nose and mouth and prevent breathing. A bag is not a toy.
 -  If using the appliance in a bathroom, make sure you unplug it after use. In fact, having it close to a water source may be dangerous even if the appliance is turned off. To ensure additional protection, we recommend installing a residual current device (RCD) with a nominal trip current of no higher than 30 mA in the electrical circuit with that supplies power to the bathroom. Ask your electrical fitter for advice.
 -  ELECTRICAL AND ELECTRONIC EQUIPMENT AT THE END OF LIFE.
In everybody's interest and to participate actively in protecting the environment:
Do not discard these products with your household waste.
Use return and collection systems available in your country.
Some materials can in this way be recycled or recovered.
- Do not use the unit if it has fallen or if it appears to be damaged.
 - Do not use any accessories other than those recommended by BaByliss.

7. Reduce the risk of injury

As with all electrical appliances, some precautions are necessary to ensure your safety when you use the Homelight epilator:

- Keep the Homelight epilator away from all sources of water!
- The Homelight epilator is an electrical appliance and should never come in contact with water.
- Do not store the Homelight epilator in a place where it could drop or fall into a bath or washbasin or any other container with liquid in it. Do not put it in water or any other liquid. This could result in severe electrical shock.

- Do not use your Homelight epilator when having a bath.
- Never use the Homelight appliance if it is damp or wet.
- Do not catch hold of the Homelight epilator if it has fallen in water.
- Unplug the Homelight appliance immediately if it has fallen in water.

Keep the Homelight epilator away from all sources of water!

- Never open your Homelight appliance!
- Do not attempt to open or to repair your Homelight appliance yourself. Opening the appliance will expose you to dangerous electrical components and to the energy of the pulsed light, which could cause serious physical injury and/or irreversible eye damage.

Do not attempt to open or repair your Homelight epilator yourself.
Only repair centres certified by BaByliss are authorised to repair any damage.

- Attempting to open the Homelight epilator could also damage the appliance and will void your warranty.
- Please contact the customer service of your Homelight epilator if your appliance is broken, damaged or needs repair;
- Only use the Homelight epilator for what it was designed to do and following the instructions of this manual.
- You appliance should never be left unsupervised while plugged in.
- Never leave your Homelight epilator unsupervised while turned on.
- Do not plug your Homelight epilator in if the cord or the plug is damaged. Keep the cord away from all heat sources.
- Do not plug your Homelight epilator in if the adapter lead or the plug is damaged. Keep the adapter away for all heat sources.
- Stop using your appliance if you notice a defect or if it appears defective or emits an abnormal noise or odour. In particular; stop using the appliance immediately if the flash lamp or the skin colour sensor appears cracked or is marked.

- Always disconnect the Homelight appliance adapter from the mains outlet immediately after use.
- Unplug your Homelight appliance before you clean it.
- Use only attachments and parts with your Homelight appliance that have been recommended by BaByliss.

8. Possible side effects with the Homelight epilator

When your Homelight epilator is used according to these instructions, side effects and complications are rare. Nonetheless, all cosmetic treatments, even those designed for use at home, involve some risks. Therefore, it is important to understand and accept these risks and possible complications when you use a pulsed light depilation system designed for use at home.

If you superimpose pulses or administer multiple pulses on the same patch:

Side Effects	Appearance of the Side Effect
Slight discomfort / annoying sensation in the skin	Light
Redness	Light
Hypersensitivity of the skin	Light
Burns or sores	Rare
Scarring	Rare
Change in pigmentation	Rare
Excessive redness of the skin and swelling	Rare
Infection	Negligible
Bruising	Negligible

■ Slight annoying sensation in the skin

Even if the home pulsed light depilation is generally tolerated well, most people who use it feel a slight annoying sensation during use, often described as a prickling sensation on the area being treated. Generally, this sensation lasts throughout the application, even for a few minutes afterwards. Any discomfort beyond this annoying sensation is abnormal and means that either you must stop the treatment, because you cannot tolerate it, or you must lower the selected power level.

- **Redness**

Your skin could become red just after use of the Homelight epilator or in the 24 hours following treatment. The redness generally disappears within 24 hours. Consult your physician if it persists longer than 2 to 3 days.

- **Hypersensitivity of the skin**

The skin is more sensitive in the area that has been treated and could dry out and scale.

- **Burns or sores**

Very rarely, the skin may be exposed to injuries or burns following the application. Several weeks may be necessary for it to heal completely and, in exceptional cases, a scar may remain.

- **Scarring**

Even if these cases are very rare, permanent scars may appear. Generally, the scars take the form of a whitish superficial lesion (hypotrophic). Nonetheless, the scar may be large and red (hypertrophic) or large and extended beyond the injured area (keloid). More advanced cosmetic treatment may be required to improve the appearance.

- **Change in pigmentation**

The Homelight appliance targets the hair roots, in particular the pigmented cells of the bulb. Nonetheless, there is a risk of temporary hyperpigmentation (increase in the brown pigmentation or discolouration) or hypopigmentation (whitening) of the surrounding skin. The risk of change in pigmentation is higher among people with darker skin.

Generally, the change in pigmentation and the discolouration are temporary. Cases of permanent hyperpigmentation or hypopigmentation are rare.

- **Excessive redness and swelling**

In some rare cases, treated skin may become very red and swollen. This reaction is more common on sensitive areas of the body. The redness and the swelling should disappear within 2 to 7 days during which you should apply ice to the affected area regularly. Gentle cleaning does not pose a problem, but you should avoid any exposure to the sun.

- **Infection**


A skin infection is extremely rare, but may occur after a burn or a lesion of the skin caused by the Homelight epilator.

- **Bruising**

Very rarely, using the Homelight epilator may cause some blue or purple bruising which could last 5 to 10 days. As the blue disappears, the skin may take on a rusty colour (hyperpigmentation) that could be permanent.

9. Long-lasting hair removal method with the Homelight epilator


The Homelight epilator is a personal long-lasting photo-epilation device. The laser or photo-epilation procedure is well known and has been mastered. It has been subjected to over 15 years of clinical testing around the world for more than 15 years and has proven its effectiveness and its safety in offering gradual but permanent hair removal.

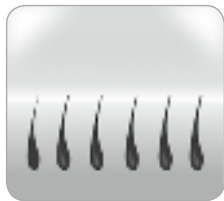
Before application
of the Homelight epilator

During the pulse emitted
by the Homelight epilator

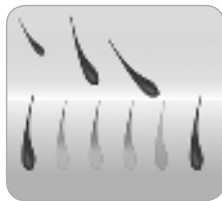
After application
of the Homelight epilator

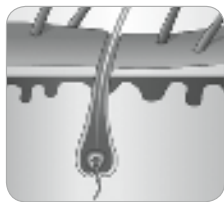Growth phase of the hair
Anagen

Transition phase
Catagen

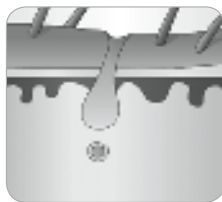Resting/dormant phase
Telogen

9.1. How does the light eliminate hair?


Photo-epilation is based on the theory of selective photothermolysis which uses optical energy to stop the growth of the hair. To obtain such thermal effects, the root of the hair must selectively absorb the photoenergy and transform it into heat. This selection is produced when most of the optical energy emitted is absorbed by the pigment of the bulb, while the surrounding skin and tissue remain cold.

Melanin is the pigment found in the hair bulb, responsible for absorbing light, generating the heat that stops the hair growth. Also, the more melanin in the hair (i.e., the darker the hair), the more light it can absorb and the more efficient photo-epilation becomes.

9.2. What impact does the growth cycle of the hair have on photo-epilation?

Every hair on our body has three growth phases: the anagen, catagen and telogen phases. These three phases play an important role in the process of photo-epilation. The anagen phase is the growth phase of the hair, while the catagen and telogen phases are both resting phases. The time the hair takes to go through its follicular cycle varies completely from one person to the next, and from one part of the body to the next. We can, however, establish an average of between 18 and 24 months.

Most hair follicles are found at any one time in one of the two resting phases. Consequently, the Homelight epilator will have no effect when they are in these resting phases. Nonetheless, it will affect the hair that is in anagen phase. It is therefore important to understand the duration of a full follicular cycle to obtain lasting results with your Homelight device.

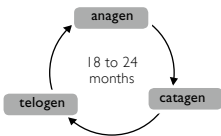9.3. Optimising the results by planning your Homelight depilation sessions

A full follicular cycle takes on average from 18 to 24 months during which many sessions with your Homelight epilator will be necessary to obtain long-lasting results. The efficiency of the appliance will vary from one person to the next and will depend on the area being treated, the colour of the hair and the way in which the appliance is used.

Schedule of your Homelight depilation sessions during a full follicular cycle:

- The first 3 or 4 sessions with the Homelight epilator should be spaced 2 weeks apart
- Sessions 5 to 7 should be spaced 4 weeks apart.
- After that, you need only use the Homelight epilator when hair regrows until epilation is long-lasting. Leave at least 2 weeks between sessions.

You can schedule your personal sessions according to your rhythm and nevertheless obtain satisfying results. However, we recommend following this schedule as it has had excellent results. Please note: treating the same zone more than once in two weeks will not improve results but it can cause undesirable effects.

You can plan your maintenance depilation sessions with your Homelight epilator after you have achieved your objectives:

Because of hormonal factors and physiological changes, the follicles at rest may become active. Occasional maintenance sessions may thus be necessary.

Expected depilation results

Skin types according to Fitzpatrick	Hair colour	Optical fluence [J / cm ²]	Surface area being treated	Average number of sessions	Efficiency (% depilation/regrowth)
I - IV	Medium brown to brown	3 - 4.5	Leg or thigh	10 to 12	60% reduction of hair
			Arm		50% reduction of hair
			Bikini line		60% reduction of hair
			Underarms		55% reduction of hair
I - IV	Brown to black	3 - 4.5	Leg or thigh	8 to 10	70% reduction of hair
			Arm		60% reduction of hair
			Bikini line		70% reduction of hair
			Underarms		65% reduction of hair

The result will vary from one person to the next and according to the kind of hair and biological factors that determine the growth cycle of the hair. Some people achieve results more or less quickly than the average number of sessions that we recommend.

Long-lasting results cannot be expected after just one or two sessions. Here again, the duration of the resting cycle of the hair follicles varies from one area of the body to the next.

10. Taking the first steps with your Homelight epilator

10.1. Installing your Homelight appliance

1. Remove the Homelight epilator appliance and all parts from the box.
2. Make sure the flash lamp is inserted firmly in the end of the epilator.
3. Plug the adapter lead into its slot in the appliance handle.
4. Connect the other end of the adapter into the mains. The 'Standby' indicator light will come on. Your Homelight epilator is now ready for use.

Wear the safety glasses included during treatment with the Homelight epilator.

10.2. The first session with your Homelight epilator

A first trial session on an arm or leg should be carried out at minimum intensity to check that your skin reacts correctly to the treatment over the following days.

You may then gradually increase the power over the course of the following sessions.

Use the lowest power level for the first 50 pulses and then increase the power gradually: from 0 to 50: level 1, from 50 to 100: level 2, from 100 to 150: level 3, from 150 to 200: level 4, ...

The power intensity level may be adjusted by pressing the + and – buttons.

Your skin must be shaved, clean, dry and free of powder or deodorant residue.

Do not epilate (either with an electric epilator or with wax) before using the Homelight epilator. Epilation (removing hair at the roots) will render the pulsed light treatment ineffective.

5. Press the On/Standby switch. The power level indicator light will come on and you will hear the sound of a fan (similar to the sound of a hair dryer).
6. About 4 seconds after pressing the switch, the 'Ready' indicator light will come on. The default level for the appliance is set at the lowest setting. It is therefore ready to release its first pulse.
7. If you would like to increase or reduce the power, press on the power intensity regulator using the '+' and '-' until the desired power intensity indicator light comes on. The power intensity indicator light that is lit corresponds to the power level administered (see details on the 'Power intensity level').

8. Press the application surface firmly against your skin so that it is in full contact. When the application surface is in full contact with the skin, the 'Ready' indicator light will begin to flash.
9. Press the release mechanism. The device will first check your skin colour:

If the skin colour is light enough for safe application, the appliance will release a pulse onto the skin and the indicator light will turn off.

You will see a flash and at the same time you will hear a light 'beep' every time it pulses, which is normal. You might feel a slight sensation of heat and prickling.

If you do not see a light pulse and the skin colour sensor indicator light begins to flash, this means that your skin is too dark for safe use. Change the patch of skin being treated or check with the Help Service for the Homelight epilator (see for more details: 'Skin colour sensor').

Your Homelight unit will begin recharging immediately after each pulse. After 4 seconds, the 'Ready' indicator light will come on again.

10. Remove the application surface from the area of treated skin.
11. Move the application surface to another area to be treated. Pressure tracks from the application surface made on your skin will help you to align the appliance correctly to administer another pulse, preventing superimposing the pulses or spacing too much. (See 'Covering the treatment area').

Warning: Never treat the same patch of skin more than once per depilation session! Repeated treatment of the same area increases the risk of side effects.

12. Repeat the application from step 5.

10.3. Epilation of the face with the Homelight epilator

As facial skin is particularly delicate and sensitive, its epilation requires very special attention.

That is why BaByliss has fitted the Homelight epilator with a precision head designed especially for the face and sensitive areas. Targeted and ultraprecise, its 2cm² window of application is designed to facilitate treatment of the face. It produces a light intensity adapted to the sensitivity of facial skin and the skin of sensitive areas. It is also equipped with a contact detector to protect the eyes such that no impulses are emitted if the device is pointing into the void.

It is imperative not to look directly at the area treated during the flash and to wear protective eye glasses when using the Homelight epilator to treat the face.

The special precision head is suited for treatment of the cheeks, upper lip and chin. Do not use for epilation of eye lashes, head hair, and eyebrows.

This product is not appropriate for full removal of male beard.

For epilation of the face with Homelight epilator a compatibility test must first be done. Your skin colour must first be tested for its compatibility with pulsed light - by flashing your arm, for example (or any part of your body, the most tanned or darkest possible: legs, back) with the special body head fitted with the skin colour detector. If an impulse is emitted, you may place the special face head on the tip end of the Homelight device and proceed to treatment of the face. Repeat this operation before each new treatment with the Homelight epilator.

A word of advice: If you have planned a complete face and body session with Homelight epilator, begin treatment with the body and end the session with the face treatment. That way you will already have had a chance to test the compatibility of your skin with pulsed light treatment using the special body head fitted with the skin colour detector.

For complete safety, use of Homelight on the face requires the aid of another person.

Your facial skin should be shaven, clean and dry with no remnants of a cream or lotion.

Please do not try to remove hair (whether with, for example, an electric epilator or with wax) before using the Homelight epilator. Epilation (removing hair by its roots) renders pulsed light treatment ineffective.

1. Apply the special precision head designed for the face and sensitive areas:
 - Press the On/Standby switch of your Homelight epilator to put it on standby.
 - Disconnect the adapter from the mains socket.
 - Grip the special body flash on both sides and pull it to remove it from its case.
 - Insert the special flash lamp for the face and sensitive areas, urging it gently into place.
2. Press the On/Stand by switch. The luminous power indicator will light up and you will hear the sound of a fan (like the noise of a hairdryer).
3. About 4 seconds after having pressed the switch the light indicator will show 'Ready'. The device automatically defaults to the minimum power setting and is now ready to emit the first pulse.
4. For the face it is strongly advised to set the device to the weakest light intensity: level 1.
5. Please wear protective eye glasses.
6. Ask someone to help you to operate the device by firmly pressing the application surface against your skin to achieve total contact. As soon as the surface is in full contact with the skin, the message 'Ready' on the light indicator will begin to blink.
7. Press the trigger. The device will now emit a pulse to the skin and the message 'Ready' will fade out.

You will see a flash of light and at the same time you will hear a faint « beep » with each pulse. This is normal. You may also feel a slight sensation of heat and a tingling. Your Homelight epilator immediately recharges itself between each pulse emission. After 4 seconds, the message 'Ready' lights up again on the indicator.

8. Move the application surface away from the area of the skin just treated.
9. Move the application surface to another area to be treated.

Avoid releasing more than one pulse on the same patch of skin!

Do not treat the same area more than once in one session!

Power level intensity

The power level intensity determines the power of the pulsed light administered into the skin during the treatment, ranging from the lowest (-) to the highest (+). The corresponding power level indicator lights on the control panel indicate the selected power level. By increasing the power, you will achieve better results, but you will also increase the risk of side effects (cf. 'Possible side effects' below).

Always use the lowest power setting (one indicator light at '-') the first time you use your Homelight epilator:

If during or after your first depilation session, you feel little or no discomfort, you can increase the power by one increment in your next session. You will see the corresponding indicator light light up and you can repeat this step in the following session if you wish.

Do not increase the power level if you felt more discomfort than usual while using or after having used the appliance (see 'Possible side effects').

Do not increase the power intensity during a session even if you have not felt any discomfort. It is possible that you won't feel discomfort during the treatment, but will do after.

Note: After having changed the replaceable flash lamp, always turn the power level one notch lower than usual. Return to your usual power level only if you experience no adverse effects during or after your treatment.

When your appliance is on, its intensity level will automatically be positioned by default at the lowest intensity level. A single power level indicator light will be lit. To adjust the intensity light to the desired level, either higher or lower, press the power intensity regulator using the '+' or '-' until the indicator light corresponding to the intensity that you want is lit.

The skin colour sensor

Depilation based on pulsed light is not appropriate for naturally dark skin and there is a risk of side effects, in particular burns, phlyctena and discolouration of the skin (hyper or hypopigmentation).

To avoid inappropriate use, the appliance comes equipped with an exclusive safety device, a sensor; that checks the skin pigmentation at the beginning of each session and from time to time throughout the session. If the sensor detects that the skin is too dark, the appliance will not be able to emit light pulses. Consequently, if you do not notice any pulsed light and the skin colour sensor indicator light starts to flash, this indicates that the skin colour measured by the sensor is too dark for safe application. Change the treatment area or contact a resource person at the BaByliss Service Centre.

Cover the area being treated

Pulses must be administered in rows, from one end to the other and proceeding systematically in successive rows. This method ensures better control of the pulses and will help you to cover the area in question better while preventing treating a same patch twice or administering pulses too close to each other in the same area.

Avoid releasing more than one pulse on the same patch of skin!

Do not treat the same area more than once in one session!

STOP the treatment IMMEDIATELY if your skin has burns or phlyctena appear!

10.4. What to expect with the Homelight epilator

Using the Homelight epilator at home could for many people be a new experience in a depilation method that uses light. The Homelight epilator is easy to handle and the depilation sessions are quick. You can expect to see, hear or feel the following during a session:

- **Fan noise** – The cooling fan of the Homelight epilator makes a noise similar to a hair dryer. This is normal.
- **Beep for each pulse** – When activating a pulse from the Homelight epilator, the unit will make a slight sound whenever the light flashes.
- **Moderate pressure on the application surface** – necessary and useful to position it correctly over the connecting area; this is an exclusive safety feature of the Homelight epilator.
- **Heat and prickling sensations** – Whenever you release a pulse, it's normal to feel some heat and prickling resulting from the energy of the light.
Do not forget that it is important to use the lowest power settings during the first sessions. The heat sensation may persist for up to an hour after your depilation session.
- **Skin slightly red or pink** – Just after use or within 24 hours of using it, it is not rare to see the skin turn slightly red or pink, especially around the hair.
If obvious redness, blisters or burns appear, stop using the Homelight epilator immediately.

10.5. After a treatment with your Homelight epilator

- At the end of your session with your Homelight epilator, turn your appliance off using the On/Stdby switch. (Make sure you remember the last level of power intensity you used, because it will not be set to that level when you turn your appliance on the next session).
- Unplug the adapter from the mains.
- We recommend cleaning your Homelight epilator after each session, particularly the application surface (see 'Cleaning your Homelight epilator').
- After cleaning your appliance, we recommend storing it in its box and away from any source of water.

Care of your skin after you depilation session

Avoid exposing the treated areas to sun. Protect your skin with a sun screen throughout the treatment and for at least two weeks after the last session with your Homelight appliance.

Side effects and complications

Some patients may notice a change in pigmentation following a treatment with the Homelight epilator. These effects, if they appear, are generally slight and temporary. **If you experience complications (see 'Possible side effects') contact your GP immediately.**

11. Maintenance of your Homelight epilator

11.1. Cleaning your Homelight epilator

We recommend you clean your Homelight epilator after every session, particularly the application surface.

- Unplug your appliance before cleaning.
- Gently clean the surface of the appliance and particularly the application surface with a clean, dry cloth and a cleaning agent formulated for cleaning electronic equipment.

Never immerse your Homelight appliance or one of its parts in water.

11.2. Changing the replaceable flash lamp

The intensity of the pulses is determined by the selected power level. The energy will not weaken during the life of a replaceable flash lamp. When it has reached 90% of its life, the replace flash lamp indicator light on the control panel will flash indicating that the bulb should be replaced soon.

Once the 5000 pulses of the replaceable flash lamp have been used, the replace flash lamp indicator light will light up indicating that no more pulses will be released. To continue your depilation, you will have to replace the replaceable flash lamp.

1000-flash replacement cartridges are available commercially under the reference code G911E.

Note: The flash lamp should also be replaced if large marks appear on the inside or if the application surface is broken.

To change the replaceable flash lamp:

1. Press the On/Standby switch of your Homelight epilator to put it in standby mode.
2. Unplug the adapter from the mains.
3. Grasp both sides of the used flash lamp, pull to remove it from the socket and dispose of it in your household waste.
4. Unwrap a new replaceable flash lamp.
5. Insert the new flash lamp gently pushing it into place.

Note: After having replaced the flash lamp, always turn the power level one notch lower than usual. Return to your usual power level only if you experience no adverse effects during or after your treatment.

12. Trouble shooting

12.1. My Homelight epilator won't start

Ensure the adapter is properly connected to your Homelight epilator and plugged into the mains.

12.2. There are no pulses when I press the release mechanism

Ensure that the application surface is evenly and firmly in contact with your skin and that the 'Ready' indicator light is flashing; this indicator light guarantees that the appliance is in full contact with your skin. For your safety, the release switch will be able to release a pulse only if the application surface is firmly pressed against your skin.

- Check the control panel. If the skin colour indicator light flashes for 3 seconds, the appliance is indicating that the skin colour sensor has assessed your skin colour as too dark for the Homelight device to be used safely. Try using the appliance on another area.
- Check the replace flash lamp indicator light. If it is lit, unplug your appliance and change the flash lamp.
- Ensure that the 'Ready' indicator light is lit.
 - If after 10 seconds, the 'Ready' indicator light is still not lit, turn your Homelight epilator off and turn it on again by pressing the On/Standby switch twice.
 - If the problem persists, contact the BaByliss customer service in your area.
- Check the appliance's appliance in use indicator light.
 - If it is lit, turn your Homelight epilator off and turn it on again using the On/Standby switch.
 - If the problem persists, contact the BaByliss customer service in your area.

If these problems persist, contact the BaByliss customer service in your area.

Do not attempt to open or repair your Homelight epilator yourself. Only certified repair centres are authorised to make repairs.

Opening your Homelight epilator could expose you to dangerous electrical parts and the energy of the pulsed light, which could cause serious physical injury and/or irreversible eye injury.

Opening your Homelight epilator could also damage the appliance and will void the warranty. Please contact BaByliss customer service if your appliance is broken or damaged and requires repair.

13. Customer Service

For more information on Homelight products, please visit www.babyliss.com

If your device is broken, damaged, needs repair or if you need any assistance, please contact BaByliss customer service.

14. FAQ

1. Is the Homelight epilator truly effective?

Yes. Clinical tests carried out by physicians have proven that the BaByliss Homelight epilator achieves long-lasting results safely.

2. On what part of the body can I use the Homelight epilator?

The Homelight epilator by BaByliss has been designed for removal of body hair. The areas most commonly treated with the Homelight epilator are the legs, underarms, arms, and the bikini line. Homelight can likewise be used for the face, but only if the precision head designed especially for the face and sensitive areas is used.

3. How long is a depilation session with the Homelight epilator?

The duration of the treatment will depend on the area being treated. Two half legs may take 40 minutes while two underarms, 10 minutes. The BaByliss Homelight epilator operates on the mains and may also be used for as long as necessary to treat the desired area.

4. Is the Homelight epilator safe?

Your safety was our first priority when we designed the Homelight epilator. Very strict standards were followed when testing and approving it for use at home. But like any electronic appliance that comes in contact with the skin, the user must follow the instructions and scrupulously follow the safety precautions.

5. Are treatments using the Homelight epilator painful?

When used properly, most Homelight epilator users report a slight heat sensation when the unit pulses. Users with darker and thicker hair sense a little more discomfort. But these sensations ease and disappear at the end of the session. For your comfort, the Homelight epilator offers 5 levels of power intensity depending on your sensitivity.

6. At what frequency should I use the Homelight epilator?

The first 3 or 4 sessions with the Homelight epilator should be spaced 2 weeks apart. Sessions 5 to 7 should be spaced 4 weeks apart. The sessions may be held when your hair grows back until you have achieved the desired results.

7. Is the Homelight epilator effective on white, grey or blond hair?

The Homelight epilator is most effective on dark hair or on those containing more melanin. Melanin, the pigment that gives colour to skin and hair, absorbs optical energy. Black and dark brown hair respond best to the treatment. If brown and light brown hair also responds, these will need a few more sessions. Red hair may also respond partially to the treatment. Generally, white, grey or blond hair does not respond to the treatment, but some users have noticed results after several depilation sessions.

8. Can I use the Homelight epilator on brown or black skin?

Do not use the appliance on naturally dark skin! The Homelight epilator eliminates hair by targeting the follicular pigment. Pigments are also found in varying quantities in the surrounding dermal tissues. The quantity of pigments in a person's dermis, visible through the colour of the skin, defines the level of risk to which he or she is exposed by using the Homelight epilator. Treating dark skin with the Homelight epilator may involve risks like burns, blisters and changes in skin colour (hyper- or hypopigmentation). The Homelight epilator was designed with an exclusive safety device, a skin colour sensor, to check the colour of the skin at the beginning of each session and from time to time during the session itself. When the appliance detects a darker area, the sensor deactivates the Homelight epilator and will not release any more pulses.

9. When will I see the first results with the BaByliss Homelight epilator?

Like all other laser or photo-epilation methods, the results are not immediate; you may even have the impression that nothing is happening. The hair will sometimes appear to continue growing after a session, but after two weeks, most of this hair will simply fall out. This is why you will need to have a several sessions before achieving the expected result.

10. Can men use the Homelight epilator ?

The Homelight epilator can be used by men. However, men's hair, especially chest hair, will require more sessions to achieve the desired results.

11. Why does hair grow back in the area that I treated last week?

It is quite common for hair to appear to continue to grow for one to two weeks after a depilation session with the Homelight epilator. This process is known as 'ejection'. After two weeks, you will nonetheless notice that this hair falls out or comes away from its follicle. That said, we recommend you do not pull the hair from the follicle – let it fall out naturally. Moreover, some hair will not be affected by the Homelight epilator either because of poor application or because the hair was in its dormant phase. This hair will be treated in the following sessions, hence the need for several sessions to achieve the expected results with the Homelight epilator.

12. I've heard that some hair grows back lighter and finer after photo-epilation. Is this true?

This phenomenon is well documented by beauticians and doctors who use laser or pulsed light appliances for depilation treatments. It is possible that some hair grows back lighter and finer after an Homelight epilator session. Generally, this hair only represents a tiny part of the original hair growth, treatments should therefore be continued to achieve the desired result.

13. Is long-term use of the Homelight epilator dangerous for my skin?

Over the last 15 years, the use of laser and of the pulsed light in aesthetic medicine has been well documented in professional journals. These journals and institutions have not reported any side effects or injuries connected with long-term use of laser or pulsed light appliances.

14. Can I use the Homelight epilator to eliminate chin or facial hair?

Homelight epilator from BaByliss allows facial hair removal (cheeks, upper lip, and chin), but only if the precision head designed especially for the face is used. However, the Homelight epilator may not be used for epilation of eyelashes, eyebrows or head hair. Please wear a pair of protective glasses when using.

15. How long should I wait before using the Homelight epilator if I am not protected from the sun's rays?

You should avoid exposure to the sun before and after treatment with the Homelight epilator no matter which part of the body has been treated. You must use suitable protection, e.g., a covering garment or a sun cream (with a protection rating above 50) which you should then use regularly for three weeks before and after treatment. Ideally, you should wait 4 weeks before using the Homelight epilator if you have been exposed to the sun without protection. If you have the least doubt about exposure to the sun consult your doctor or contact the BaByliss customer service.

16. What should I do before I use the Homelight epilator?

Before each session with the Homelight epilator, it is important that the area to be treated has not been exposed to the sun for at least four weeks. A sun screen with a high protection factor (protection factor of 50+) could be some help, as well as clothes that cover the area to be treated. Moreover, the area to be treated should be washed with soap and water beforehand and then shaved clean of hair.

17. How do I take care of the treated area after using the Homelight epilator ?

The area treated with the Homelight epilator may be cleaned and cared for with well-known skin care products. Specific precautions should be taken to avoid unprotected exposure to the sun. Sun screens with a higher protection factor (protection factor of 50+) and long clothing are good protection against the sun.

18. Should I remove hair before, between and after treatments?

No, do not remove hair before, during or after treatments. The hair will fall out on its own naturally and gradually. This stage may last up to two weeks.

15. Specifications

Point of impact, Specs	1,5 cm x 3 cm [3,9 cm ²]	
Speed	1 pulse every 4 seconds: 1,1 cm ² / sec	
Technology	Home Pulsed Light™	
Light intensity	2,5 – 4,5 J/cm ²	
Wavelength	475 - 1200 nm	
Charge time / Energy source	Continuous use	
Required electrical installation	100 - 240 VAC, 0,9A	
Time needed for half-leg treatment	40 minutes	
Handling and safety	New: a colour detector sensor ensures safe use.	
Temperature	Use	10 to 40°C
	Storage	-40 to +70°C
Relative Humidity	Use	30 to 75%
	Storage	10 to 90%
Atmospheric pressure	Use	700 to 1060hPa
	Storage	500 to 1060hPa

16. Packaging - Labelling

Level of protection against electrical discharge: type BF appliance

Waste electrical and electronic equipment 

CE Marking

Following the usage guidelines

CSA Marking for the United States and Canada

Moisture protection: ordinary

This appliance is not intended to be used near a flammable mixture of anaesthetic and air or with oxygen or nitrous oxide.

17. Quick start-up

Installation of the Homelight appliance

1. Remove the Homelight epilator appliance and all the parts from the box.
2. Ensure that the flash lamp is firmly inserted in the end of the epilator.
3. Plug the adapter lead into its slot in the appliance handle.
4. Plug the other end of the adapter into the mains. The 'Standby' indicator light will come on. Your Homelight epilator is ready to use.

Use the safety glasses included when using the Homelight epilator.

Your first epilation session with your Homelight appliance

Your skin should be shaved, clean, dry and free of residue from powders or deodorants.

5. Press the On/Standby switch. The power indicator light will come on and you will hear the noise of a fan (similar to the noise made by a hair dryer).
6. About 4 seconds after having pressed the switch, the 'Ready' indicator light will come on. The default level for the appliance is set at the lowest setting. It is therefore ready to release its first pulse.
7. Turn the light intensity down to the lowest setting for the first use. If you would like to increase or reduce the power, press on the power selector using '+' and '-' until the indicator light is at the desired power level. The power level indicator light that is lit corresponds to the power level administered (see details on the 'Power intensity level').
8. Press the application surface firmly against your skin so that it is in full contact. When the application surface is in full contact with your skin, the 'Ready' indicator light will begin to flash.
9. Press the release. The mechanism will first check skin colour. If the skin colour is light enough for safe application, the appliance will release a pulse onto the skin and the 'Ready' indicator light will switch off. You will see a flash of light and at the same time you will hear a slight 'beep' with each pulse, which is normal. You may also feel a slight heat and tingling sensation. Your Homelight appliance recharges immediately after each pulse. After 4 seconds, the 'Ready' indicator light will again light up. If you don't see pulsed light and the skin colour indicator light starts to flash, this means the your skin is too dark for safe application. Change the area of treatment or check with the Homelight epilator Customer Service (see details: 'Skin colour sensor').
10. Remove application surface from the treated area of skin.
11. Move the application surface to another area to be treated. Pressure tracks from the application surface made on your skin will help you to align the appliance correctly to administer another pulse, preventing superimposing the pulses or spacing too much. (See 'Covering the treatment area').
12. Repeat the application from step 5.

For epilation of the face, please refer to the instructions on page 150.

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