Hello Kitty  
GESKE

Quickstart Guide

[www.geske.com](http://www.geske.com)

Sanrio  
LICENSE

©2022 SANRIO CO., LTD.


GESKE Beauty Tech GmbH,  
Leipziger Platz 18, 10117 Berlin, Germany  
[info@geske.com](mailto:info@geske.com), [www.geske.com](http://www.geske.com)

SmartAppGuided™  
MicroCurrent Face-Lift Pen | 6 in 1

# Short Introduction

Scan for the full instructions  
in your language

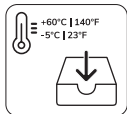GESKE German Beauty Tech App  
Available in 45 Languages

Apple and the Apple logo are trademarks of Apple Inc.  
Google Play and the Google Play logo are trademarks of Google LLC.

Read the safety information before first use.  
Keep both instructions for use.

## YOUR PRODUCT AT A GLANCE

Your product is a multi-functional MicroCurrent device. It is safe to be used up to 3 times per week but not more than once per day. Please do not use more frequently. For maximal safety and ideal results, always follow your personalized Skin Care Routine created by the GESKE App.

**Download the GESKE App and get your personalized routine.**

**Intended Use:** The Geske MicroCurrent Face-Lift Pen & Hello Kitty MicroCurrent Face-Lift Pen Devices are intended for the stimulation of facial skin and are indicated for over-the-counter cosmetic use.

## SPECIFICATIONS:

**Materials:** ABS, silicone, zinc alloy | **Storage temperature:** -5 - +60°C/ 23 - 140°F | **Battery:** Li-Ion 3,7V | **Battery Capacity:** 220 mAh | **Maximum Charging Current:** 220 mA | **Charging Voltage:** 5V | **Charging Time:** approx. 2,5 hours | **Operating Time:** Operating time differs depending on activated technologies & intensities, | **Interface:** USB | **Programs:** Vibration: 3 speed levels, MicroCurrent: 3 intensities

## IMPORTANT SAFETY NOTES! READ BEFORE USAGE!

The product is not suitable for use under the following conditions:

- If you have an implanted medical device or other electronic device such as a pacemaker
  - The safety and effectiveness of the product has never been assessed in combination with non-GESKE skin care products and masks. We recommend using the product only with GESKE products. Use with other products at your own risk
  - After laser treatment or chemical peeling
  - After facial surgery, both cosmetically and medically indicated
- In the following cases, a consultation with a physician is necessary as to whether the product is suitable for you to use:**
- If you have a heart disease
  - In case of pregnancy. The safety of use during pregnancy has not been established

## Further information on use:

- Before using the MicroCurrent Face-Lift Technology, the product must be completely dry
- Use of the MicroCurrent Face-Lift Technology only in combination with the GESKE serum
- Do not use the microcurrent function on the chest area, around the eyes,

midline of the neck, genital/groin area, or on red and inflamed skin, over birthmarks and enlarged blood vessels

- Nothing is known about the long-term effects of microcurrent
- As a result of stimulation of the optic nerve, a flashing light may be perceived during use. If this persists when not using the product, seek medical advice
- A slight tingling sensation may occur during application. Decrease the intensity of the microcurrent to reduce or eliminate tingling
- Electronic monitoring devices may not function properly while the device is in use
- Never use the product while driving or operating heavy machinery
- Do not poke any sharp objects into the product's openings
- The product may only be used with a safety extra-low voltage (SELV) adapter
- We recommend using standard IEC60335 power supplies to charge the device
- Start with the lowest intensity and slowly increase as needed
- Do not use the product on irritated or injured skin or open wounds (after any kind of injury), in the case of bacterial or viral infections of the skin (either in the usage area and elsewhere)
- This product is not intended to treat medical skin conditions or diseases. If you have medical skin conditions or diseases, consult your dermatologist before using this product
- Do not use the product if it is 1) untypically warm, 2) mechanically damaged or deformed, 3) discolored
- If pain and/or discomfort occurs during use, discontinue use and seek medical advice
- The eye area is very sensitive: Be extra careful in this area and avoid any contact of the product with the eyelids or the eye itself
- Only use the product as recommended in your personalized routine
- For hygiene reasons, we recommend you do not share your product with anyone else
- The product is intended for cosmetic use on the face. GESKE assumes no liability for harmful consequences arising from the use on other body areas or other improper use, connection to unsuitable voltage sources, contaminated conducting liquid or balls
- Always make sure the plug and socket are completely dry before charging. Otherwise, electric shocks, a short circuit or fire may result. The product must not be used while connected to the power outlet

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