GESKE

GERMAN BEAUTY TECH

Quickstart Guide


www.geske.com

GESKE Beauty Tech GmbH,
Leipziger Platz 18, 10117 Berlin, Germany
info@geske.com, www.geske.com

SmartAppGuided™
Sonic Facial & Body Roller | 4 in 1

Short Introduction

Scan for the full instructions
in your language

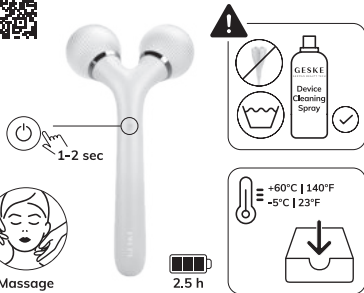GESKE German Beauty Tech App
Available in 36 Languages

Apple and the Apple logo are trademarks of Apple Inc.
Google Play and the Google Play logo are trademarks of Google LLC.

EN Read the safety information before first use.
Keep both instructions for use.

YOUR PRODUCT AT A GLANCE

Your product is a multifunctional derma roller. It is safe to be used up to 3 times a day, please do not use more frequently. For maximal safety and ideal results, always follow your personalized Skin Care Routine created by the GESKE App.

Download the GESKE App and get your personalized routine.

SPECIFICATIONS:

Materials: ABS, silicone | **Storage temperature:** -5 - +60 °C/23 - 140 °F | **Battery:** Li-Ion 3.7V | **Battery Capacity:** 140 mAh | **Maximum Charging Current:** 140 mA | **Charging Voltage:** 5V | **Charging Time:** approx. 2.5 hours | **Operating Time:** Operating time differs depending on activated technologies & intensities. | **Interface:** USB | **Programs:** Vibration, 1 intensity

IMPORTANT SAFETY NOTES! READ BEFORE USAGE!

In the following cases, a consultation with a physician is necessary as to whether the product is suitable for you to use:

- In case of pregnancy. The safety of use during pregnancy has not been established
 - Do not use the product on irritated or injured skin or open wounds (after any kind of injury), in the case of bacterial or viral infections of the skin (either in the usage area and elsewhere)
 - This product is not intended to treat medical skin conditions or diseases. If you have medical skin conditions or diseases, consult your dermatologist before using this product
 - Do not use the product if it is 1) untypically warm, 2) mechanically damaged or deformed, 3) discolored
 - If pain and/or discomfort occurs during use, discontinue use and seek medical advice
 - The eye area is very sensitive: Be extra careful in this area and avoid any contact of the product with the eyelids or the eye itself
 - Only use the product as recommended in your personalized routine
- For hygiene reasons, we recommend you do not share your product with anyone else
- The product is intended for cosmetic use on the face. GESKE assumes no liability for harmful consequences arising from the use on other body areas or other improper use, connection to unsuitable voltage sources, contaminated conducting liquid or balls
 - Always make sure the plug and socket are completely dry before charging. Otherwise, electric shocks, a short circuit or fire may result. The product must not be used while connected to the power outlet

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