
LIROMA

RED LIGHT THERAPY

Manual

FORCE

RL60 / RL120 / RL240 / RL300Max /
RL600Max

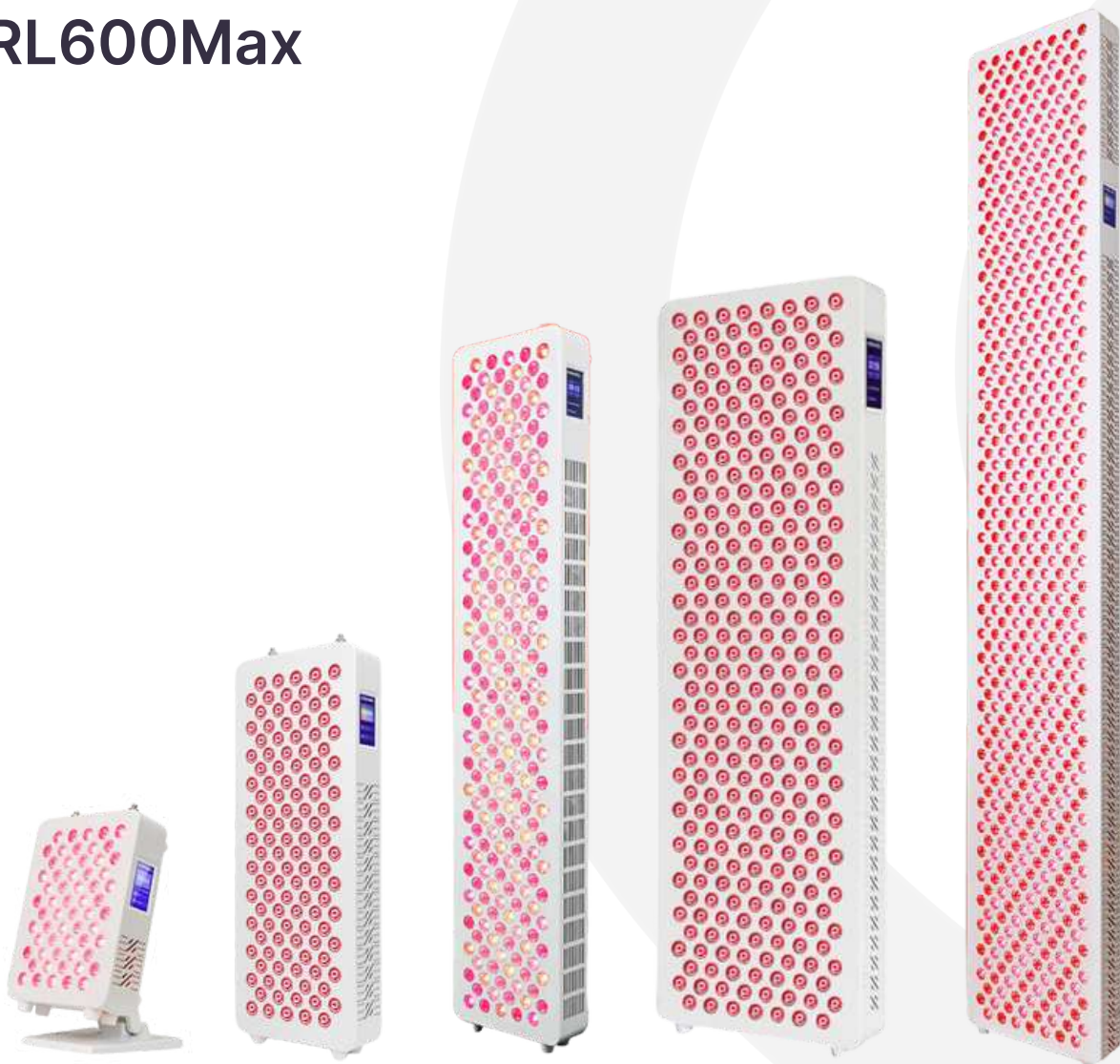- Blue light, Red light, Infrared light (NIR)
- Seven Wavelengths

Table of content

- 1.Thank you for your purchase. (3)
- 2.Red light therapy. (3)
- 3.Specifications. (4-6)
- 4.Benefits Liroma red light therapy panels. (7)
- 5.Instructions for use. (8)
- 6.Control panel. (9-10)
- 7.Which setting should I choose. (11)
- 8.PULSE function. (11)
- 9.Dimmer function. (12)
- 10.Standard - Suspension system (13)
- 11.Accessories (14)
- 12.Warnings (15-16)
- 13.Important notes. (16)
- 14.Liroma. (17)

1. Thank you for your purchase

Thank you for choosing Liroma Red Light Therapy!

Liroma wants its customers to feel healthy, fit and relaxed.

To achieve this, Liroma has a wide range of high quality red light therapy products. Our products help with complaints such as: fatigue, poor circulation, pain in muscles and joints, pain due to underlying diseases, poor sleep, skin problems, stress and many other ailments. Our products target both mental and physical ailments.

Note: Red light therapy lamps do not get hot! During treatment, however, you may feel heat and sweat inside.

2. Red light therapy

Red light therapy uses a combination of red light and near-infrared light (NIR) to stimulate and repair the body at a cellular level. This therapy works by allowing the light to penetrate deep into the skin, with red light targeting the upper skin layers and NIR light reaching deeper into muscles and joints. A key effect of this light therapy is the stimulation of mitochondria—the energy factories of our cells—helping them produce energy more efficiently in the form of ATP (adenosine triphosphate).

With this increase in cellular energy, cells can recover more quickly, leading to reduced inflammation and improved blood circulation. Red light also promotes collagen production, which can contribute to firmer, more youthful skin and improved skin texture. Thanks to the combination of red light and NIR light, users often experience benefits such as reduced muscle and joint pain, enhanced skin health, faster muscle recovery, and increased energy levels. These effects make red light therapy a popular choice for those looking to support their well-being naturally.


Compared to other providers, Liroma's panels also feature blue light (480 nm). Blue light is known for its powerful antibacterial properties and is often used in the treatment of acne and skin impurities. It penetrates the superficial skin layers and helps combat bacteria that cause inflammation. Additionally, blue light can contribute to improved skin texture and a more even complexion. Combined with red and NIR light, this new functionality offers an even more comprehensive approach to skin health and overall well-being.

3. Specifications

Model	Force RL60
Number of LED's	60
Wavelength (nanometer)	480/630/660/ 810/830/850/1060
Light spectrum	Blue light/ Red light/ NIR-infrared
Size	28,5 x 20,9 x 6,4 cm
Weight	2,4kg
Beam angle	30 degrees
Current	182W±10%
EMF	-0,0 UT
Radiation	219 mw/cm2 on the surface 204 mw/cm2 on 7,5cm
Lifespan	100.000 hour
Input voltage	100V - 240V 50/60Hz

Model	Force RL120
Number of LED's	120
Wavelength (nanometer)	480/630/660/ 810/830/850/1060
Light spectrum	Blue light/ Red light/ NIR-infrared
Size	54 x 20,8 x 6,4 cm
Weight	4kg
Beam angle	30 degrees
Current	278W±10%
EMF	-0,0 UT
Radiation	219 mw/cm2 on the surface 204 mw/cm2 on 7,5cm
Lifespan	100.000 hour
Input voltage	100V - 240V 50/60Hz

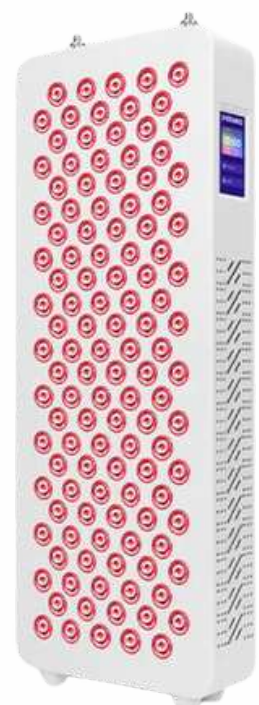Model	Force RL240
Number of LED's	240
Wavelength (nanometer)	480/630/660/ 810/830/850/1060
Light spectrum	Blue light/ Red light/ NIR-infrared
Size	99,8 x 20,8 x 6,4 cm
Weight	7,2kg
Beam angle	30 degrees
Current	374W±10%
EMF	-0,0 UT
Radiation	219 mw/cm2 on the surface 215 mw/cm2 on 7,5cm
Lifespan	100.000 hour
Input voltage	100V - 240V 50/60Hz

Model	Force RL300Max
Number of LED's	300
Wavelength (nanometer)	480/630/660/ 810/830/850/1060
Light spectrum	Blue light/ Red light/ NIR-infrared
Size	92 x 30 x 6,4 cm
Weight	10,5kg
Beam angle	30 degrees
Current	460W±10%
EMF	-0,0 UT
Radiation	230 mw/cm2 on the surface 220 mw/cm2 on 7,5cm
Lifespan	100.000 hour
Input voltage	100V - 240V 50/60Hz

Note: NIR infrared light is not visible to the naked eye. You may see a faint glow. Additionally, red light therapy lamps do not get hot! During the treatment, you may feel warmth from within and start sweating.

Model	Force RL600MAX
Number of LED's	
Wavelength (nanometer)	480/630/660/ 810/830/850/1060
Light spectrum	Blue light/ Red light/ NIR-infrared
Size	182 x 29,9 x 6,4 cm
Weight	18kg
Beam angle	30 degrees
Current	374W±10%
EMF	-0,0 UT
Radiation	219 mw/cm2 on the surface 215 mw/cm2 on 7,5cm
Lifespan	100.000 hour
Input voltage	100V - 240V 50/60Hz

4. Benefits of Liroma red light therapy panels

- Penetrates deep into the skin (up to the bones)
- Combination of blue, red, and NIR infrared light
- 7 wavelengths > 480 / 630 / 660 / 810 / 830 / 850 / 1060nm
- User-friendly display
- Touchscreen
- Various programs
- Timer
- Pulse mode
- Can be used with a rolling foot stand and a large rolling stand (see the website for this product)
- LED lamps (energy-efficient)
- Reduction of muscle and joint pain
- Acceleration of wound healing
- Improvement of skin condition
- Anti-inflammatory effect
- Stimulation of blood circulation
- Improvement of sleep quality
- Hair restoration
- Relief from stress and mood improvement
- Better muscle recovery after exercise
- Enhanced immune function
- Relief from seasonal depression
- Support for the nervous system
- Improvement of mental clarity and concentration
- Reduction of acne and other skin issues
- Increase in cellular energy production
- Improvement of hormonal balance

5. Instructions for use

- Ideal usage frequency: 3-7 times per week, with a maximum of once per day per body part, at a distance of 30-50 cm.
- Maximum distance: 90 cm.
- Recommended treatment example: 12 minutes at a distance of 30 cm.
- Use the lamp for a maximum of 20 minutes. Be mindful of the distance—if you wish to use the lamp for 20 minutes, increase the distance.
- Do not look directly into the lamp and wear protective goggles.
- Ensure that the body part being treated is fully exposed to the light from the lamp and is not covered by clothing.
- Adjust the lamp and then put on the protective goggles. Try to avoid looking directly into the lamp.
- Always consult a doctor before use in case of contraindications, medical conditions, or if you are unsure.
- Note: The lamp does not emit heat! It is designed solely for LED therapy.

6. Control panel

Turn on the device using the red button at the back of the unit.

You will now be on the home screen.

Below, we explain the different buttons:


- Scene
- Settings
- Start/Pause

Scene

With this button, you can choose from different programs. We always recommend activating both red light and NIR infrared light for any purpose.

- The difference between the programs mainly lies in the brightness and duration. You can scroll between the different programs.
- Once you select a program, it will start immediately. You can pause it using the Start/Pause button.

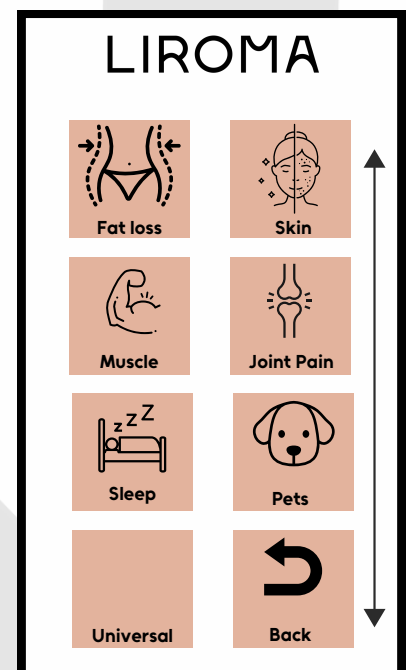Settings

- **Timer**

Press this button to adjust the timer. You can set the time by rotating the circle.

- **Mode**

Press this button and choose from RED/NIR, RED, or NIR. Use red light (RED) and near-infrared light (NIR) together or activate them separately. We recommend always using both lights simultaneously for the best results.

- **Brightness**

Adjust the brightness of the lights with this button. Many people start with a lower brightness if the light feels too intense at first.

- **Pulse**

Control the pulse function here. This is displayed in Hertz (Hz). The higher the Hz, the faster the light flickers. We recommend experimenting with pulse settings between 1 Hz and 10 Hz to find what feels most comfortable for you.


- **Blue Light**

Turn the blue light on or off here.

- **Group**

You can link two panels together by connecting an additional adapter to both devices. In the settings menu, go to "Group" to configure this feature.

On one panel, select "Primary," and on the other panel, select "Secondary." You can then control both panels from the panel set to "Primary."

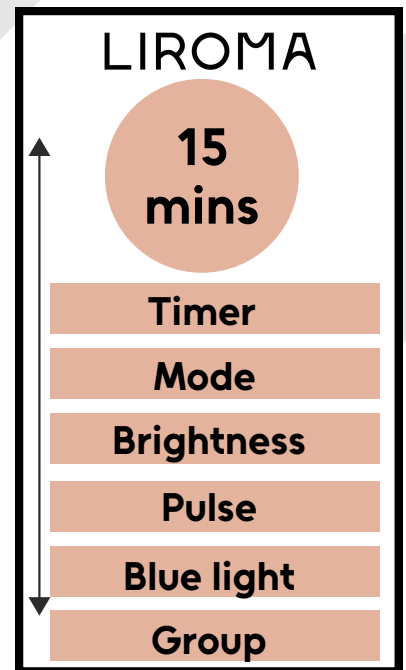Start/Pause

With this button, you can turn the lamps on and pause them.

7. Which setting should I choose?

Blue light (480nm)

- Reduces acne
- Reduces redness
- Regulates sebum production
- Helps stay awake

Red light (630 & 660nm)

- Improves skin health
- Stimulates blood circulation
- Stimulates collagen and elastin production

NIR (near-infrared light) (810, 830, 850, 1060nm)

- Reduces inflammation
- Relieves muscle and joint pain
- Restores deep tissues, such as muscles, joints, and organs

We always recommend turning on all the lights at the same time for the best results.

8. PULSE-function

Pulsed light in red light therapy penetrates deeper into the tissue, stimulates cell repair, and increases ATP production. It enhances circulation, reduces heat stress, and can activate specific biological responses at different frequencies. Additionally, research suggests it may offer neurological benefits, such as in Alzheimer's and Parkinson's diseases. A combination of continuous and pulsed light may yield the best results.

We recommend experimenting with pulsed settings between 1 Hz and 10 Hz to find what feels most comfortable for you.

9. Brightness

A dimmer function on a red light therapy panel offers users the flexibility to adjust the light intensity. This provides greater comfort, especially for sensitive skin types, and allows for customized therapies. Additionally, a lower intensity helps during longer sessions and can reduce energy consumption, which in turn extends the lifespan of the device.

10. Standard - Hanging System

Standard:

Assemble the foot stand and ensure that the bottom of the lamp is placed into the stand. The stand is delivered fully assembled.

Hanging system:

- Step 1

Screw the silver hanging system onto the top of the device.

- Step 2

Place the door hanger at the top of the door.

- Step 3


Connect the pulley to the pulley of the silver hanging system. Hang the pulley onto the door hanger.

- Step 4


To raise the lamp, pull the loose black cable down while supporting the lamp with your hand.

- Step 5

Push the lever on the height adjustment pulley upwards to lower the lamp. Support the lamp again with your hand.

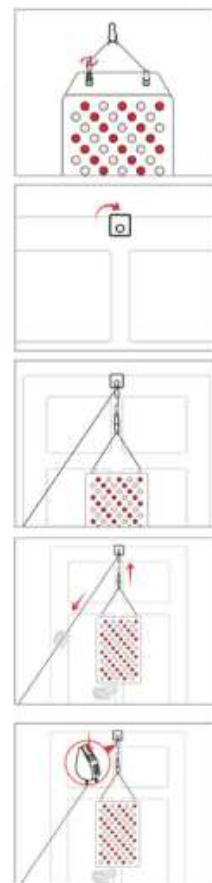11. Accessories

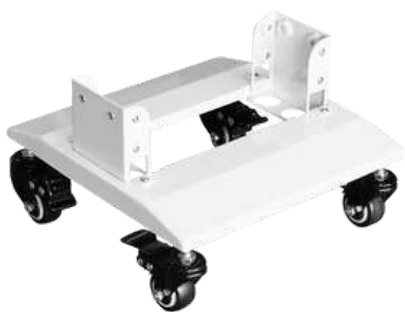Footstand on Wheels (RL120, RL240, and RL300Max)

The footstand is compact and adjustable, making it ideal for use while sitting or standing. The wheels allow you to easily move the panel to the desired location, ensuring both flexibility and stability.

Compatible with the RL120, RL240, and RL300Max.

Large Stand on Wheels (Horizontal Use)

This large stand enables horizontal use, making it ideal for lying down during treatment. With its wheels, it is easy to move and perfect for treating larger body areas, such as the back or legs.

Compatible with the RL120, RL240, and RL300Max.

Footstand on Wheels (RL600)

The footstand for the RL600 is made of sturdy material. The wheels allow you to move the panel smoothly and easily, making it ideal for flexible use in different spaces.

Compatible with the RL600.

Large Stand on Wheels (RL600, Horizontal Use)

With the large stand, you can use the RL600 both horizontally and vertically, making it perfect for clinics and salons. The sturdy construction and wheels allow for easy movement and positioning.

Compatible with the RL600.

12. Warnings

Warnings and Precautions

- Do not use this device on babies.
- Do not allow sleeping, disabled, or unconscious individuals to use this device.
- Do not use the device for more than 20 minutes per therapy session.
- Do not use the device on open wounds.
- Do not use the device in combination with liniments, ointments, or balms that contain heat-producing ingredients, as this may cause burns.
- Do not attempt to modify or repair the device yourself. If it is damaged or not functioning properly, please contact us as soon as possible at info@liroma.nl.

This device is electrical. To prevent electric shock and other hazards, please follow these instructions:

- Do not wash the product with water or other liquids.
- Do not place the device in water or in areas where it may be exposed to water or other liquids.
- Do not expose the device to running water.
- If the device comes into contact with or is submerged in water or any other liquid, immediately unplug it from the socket. If using the device near water, we recommend using a socket with a ground fault circuit interrupter (GFCI).
- Do not touch the device or its plugs with wet hands.
- Do not use or touch the device with wet or damp skin.
- Do not lift the device by the power cord.
- Do not use the device if the power cord is damaged.
- Do not use the device on or near heated surfaces.
- Do not use the device in areas where it may be exposed to flammable or combustible products or fumes, as this could cause a fire or explosion.
- Do not store the device at temperatures above 140°C (284°F).

Additional Warnings

- Do not use the device if it is damaged, as continued use of a damaged device may cause electric shock or injury.
- Do not use extension cords unless they are designed to handle the total power of the device.
- Do not unplug the device by pulling on the cord. Always grasp the plug itself when removing it from the outlet.
- This device is not suitable for individuals with physical, neuro-sensory, or mental impairments, or those lacking experience and/or knowledge.
- Do not use the lamp in combination with implanted electronic devices, such as artificial hearts or artificial lungs, without first consulting your doctor.
- Do not use the lamp on pregnant women, individuals with severe heart conditions, critically ill patients, or emergency cases without consulting your doctor.
- Discontinue treatment immediately if you experience hypersensitivity or allergic skin reactions due to medication, cosmetics, or food.
- Under no circumstances should you open the device or attempt to repair it yourself, as this may cause injury and void your warranty. If the device is damaged or not functioning properly, contact the retailer or info@liroma.nl as soon as possible.

Do not use this product

Without first consulting your doctor if you suffer from:

- Severe illnesses such as heart disease, acute inflammation, thrombosis, blood clotting disorders, adrenal suppression, systemic lupus erythematosus, or malignant conditions.
- Hypersensitivity to infrared light or have edema.
- A condition related to nerve demyelination, such as carpal tunnel syndrome, multiple sclerosis, or ulnar neuropathy.
- Pregnancy.
- Certain medications or plant extracts that cause photosensitivity are not compatible with light therapy. (Examples: antibiotics, antidepressants, Roaccutane, MTX, antifungal medication, and St. John's Wort).


13. Important notes

- Half of the LED lights in the panels emit red light, while the other half emits NIR infrared light (Near Infrared).
- Red light is visible, while NIR infrared light is invisible to the naked eye. However, with NIR infrared light, you may perceive a slight glow or dimming, which is completely normal.
- Always use it on exposed skin and wear the included glasses for added protection.

LIROMA

Red Light Therapy

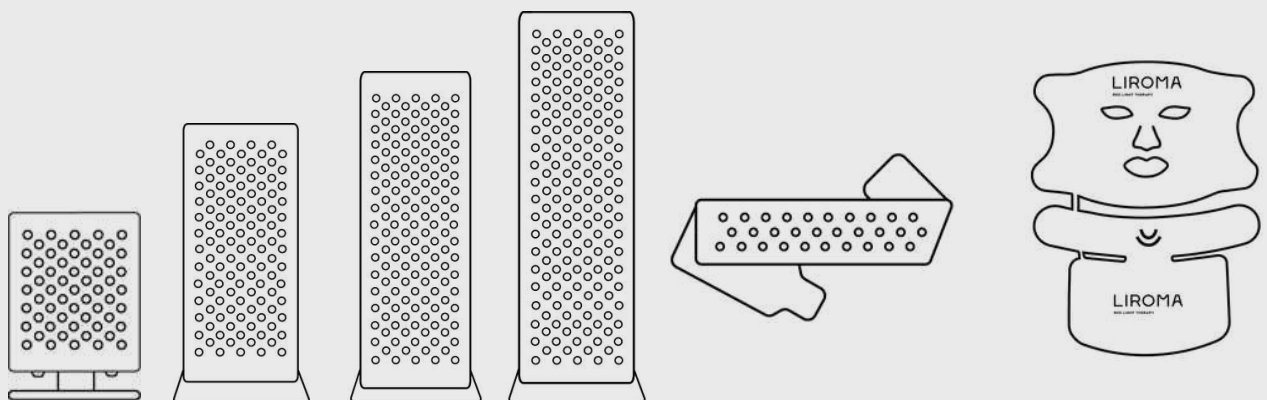Liroma has been a specialist in red light therapy for years. We offer lamps that provide full-body exposure, as well as smaller devices that allow for targeted treatment on specific body parts. Additionally, we have red light devices specifically designed for the back and face. Visit our website to explore all our products.

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