
Manual version 1.0

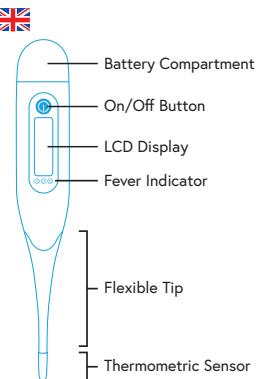This device is intended for measurement of human body temperature.

⚠ Warning and safety


- This product is not intended to diagnose, treat, cure, or prevent any disease. Consult your doctor for diagnosis and treatment.
- Do not immerse the thermometer under water deeper than 15cm for longer than 5 minutes.
- Do not drop the thermometer. It is neither shock-proof nor impact-resistant.
- Do not clean with thinners, petrol or benzene. Only clean with water or disinfectant.
- The thermometer contains small parts. Do not leave the thermometer unattended in the hands of children.
- Do not allow the thermometer to come into contact with hot water or fire.

Method of Correct Measurement

The measurement site must always be specified in order to ensure that a correct temperature reading is recorded.

Rectal Measurement

From a medical point of view, this is the most accurate method because the temperature at this measurement position is closest to the core body temperature. The thermometer tip is inserted carefully into the rectum for a maximum of 2cm (1inch). The usual measuring time is approximately 8 to 15 seconds.

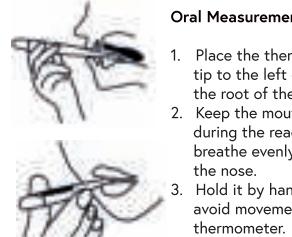Axillary Measurement is not recommended

The following situations will cause incorrect temperature

- After activities, such as bathing or exercising
- Drinking or eating before oral measurement
- Breathing through the mouth during oral measurement

Usage

- Press the power button. A beep sounds and the thermometer runs a check after which the last taken temperature is shown for 2 seconds.
- As soon as - - - appears on the display, the thermometer is ready for use.
- During the reading, the current temperature is displayed continuously and the °C or °F symbol flashes. The measurement is completed when a constant temperature value has been reached. A beep will sound four times, and the °C or °F symbol will stop flashing.

Fever alarm

If the detected body temperature exceeds 37.8°C (100.0°F) the thermometer will sound an alarm: "bee-bee-bee — bee-bee-bee —".

Fever indication

- Below 37.2°C (99.0°F); healthy temperature
- 37.3°C - 37.7°C (99.1°F - 99.9°F); slightly elevated temperature
- Over 37.8°C (100.0°F); elevated temperature or fever

Switch between Celsius and Fahrenheit

- Make sure the thermometer is turned off.
- Press and hold the power button until °C or °F appears on the screen.
- Keep holding the button and release only after it automatically switches to °C or °F.

Auto-off function

The thermometer will shut off automatically after 10 minutes of not using the thermometer.

Memory function

When the thermometer is switched on, the last taken temperature will be shown for 2 seconds.

Cleaning

The best way to clean the thermometer is by applying disinfectant (e.g. 70% ethyl alcohol) with a damp cloth.

Battery replacement

The battery is empty and needs replacing when the battery symbol █ appears on the right of the LCD.

Battery: Alkaline button battery LR41

Technical data

Type: maximum thermometer
Measurement range: 32.0-43.0°C (89.6-109.4°F)

Measurement accuracy:

±0.1°C at 35.5-42.0°C at ambient temperature 18-28°C

±0.2°C at 32.0-35.5°C at other ambient temperatures. (±0.2°F at 95.9-707.6°F at ambient temp. 64.4-82.4°F

±0.4°F at 89.6-95.9°F at other ambient temperatures.)

Storage/transport conditions: -20-55°C, 75%-95% RH

Ambient operating conditions: 5-35 °C, 75%-95% RH

Min Scale: 0.1°C/0.1°F

Battery type: Alkaline manganese battery, type LR41,1.5V.

Weight: Approximately 12g

This thermometer is latex free.

Explanation of LCD symbols

The battery is empty

Lo°C/ Lo°F The temperature is under 32 °C or 89.6 °F

Hi°C/ Hi°F The temperature is over 43 °C or 109.4 °F

Explanation of symbols on the label

The battery in this product complies with the requirements stated in European Directives 2006/66/EEC.

Type BF equipment.


Consumed electronic products must not be mixed with conventional domestic waste. Take this product to an appropriate waste collection point for correct handling and recycling.

Warranty

This limited warranty covers any defects in material or workmanship under normal use during the warranty period of 365 days after purchase. iProven will either replace the product or repair the product at no charge, using new or refurbished replacement parts. A replacement product or product part assumes the remaining warranty of the original product purchase. This Limited Warranty does not cover batteries and packaging. It also does not include any problem that is caused by: conditions, malfunctions or damage not resulting from defects in material or workmanship.

Customer support

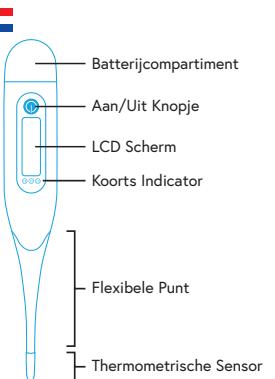For obtaining warranty service and for any other inquiries or suggestions please contact iProven customer support at www.iproven.com

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