MEDISANA®

- DE** Haarentfernungssystem IPL 850
GB Hair Removal System IPL 850
FR Système d'épilation IPL 850
IT Sistema di epilazione IPL 850
ES Sistema de depilación IPL 850
NL Ontharingsysteem IPL 850

Art. 88588

Gebrauchsanweisung
Instruction Manual
Mode d'emploi

Bitte sorgfältig lesen!
Please read carefully!
A lire attentivement s.v.p.!

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IMPORTANT INFORMATION! RETAIN FOR FUTURE USE!

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use.

Should you give this device to another person, it is vital that you also pass on these instructions for use.

Explanation of symbols

This instruction manual belongs to this device. It contains important information about starting up and operation. Read the instruction manual thoroughly. Non-observance of these instructions can result in serious injury or damage to the device.

WARNING

These warning notes must be observed to prevent any injury to the user resp. damage to the device.

LOT

LOT number

Manufacturer


**Protect against moisture!
Store in dry location!**

**Warning - Optical radiation
Do not look into the light source!**

Only for indoor use!

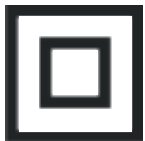Protection category II

WARNING

Possible eye injury (potentially leading to loss of vision) or skin injury if instructions are not followed. Read and follow the instructions. Protect the eye from exposure. Do not use the appliance over the eye lids or close to the eye.

- Only use the device for the specific purpose described in the instruction manual. Any misuse will void the warranty.
- Never use the device if you have skin type VIII/IX (label: rarely sunburnt, very dark suntan) or very dark brown skin. Observe the “light intensities” table in these operating instructions!
- This device is not suitable for use on white or grey hairs. In this case contact your doctor before use.
- Never use the device after sunbathing or artificial tanning, in order to avoid damage to the skin. Leave at least 4 weeks between skin tanning and use of the device.
- After use for at least 2 weeks no solar radiation or artificial tanning may take place for more than 15 minutes. Use a sun protection cream with minimum protection factor 30 to avoid complications.
- For some people treatment with this device may cause the opposite effect (stronger hair growth). Consult your doctor.
- In any case contact your doctor before using the device, to find out if the use of the device could be questionable for your health.
- Never use the device while pregnant or while breastfeeding.

- Never use the device if you suffer from the following diseases or problems:
 - Eczema, psoriasis, lesions or open wounds
 - Infections, scars, abrasions or skin burns
 - Sunburn, redness or blistering
 - Diabetes or other metabolic diseases
 - Skin cancer, potential malignoma, epilepsy, collagen disorders, Herpes or blood clotting disorders
 - Immunodeficiency (e.g. HIV/Aids)
- For skin wounds that are healing wait until the skin has completely healed before you use the device.
- The device may not be used after radiation or chemotherapy within the last 3 months.
- Use is also to be refrained from in case of photosensitivity or other light sensitivity (e.g. polymorphic light dermatosis, light urticaria, Lupus etc).
- Never use the device if you are taking one or more of the medications listed below:
 - Medications which make the skin sensitive to light, including non-steroidal anti-inflammatory drugs (e.g. aspirin, ibuprofen, paracetamol), Tetracycline, Phenothiazine, thiazide diuretics, sulphonylureas, sulphonamides, DTIC, fluorouracil, vinblastine, Griseofulvin, alpha hydroxy acids (AHAs), beta hydroxy acids (BHAs), Retin-A®, Accutane® and/or topical retinoids.

- Also do not use the device if:
 - if your skin is currently or has in the past been treated with alpha-hydroxy acids (AHAs), beta-hydroxy acids (BHAs), locally applied isotretinoin and azelaic acid.
 - if you have been treated within the last 6 months with Accutane®.
 - If in the last 3 months you have had a steroidal course of treatment
 - if in the last 8 weeks you have undergone skin peeling or other skin smoothing methods
 - in case of active implants such as a heart pacemaker, a device for incontinence, an insulin pump etc.

If you are unsure if the application of the device is appropriate contact your doctor. The points listed in these operating instructions are in no way intended to represent a comprehensive list!

- This device is only suitable for persons over 18 years of age.
- Never look directly into the light source!
- Never use the device close to wet environment such as bath tubs, showers, washbasins etc - risk of electric shock!
- Keep the device away from moisture and liquids! A wet device may not be operated.
- Do not push anything into the device.

- Do not use the device if you find wear, damage or signs of improper use.
- If a fault occurs, do not attempt to repair the unit yourself. Repairs must only be carried out by an authorized specialist dealer or other suitably qualified personnel.
- This device is not designed to be used by persons with limited physical, sensory or mental abilities, or by persons with insufficient experience and/or knowledge, unless under observation by a person responsible for their safety, or unless they have been instructed in the use of the device.
- Children must be supervised to ensure that they do not play with the device.
- Keep the device out of the reach of children.
- The device and the power supply are only intended for use in interior rooms in a private area, and not for commercial use.
- Never use the device on the following parts of the body:
 - close to the eyes or the eyebrows
 - on nipples, areola, inside the labia, vagina, anus, on mucous membranes and in the nose or ears
 - for men on the face or on the scrotum
 - over piercings or other metal objects, such as earrings or jewellery.
 - on tattoos, make-up, warts, paint, freckles, dark areas of the skin etc.

- on damaged or cracked areas of skin, or on areas, that could be damaged
- on parts of the body on which you may later wish to have hair
- If you find redness, blistering or burns on the skin, stop use immediately!
- If the light area ⑤ or the area of skin itself should become hot, end the application.
- Never try to emit a pulse of light into the air. The light area ⑤ must always be in full contact with the skin surface.
- Never use the device several times in succession on the same area of skin.
- Avoid overlapping light pulses.
- Cover body parts with birthmarks, tattoos, dark brown or black spots, freckles, warts etc with a cloth before use.
- Never use flammable liquids (e.g. perfumes, disinfection agents or alcoholic substances) or acetone for skin cleaning before the application.
- Use of the device can lead to pigmentation changes in the skin.

- Never alter or dismantle the device.
- Never reach into the device, if it has fallen into water. Immediately disconnect the power supply and contact customer service.
- Never leave the device unattended when it's connected to the mains.
- The device will be switched off automatically after approx. 3 minutes.
- Keep the mains cable away from sources of heat. Never use the device if the cable or power supply is defective.
- Never use the device if it emits smoke or the fan doesn't work.
- Always disconnect the device from the power grid after use and before any cleaning.
- Never use spare parts or attachments that are not recommended by the manufacturer.
- Never operate the device with any power supply other than the original one.
- This device is not suitable for use in environments with flammable, anaesthetic air mixtures, oxygen or laughing gas.

**Possible side effects:**

If this device is used as intended then side effects and complications are very rare. Nevertheless any cosmetic application involves a certain risk potential. It is therefore important that you are familiar with the possible dangers in conjunction with the hair removal by IPL technology (IPL=intense pulsed light) in home use:

- It is possible that the skin is slightly reddened after treatment, a higher amount of skin sensitivity occurs or that you experience a general sense of discomfort. These reactions are mostly harmless and quickly ease. However, if these reactions are still present after 24 hours contact your doctor.
- After treatment skin reactions may occur similar to sunburn or scarring or pigmentation changes may also occur. Swellings or a strong degree of redness may also occur. These symptoms should subside within 2 to 7 days. Contact your doctor to rule out further complications.
- In extremely rare cases infections or bruises / contusions may occur. In such cases contact your doctor immediately.
- The use of the device can in very rare cases lead to blue or pink coloured haematomas which can last 5 to 10 days. When the haematomas fade, this can lead to a rusty brown discolouration of the skin (hyperpigmentation) which may be permanent.

- Increasing the light intensity (can be set in levels 1 to 4) on the one hand means improving the effectiveness of the application, but on the other it means a higher risk likelihood for the occurrence of side effects. Therefore always start an application with the lowest level (level 1). Only if after several applications over a number of days you observe no side effects, you should successively increase the light intensity (depending on your skin and hair type, up to the values given in the table).
- In combination with shaving itching or dry skin may occur. These symptoms are generally harmless and disappear after a few hours.
- Strongly felt pain may occur during or after treatment in the following instances:
 - if the light intensity is too strong for your skin tone;
 - if you apply the device to unshaven skin;
 - if the same skin area was treated several times with light pulses;
 - if you use the device on open wounds, infections, tattoos, burns etc.
- Cleaning and user maintenance shall not be made by children without supervision.

**Thank you
very much**

Thank you for your confidence in us and congratulations!

With the **IPL 850** hair removal system you have purchased a quality product from **MEDISANA**. The device is designed for the removal of unwanted body hair as well as for the reduction of hair regrowth in women and men. This device is exclusively suitable for use on the legs, armpits, bikini area, breast, stomach, back and face below the cheek bones (on face only for women).

So that you achieve the desired results and have a long period of enjoyment with your **MEDISANA IPL 850** we recommend that you carefully read the following instructions on usage and care.

**3.1
Scope of
delivery
and
packaging**

Please check first of all that the device is complete and is not damaged in any way. In case of doubt, do not use the appliance and contact your supplier or your service centre.

The delivery scope includes:

- **1 MEDISANA** Hair Removal System **IPL 850**
- 1 Cartridge with 200,000 light pulses
- 1 Power supply
- 1 Brief guide
- 1 Operating instructions
- 1 Storage bag

The packaging can be reused or recycled. Please dispose properly of any packaging material no longer required. If you notice any transport damage during unpacking, please contact your supplier without delay.

**WARNING**

Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!

3.2

Functionality

The IPL (intense pulsed light) technology sends light pulses to the skin, which are absorbed by the hair roots. These light pulses lead to a selective heating of the hair shaft, that finally stimulates the hair follicle to go into a rest phase. Especially the pigment melanin found in hair or the red blood cells receive the light energy well - that is also the reason why for example thick black hair responds better to IPL treatment than fine, blonde hair. Through the treatment existing hair is shed and the subsequent growth of hair reduced. The growth cycle of hair consists though of different phases - in order to work optimally, the application of IPL technology must be present during the growth phase. As not all hair of the human body are in a growth phase at the same time, the application of the device is required over a longer period of time.

A complete hair growth cycle depends on individual factors and can - depending on type, hairiness and part of the body - last from 18 to 24 months. So permanent hair removal is only possible if the device is applied over the entire period of time.

For your protection the device has a built-in skin contact sensors in the skin contact area **4**. The sensor then only allows the emission of a light pulse when the device is in full contact with the skin. So light pulse emissions into the air are prevented.

The device has a built-in UV filter (integrated into the light area **5**), by which harmful UV rays are blocked.

Skin tone sensor: The device has a skin tone sensor **6**, which is integrated in the cartridge. The sensor measures skin tone before and during application. If the skin tone within treated area is too dark, light impulses are prevented.

3.3 Light intensities

The device allow setting the light intensity levels 1 to 4. To start with you should initially begin with Level 1 and check whether the application runs without complications (see chapter 4 “Operating”). Then slowly you can subsequently (depending on your hair and skin colour) switch to higher levels, if necessary. Observe the last page of these operating instructions for this. The table shows which light intensities are best suited in conjunction with your hair and skin colour or whether these methods are suitable for you at all (if not, this is shown by an ‘X’ on the table).

4 Operating

4.1 Preparation

The **IPL 850 Hair Removal System** is only suitable for skin types I to VII. Skin type VIII (darker skin type) as well as heavily browned skin are not suitable for treatment with this device. Before starting the application, the following points must be observed:

1. Observe all safety instructions and the possible side effects, which are described in this guide.
2. Using the light intensity table (you will find it on the rear side of these operating instructions) which intensity level is suitable for you (depending on skin and hair colour). An “X” means that the device is not suitable for you.
3. The areas of skin which should be treated with the device should be as free of hair on the skin surface as possible. You can shave the areas or shorten the hairs situated there. But do not use a hair removal cream.
4. The areas to be treated must be hairless, completely dry and free from oily or creamy substances, powder, make-up or deodorant etc.
5. Before the very first use of the **IPL 850** first a skin test is to be carried out (see next chapter).

4.2 Skin test

The skin test serves to check tolerance and for a better estimation as to how your skin will react to the treatment. Perform this at least 48 hours before the first application on the area of skin to be treated. Proceed as follows:

1. Plug the end of the device plug of the power supply cable into the device **2** and the power supply plug into a socket. The readiness indicator **9** flashes green.
2. Press the ON/OFF button  **1**. The intensity display **7** shows level 1, (blinking) and the device fan switches itself on. The device is now ready for positioning on the skin.
3. Place the light area **5** onto the skin. The intensity display **7** “Energy” stops blinking. (Display still blinking? -> chapter 5.3)
4. Press the ON/OFF button  **1**. The intensity display **7** shows level 1 again.
5. Press the trigger button **8** to trigger a light pulse at intensity level 1. The skin contact sensors **4** must lie fully on the skin, so that a light pulse can be emitted by the device. When the pulse is emitted, you will see a bright flash and may hear a clicking sound.
6. Now place the device on another area of the skin while it charges for the next light pulse.
7. Then emit a new light pulse on the other skin area.
8. Repeat this procedure for further areas of skin, for which a treatment is planned.
9. End the skin tolerance test (switch off the device by pressing and holding the ON/OFF button  **1**) and wait at least 48 hours. Check the skin areas that have been treated. If they appear normal (no reactions or only slight redness), the areas of skin tested can be treated. However if the areas of skin display strong, painful or other unusual reactions, you may not carry out any treatment on the corresponding skin areas. In this case contact your doctor.

4.2 Skin test

10. Perform the skin tolerance test again, if you intend to increase the intensity (then test with the new planned intensity) or if you wish to treat new skin areas.

WARNING

If pain, swelling or other unusual reactions occur during the skin test, immediately stop the test and do not carry on using the device. Contact your doctor.

4.3 Application plan

So that optimum results are achieved with your **IPL 850 Hair Removal System**, several applications should be carried out at certain intervals and repeated over a longer time scale of at least 18 to 24 months. Here should be differentiated between body applications and facial applications (facial applications only for women below the cheekbones).

For bodily applications:

- The first 3 to 4 hair removal applications should each be 2 weeks apart.
- The 5th to 7th hair removal applications should each be 4 weeks apart.
- From the 8th application an application should take place from time to time (about every 4 to 8 weeks) until long-term results have been achieved.

For facial applications:

- The first 6 hair removal applications should each be 2 weeks apart.
- The 7th to 12th hair removal applications should each be 4 weeks apart.
- From the 13th application an application should take place from time to time (about every 4 to 8 weeks) until long-term results have been achieved.

This application plan is our recommendation for the best possible results. Depending on personal skin type, hair colour, application area and other biological factors for various people this may result in different treatment results.

4.4 Starting application

Make sure that the skin areas to be treated are shaved, dry, clean and free of creams, deodorants etc. Proceed as follows:

1. Plug the end of the device plug of the power supply cable into the device **2** and the power supply plug into a socket. The readiness indicator **9** flashes green.
2. Press the ON/OFF button  **1**. The intensity display **7** shows level 1, (blinking) and the device fan switches itself on. The device is now ready for positioning on the skin.
3. There are two ways for choosing the intensity:
 - a. Set the intensity that you have tested and which is suitable for your skin, by pressing the ON/OFF button  **1** briefly in succession so long, until the desired intensity (Level 1 to 4) is lit on the intensity display **7**.
 - b. or: Place the light area **5** on the skin directly. By taking use of the build-in skin tone sensor, the device can automatically choose the safest intensity level. The readiness indicator **9** illuminates green permanently and the (automatically chosen) level is shown in the intensity display **7** (“energy”). If intensity display **7** is still blinking, skin tone of treated area is too dark. Choose another area, if possible.
4. The light area **5** must lie fully on the skin, so that a light pulse can be emitted by the device. Press the trigger button **8** to trigger a light pulse. When the pulse is emitted, you will see a bright flash and may hear a clicking sound.
5. Now place the device on another area of the skin while it charges for the next light pulse.
6. Then emit a new light pulse on the other skin area.
7. Repeat this procedure for other planned skin areas.
8. Switch the device off by pressing and holding the ON/OFF button  **1**. Readiness indicator **9** will flash green and cooling fan will stop immediately.

**WARNING**

- **Never emit several light pulses in succession onto the same area of skin!**
- **If pain, swelling or other unusual reactions occur during the application, immediately stop the test and do not carry on using the device. Contact your doctor.**

4.5 Display / Intensity display

The display of **IPL 850** features two lines. First/upper line (“Energy”) shows intensity and readiness (flashing or non flashing number). Second line (“Flash”) shows the number of emitted flashes. It will refresh every 1,000 flashes. Here a value of “001” stands for 1000 flashes. If the number of flashes passes “180” (180,000) the display will show 

That means: The cartridge will soon be depleted and needs to be replaced. (-> chapter 5.2 Replacing the cartridge)

While transferring the device to other skin areas you can proceed in two different ways:

Method 1: “Gliding and flashing”

Keep the trigger button **8** pressed down. The device emits a light pulse in succession as long as the light area **5** is in full contact with the skin. After each light pulse slide the device to the next area. Keep the light area **5** pressed onto the skin. This mode is especially suited to treating larger areas, such as the legs, for example.

Method 2: “Stopping and flashing”

Release the trigger button **8** after each light pulse. For the next light pulse position the device on the area directly next to the area that has just been treated. Make sure that the light area **5** and the skin contact sensors **4** lie fully on the skin. This mode guarantees precise treatment for example of the knees or ankles.

**WARNING**

If a device fan fault should occur during the application, immediately stop use of the device, otherwise skin injuries may be caused by the high heat. Contact the customer service.

After switching the device off take the power supply plug out of the socket and disconnect the device plug end of the power supply cable from the device. Then clean the device in accordance with section 5.1 “*Cleaning, care and storage*” in these operating instructions.

NOTE

The device will be switched off automatically after approx. 3 minutes. To continue the application, press the ON-/OFF-button  ❶ and select the desired intensity.

**4.6
Instructions
for use on
face**

This device may not be used on men’s faces and in women only below the cheekbones. Do not use the device on the mucous membranes in the ear and nose region. As the facial skin is very sensitive, we recommend only use Level 2 as a maximum intensity here. When applying to the upper lip, cheek, jaw, chin or neck use the utmost care. You could for example move your lips inwards during the application or press them together, to tauten the treatment area. If necessary you could apply white eyeliner pencil to the lips and lip edges - this will ensure that the energy of the flash is reflected and does not reach your lips, if you accidentally aim a pulse at the lips.

**4.7
Post-
treatment**

After the treatment the skin could possibly be slightly red, tingle or feel warm. This is however completely harmless and should abate in a short time. Nevertheless you can still take simple steps to relax your skin after the application of the **IPL 850** and to accelerate the regeneration.

You can apply a soothing or moisturising and non-perfumed lotion, e.g. aloe vera, to the treated areas. You can also cool down the area by applying a wet wash cloth or an icepack for a few minutes. However, avoid the application of peels, bleaching creams, deodorants etc for at least 48 hours after the treatment. Also direct sunlight and artificial tanning should be avoided for at least 48 hours. After this time the treated skin should be treated with a sun protection cream light protection factor 30, if the skin will be exposed to direct sunlight within the following 2 weeks. Avoid removing hair from the treated areas with wax or by plucking.

5.1 Cleaning, care and storage

The device must be cleaned after each use. Do not use aggressive cleaning agents or strong brushes so that the surface is not damaged. Never immerse the device or parts of it in water and never rinse them under running water! To clean properly proceed as follows:

- Make sure that the device has been disconnected from the power grid.
- Clean the housing with a clean, lint-free and dry cloth.
- Moisten a soft, lint-free cloth with a few drops of water and clean all sensors and the light area **5** thoroughly.
- Store the device in a dry and dust-free location, ideally in the original packaging according to the storage conditions in chapter 5.5 “*Specifications*”.

5.2 Replacing the cartridge

1 cartridge is included in the delivery scope with a lifetime of up to 200,000 light pulses. In case of normal use in accordance with these instructions one cartridge is enough for an individual user for several years of use. If spots are visible on the cartridge, no or hardly any light emissions take place or the lens is broken, the cartridge should be replaced. Switch the device off and take the power supply plug out of the socket. Pull out the old cartridge and insert the new one.

5.3 Trouble- shooting

Here you will find a list of possible faults and how to fix them.
If this list should not help you please contact the customer service.

Problem: Device cannot be switched on.

Fix: Make sure that the device plug end of the power supply cable is connected to the device and the power supply plug is plugged into a socket. Press the ON/OFF button  ①.

Problem: No light pulse emission takes place when actuating the trigger button ⑧.

Fix: The readiness indicator ⑨ should light up in green permanently before the trigger button ⑧ is pressed. This only happens if the skin contact with the light area ⑤ is sufficient. If intensity display ⑦ is blinking, skin tone of treated area is too dark. Choose another area, if possible.

Problem: There is sufficient skin contact and skin type is I to VII but still the readiness indicator does not flash green and it is not possible to trigger the light pulse.

Fix: If the readiness indicator ⑨ flashes or is lit in another colour, possibly there may be a system fault. Switch the device off and disconnect it from the power supply. After about 1 minute restore the power supply and switch the device back on. If it still does not work contact customer service.

5.4 Disposal

This product must not be disposed together with the domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

5.5	Name and model	: MEDISANA Hair Removal System IPL 850
Specifica- tions	Power supply	: 100-240V~ 50/60Hz 1 A
	Adaptor	: SW24C-12002000-EU
	Technology	: IPL (Intense Pulsed Light for home use)
	Max. Energy level	: max. 5 joules/cm ²
	Light area	: 2.45 cm ²
	Wavelength	: 470 - 1200 nm
	Selectable	
	intensity levels	: 4 (6 in automatic mode)
	Flash duration	: 6400 - 7600 µsec
	Operating	
	conditions	: +10°C to +35°C; rel. humidity: 3 % to 75 %; Pressure: 700 to 1060 hPa
	Storage conditions	: -35°C to +70°C; rel. humidity: 10 % to 90 %; Pressure: 500 to 1060 hPa
	Dimensions	: approx. 14 x 8.5 x 4.5 cm
	Weight	: approx. 260 g
Item number	: 88588	
EAN number	: 4015588 88588 4	
Accessories	: Spare cartridge (capacity: approx. 200,000 light pulses) item no. 88589, EAN 4015588 88589 1	

In accordance with our policy of continual product improvement, we reserve the right for technical and design changes without notice.

The current version of this instruction manual can be found under
www.medisana.com

6.1 Warranty and repair terms

Please contact your supplier or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

1. The warranty period for **MEDISANA** products is 3 years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be corrected free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

MEDISANA AG

Jagenbergstraße 19, 41468 NEUSS
GERMANY

E-Mail: info@medisana.de

Internet: www.medisana.com

The service centre address is
shown on the attached leaflet.

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