# Silk'n®

Beautiful Technology

## SKINVIVID

### TB-1389

EN USER MANUAL

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IT MANUALE D'USO

EN

DE

FR

NL

ES

IT

## **Legal notice**

Home Skinovations Ltd. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Home Skinovations Ltd. is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd. Data is subject to change without notification.

Home Skinovations Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Home Skinovations Ltd. Specifications are subject to change without notice.

Silk'n and the Silk'n logo are registered trademarks of Home Skinovations Ltd.

Invention Works B.V.


Donk IB, 2991 LE Barendrecht, The Netherlands

[www.silkn.eu](http://www.silkn.eu) – [info@silkn.eu](mailto:info@silkn.eu)

# Table of contents

|            |  |           |
|------------|--|-----------|
| <b>1.</b>  | <b>Overview</b>                                  | <b>4</b>  |
| <b>2.</b>  | <b>Package contents/device parts</b>             | <b>5</b>  |
| <b>3.</b>  | <b>General information</b>                       | <b>5</b>  |
| 3.1        | Reading and storing the user manual              | 5         |
| 3.2        | Explanation of symbols                           | 6         |
| <b>4.</b>  | <b>Safety</b>                                    | <b>7</b>  |
| 4.1        | Proper use                                       | 7         |
| 4.2        | Safety instructions                              | 8         |
| <b>5.</b>  | <b>First use</b>                                 | <b>15</b> |
| 5.1        | Checking the SkinVivid and package contents      | 15        |
| 5.2        | Charging the SkinVivid                           | 16        |
| <b>6.</b>  | <b>Preparation</b>                               | <b>16</b> |
| <b>7.</b>  | <b>Massage mode (hot, cold, sonic vibration)</b> | <b>17</b> |
| 7.1        | Using the hot massage                            | 17        |
| 7.2        | Using the cold massage                           | 17        |
| <b>8.</b>  | <b>Operation</b>                                 | <b>17</b> |
| <b>9.</b>  | <b>Tips for daily use</b>                        | <b>19</b> |
| 9.1        | Morning ritual                                   | 19        |
| 9.2        | Evening ritual                                   | 19        |
| <b>10.</b> | <b>Cleaning</b>                                  | <b>19</b> |
| <b>11.</b> | <b>Troubleshooting</b>                           | <b>20</b> |
| <b>12.</b> | <b>Technical data</b>                            | <b>21</b> |
| <b>13.</b> | <b>Declaration of conformity</b>                 | <b>21</b> |
| <b>14.</b> | <b>Disposal</b>                                  | <b>21</b> |
| 14.1       | Disposing of the packaging                       | 21        |
| 14.2       | Disposing of the SkinVivid                       | 22        |
| <b>15.</b> | <b>Warranty</b>                                  | <b>22</b> |
| <b>16.</b> | <b>Customer service</b>                          | <b>23</b> |

## Overview
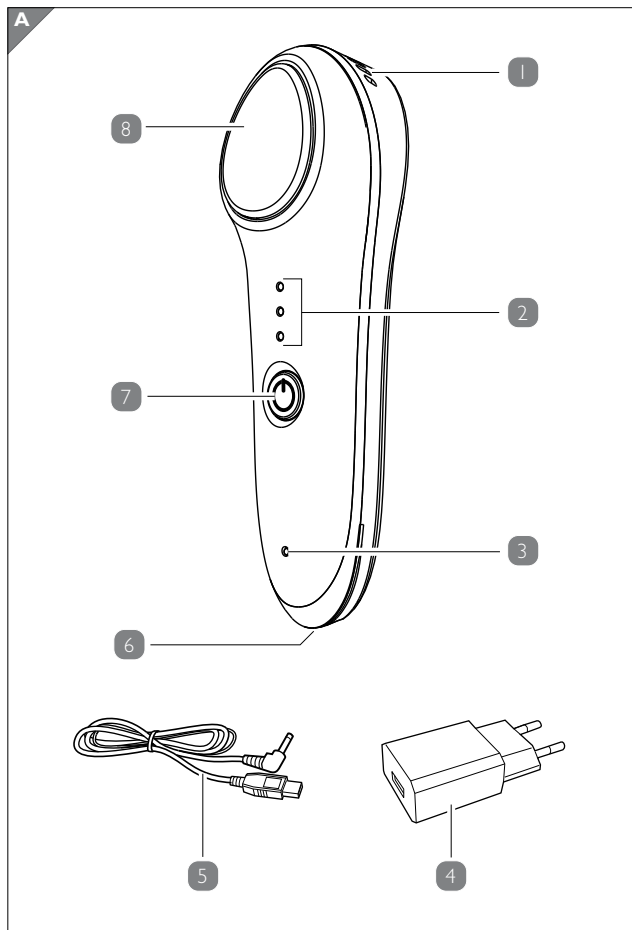
## 2. Package contents/device parts

- 1 Vent
- 2 Indicator for the massage mode
  - Indicator for hot massage
  - Indicator for cold massage
  - Indicator for sonic vibration
- 3 Charging indicator
- 4 USB mains adapter
- 5 USB charging cable
- 6 Charge connector
- 7 Control knob
- 8 Massage head

## 3. General information

### 3.1 Reading and storing the user manual

This user manual accompanies the SkinVivid facial massager (hereinafter referred to only as the "SkinVivid"). It contains important information about start-up and handling.

Before using the SkinVivid, read the user manual carefully. This particularly applies to the safety instructions. Failure to comply with this user manual may result in severe injury or damage to the SkinVivid.

The user manual is based on the standards and rules in force in the European Union. When abroad, you must also observe country-specific guidelines and laws.

Store the user manual for future use. If you pass the SkinVivid on to third parties, please be absolutely sure to include this user manual.

## 3.2 Explanation of symbols

The following symbols and signal words are used in this user manual, on the SkinVivid or on the packaging.

### **WARNING!**

This signal symbol/word designates a hazard with moderate degree of risk which may lead to death or severe injury if not avoided.

### **CAUTION!**

This signal symbol/word designates a hazard with low risk that, if not avoided, may result in minor or moderate injury.

### **NOTICE!**

This signal word warns against potential damages to property.

This symbol provides you with useful supplementary information about use.

Declaration of conformity (see chapter “Declaration of conformity”): Products marked with this symbol meet all applicable Community regulations for the European Economic Area.

Do not use the SkinVivid near bath tubs, shower trays or other vessels containing water.

This symbol identifies electrical devices that belong to the protection class II (mains adapter).

Electrical devices marked with this symbol are only suitable for use in dry indoor rooms.

This symbol identifies voltage and the current rating for direct current.

The word "GS" stands for tested safety. Products marked with this symbol meet the requirements of the German Product Safety Act (ProdSG).

**IP20**

The USB wall adapter is protected against solid objects up to 12 mm.

The international efficiency mark indicates that the power supply meets the level VI requirements (USB wall adapter).

The USB wall adapter is equipped with a short-circuit proof safety transformer.

This symbol indicates the polarity of the pins in the USB wall adapter.

The main adapter has a switching power supply.

**ta**

Ambient temperature (USB wall adapter)

## 4. Safety

### 4.1 Proper use

The SkinVivid is only intended for smoothing, strengthening and improving the elasticity of the skin using hot and cold massage. The SkinVivid is suitable for all skin types. Pregnant persons as well as persons with burns, acne, allergies or skin that is otherwise irritated must not use the SkinVivid.

The SkinVivid is intended for private use only and is not suitable for commercial purposes, e.g. for cosmetic salons. Do not expose the SkinVivid to weather conditions and do not use it in wet rooms. The SkinVivid is not a children's toy.

Only use the SkinVivid as described in this user manual. Any other use is considered improper and may result in damage to property.

Silk'n accepts no liability for damage caused by improper or incorrect use.

## 4.2 Safety instructions

### **WARNING!**

#### **Risk of electric shock!**

A faulty electrical installation or excessive mains voltage may result in an electric shock.

- Only connect the SkinVivid if the mains voltage of the socket corresponds to the rating plate specification.
- Only connect the SkinVivid to an easily accessible socket so that you can quickly disconnect it from the mains supply in the event of a fault.
- Do not operate the SkinVivid if it is visibly damaged or the USB charging cable or the USB mains adapter is defective.

- If the charging cable of the SkinVivid is damaged, it must be replaced with a USB charging cable with the same specifications to avoid risks. Please contact the specialist shop or customer service.
- Do not open the housing; instead, have a qualified professional perform the repairs. Contact a specialist workshop or customer service.
- Liability and warranty claims are waived in the event of repairs performed by the user, improper connection of the device or incorrect operation.
- Only parts that comply with the original device information may be used for repairs. Electrical and mechanical parts, which are essential for providing protection against sources of danger, are located in the SkinVivid.
- Do not immerse the SkinVivid or the USB charging cable or USB mains adapter in water or other liquids.
- Never touch the USB mains adapter with damp hands.

- Never pull the SkinVivid out of the socket by the USB charging cable; instead, always do so by taking hold of the USB mains adapter.
- Keep the SkinVivid, the USB mains adapter and the USB charging cable away from open flames and hot surfaces.
- Do not kink the charging cable and do not lay it over sharp edges.
- Only use the SkinVivid indoors.
- Do not store the SkinVivid in wet rooms and do not in any case expose it to the elements such as rain.
- Never store the SkinVivid such that there is a risk of it falling into a tub or a basin.
- Never reach to retrieve an electrical device connected to the mains supply if it has fallen in water. In such a case, immediately disconnect the USB mains adapter from the socket.
- Never touch the USB mains adapter with damp hands.
- Once the SkinVivid is completely charged or if you would like to clean it, leave the Skin-

Vivid switched off and pull the mains adapter out of the socket.

- Do not go too close to strong magnetic fields with the SkinVivid.
- Do not use the SkinVivid while charging.

## **WARNING!**

### **Risk of explosion and fire!**

Improper handling of the rechargeable battery poses an increased risk of explosion and fire.

- Never expose the SkinVivid to heat such as direct sunlight, open flames or similar.
- Once the SkinVivid is fully charged, pull the USB mains adapter out of the socket.

## **WARNING!**

**Danger for children and persons with impaired physical, sensory or mental capacities (e.g. partially disabled persons, older persons with reduced physical and mental capacities) or lack of experience and knowledge (e.g. older children).**

The following groups of persons must not use the SkinVivid: Children, persons under 16 years of age, adults with impaired physical, sensory or mental abilities or persons lacking experience handling the device.

- Cleaning and user maintenance must not be performed by unsupervised children.
- Keep the SkinVivid and its accessories away from children under 16 years of age.
- The package contents include the packaging wrapper. Keep it away from children as there is a risk of choking if swallowed.

### **CAUTION! Risk of chemical burns!**

Leaked fluid from the rechargeable battery could cause injuries if it comes into contact with skin or other body parts.

- If the rechargeable battery in the SkinVivid has leaked, put on protective gloves and ensure that the fluid from the rechargeable battery does not come into contact with your skin, eyes or mucous membranes.

- In the event of contact with acid from the rechargeable battery, immediately flush the applicable areas with plenty of clean water and consult a physician immediately.

### **CAUTION! Risk of injury!**

Handling the SkinVivid improperly may result in injury.

- Do not use the SkinVivid if you:
  - are pregnant,
  - suffer from neurological diseases,
  - have undergone cosmetic surgery in the area of the face,
  - are allergic to cold, heat or soft rubber,
  - have metal, silicon or plastic implants due to fractures or other plastic surgery in the area of the face,
  - are unable to feel heat or are not sensitive to warmth.
- Do not use the SkinVivid if you suffer from the following skin irritation or diseases:

- allergies,
  - acne,
  - eczema,
  - irritated skin, reddening or irritation (including a sun burn),
  - swelling in the area of the face,
  - wounds around the face,
  - psoriasis,
  - neurodermatitis.
- Do not use the SkinVivid near open wounds or to remove warts.
  - Immediately discontinue treatment if you experience pain or skin irritation.
  - Do not continuously treat your skin with the SkinVivid for more than two to three seconds.
  - If you are not sure whether SkinVivid is suitable for you, ask your physician before using it for the first time.

## **NOTICE! Risk of damage!**

Improper handling of the SkinVivid may result in damage to the SkinVivid.

- Protect the SkinVivid from shocks, extreme temperatures and direct sunlight.
- Only use the SkinVivid with the original accessories.
- Keep the SkinVivid away from corrosive liquids.

## **5. First use**

### **5.1 Checking the SkinVivid and package contents**

## **NOTICE! Risk of damage!**

If you are not cautious when opening the packaging with a sharp knife or other pointed objects, you may quickly damage the SkinVivid.

- Be very careful when opening the packaging.
  1. Take the SkinVivid out of the packaging.
  2. Check to make sure that the delivery is complete (see **Fig. A**).
  3. Check whether the SkinVivid or the individual parts exhibit damage. If this is the case, do not use the SkinVivid. Please contact the local Silk'n service centre (see chapter "Warranty").

## EN 5.2 Charging the SkinVivid

It is normal for the SkinVivid to become warm while charging.

1. Make sure that the SkinVivid is switched off.
2. Completely unwind the USB charging cable **5**.
3. Insert the jack plug of the USB charging cable in the charge connector **6** of the SkinVivid.
4. Insert the USB plug of the USB charging cable in the USB jack of the USB mains adapter **4**.
5. Insert the USB mains adapter in a properly installed socket.

The charging indicator **3** is illuminated in red during the charging process. The charging time for fully charging the rechargeable battery of the SkinVivid amounts to approx. 3 hours.

6. Once the rechargeable battery is fully charged, the charging indicator will illuminate in blue. Disconnect the SkinVivid from the mains supply.

The rechargeable battery can now be operated for 40 minutes.

## 6. Preparation

1. Before using the SkinVivid, clean your face and the parts of your skin you intend to treat with your usual cleaning products.
2. If necessary, groom your face and the parts of your skin you intend to treat with the usual care lotions and moisture creams.
3. If you would like these skin care products to absorb faster, use the hot mode with sonic vibration after applying the cream (see chapter "Using the hot massage" and "Operation").

## 7. **Massage mode (hot, cold, sonic vibration)**

Every mode you select runs automatically for 2 ½ minutes. The SkinVivid switches off afterwards.

### 7.1 **Using the hot massage**

The SkinVivid's hot massage at 42 °C promotes circulation in the skin and allows skin care products absorb faster for a smoother complexion and improved moisture balance. Also switching sonic vibration on will enhance the effect of the hot massage.

Ideally, you should use this mode after your daily facial cleanse.

### 7.2 **Using the cold massage**

The cold massage firms up your skin and reduces pore size at 6 °C. In this mode, the SkinVivid also relieves tense muscles, reduces puffy eyes and wrinkles. Also switching on sonic vibration will enhance the effect of the cold massage.

Ideally, you should use this mode after your daily facial cleanse and complete absorption of your skin care products and before you apply your make-up.

For smoother, more radiant and more supple skin, we recommend alternating between the hot and cold massage.

## 8. **Operation**

The SkinVivid cannot be used while charging.

1. To switch the SkinVivid on, hold the control knob **7** down for 3 seconds.

A beep sound will be emitted, the top light of the indicator for

the massage mode **2** will illuminate in red and indicate that the SkinVivid is in the mode for the hot massage.

2. Push the control knob again to additionally switch on the sound vibration.

A beep tone will sound and, in addition to the top light, the bottom light of the indicator for the massage mode will illuminate in yellow at the same time. The SkinVivid is in the mode for hot massage with sonic vibration.

3. To switch to the mode for cold massage, push the control knob again.

A beep sound will be emitted, the light in the middle of the indicator for the massage mode will illuminate in blue and indicate that the SkinVivid is in the mode for the cold massage.

4. Push the control knob again to additionally switch on the sound vibration.

A beep tone will sound and, in addition to the middle light, the bottom light of the indicator for the massage mode will illuminate in yellow at the same time. The SkinVivid is in the mode for cold massage with sonic vibration.

5. Move the SkinVivid up to your face and place the massage head **8** on the desired part of your face.
6. Depending on the part of your skin you are treating, allow the SkinVivid to glide in circular movements from the inside to the outside or from top to bottom.
7. To switch the SkinVivid off, hold the control knob down for 3 seconds.

A beep tone will be emitted and the lights in the indicator for the massage mode will go out.

It is recommended to fully recharge the SkinVivid after 5 uses for maximum power.

## 9. Tips for daily use

Regardless whether you have dry, oily or sensitive skin, Silk'n SkinVivid is suitable for any type of skin.

### 9.1 Morning ritual

You can use the hot massage e.g. in the morning after cleaning your face. Then apply your daily care products. Once they have absorbed well, use the cold massage, which e.g. is perfect for relieving puffy eyes and is suitable for creating a smooth base to apply make-up to.

### 9.2 Evening ritual

Or opt for a complete evening beauty ritual after a stressful day. After cleaning your skin, first use the hot massage to open your pores. Then apply eye gel, a nurturing mask or moisturiser cream. Once your care products have been absorbed, use the cold massage to close your pores and rejuvenate your skin.

## 10. Cleaning

### **NOTICE! Risk of damage!**

Using the wrong cleaners and cleaning utensils could damage the SkinVivid.

- Do not use any hard, scratching or scouring cleaners and cleaning utensils.

- Do not immerse the SkinVivid or the USB charging cable or USB mains adapter in water or other liquids and never put any of the parts in the dishwasher.
- Clean the massage head  and the SkinVivid with the enclosed cleaning cloth after each use.
- Do not clean the SkinVivid with water.

## II. Troubleshooting

In the event of an error, see if the following steps resolve the problem:

| Fault                         | Possible causes and solution  |
|-------------------------------|---|
| The SkinVivid doesn't work.   | <p>The rechargeable battery may be empty or not adequately charged.</p> <ul style="list-style-type: none"> <li>– Fully charge the SkinVivid before using it (see chapter "Charging the SkinVivid").</li> </ul>          |
| There are no visible effects. | <p>You have forgotten to clean or groom the skin of your face before use.</p> <ul style="list-style-type: none"> <li>– Clean and groom your skin according to the instructions in the chapter "Preparation".</li> </ul> |
|                               | <p>You have not been using the SkinVivid for very long or not on a regular basis.</p> <ul style="list-style-type: none"> <li>– Use the SkinVivid regularly.</li> </ul>  |

If these steps are not successful or if your SkinVivid has to be repaired, contact our customer service team.

## 12. Technical data

Model: TB-I389

### Adapter

Model: PS06H050K1000EU  
 Input voltage: 100–240 V~, 50/60 Hz, 0.25 A  
 Output voltage: 5 V  $\overline{=}$  1000 mA  
 Protection class: II

### Hand piece

Rechargeable battery: 5 V  $\overline{=}$  800 mAh  
 Charging time: approx. 3 h  
 Operating time: approx. 40 minutes (on a full charge)  
 Dimensions (L x W x D): approx. 15.2 x 4.8 x 4.3 cm  
 Temperature massage head: hot massage 42 °C  $\pm$  3 °C  
 cold massage 6 °C  $\pm$  3 °C  
 Automatic deactivation: after 2.5 minutes in any mode  
 Weight: approx. 109 g

## 13. Declaration of conformity

The EU declaration of conformity can be requested from the address specified on the back.

## 14. Disposal

### 14.1 Disposing of the packaging

Dispose of the packaging separated into single type materials. Dispose of cardboard and carton as waste paper and film via the recyclable material collection service.

## EN 14.2 Disposing of the SkinVivid

(Applicable in the European Union and other European states with systems for the separate collection of reusable waste materials)

### **Old devices must not be disposed of with household waste!**

If the SkinVivid is no longer operational, every consumer is **required by law to dispose of used devices separately from household waste**, e.g. at a collection point run by the communal authority or borough. This ensures that old devices are recycled in a professional manner and also rules out negative consequences for the environment. For this reason, electrical equipment is marked with the symbol shown here.

### **Batteries and rechargeable batteries do not belong in the household waste!**

As the end user you are required by law to bring all batteries and rechargeable batteries, regardless whether they contain harmful substances\* or not, to a collection point run by the communal authority or borough or to a retailer, so that they can be disposed of in an environmentally friendly manner.

\*labelled with: Cd = cadmium, Hg = mercury, Pb = lead

## 15. Warranty

This product is covered by a 2-year warranty according to European regulations and law. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may

be able to solve your problem without needing to return the product to the store or our service center. Our customer service is always pleased to help you!

## **16. Customer service**

For more information about the products of Silk'n visit your regional Silk'n website: [www.silkn.eu](http://www.silkn.eu). If the device is damaged, defective, in need of repair or you need our assistance, please contact your nearest Silk'n service center. This manual is also available as a PDF download from: [www.silkn.eu](http://www.silkn.eu).

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Copyright 2020 © Home Skinovations Ltd. Tutti i diritti riservati.
All rights reserved. SILK'N and the SILK'N logo are registered trademarks of Home Skinovations Ltd.
InnoEssentials Deutschland GmbH
© 2020 Home Skinovations Ltd.