Smart Personenweegschaal Smart Body Fat Scale

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
NL

p.2-10

User Manual

ENG

p.11-18

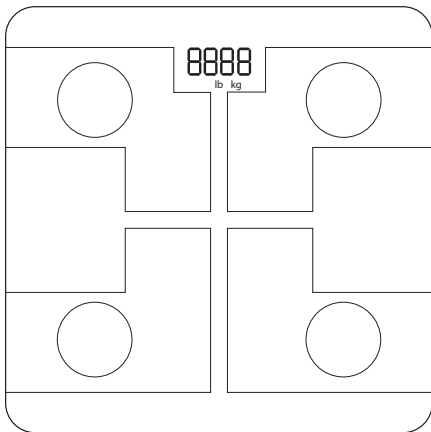RoHS  CE

MADE IN CHINA

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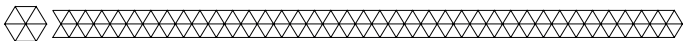
Introduction

- Thank you for choosing our LifeGoods product! To ensure proper and safe functioning of the product, please read the user manual carefully before use. Please always follow the instructions for safe usage. If you have any questions about the product or experience any problems, please contact us by email: service@lifegoods.nl
- The LifeGoods Smart Body Fat Scale is equipped with BIA technology and capable of 11 different body measurements. Turn on your bluetooth, download the Feelfit app, create an account and keep track of your health! The scale is compatible with iOS 7.0/Android 4.3 or newer devices. With a max. load bearing capacity of 180 kg and including batteries.

Warnings and Safety Instructions

Please ensure that the packaging is kept out of the reach of children!
Risk of suffocation!

Make sure to keep this user manual. Before using your Smart Body Fat Scale, read the following:

General

- If you have not used the LifeGoods Smart Body Fat Scale before, it is important that you read and follow the instructions as stated in the user manual. Be aware of the risks and consequences resulting from misuse of the Smart Body Fat Scale.
- Misuse of the Smart Body Fat Scale can result in personal injury, or destruction of the instrument/appliance/device and of the installations.
- Use the Smart Body Fat Scale only for purposes described in this user manual.
- If the Smart Body Fat Scale malfunctions, is dropped or (appears to be) damaged in any manner, do NOT use the device and contact Customer Service.
- The Smart Body Fat Scale is designed and intended for private, non-commercial use only.
- The Feelfit App is a third party software and is property of en developed by Shenzhen QingNiu Software Co. Ltd. LifeGoods bv is not responsible for the services of the app and its privacy and cookie policy.
- Always weighs yourself at the same time of day and under similar circumstances to ensure that the results are comparable.

- The Smart Body Fat Scale has a max. load bearing capacity of 180 kg (400 lbs.). Exceeding this weight limit can cause damages to the device and personal injuries.

Persons and use

- This device is equipped with BIA (bioelectrical impedance analysis) to measure the body fat. persons (including children under the age of 18) with physical, sensory or mental disabilities, or persons (including children) with lack of experience with and knowledge of the device should NOT measure body fat. Measuring weight is allowed ONLY if they are supervised or have been given instruction concerning the use by a person responsible for their safety.
- Do NOT use the device or consult your physician: If you are pregnant; If you have medical implants (such as a pacemaker). The Smart Body Fat Scale can disrupt the working of the implant; During or immediately after strenuous activities; Immediately after taking a shower, a bath or going to the sauna; When you have consumed alcohol; During or immediately after eating; When you are running a fever; If you have any unexplained pain, if you are undergoing medical treatment or using medical equipment.
- The measurement of body fat can be inaccurate when you have diabetes or another medical or physical condition.
- The measurement of body fat can be inaccurate when you are an athlete
- This device is NOT a substitute for treatments administered by doctors or physiotherapists.

Electronics

- Do NOT reconstruct, repair or disassemble the Smart Body Fat Scale. This can damage the device. For safety reasons, these parts may only be replaced by authorized service personnel. In order to prevent hazards, always send the Smart Body Fat Scale to the service centre for repairs.
- Do NOT use the Smart Body Fat Scale in or near a bath, shower, swimming pool, or over a basin filled with water.
- Make sure to prevent any contact with water. Do NOT submerge the Smart Body Fat Scale in water. If water does get inside the casing, turn the device off immediately and contact our customer service.
- Do NOT operate the device with damp or wet hands
- Keep the Smart Body Fat Scale away from moisture, heat and in or close to areas where the surrounding temperature is relatively hot (e.g. stoves or other sources of heat).
- Do NOT place the device in direct sunlight.

- Leaking batteries can damage the device. Remove the batteries if you are not using the device for an extended period of time.
- Use gloves to remove leaking batteries. Clean the battery compartment with a dry cloth.
- Do NOT expose the batteries to high temperatures.
- Batteries prove fatal when swallowed. Keep the batteries away from children.

Other

- Do NOT stand on the outer edges of the scale. This can cause the scale to tilt over.
- Always place the scale on a flat and hard surface, and only use the device indoors.
- The surface of the scale can be slippery.
- Protect the scale against impacts and shocks.
- Do NOT clean the scale with aggressive, abrasive or foaming cleaning agents.

▶ Receiving your Shipment


- When receiving your shipment, make sure that the contents are consistent with the packing list. Notify your distributor of any missing items. If the equipment appears to be damaged, file a claim immediately with the carrier and notify your distributor at once, giving a detailed description of any damage. Save the damaged packing container to substantiate your claim.

▶ Contents of Package

- LifeGoods Smart body Fat scale - Model LG1004 - EAN 8720195250623

Included Accessories:

- 1x User Manual
- 3x AAA batteries

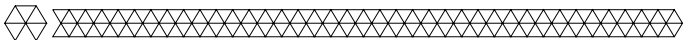▶ Product Features

- Technical Data

Size	26,00 x 26,00 cm
Weight	1200 gr
Colour	Black

Material	ABS - Glass - Plastic
Diameter Scale	280 - 300 mm
Max. Load Bearing Capacity	180kg (400 lbs.)
Measurement Resolution	0,1 kg (0.2 lbs).
Bluetooth	4
Name Bluetooth	#1004
Name App	Feelfit
Compatible with	iOS 7.0/Android 4.3 or newer devices
Synchronized with	Apple Health, Google Fit, FitBit
Input - Output	50mA - 50mA
Power	3x AAA batteries

Other Features

- LED-display
- Displays 11 different body measurements
 - Body Fat
 - BMI
 - Visceral Fat
 - Metabolic Age
- Weight
- BMR
- Protein
- Body Water
- Bone Mass
- Muscle Mass
- Skeletal Muscle


Operating Instructions

- Overview

6800 x kg	No display	6800 kg
Normal	Change Battery	No Connection
Lo	EEr	CAL
Change Battery	Overload	Calibrate

- Power on and Measuring

1. Put the AAA batteries (included) in the battery compartment and turn on the Smart Body Fat Scale.
2. Check if the scale displays 0.0lb before stepping on the scale.
3. The units of measurement can be changed from pounds (lb) to kilograms (kg), or the other way around. The height measurement units (meter vs. inch/foot) will change accordingly.
4. Measure body fat and other body composition parameters correctly: have both feet uniformly touch the electrodes as shown in the figure below:

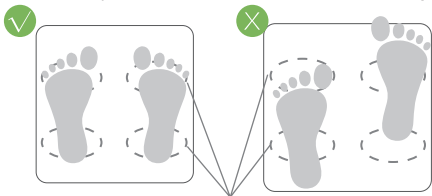Conductive areas

5. Your weight will be measured. Download the app to measure other parameters.

- Installing App




The Smart Body Fat Scale can connect wirelessly to mobile phones or tablets when you install the Feelfit app.

1. Find and download the Feelfit app in the App Store ((iPhone/iPad) or the Play Store (Android).
2. Make sure to turn on your bluetooth before using the app and to put the batteries in the scale. When the bluetooth is turned on, connect with the name #1004 to connect to the Smart Body Fat Scale.

- Register with Feelfit

1. Choose your preferred unit of measurement (kg/lb/st - cm/inch).
2. Create an account.
Use your email address to register. The app can realize synchronous cloud storage of data and synchronizes with other devices when logging in.
3. Fill in your personal information

Note: See next page for visuals

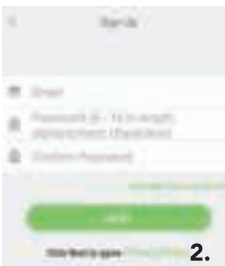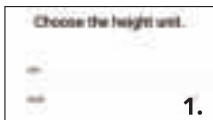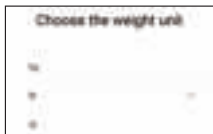
1.

2.

3.

- Add a device

1. Click on the "+" symbol in the top right corner to access Device Management.
2. Step on the scale when your device is paired and the scale will power on. Click on the CS20M symbol on the interface to add the body fat scale.

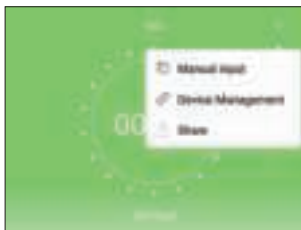1. All Data

Step on the scale when you've added your device. Make sure to:

- be barefoot with clean feet and stand up straight;
- NOT stand on the outer edges of the scale.

Once the transmission of data is completed all the measured data will be displayed. Click on the data to obtain results.


2. Historic Data

- Click on the Trend icon at the bottom of the screen to view your historical data. You can select a week, month or year.
- Select the type of data you want to view (such as weight or BMI) to view changes in your data.

- Other Options

Click on "My Account" (bottom of the page) to change your settings, view FAQ and add users.

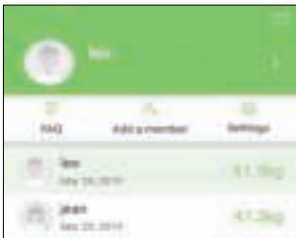1. Settings

Enter the settings interface to change your settings. Android users can connect with Google Fit and FitBit accounts. iPhone and iPad users can connect with their Apple Health Kit account.

2. FAQ


Open FAQ to view frequently asked questions and their answers.

3. Add Users

You can add other users or Feelfit friends by clicking on Add Measuring User. Add the desired personal information and click on Confirm to add the new user. To independently measure the data of each user, select the desired user before measuring. The new user can also download the Feelfit app on their own phone after registering.

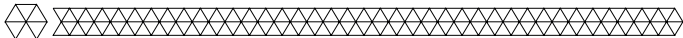Maintenance and Cleaning

- Clean the Smart Body Fat Scale with a moist cloth. Do NOT use aggressive, abrasive, or foaming cleaning agents.
- **Storage**
Store the Smart Body Fat Scale in a clean and dry place. Do not store the device in extreme high or extreme low temperatures.


Disposal - Recycle

- Contribute to a cleaner living environment! Do not dispose of the LifeGoods Body Fat Scale with the household waste. For proper recycling, please take these products to a designated collection point. Please contact your local authority or your household waste disposal service for further details of your nearest designated collection point.

Service and Warranty

- If you need service or information regarding your product, please contact the LifeGoods Customer Support at service@lifegoods.nl
- LifeGoods provides a 2 year warranty on its products. To obtain service during the warranty period, the product needs to be returned with a proof of purchase. Product defects have to be reported within 2 years from the date of purchase. This warranty does not cover: damages caused by misuse or incorrect repair; parts subject to wear; defects that the customer was aware of at the time of purchase; damages or defects caused by customer neglect; damages or defects caused by third parties.

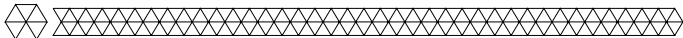Disclaimer

- Changes reserved; specifications are subject to change without stating grounds for doing so.

LifeGoods B.V.

Wisselweg 33
1314CB Almere
Netherlands
06/2020 - V.1