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User Manual  
Light Therapy Lamp

Dear customer,

Thank you for choosing Light Therapy Lamp. This introduction will help you to properly install and use your new lamp. Please read these carefully before using and keep them for future reference.

PACKAGE CONTENTS

- 1x Light Therapy Lamp
- 1x Adapter
- 1x User Manual

SAFETY INSTRUCTIONS

1. Please read the instruction manual before using the lamp.
2. To ensure the longer life, please do not install near heat, corrosive or damp place.
3. Installation at the point that no shock, momentum and fire. Do not let it fall from high place and avoid a crash.
4. Before cleaning, make sure that the lamp is disconnected from the power supply.
5. Before using the lamp, please check if the lamp is damaged during transport. If the lamp has shown any damage, do not use or install the lamp.
6. Please do not disassemble or repair the lamp yourself.


INSTALLATION

Installation Instructions

1. Please check the accessories of the lamp carefully before installation.
2. To ensure the normal work of the lamp, please use the appropriate voltage.


Installation steps

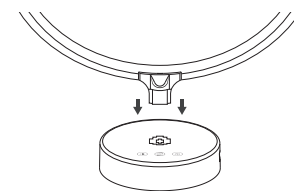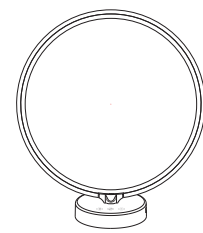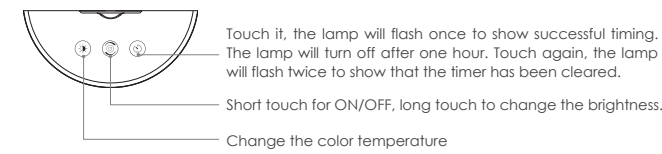
Insert the lamp ring part into the base. Then connect the power supply.

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