




# ULTIMA M7-W

## User Manual



### Parts Name

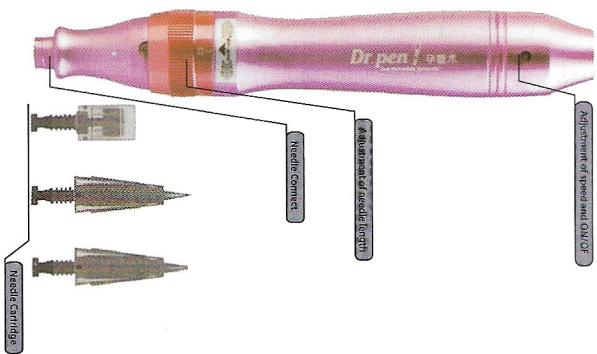

### How to Operate

- Adapler Mode**  
 (1) Plug in DC adaptor to the jack.  
 (2) Plug in DC adaptor to AC outlet.

- Battery Mode**  
 1) Turn on Dr. pen directly with the fully charged battery.  
 2) To charge the battery: Plug the adaptor into an outlet and connect it to charger. The battery indicator light will illuminate in red. Indicator light will turn green when fully charged



- 3) Battery fully charges in about 4 hours.  
 4) Battery can be charged while working.



### How to Operate

- 2. Preparation**  
 (1) Take out needle cartridge in the sterilized packing.  
 (2) Fitting round slots 1 into 2 of machine, push needle cartridge horizontally up to the end.


- (3) Lock the cartridge by turning it to the counter clockwise.  
 (4) Adjustment of needle length should be done after turn on the machine.  
 (5) If you want to remove the needle cartridge, turn it to the clockwise and push out horizontally, fitting round slots(slot1 slot2) together.

- (3) Adjustment of needle length**  
 You can adjust the length of needle by turning the adjustment ring clockwise or counter clockwise. Turning it clockwise, length of it will be longer. Turning it counter clockwise, length of it will be shorter.



### How to Operate

- 3. Adjustment of speed and needle length**  
 (1) Adjustment of Speed and ON/OFF

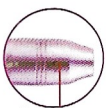

**Power switch**  
 Press 2 seconds to open the device

- (2) 5 speed level control**

Pressing the switch 2 seconds to open the device. It starts from the lowest speed to highest speed when you turn on the power firstly. The speed cycles low to high. Please press the switch 2 seconds when you need to stop the device.  
 The flicker frequency of light faster, the speed is higher

# ULTIMA M7-W

## User Manual



### Maintenance

- The period when the parts of Dr pen is 6 months after stopping of the production of the parts.
  - We offer free repair to our products within 600 hours basis(warranty) from the day they are purchased. If any of our products has damages or natural trouble, we offer repairs or replacements.
  - If they are included in the following problems, you have to pay for the repair even during warranty period
    - Troubles due to a user's careless treatment.
    - Voluntary repair or alteration
    - Trouble due to the repairs by a man other than the technician from a dealer shop or the customer center appointed by manufacturer
    - Troubles or damages due to shocks while in movement
    - Trouble due to the accessories or options not appointed by manufacturer
- 1. General maintenance**
- Keep it with its power off. Wipe its surface out with a dry cloth.
  - Keep it away from direct sunlight or high temperature and humidity.
  - Do not leave the measurement, repair or alter of the device to non-professional staff.
- 2. Routine inspection**
- Make sure to inspect the device and adapter before use
  - When you use the device which has not been used for a long time, make sure to check if it works properly with safety.



### Warranty

<input checked="" type="checkbox"/>	Please be sure to read instructions before use and follow as instructed.
<input checked="" type="checkbox"/>	Do not apply it to the skin as follows: <ol style="list-style-type: none"> <li>When the skin is dry or peeling.</li> <li>When the skin is red or itchy.</li> <li>When the skin is sensitive or itchy.</li> <li>When the skin is damaged.</li> <li>When the skin is infected.</li> <li>When the skin is 12 years old.</li> </ol>
<input checked="" type="checkbox"/>	Do not use it on the power chair than AC 200V.
<input checked="" type="checkbox"/>	Do not use it near the place with high temperature, dust or humidity.
<input checked="" type="checkbox"/>	Do not use it near the flammable or highly volatile materials.
<input checked="" type="checkbox"/>	Do not use it after the controller near the body.
<input checked="" type="checkbox"/>	Do not touch the device directly by your hands, feet and hair.
<input checked="" type="checkbox"/>	Do not use the device when it is wet.
<input checked="" type="checkbox"/>	Do not use the device when it is exposed and it may cause electric shock or fire.
<input checked="" type="checkbox"/>	Do not use the device when it is damaged or it may cause electric shock or fire.
<input checked="" type="checkbox"/>	When you connect a plug, always check if there is the correct plug.
<input checked="" type="checkbox"/>	Do not use the device when it is not in the correct position.
<input checked="" type="checkbox"/>	Do not use the device when it is not in the correct position.

### Packing list

Host	1 pcs
Cartridges	2 pcs 12 pins
Manual	1 pcs
Adapter	1 pcs



### Specification

- 1.Object**  
Apparatus for applying a Micro-needle Therapy by vibration.
- 2. Performance**  
A device for puncturing skin-penetrating drugs or cosmetics supports with a micro-needle guided between extended and retracted positions.
- 3. Classifications**  
-Type of protection against electric shock: Class II  
-Degree of protection against electric shock: Floating  
-Degree of protection against the ingress of water: Ordinary  
-Mode of operation: Continuous operation
- 4. Specification**  
1) Product name: Auto-stamp Motorized Meso Machine  
2) Model name: Ultima M7  
3) Adapter: SVJ-1000MA  
4) Weight: Main body / 78g  
150\*23 (mm)  
5) Size:



### Reference Chart

REGION	THIN SKIN	THICK SKIN
Forehead	0.25-0.5mm	0.5-0.75mm
Between Eyebrow	0.25-0.5mm	0.5-1.0mm
Nose	0.25mm	0.5mm
Eye Area	0.25mm	0.25-0.5mm
Cheek Bone	0.5-1.0mm	0.5-1.0mm
Cheek and Chin	0.25mm	1.0-2.0mm
Upper Lip	1.0-1.25mm	0.25-0.75mm
Facial Scarring	1.5-2.0mm	1.5-2.0mm
Stretch Mark	1.5-2.0mm	1.5-2.0mm
Surgical & Burn Scars		2.0-3.0mm



### Precaution

	Alternating current		Serial No
	Direct current		EU Representative
	Class II equipment		CE mark
	Read the manual		Manufacturer
	Date of manufacture		Finishing applied part