FYSIO PRO MASSAGE GUN

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READ THIS BEFORE USE

WARNING

- Never use this device near water
- Keep the ventilation grille free from dust and dirt
- Do not remove any screws and never try to open the device yourself
- Do not use this appliance for more than one hour at a time rest 30 minutes before using again
- Never charge the device overnight and only under supervision
- Disconnect the charger from the battery before use
- Suitable for adults only. Do not use this device on injuries. Consult a doctor before using this device

WARNING

DANGER: To reduce the risk of electric shock, injury to limit fire, this appliance may only be used in accordance with instructions below

- Only suitable for adults.
- Never leave the FYSIO PRO massage gun unattended while charging.
- Never use this massage gun in any other way than indicated.
- Check the appliance carefully for defects before each use.
- Charge only with the supplied FYSIO PRO charger.
- Never use water near the FYSIO PRO massage gun.
- Keep hair, fingers and other body parts away from the shaft of the massage heads because these can get pinched.
- Bruises may occur when using a massage gun.

Check treated parts of the body regularly and stop immediately in case of pain or discomfort.

- Use the FYSIO PRO massage gun only on soft body tissues, Without discomfort or pain. Don't ever use these on head and bone.
- Use the FYSIO PRO massage gun only on dry skin or clothes in an up and downward movement for about 60 to 90 seconds.

DISCLAIMER READ THIS BEFORE USE

Never use the FYSIO PRO massage gun without approval from your doctor in the cases below:

Pregnancy, diabetes, use of pacemakers, recent surgeries, epilepsy, migraines, hernias, spondylolisthesis, spondylol-ysis, spondylosis, joint replacements, metal plates in the body or other concerned present inc. your physical health. These cases exclude the use of the FYSIO PRO massage is not an option, but advice from a medical expert is needed.


GUARANTEE


If a manufacturing error or defect occurs within two years of purchase then FYSIO PRO will repair or replace the massage gun free of charge. In case of defect please contact reseller.

SPECIFICATIONS

IMPORTANT: The massage gun shuts down after 10 minutes of use a safety measure

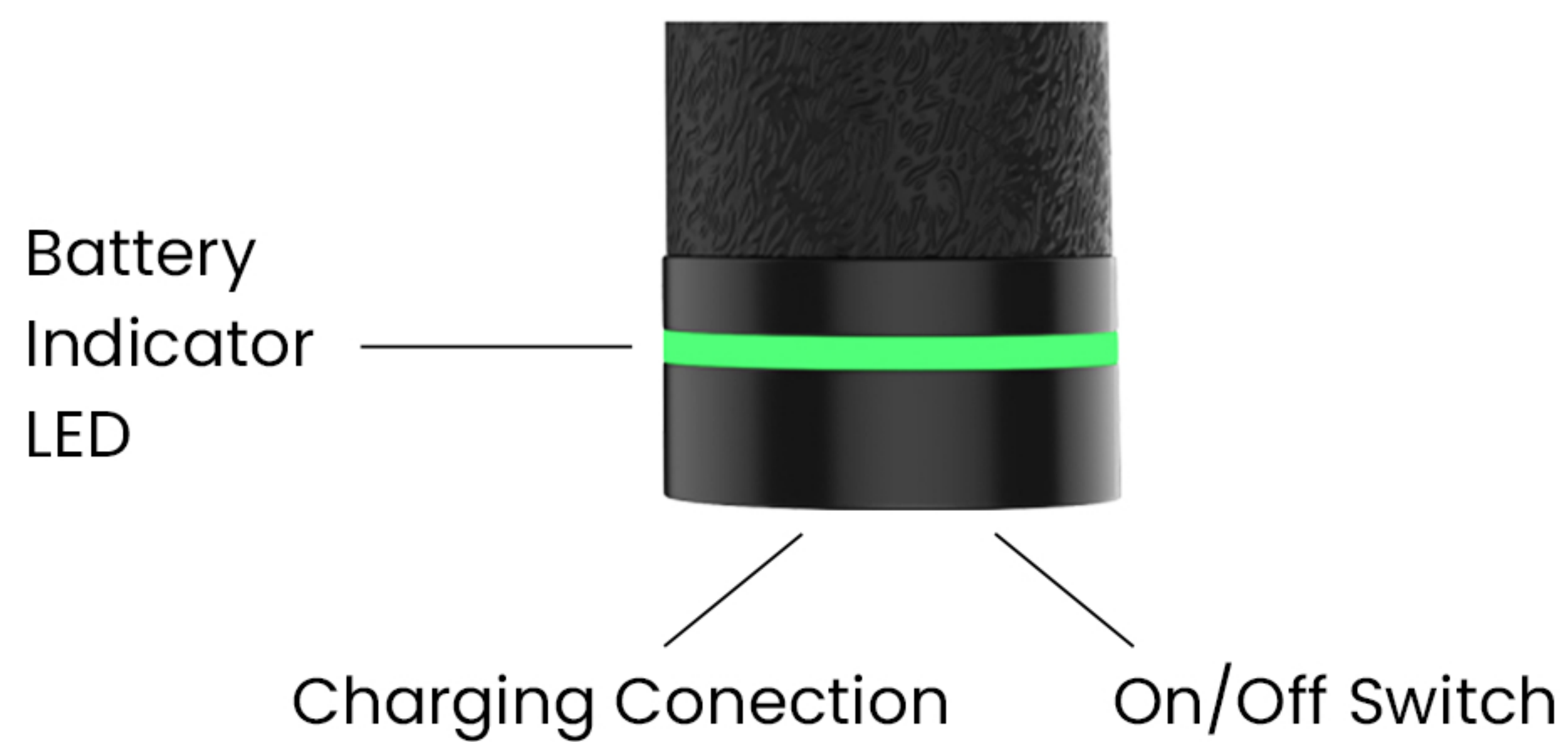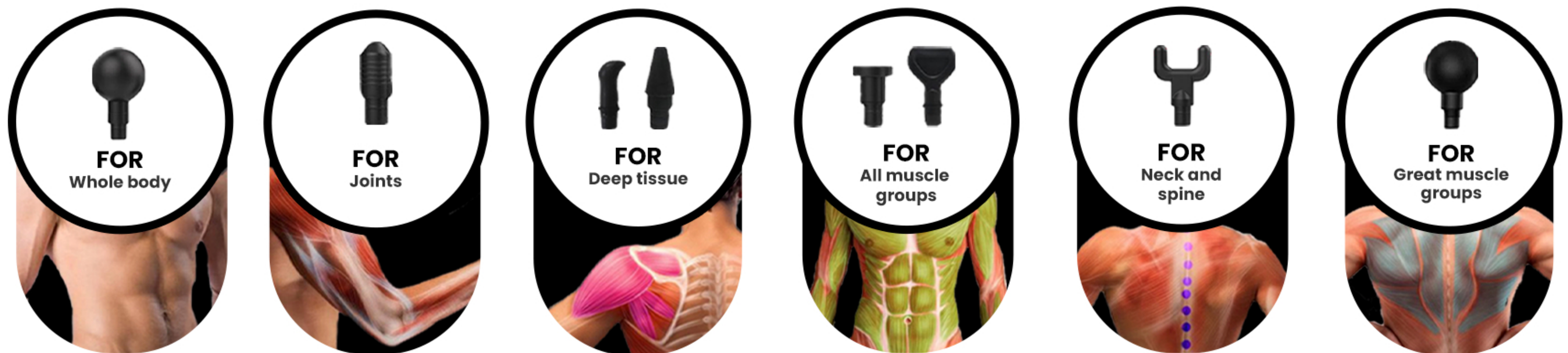MANUAL

- PAY ATTENTION! The massage gun switches itself off after 10 minutes for security
- Charge the battery for up to 2 hours before first use
- The average working time is up to +- 2 hours, depending on the vibration position and the applied pressure
- Attach a massage head before turning on the device. You can push in the head and pull it out after use
- Switch the device on and off by pressing the On/Off button.

TECHNICAL DATA

Brand: FYSIO PRO
Battery capacity: 2500 mAh
Power: 24W
Voltage: 14V
Weight: 2KG
Charging time: 2 hours
Vibration levels: 30

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