# Xclusive

**LIFE STYLE**

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XLS-BLD-307


# Product Overview

Massage head  
(Replaceable)





Gear light

Red breathing light

Power button
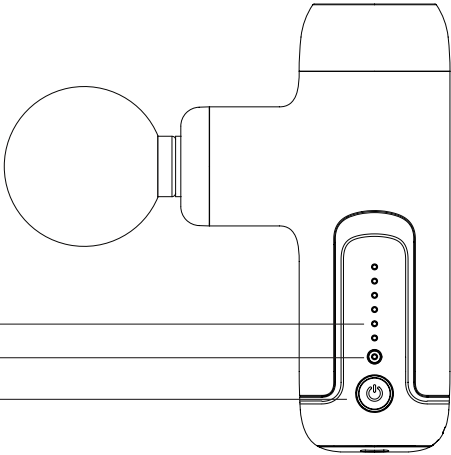
## Replace the massage head

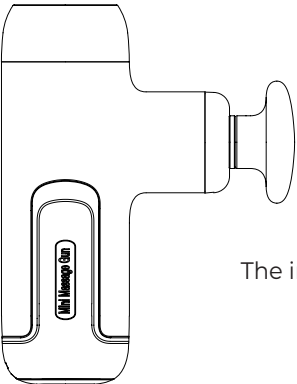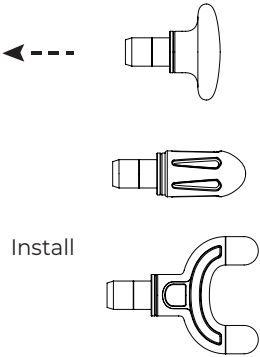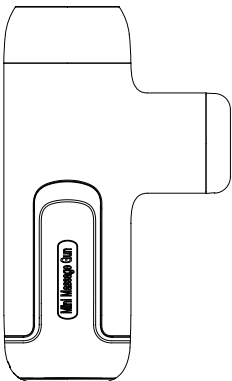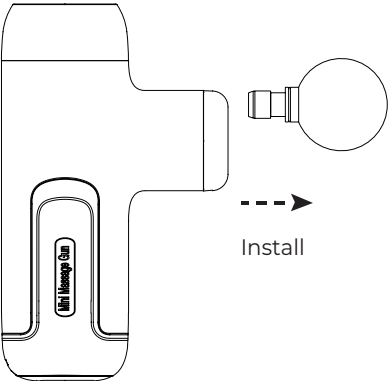
The installation is complete

# Product Use

## ON/OFF

1. Press and hold the power button for three seconds to power ON/OFF.
2. In the power-on state, short press the button once to increase the gear.
3. After the product is turned on, every 10 minutes of operation, It will stop automatically, if you want to continue using it, please just press the switch.
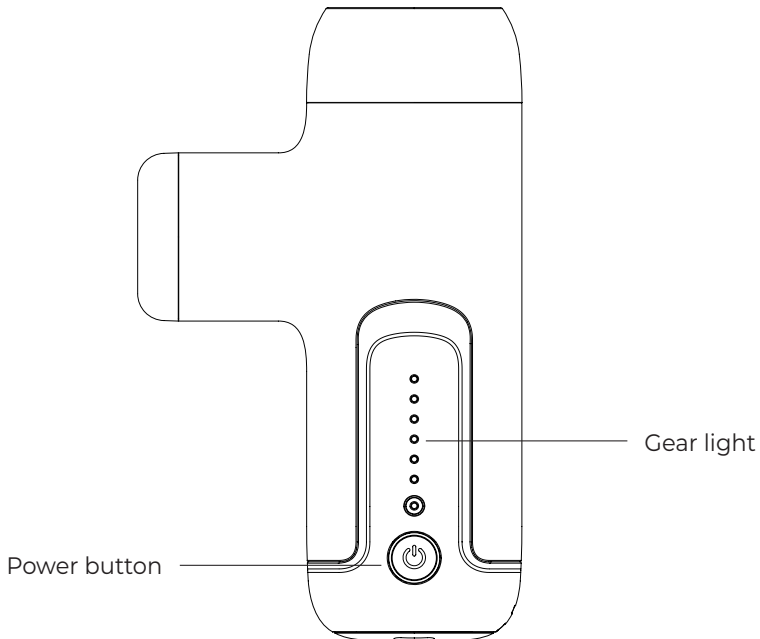
## Strength adjustment


1. The initial gear when power on is the first gear. Click the power button once to increase by one step.
2. The maximum gear is the sixth gear. When the sixth gear is clicked, the power button will change to the first gear.

**Note:** The initial gear when power on is the first gear. Click the power button once to increase by one step.


The maximum gear is the sixth gear. When the sixth gear is clicked, the power button will change to the first gear.

## Gear light, red breathing light





- The gear light indicates the gear size of the fascia gun.
- When one light is on, it is the first gear, when two lights are on, it is the second gear.
- When three lights are on, it is the third gear, when four lights are on, it is the fourth gear. And so on.

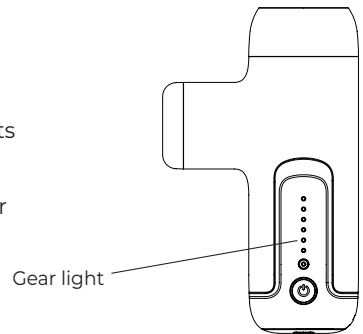- When charging, the red breathing light turns red,
- Do not use and switch the machine while charging!
- Fully charged, the red breathing light turns green,
- Can be unplugged to start using.
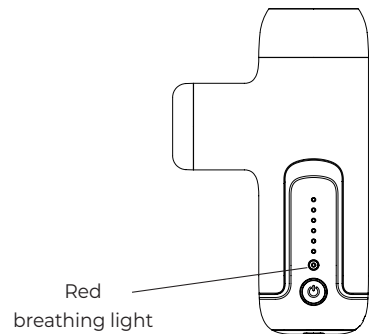
## Instructions for charging

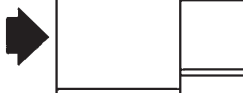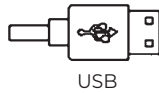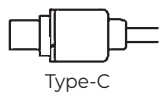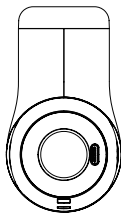USB (bring your own)  
Note: The output is 5V-2A or 5V-1A  
mobile phone charging head

1. It cannot be turned on during charging, and the USB charging cable must be unplugged after charging.
2. It takes about 3 hours to fully charge. After full charge, the first gear test can be used for about 4 hours (laboratory data, there is a certain error in actual use). The specific use time will decrease depending on the gear used and the pressing force.
3. Plug one end of the USB cable into the bottom of the host, and the other end can be connected to the USB port of a charger, power bank or desktop computer for charging. Please note that if you use a USB power supply for charging, please do not use a power adapter with an output specification exceeding DC 5V. Otherwise, the host will be damaged.

**\*In order to ensure the safety of use, the massage gun cannot be used when charging.**



## ♥ Indications

- Pain and cramps from muscle contusions, sprains, strains
- Loose and thickened of connective tissue and fascia
- Reduces the accumulation of lactic acid in the muscles
- Increase joint mobility
- - Remove muscle fatigue

## ! Taboo

- Aneurysm, bleeding, use of blood thinners.
- Heart disease, have a pacemaker or defibrillator.
- Pregnancy, cancer.
- Within 90 days of internal fixation, the device cannot be used within 3 inches of internal fixation.
- Sensitive areas head, face, cervical spine, spine vertebrae.
- Positions with less muscle coverage, such as shinbone, dorsum of the foot, dorsum of the hand, etc., Where muscle coverage is less close to bones.

## List
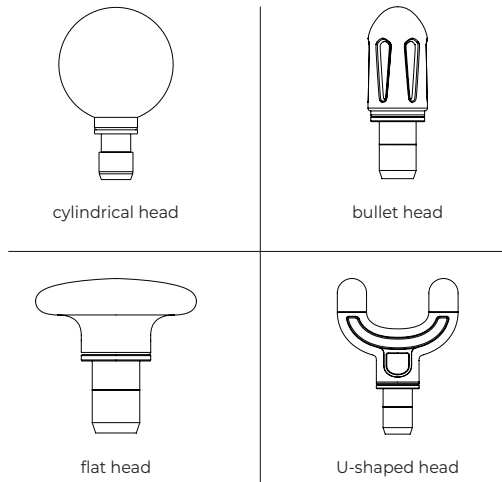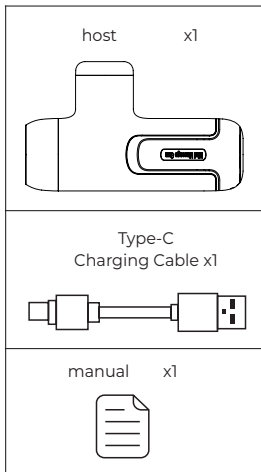
## Product Parameters

Model: BLD-307	Input parameters: 5V—1A
Rated voltage: 7.4V	maximum speed: 3300r/min
Rated current: 1A	Battery capacity: 1800mAh
	Executive standard: GB4706.1-2005 GB4706.10-2008

# Thank You

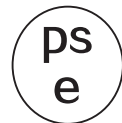
Thank you for your purchase and your trust in us. We select our products with the greatest care. If you have any questions about our product, please send us an email to [bol@xclusive-lifestyle.nl](mailto:bol@xclusive-lifestyle.nl)

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