


The instructions








Product structure

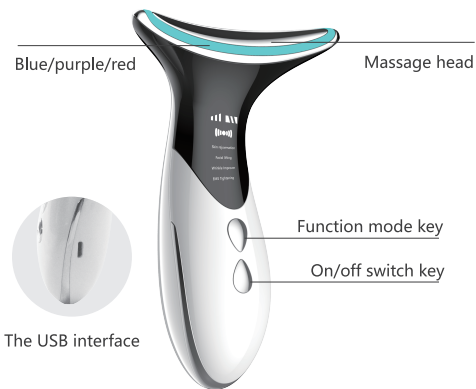

Name: Beauty neck instrument	Rated voltage :5V
Product model :2032	Material :ABS+PC material
Product size :92* 126* 37.4mm	Rated input power :2.5W
Charging time :3.5 hours	Battery capacity :750mAh
Main functions: lifting, tightening, relaxing and activating muscle	

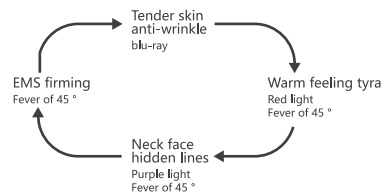Model explanation

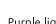

- 
Blue light Blue light can promote the synthesis of protein and collagen, activate skin, tighten loose skin, shrink pores, and greatly improve oily skin and severely sensitive skin.
- 
Red light Red light can enhance cell vitality, accelerate blood circulation, promote the growth of fiber cells and collagen, smooth fine lines, wrinkles, tighten the skin, lighten the pigment, dilute scar and other effects.
- 
Purple light Purple light can improve the alternate function of cell oxygen utilization, promote the microcirculation on the surface of the body, dredge lymph and drainage edema, at the same time have the effect of sedation skin, balance and stability.
- 
EMS EMS microcurrent releases intermittent pulse current, stimulates ATR production, improves facial muscle fiber elasticity, activates collagen protein activity, and tightens skin.

Product characteristics

- The streamline massage head is specially designed for the neck. The streamline design conforms to the ergonomic theory and transmits energy to every inch of the neck skin.
- Massage head acoustic vibration with four modes, there are three color light blue, red, purple, EMS massage, promote blood circulation, stimulate the points on the face and neck, dredge meridian, improve circulation to promote metabolism, promote waste discharge toxins accumulated for a long time, slow down the aging epidermis cells, to downplay the neck lines, minus the dropsy, relax the neck, acupuncture therapy, quickly eliminate fatigue, restore physical function.
- The four modes are: anti-wrinkle skin tender (blue light), temperature sensing pull (red light), neck and face cryptogram (purple light), EMS firming, long press the on key to start the machine, each press the mode conversion, one cycle, long press the on key to shut down.
- Each mode has three gears, and each mode has vibration
- Skin anti-wrinkle pattern dot long bright, small brackets from the inside out
- The heating mode has a heating temperature of 45 degrees, the symbol dot flashes and stops, and the bracket light does not light up
- The neritic pattern has a fever of 45 degrees, the symbol dots do not shine, and the brackets flash and stop
- EMS fever 45 degrees, long bright symbol dot, bracket from the outside to the inside horse
- USB cable charging, safe and environmentally friendly.
- Shape mini portable, anytime and anywhere use
- The new introduction of bionics design concept, the use of dolphin tail modeling, exquisite novel.

- Warm tips:
1. Please charge to saturation before using the product.
 2. If the product is not used for a long time (about 3 months), please take out the product and store it after full charge.

Procedures for neck use

1. Evenly apply skin care products to the neck.
2. Choose the mode according to personal needs. Massage head in the neck from the bottom to the top of the support pull gauge.
3. after the ministry also need to massage, the use of points, bridge, massage, push and other ways to massage.

Steps to use face

1. Apply skin care products evenly on the face.
2. Choose the mode according to personal needs.
3. Point pressure to relieve muscle tension.
4. Fit facial skin, pull lift from bottom to top, three-dimensional shape.

Cleaning steps

1. Long press the on key for 1 second to close the product.
2. Wipe the product with a clean flannelette without using detergent.
3. Place in a dry and ventilated place.

Matters needing attention

1. Before using the product, please read the instruction carefully and keep it for future reference.
2. The product is charged by USB cable. When the battery is low, it takes about 3.5 hours to reach saturation.
3. Flash red light when the battery is low, red light and green light when charging.
4. The product has the function of memory mode. When restarting, it will automatically enter the mode in which it was shut down last time.
5. Do not use this product when charging.
6. The product has the function of heating, about 45°C. There is a little difference in the temperature perception of each skin, which may cause it to feel hot.
7. If you feel uncomfortable or have abnormal skin during use, please stop using it immediately and consult your doctor.
8. The product has no waterproof function, so please do not place it in water.
9. Please place the product in a dry and ventilated place, avoid direct sunlight, and there are no corrosive chemicals around.

10. When nursing, please take an appropriate amount of skin care products. If excessive use of skin care products, the skin can not be completely absorbed resulting in waste, at the same time there is the possibility of residual liquid into the product internal, damage the product service life. Please wipe it with a clean cloth in time after use.
11. Do not let children touch this product to avoid danger.

Common Use Q and A

The problem	The solution
The product suddenly didn't vibrate. What's the reason for not working?	Increase the amount of maintenance products used. Please check whether the power of the product is insufficient. Is it turned on correctly? Restart the product after shutdown. Please contact the brand dealer consultation, or dial After-sales service hotline consultation.
I can use Ben every day Products?	Yes, the vibration intensity and temperature of this product are adjusted All are designed by professional engineering teams, but one Use it many times a day.
I am allergic skin, Can I use it?	Special selection of environmental protection, in line with international standards of materials, Safe, healthy, non-irritating to skin Please feel free to use it daily.
What time of day is it in use?	It can be used at any time. It is recommended to use 5- each time Ten minutes is better.

CERTIFICATE

DATE OF PRODUCTION: _____

EXAMINATION CLERK: _____

IMPLEMENTATION CRITERIA: GB4706.1-2005
GB4706.1-2005

The neck instrument

Please read the instruction carefully before use