ELECTRIC **BETTER** **TOOTHBRUSH**

with 5 BRUSHING MODES

IMPORTANT SAFETY INSTRUCTIONS

| | |
|----------------------|-------|
| English (EN)..... | P.1 |
| 中文繁體 (CN)..... | P.11 |
| 中文簡體 (SCN)..... | P.16 |
| čeština (CZ)..... | P.21 |
| dansk (DK)..... | P.32 |
| Deutsche (DU)..... | P.43 |
| Español (SP)..... | P.55 |
| Français (FR)..... | P.66 |
| Italien (IT)..... | P.74 |
| Nederlands (NI)..... | P.85 |
| Polskie (PI)..... | P.97 |
| svenska (ZW)..... | P.108 |

Read all instruction before using this product.

IMPORTANT SAFEGUARDS

This Electric Toothbrush has been carefully designed to offer you and your family a new and unique brushing experience, that is both safe and effective. However, when using this electrical product, especially when children are present, there are some basic safety precautions that should be followed.

WARNING

- Use this appliance only for its intended purpose described in this manual. Use only manufacturer recommended attachments.
- Not following any of the below mentioned instructions may cause malfunction, fire, burns, explosions, short circuiting, electric shock or any other injury.
- Keep this appliance out of the reach of children or infants as this is not a toy.
- This appliance is not intended for use by persons (including children and infants) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given proper supervision and instruction concerning the use of this appliance by a person responsible for their safety.
- Unplug the appliance immediately if you have any problem while using it.
- If the plug does not fit fully into the electric socket, reverse the plug. If it still does not fit, contact a qualified electrician.
- Always inspect the appliance before use, pay particular attention to the condition of the cord and plug.
- Never use while sleeping or drowsy.
- Connect this product to a properly grounded outlet only.
- Do not use the appliance if it has a damaged cord, plug or connector.

- Do not pull or carry the product charging unit by the cord and keep power cords away from heated surfaces.
- Do not unplug by pulling on the cord, instead unplug by pulling the head of the plug.
- Do not use the appliance if it is not working properly. Stop using the appliance immediately.
- Do not wind/wrap the main cord around the appliance.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a bathtub or sink.
- Do not place in or drop the charging unit into water or other liquid.
- Do not reach for the charging unit if it has fallen into water. Unplug immediately.
- We do not accept any responsibility for damage to the appliance, person(s), or other items caused by misuse, abuse or non-compliance with the instruction mentioned in this User Manual.

PLEASE READ CAREFULLY BEFORE USE

- This product has no user serviceable parts.
- Store the charging unit in a dry and safe place when it is not in use.
- Untwist the power cord before use.
- In the unlikely case this product does not function properly, you should not dismantle, modify or manipulate any part of it otherwise there is a risk of electrical shock. You should contract the manufacturer.
- If the power supply cord is damaged, it must be replaced by the manufacturer in order to avoid any risk of electrical shock.

- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.
- This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, turn the plug head. If it still does not fit, contact a qualified electrician. Do not attempt to force the electrical plug into the outlet.

CARE AND MAINTENANCE

- Never use any abrasive materials that may scratch or damage the appliance. Regular cleaning will maintain your electric toothbrush's optimum performance level. Proper use and storage of the power supply cord is essential to the safety and function of the appliance. The cord should never be twisted, pulled or bent. Do not wrap the cord around the appliance. Routinely examine cord for damage including tears in the cord or exposed wires.
- After brushing, rinse you brush head thoroughly under running water.
- Remove the brush head from the handle and clean both parts separately under running water. Wipe them dry before you put the handle onto the charging unit.
- Before cleaning your charging unit, make sure you unplug the adapter.
- Only use a dry or damp cloth to wipe the surface of the charging unit.

IMPORTANT

- Do not clean the charging unit by immersing it under running or still water.
- Place the charging unit onto a dry and safe place.
- Do not plug in the charging unit until it is completely dry.

- This product contains rechargeable batteries in the interest of protecting the environment.

DECLARATION

This toothbrush conforms to the applicable provisions of US Federal requirements, including 21 CFR, Chapter 1, Subchapter J, Radiological Health. This product complies with applicable EU-regulations for household products and standards including electromagnetic fields (EMF). If handled properly and according to the instructions in this user manual, the product is safe to use based on scientific evidence available today.

LIMITED WARRANTY

USA

World Wide Daily will warranty this electric toothbrush for one (1) year on all manufacturer defects in material and workmanship. This warranty is only valid with a receipt from authorized reseller and registration of the appliance at www.abetterdental.com. World Wide Daily will not warranty any electric toothbrush that is damaged or abused due to user error (dropping, wrong voltage supply, toothpaste build up, etc)

Asia

Hong Kong:

World Wide Daily will warrant this electric toothbrush for one (1) year on all manufacturer defects in material and workmanship. This limited warranty does not cover any normal wear and tear parts including brush head and degradation of battery. It does not warrant any electric toothbrush that is damaged or abused due to user error (dropping, wrong voltage supply, toothpaste build up, etc). This warranty is void immediately also, however, if the electric toothbrush is performed any of modifying or repairing jobs by any unauthorized person or third party.

China:

深圳市美活貿易有限公司 will warrant this electric toothbrush for two (2) years on all manufacturer defects in material and workmanship. This limited warranty does not cover any normal wear and tear parts including brush head and degradation of battery. It does not warrant any electric toothbrush that is damaged or abused due to user error (dropping, wrong voltage supply, toothpaste build up, etc). This warranty is void immediately also, however, if the electric toothbrush is performed any of modifying or repairing jobs by any unauthorized person or third party.

Other Asian countries:


For the terms and conditions of warranty service in all the other Asian countries / regions other than China and Hong Kong, please contact / consult your local authorized resellers or distributors, visit web page www.worldwide-daily.com or email to sales@worldwide-daily.com to get details and arrangement.

EU & Africa

Perfect Care BV will warrant this electric toothbrush for the legal period valid in the European Union on all manufacturer defects in material and workmanship. This warranty is only valid with a receipt from authorized reseller. Perfect Care BV will not warrant any electric toothbrush that is damaged or abused due to user error (dropping, wrong voltage supply, toothpaste build up, etc). Please check the European Union website for all obligations and rules for the legal warranty period.

LIST OF COMPONENTS

- a** Brush head
- b** Toothbrush Handle Grips
- c** On/Off & Brushing Mode Button
- d** Brush Mode Indicators
- e** Power & Charge Indicator
- f** Charger Unit
- g** Handle Stander
- h** Brush head Stander (x2)

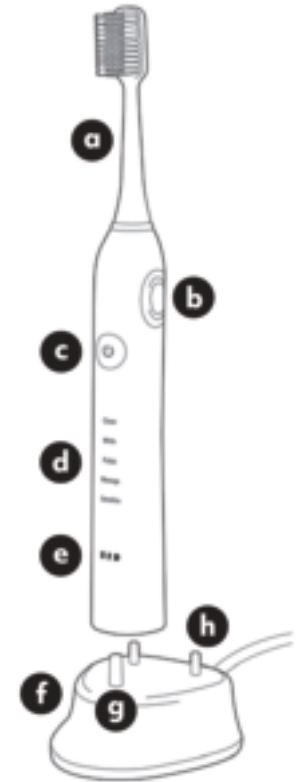CHARGING

This Electric Toothbrush has a waterproof handle. It is electrically safe and is designed for use in the bathroom. But it is not meant to be used while fully-submerged in water.

When charging the appliance, the green charging indicator will be flashing, once the appliance is fully-charged the green indicator light will become continuous.

When your electric toothbrush is not in use, you can leave it on the charging unit as the battery will not overcharge.

It is always recommended that if the charging base is not in use, to unplug the unit from the electrical outlet.

HOW TO USE THE ELECTRIC BETTER TOOTHBRUSH

Attach a brush head to the handle unit. The attached brush head should snap in and fit snug to the handle.

To avoid splashing of toothpaste or water, guide the brush-head to your teeth before switching the appliance on.

BRUSHING TECHNIQUE

1. Guide the brush head slowly from tooth to tooth.
2. Hold the brush head in place for a few seconds before moving on to the next tooth.
3. Brush the gums as well as the teeth. Start with the outside of the teeth, then the insides ending with the chewing surfaces.
4. Do not scrub or press hard, simply let the brush do all the work.

BRUSHING MODES

This Electric Toothbrush offers five effective brushing modes for varying oral care needs. It features an auto shut-off after two minutes of brushing and a dentist-recommended auto-alert when brushing the same area for more than 30 seconds.

CLEAN: With up to 31,000 strokes per minute, exceptional mouth cleaning—great for everyday use.

WHITE: With up to 41,000 strokes per minute, removes microscopic plaque and bacteria left behind by other toothbrushes.

POLISH: With up to 48,000 strokes per minute, will effectively scrubbed away up to 2x more plaque than other toothbrushes.

MASSAGE: With up to 31,000 strokes per minute, pulsates between strong to soft and then from soft to strong in a circular cycle which is great for gum stimulation.

SENSITIVE: With low pulse 31,000 strokes per minute, great for gentle cleaning for those with sensitive teeth and gums.

TWO MINUTE AUTO-TIMER

This Electric Toothbrush is designed to reach the proven and best brushing timing of 2 minutes. After 2 minutes of continuous brushing, the toothbrush will automatically stop.

HOW TO SWITCH BETWEEN BRUSHING MODES

1. Power On/Off and brushing mode button (See list of Components).
2. Keep pressing the mode button to cycle through the operating sequence from power on, to Clean mode/White mode/Polish mode/Massage Mode and Sensitive Mode.
3. The device will turn off if the mode button is depressed after Sensitive mode or by pressing and holding the mode button for approximately 5 seconds.
4. The Electric Toothbrush automatically starts in the “Clean” mode when pressing the brushing mode button once.
5. By hitting the mode button the appliance will start in Clean mode, then by pressing the mode button the modes will switch in sequence to White, Polish, Massage and lastly Sensitive.
6. The device will turn-off after Sensitive mode.
7. Attach a brush head to the handle, dispense your favorite toothpaste onto the brush head, to avoid splashing guide the brush head to your teeth before switching the appliance on.
8. Place the bristles against the teeth at a slight angle towards the gum line.
9. Press the mode button to turn on the toothbrush.

To replace a brush head, hold the neck of the brush head (“a”) and gently pull up to remove. To re-attach a brush head, align the brush head replacement to the center of the toothbrush, gently press down the brush head to the handle until it locks into place. There will be a small (0.5mm) gap between the brush head and handle.

After brushing, rinse both the brush head and handle under running water, then wipe dry and dock the toothbrush on the charging unit.

WHAT TO DO IF YOU HAVE A CLAIM

USA

If a product fails to operate properly under normal use because of a covered defect and it is WITHIN the Warranty Period, then please:

1. Visit www.abetterdental.com to obtain a Return Authorization Number and return instructions or email us at support@abetterdental.com.
2. To return the defective product with the Return Authorization Number, shipping, handling and replacement fee of \$35.00 ("SHR Fee"). The SHR Fee must be paid in the form of money order or check drawn on a US bank and made payable to "JP Products International LLC" (PLEASE DO NOT SEND CASH). We WILL NOT process any returns without a Return Authorization Number or the SHR Fee. We are not responsible for products lost or damaged in transit.

Asia

Hong Kong

If a product fails to operate properly under normal use because of a covered defect and it is WITHIN the Warranty Period, then please:

1. Please call +852 2426 4921 or email to sales@worldwide-daily.com to describe your case. Then, if your defective unit is confirmed to be qualified for this limited warranty service, please send the defective unit to our service center together with the original purchase invoice for the warranty validation. Please be aware that freight costs that occur are for the consumers account. We WILL NOT provide any warranty service without the present of the original purchase invoice.

China

If a product fails to operate properly under normal use because of a covered defect and it is WITHIN the Warranty Period, then please:

1. Please call +86 755 8325 6458 or email to sales@worldwide-daily.com to describe your case. Then, if your defective unit is confirmed to be qualified for this limited warranty service, please send the defective unit to our service center together with the original purchase invoice for the warranty validation. Please be aware that freight costs that occur are for

the consumers account. We WILL NOT provide any warranty service without the present of the original purchase invoice.

Other Asian countries:

For the warranty validation of terms and conditions in all the other Asian countries / regions other than China and Hong Kong, please contact / consult your local authorized resellers or distributors, visit web page www.worldwide-daily.com or email to sales@worldwide-daily.com to get details and arrangement.

EU & Africa

If a product fails to operate properly under normal use because of a covered defect and it is WITHIN the Warranty Period, then please:

Download a Return Form at <https://shop.bettertoothbrush.com/wp-content/uploads/Retourformulier-EN.pdf> and see further instructions. You may also call at +31 85 0643011 for further help. To return the defective product with Return Form, please be aware that freight costs that occur are for the consumers account. We WILL NOT process any returns without a valid invoice and a duly filled out Return Form. We are not responsible for products lost or damaged in transit.

ENVIRONMENT

Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. Dispose of the rechargeable battery at an official collection point for batteries. By doing this, you help to preserve the environment. Only remove the rechargeable battery at the end of the product's life. This process is not reversible! Make sure the rechargeable battery is empty. Observe basic safety precautions.

To remove the rechargeable battery from the handle:

Place the handle on a solid surface and lightly tap the housing about 2.5cm from the bottom edge with a hammer. Turn the handle 90 degrees and repeat until the bottom cap separates from the handle.

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