# BEAUTICAL SONIC EYE CARE BOOSTER

### EN: Instructions for use

Thank you for purchasing this product. Before using this device, please read this user manual carefully and keep it for future reference.

### Precautions

- This device must not be used by children, pregnant women, patients suffering from epilepsy, cancer or other serious conditions, patients with a pacemaker or internal defibrillator, patients undergoing radiotherapy or chemotherapy.
- Please follow the instructions. Do not use this device if you suffer from eye inflammation such as conjunctivitis or any active condition in the treatment area such as sores, eczema or rash.
- Remove all metals such as earrings, necklaces, bracelets etc. before use.
- Do not use the unit while charging.
- Do not repair, disassemble or modify this device without authorization.

### Product Description


This highly efficient eye massager features 3 different functions:

- high-frequency vibrations – to refresh and revitalize the eye contour and to reduce the appearance of wrinkles, dark circles and puffiness
- thermotherapy – to open pores and stimulate blood circulation
- red-light treatment – to smooth and firm the delicate skin around the eyes

When used with your favourite eye cream or eye serum, the device increases the absorption of active ingredients and therefore maximizes their efficacy for rejuvenated, firmer and more radiant skin.

### 3 working modes:

- Heat + vibration + light
- Heat + light
- Heat + vibration

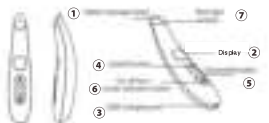- 1- Metal massage head
- 2- Display
- 3- USB charging port
- 4- Upshift button
- 5- Downshift button
- 6- On/off key / mode selection button
- 7- Red light window

### Instructions for use

- Take an appropriate amount of eye cream/serum and spread it evenly on the eye area.
- Long press **🔴** to turn on the device, then short press **⏮** to select the mode and press "+/-" to adjust the temperature as needed.
- Holding the device by the upper part, use the massage head to massage **each treatment zone for 3 minutes.**
- After using, long press **🔴** to turn off the device.

### Charging Method

- When the device's power is too low, turn off the device and plug in the charger.
- The battery symbol on the display flashes during the charging process and stops flashing when charging is complete.
- When charging is completed, please unplug the charger.
- Please charge occasionally to avoid ageing of the device after a long period of non-use.

### Maintenance

- Please distribute the skin care products evenly before each use in order to avoid the accumulation of an excessive amount of product on the massage head.
- Please wipe the product with a dry towel after each treatment.
- Do not pour water on the device or wash the device with water.
- Do not use corrosive products to clean the unit.
- Store it in a cool, dry place.
- Do not use the device if it is damaged.
- Avoid collisions or violent shocks.

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