
Handleiding - Infrarood rugband
Manual - Infrared back wrap

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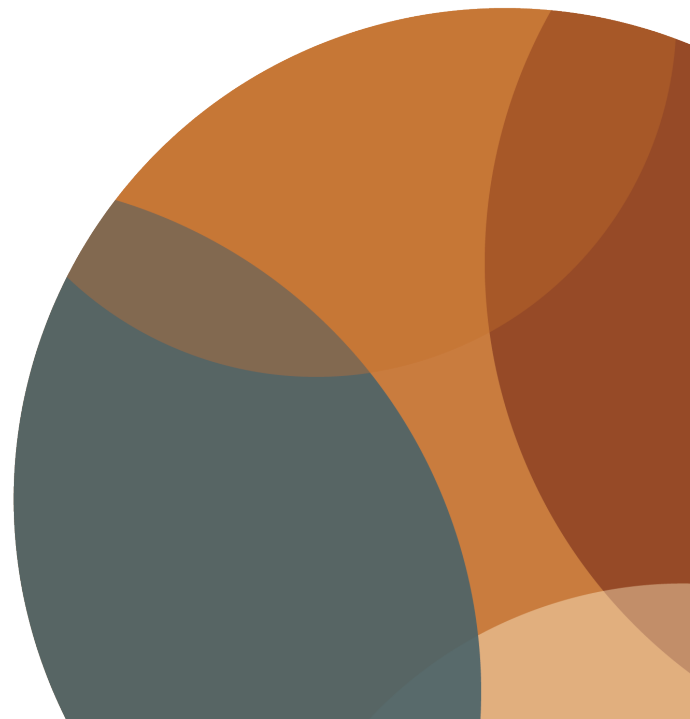Thank you for choosing one of Liroma's products!

The Liroma infrared belt offers an extensive list of possible benefits, including detoxifying the body, burning calories, skin rejuvenation, strengthening the immune system and relieving pain.

The infrared belt has many benefits

- Rejuvenates the skin
- Improves sleep quality
- Faster recovery
- Burns calories
- Detoxifies the body
- Stimulates the happiness hormone
- Reduces stress
- Increases collagen production
- Reduces swelling
- Promotes blood circulation
- Strengthens immunity

READ BEFORE USE

IMPORTANT:

1. Before use, make sure your blood pressure is normal and remove metal objects from your body.
2. Before, during and after use, it is recommended to drink 200-300 ml of water each time (with a little salt if necessary), to provide fluids to the body and flush toxins from the body.
3. For the first time, use the low or medium temperature setting. The recommended user time is 30-45 minutes.
4. If possible, use the device after taking a bath.
5. Do not use the device one hour before and after meals. Do not use it when extremely hungry or full. Do not use it if you are extremely debilitated.
6. Do not use air conditioning after use or take a bath within 1 hour.
7. It is prohibited for use by humans:
 - With high blood pressure
 - By people who are sensitive to heat
 - Who have a pacemaker
 - Who cannot take care of themselves.
 - Pregnant women, babies and children should not use the product.
8. Start using the belt at around 60°C. If you want a sauna-like experience, you can gradually increase the temperature with caution. Always place your belt on a heat-resistant surface.
9. Always disconnect the belt from the power outlet when not in use. The belt is designed for use with regular power. Do not use it with generators, power converters or inverters.
10. Avoid exposing the controller or connector to water, as they are not waterproof.
11. If you have medical conditions, we strongly advise against using the belt without first consulting your doctor.
12. When using the belt. Please do not put the remote control on the belt.
13. Do not fold the belt while it is hot and cooling down. This may cause damage to the heating wires and internal material of the blanket.
14. Let it cool down completely before folding it.

Operating instructions

1. Wear the belt around your waist and adjust it to a comfortable position (**Note:** you can also use the belt over the shoulder or around your leg)
2. Plug in the power cord.
3. Press the power button.
4. Press the temperature + and - buttons to adjust the temperature.
5. Press the time + and - buttons to adjust the timer.
6. Turn off the power after use.

Care and storage

1. Do not use sharp objects with the product.
2. After use, clean the product and store it out of direct sunlight.
3. Keep the product clean by using a soft cloth.
4. Do not use water, flammable or chemical solutions to clean the product. Also, do not store it in places with high humidity.
5. Do not place heavy objects on this product.
6. Do not immerse this product in water.

If you have any questions about the product please contact the seller or at info@liroma.nl