BABY SOOTHER THERMOMETER 110

**USER GUIDE / HANDLEIDING / MANUEL D'UTILISATEUR /
BEDIENUNGSANLEITUNG / MANUAL DE USUARIO /
BRUKSANVISNING / BRUGERVEJLEDNING / ANVÄNDARHANDBOK /
KÄYTTÖOHJE / MANUALE D'USO / MANUAL DO UTILIZADOR /
РЪКОВОДСТВО НА ПОТРЕБИТЕЛЯ / UŽVATELSKÁ PŘÍRUČKA /
ΟΔΗΓΙΕΣ ΧΡΗΣΗΣ / HASZNÁLATI ÚTMUTATÓ / INSTRUKCJA
OBSŁUGI / MANUAL DE INSTALARE SI UTILIZARE / РУКОВОДСТВО
ПОЛЬЗОВАТЕЛЯ / UŽÍVATEĽSKÝ MANUÁL / KULLANICI KILAVUZU**

V1.0 - 04/09

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
1 SAFETY ADVICE

- Read this user's manual carefully before use.
- Keep this user's Manual for future reference.
- This product is for household use only. It is not intended as a substitute for consultation with your doctor.
- Examine the pacifier bulb before each use for cracks, tears, tackiness or other signs of wear. Discard the worn unit immediately.
- Do not expose to direct sunlight, high temperature, moisture or dirt.
- Do not attempt to pull or twist the nipple apart from the body as this will damage the thermometer's sensitivity probe and electronics.
- Never allow a child to play with this product. It is not a toy. It is a sensitive electronic device and may be damaged if not used properly under parental supervision.
- Use this thermometer for oral temperature measurements only.
- The battery is not replaceable. Do not attempt to remove the battery.

2 DESCRIPTION

1. Temperature sensor
2. On/Off button
3. Display

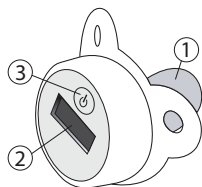3 HOW TO USE

- Remove the protection cover carefully.
- Press the On/Off button briefly. The display indicates first 36.5°C and afterwards 'L°C'. When the '°C' symbol starts flashing, the thermometer is ready to use.
- If the ambient temperature is higher than 44°C, the display indicates 'H°C'.
- Position the pacifier thermometer bulb in your child's mouth. Allow your child to suck at least 3 minutes (minimum). To have a more accurate temperature, keep the thermometer in the child's mouth for 6 to 8 minutes.
- When the peak temperature is stable for 32 seconds, the alarm tone will be heard.
- To turn the unit off, press the On/Off button briefly.
- Re-install the protection cover immediately after use.

4 MEMORY

The last 5 measurement results are stored in the memory.

To recall the memory:

- Press and hold the On/Off button. The display will show M1, followed by the last measurement result in °C.
- When the result is shown, press the On/Off button to recall the previous result.

Note: If you wait too long, 36,5°C followed by 'L°C' will be displayed and the thermometer will initialize for starting a new measurement.

- Press the On/Off button repeatedly to scroll through the 5 memories.

This product is in compliance with the essential requirements and other relevant provisions of the directive 93/42/EEC. This is confirmed by the CE-mark.

The Declaration of conformity can be found on:

<http://www.topcom.net/cedeclarations.php>

5 CLEANING AND MAINTENANCE

The thermometer may be sterilized immersing the silicon rubber nipple ONLY in a sterilizing solution such as alcohol. Immerse the nipple approximately 1 second in the solution. The protective clip-on plastic supplied may be used as container for this purpose.

The Pacifier Thermometer should never be completely immersed in sterilizing solution or sterilized in boiling water as damage will occur to the electronic components.

6 DISPOSAL OF THE DEVICE (ENVIRONMENT)

At the end of the product lifecycle, you should not throw this product into the normal household garbage but bring the product to a collection point for the recycling of electrical and electronic equipments. The symbol on the product, user guide and/or box indicate this.

Some of the product materials can be re-used if you bring them to a recycling point.

By re-using some parts or raw materials from used products you make an important contribution to the protection of the environment. Please contact your local authorities in case you need more information on the collection points in your area.

7 TECHNICAL SPECIFICATIONS

Model number	VT-901 W5
Temperature range	+32,0°C to +44°C
Accuracy	+/- 0,1°C (+35°C to +41°C) +/- 0,2°C (rest)
Resolution	0,1°C
Working conditions	Temp: 0 to +50°C - Hum: 10-83% Temp: -20 to +60°C - Hum: 10-83%
Battery	1,55V DC (Alkaline) - SR41/LR41/UCC392/ UCC192 (non replaceable)
Battery lifetime	Approx. 100 hours continuous use Approx. 18 months while used 10min/day
Power consumption	0,15mW
Sensor type	Thermistor
Display	Liquid Crystal Display 3-1/2 digits
Memories	5 last results
Auto-shut-off	10 minutes after last temperature increase
Dimensions	55x50x50mm
Weight	+/- 12gr (Including batteries)
Classification 	<ul style="list-style-type: none"> • Internally powered equipment • Classification of type BF • IP44 • Not suitable for use in presence of flammable anaesthetic mixture with air or with oxygen or nitrous oxide • Continuous operation with short-time loading

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