Attentie

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
Waterproof: IPX

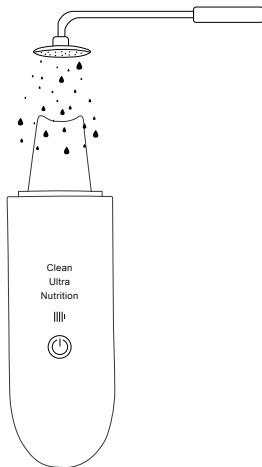Attention

1. Do not exceed the scope of use specified in this user manual.
2. Pregnant women should follow a doctor's advice.
3. Do not use this near an area with metal prosthesis, such as eye balls, inner ear, bone areas and joints.
4. Prohibit use in vein or arterial pathology areas or patients with thrombophlebitis or blood coagulation precedents.
5. Avoid physical impact.
6. Only use according to the user manual. Prohibit disassemble, transforming, maintenance or retrofit the product.
7. Please store the device in a clean and stable place. Do not store it in a place that is overheated or over cold. Keep it away from fire, such as a stove.
8. Charging base is not waterproof. Keep the charging base dry and clean.
9. Unplug the device after it finishes charging.
10. Keep the device away from children.

Introduction

Thank you very much for purchasing our product and your confidence in us. We hope our product can bring you a perfect experience. Please read the user manual and all related information carefully and keep it for future reference.

Product specification

NW: 205G

Measurement: 21,5 * 9,8 * 8 cm

Working environment: 0-35 ° C

Input: AC 100-240V, 50 / 60Hz

Output: DC-5V, 1A (rechargeable)

Power: 2. 2W

Battery volume: 800 mAh

Ultrasonic frequency: 24 kHz

Blue light wavelength: 465 ± 3 nm

Waterproof level: IPX5 (device only, charging base is not waterproof)

Accessories

1. Main device
2. 5-pin USB charging cable
3. Charging base
4. User Manual
5. Gift box
6. Water transfer bottle

Features

1. Standby mode: Press the button one time to turn on the device and to enter standby mode. In this mode you can see the lights for 'clean', 'ultra' and 'nutrition' flashing.

2. Clean mode: Press the button three times to enter this mode. You can hear the sonic vibrations and you see the mist spray. The 'clean' indicator light is always on.

3. High level deep cleaning mode: Press the button four times to enter this mode. You can hear the sonic vibrations The 'ultra' indicator light is always on.

4. Sonophoresis mode (ultrasonic-treatment): Press the button five times to enter this mode. You can hear the sonic vibrations and can feel ions. The 'nutrition' indicator light is always on.

5. Standby mode again: Press the button six times to go back to standby mode. No functions are active now. Press the functions circle to activate.

6. Turn off the device: Long press the button for one second to turn off the device.

7. The device is set to run for ten minutes. After ten minutes the device will turn off automatically.

Charging

1. When the battery is low, the red indicator light lights up and the device turns off automatically.

2. Put the device on the charging station and make sure it is well connected. Then charge it with a USB cable. While charging, the red light flashes. When it's completely charged, the light turns green and stays on all the time.

*While charging,
the red light flashes.*

*Done charging,
green light stays on.*


Functions

1. **Ultrasonic deep cleaning technology:** The ultrasonic skin scrubber produces 24,000 Hz per second of ultrasonic waves. These waves effectively releases clogged pores and exfoliates the skin to kill and remove dead skin cells painlessly.
2. **Firms and lifts skin-micro-current:** Pulse technology intermittently on ultrasound helps to stimulate cell collagen regeneration and activates cell vitality, keeping your skin hydrated and stays healthy.
3. **Advanced ion technology:** The positive ion helps to remove dirt and deeply cleanses the skin. It accelerates the speed of the skin cell metabolism, softens old keratin, restores the elasticity of the skin and helps the skin to tighten and be smooth. The leading function of the negative ion helps to break down and refine larger molecules like collagen. It helps to massage and support the skin. It also absorbs nutrients easily.
4. **Bio Blue Light:** Wavelength $465 \pm 3\text{nm}$ - Reduce inflammation, kills bacteria and propionibacterium acnes. It removes acne and scars.

Tips for effective cleaning

Deep cleaning

1. Apply skin cleaning serum or oil on your face and keep your skin moist at all time.
2. Switch on the device. Choose the first 'continuous wave deep cleaning mode' or the second 'high level deep cleaning mode'. Hold the appliance with the control panel away from the face. Clean the face according to the diagram below. The treatment time is usually three to five minutes each time, no more than five times a week. It is recommended using 'continuous wave deep cleaning mode' for first time users.
3. Wash your face with water after the treatment.

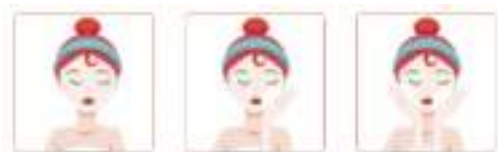Cleaning

It is recommended to use the device according to the diagram, because it supports the body's natural lymph drainage.

1. Sweep the machine upwards onto the face in horizontal strokes. Repeat on both sides of the face.
2. Work methodically over all areas of the face as shown in the diagram below.
3. Change the direction of the handset and glide in vertical strokes. Work from jaw to hairline.
4. Clean around the curve of the nostril using the curved inner edge of the blade.

Use tips

Use the scooped-out middle section of the blade to clean the chin area under the lip with vertical strokes. Select the first 'ultrasonic deep cleaning mode' to work around the eyebrows. Avoid the sensitive areas below the eyes when using the second 'high level deep cleaning mode'".

lift up & sonophoresis

1. Apply your favorite serum or moisturizer onto your face. Hold the appliance so that the blade is directly against the skin.
2. Turn on the device and turn on the fifth 'sonophoresis mode'.
3. Slide the blade over the surface of the skin with light pressure using backward movements. For best results, work the knife slowly at all times.
4. Three to five minutes is sufficient for each treatment. Use three to five times weekly.

Tips lift up & sonophoresis

Let the blade slide over the surface of the skin, with a backward movement. Work from the center of the face outwards.

Continue upwards from the jawline to the hairline. Sweep over the face until all the product has been absorbed.

Maintenance

1. After each use, clean the device with a wet cloth and than keep it in a dry and cool place.
2. Before cleaning the device, make sure that it is turned off and not connected to power.
3. Personal use only. Do not share the device with others.
4. Keep this user manual for future reference.

Failure

If the device does not work normally, please check the following:

1. Make sure the device is fully charged.
2. Be sure to press the 'power' button and that the indicator lights are on.
3. If there is no response after turning on the device, maybe the battery is low. Charge the device and try it again.
4. If no action is taken for thirty seconds, the device automatically shuts down to protect the circuit.


Adding water

Caution! Use the VRL Skin Scrubber with demineralized water for better misting results and to prevent limescale.

Waterproof

Waterproof: IPX

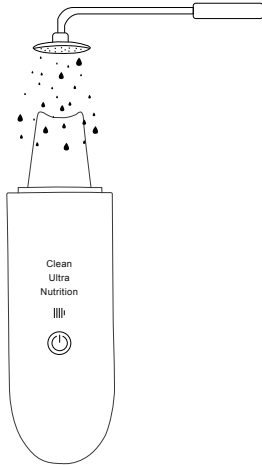VRL *Beauty*

Fabrikant
VRL Shop
Postbus 4189,3130 KD, Vlaardingen, Nederland
info@vrl-shop.nl | www.vrl-shop.nl