BLUE LION
LIFESTYLE

BLUE LION
LIFESTYLE

MARULA

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Congratulations!

You must be proud that you have invested in yourself and your physical health. When you are working hard every day, you also need to take good care of your body. In this user guide you can read how you can safely give your back and neck a nice and loose massage by your new personal masseur, so that you can get after it again!

About your new personal back masseur: Marula

Your new back and neck massager, called Marula, solves a problem that many people face these days. More and more people are experiencing problems with their backs and sitting a lot does not help. Until now. It is Marula's mission to help you improve your physical health. Put him in your favorite seat and experience a deep massage that loosen up your back muscles in three different zones:

1. Full back
2. Upper back
3. Lower back

At the same time, you can also enjoy a firm neck massage.

In addition to the kneading 'shiatsu' massage, Marula can also give you a firm tapping massage.

Do you have a specific spot where you want to be massaged? If you press the 'spot' button on the remote control, he will focus his massage on a specific place and you can manually guide him up or down with the arrows. You can determine the massage intensity with the remote control and by using your massage cushion with or without a protective cover. You can also put an extra towel over it or put on a thick sweater if you want to make it even less intense. The seat has a vibration option with two different intensities so you can enjoy those 'good vibes'.

Use your massage cushion on the couch, in the chair or in the car. Loosen your muscles before a workout, promote recovery after an injury, stimulate your circulation while working in the office, or enjoy a relaxing massage after a long day at work. However you choose to use it, we hope you're as excited about it as we are!

Before you start...

Read the safety instructions on the next page. We understand that you are excited and want to use your massage cushion right away, so we have kept it as short as possible to save your time and also save paper.

Warning

Careless use of the massage cushion can potentially cause personal injury and damage to the device or the environment. Read the safety instructions carefully to continue to enjoy your back massage cushion in the long term.

Safety instructions

1. This device is not suitable for animals or use while driving (only use in a stationary car)
2. Only use your your massager after consulting a doctor or specialist if:
 - you are pregnant;
 - you have a pacemaker;
 - you have artificial joints or implants;
 - you experience intense pain during the massage or have physical health conditions such as: heat sensitivity, blood circulation problems, varicose veins, bruises, muscle or tissue inflammation, injuries and open wounds.
3. Unplug the massager after every use.
4. Grab and carry your massager by itself: do not lift or pull on the cable.
5. Avoid your massager from getting wet and only handle it with dry hands.
6. Only use the original adapters and power cable supplied by the manufacturer.
7. Keep your massager and its power supplies away from heat sources, highly inflammable substances, vapour, liquids or sharp objects.
8. Before use, check if the cables are free from knots, twists or damage.
9. Provide adult supervision when a child, elderly or disabled person uses the massager.
10. Clean with a soft absorbent cloth without detergents (see p.8).
11. For your own safety, parts should only be replaced by professional / authorized service personnel.
12. Send the device to a professional service center for repair and maintenance.
13. Lay the cable carefully so that no one can trip over it.
14. Your massager is intended for indoor use.
15. Your massager is intended for private use.
16. Keep all parts and packaging materials out of the reach of children.

Important to know:

The warranty is void if these safety instructions are not observed or if the back massager is used for purposes other than those described above. The supplier cannot be held liable for damage in these cases.

What's in the box?

1x Blue Lion Marula-1 Back & Neck Massager

1x Remote control

1x 12 V Home adapter

1x 12 V Car adapter

1x User guide

Technical details

Brand	Blue Lion
Model	Marula-1
Voltage Supply	100-240V - 50/60Hz
Power Supply Unit	12V - 3,5A
Power Consumption	42W
Timer	5 - 10 - 15 minutes
Dimensions	32 x 6.9 x 17.3 inches
Weight	16 lbs
Power cable length	118 inch
Car adapter cable length	71 inch
Test reports and certifications	CE, RoHS, LVD, EMC, ERP

RoHS

Did you know...


The massage cushion is named after the Marula tree in Africa. Oil is extracted from the seeds and peels of the fruits of the Marula tree, which is used, among other things, for massage oil.


8 Steps to your perfect massage experience

- 1 Place your Blue Lion back and neck massager on your favorite chair or sofa.
- 2 Use the elastic straps on the back of your massager to attach it to the chair for extra support.
- 3 Connect the adapter cable to the cable attached to the cushion as shown in the picture below:

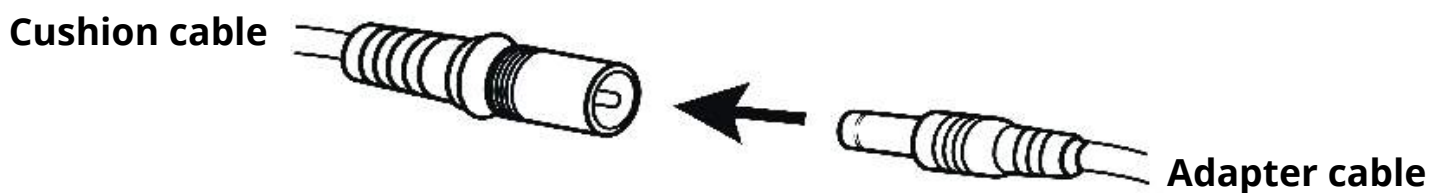- 4 Plug the home adapter into an electric outlet.
- 5 Turn on your back and neck massager with the remote control and choose your favorite zone and mode (see controller instructions on next page).
- 6 Choose your preferred massage intensity.
 - Too intense? Use the detachable protection cover for a gentler massage.
 - Not intense enough? Remove the detachable protection cover for a more intense massage.

Is the massage at the lowest intensity with the protection cover still too intense? You can use an extra towel or blanket or put on a nice thick hoodie to make it more soothing.

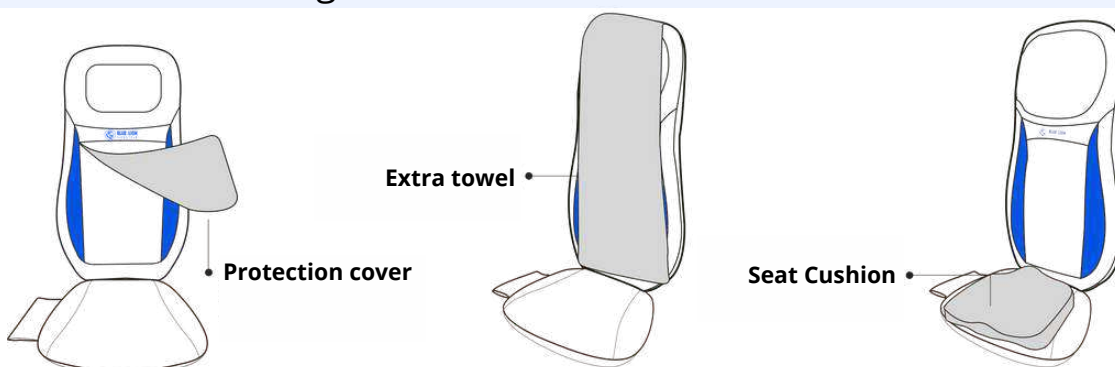- 7 If necessary, use a seat cushion to adjust the height of your back. Order your seat cushion at www.bluelionlifestyle.nl
- 8 Put the remote control in the pouch on the right side of the seat and enjoy your personalized massage!

Controller instructions

Neck

Activate a nice neck-massage. Press again to reverse the movement. Adjust the height with the arrows.

Tapping message

Activate a firm tapping massage for your back in two modes. Continuous and pulsating.

Massage zones

Where do you want to be massaged?

- Full back
- Upper back
- Lower back

Select your favorite massage zone.

Vibration

Activate vibration in the seat in two different intensities

Timer

Choose a massage time of 5, 10 or 15 minutes

Demo

Quick demonstration of all different functions. The flashing lights indicate the active functions. The back massager automatically turns off after completion.

Intensity tapping message

Change the intensity of the tapping massage.

Spot message


Press to keep Marula massaging in one place. Adjust the position up and down with the arrow buttons.

Heat

Activate the heating in the massage heads. (perfect for after a hard workout).

Power

Switch your back massager on or off.

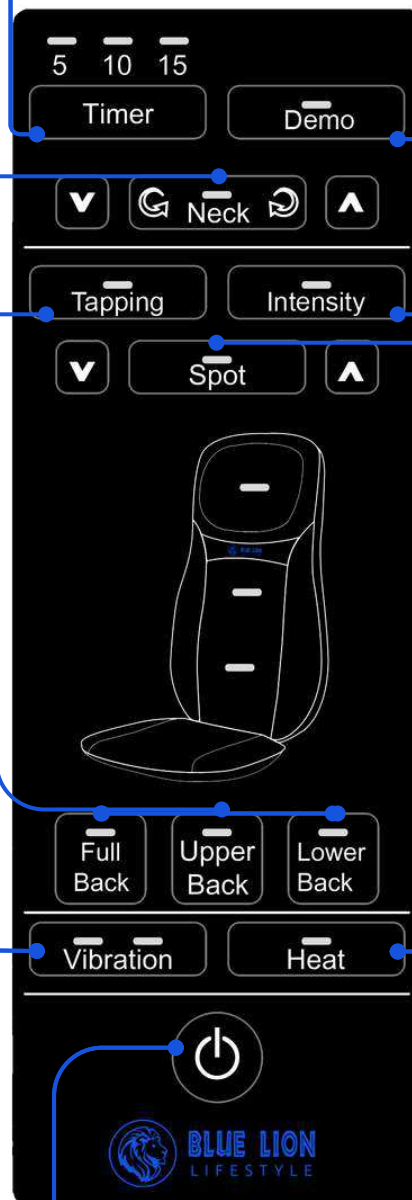Marula will automatically turn off after the time you have selected has passed. Do you want Marula to continue massaging? Press the power button again.

Cleaning

Take good care of Marula and he will take good care of you. Careless maintenance can shorten its life and void the warranty. This is how to keep Marula clean:

1. Unplug the massage cushion.
2. Use a soft absorbent (slightly damp) cloth.
3. Immediately dry your back-massager cushion afterwards.
4. Only use your massager again until he is completely dry.

Service & Warranty

We believe in the quality of our Marula-1 back-massager with great confidence, so we offer you a full year warranty against defects due to material or manufacturing defects. Please keep your proof of purchase and this manual safe in case you need it.

Note: The warranty only applies if the user has carefully followed the instructions in this manual. The warranty does not apply to damage caused by improper/careless use, or wear and tear. Read the extended warranty conditions on our website for more information.

Contact

Do you have any questions, comments or just want to share something with us? You can chat or mail us and send us a WhatsApp message too. Whatever it is, we are here for you.

help@bluelionlifestyle.nl

085-0604658

www.bluelionlifestyle.nl

Help nature grow

We hope you enjoy Marula's massages for many years to come. Has your massager reached the end of its life after a long enjoyable time?

Don't just take it to the dump, but hand it in at a recycle collection point and be a hero! By recycling the cardboard, plastic and equipment, together we can contribute to the growth of nature and a cleaner living environment.

Share your photos, video's and stories about your personal masseur with us!

Tag us @bluelionlifestyle and we'll share your post with our followers so you can help and inspire each other.

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Download your digital version of this manual here

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Did you know that we are climate positive because we plant trees all around the world, that absorb more CO₂ than we emit with our shipments?

BLUE LION
LIFESTYLE

BLUE LION
LIFESTYLE

Blue Lion Lifestyle
De Smidse 21
6931KJ Westervoort
The Netherlands
www.bluelionlifestyle.nl

MARULA

© 2024 Blue Lion Lifestyle