ENGLISH
Your WildOnes® toothbrush has been carefully designed to offer your child a fun and unique brushing experience that is both safe and effective.

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USING THE APPLIANCE

- Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or a lack of experience and knowledge if they have been given supervision or instruction concerning use of appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- The toothbrush must only be charged using the supplied charger cable plugged into a list/certified power supply or other approved sources incorporating a standard 5V DC USB 2.0 charging socket.
- Periodically check the USB cord for damage. If cord is damaged please dispose of and replace. A damaged or non-functioning unit should no longer be used. Do not modify or repair the product. This may cause fire, electric shock or injury.
- Children should not play with the appliance.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.

DANGER

- Do not use as a toy.
- Always unplug this product immediately after using.
- Do not use whilst bathing.
- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Do not reach for a product that has fallen into water. Unplug immediately.


WARNING

- A product should never be left unattended when plugged in.
- Close supervision is necessary when this product is used by, or near children or persons with reduced physical, sensory or mental capabilities.
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged or dropped into water. Return the product to a service centre for examination and repair.
- Never use the product if the cord is frayed, frayed, or damaged.
- Keep the cord away from heated surfaces.
- Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair and the like.
- Never use whilst sleeping or drowsy.
- Never drop or insert any object into any opening or hose.
- Do not use outdoors or operate where aerosol (spray) products and being used or where oxygen is being administered.
- Follow all charging instructions and do not charge the battery pack or appliance outside of the temperature range specified in the instructions. Charging improperly or at temperatures outside of the specified range may damage the battery and increase the risk of fire.
- Do not incinerate this appliance even if its severely damaged. The batteries can explode in a fire.
- Only use appliance and battery in temperatures between 5-40°C. Store in temperatures between -10-20°C (up to 1 year) or -10-45°C (up to 3 months). Recommended ambient temperature range for charging system during charging is between 0-40°C.

SAVE THESE INSTRUCTIONS

- If the product is dropped, the brush head should be replaced before the next use even if no visible damage.
- Do not place the USB charger in water, liquid or store it where it can fall or be pulled into a tub or sink. Do not reach for it if it has fallen into water. Unplug immediately.
- This appliance contains batteries that are non-replaceable. Do not disassemble the product except when disposing of the battery. When taking out the battery for disposal of the unit, use caution not to short the positive (+) and negative (-) terminals.
- When unplugging, always hold the USB power plug instead of the cord. Do not touch the USB cable or plug in with wet hands. This can cause an electric shock.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- This toothbrush is a personal care device and not intended for use on multiple patients in dental practice or institution.

Description Connecting and Charging

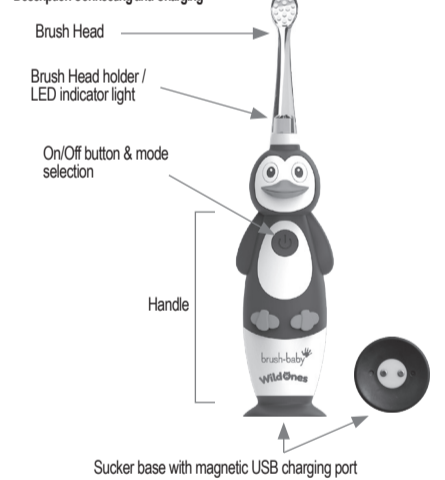The WildOnes® toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Connect the USB charging cable to a USB Power Plug and place USB magnetic coupling device to the receiver on toothbrush handle.
- The LED light flashes white indicating charging in process. Once it is fully charged the light turns off. A full charge can take up to 4 hours and enables up to 30 days of regular brushing (twice a day, 2 minutes). Note: After deep discharge, the white light might not flash immediately; it can take up to 10-15 minutes. Recharging is impossible.
- If the rechargeable battery is low the LED light will flash white for a few seconds when turning your toothbrush on/off.
- For everyday use, the toothbrush handle can be stored without connecting to the USB charging cable.
- To maintain the maximum capacity of the rechargeable battery, fully discharge the handle with regular use at least every 6 months.

Brushing modes. This produce has two brushing modes – Gentle and Moderate. Successively press the on/off button to move from the on mode through gentle mode, moderate mode and to off mode. The toothbrush pulses every 30seconds (to encourage brushing a different quarter of the mouth) and switches off after 2 minutes. Brushing techniques. Children must be supervised during use of this toothbrush. Wet brush head and apply Brush-baby toothpaste. To avoid splashing, guide the brush head to your teeth before switching on the appliance. Guide the toothbrush slowly from tooth to tooth, spending a few seconds on each tooth. Brushing in small circles as you would do with a manual toothbrush. Brush all surfaces of the tooth, move to a different quadrant of the mouth when the brush pulses (every 30 seconds). Do not brush too hard, simply let the toothbrush do all the work.

- Brush heads. Use only Brush-baby heads with the product.
- 3+ Brush head: With deep clean bristles specially designed for kids developing teeth. Good for gaps, missing teeth, mixed teeth, kids who won't floss and braces.

For best brushing results change the brush heads when the bristles become worn, played or if your child has a cold, is unwell or every 3 months for hygienic purposes. Replacements brush heads are available to purchase. Cleaning recommendations. After brushing: Rinse your brush head under running water for a few seconds with the brush switched on. Switch off brush, stand upright and allow to air dry. To clean the handle: Clean under running water and wipe dry. Never place in the dishwasher or water, clean with a damp cloth only.

Environmental notice: This product contains batteries and recyclable electric waste. Take electric waste to recycling collection points.

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