


User manual of hair removal instrument

Catalogue

- I. Product summary 2
- 1.1 Button description 2
- II. Principle exposition 2
- III. Use taboo 2
- IV. Security matters 4
- V. Introduction to use 4
- 5.1 Use preparation 4
- 5.2 Steps for usage 5
- 5.3 Skin detection 7
- 5.4 End of use 7
- 5.5 Clean and save 7
- VI. Troubleshooting 8
- VII. Product specifications 9
- VIII. Disclaimer 10
- XI. After-sales service card 11

I. Product summary

1.1 Button description

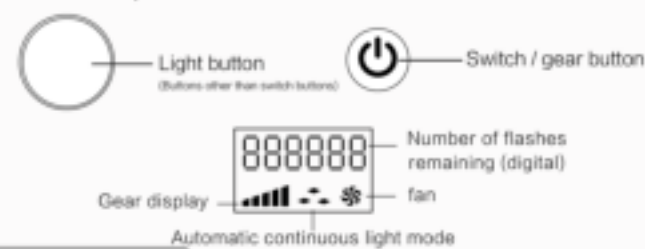

II. Principle exposition

The instrument adopts the international advanced photon cosmetic IPL technology, which directly irradiates the skin with strong pulsed light to produce the photochemistry effect. On the one hand can make the collagen fiber and the elastic fiber in the skin regenerate, restore the skin elasticity, and improve the facial microcirculation and eliminate or lighten the wrinkles. On the other hand, strong pulsed light can penetrate the skin and be absorbed selectively by the pigmented mass and capillaries in the tissue. On the premise of not destroying the normal skin tissue, heat energy makes the blood coagulation and pigmented mass destroyed by thermal decomposition, to treat capillary dilatation, remove color spots, and remove excess hair and so on.

III. Use taboo

Important safety information—please read before use

Important safety information—please read before use:

— Refer to the skin type on page 2 [skin color control picture]. If you are VI, do not use this instrument, otherwise it is likely to cause adverse reactions to the skin, such as pigmentation and depigmentation, severe redness and swelling or burning.



IV. Safety precautions

Special attention

- Do not apply any gels or care products (also veneers, etc.) before use.
- Be sure to assemble and use the product in accordance with the regulations contained in the instruction manual.
- Please read the "Instruction Manual" carefully and keep it for future reference.
- If you have any questions during use, please contact the online customer service in the required language.
- **DANGER: Contents related to personal injury due to incorrect operation**
- **WARNING: Dangerous content due to incorrect operation**
- **CAUTION: Contents related to personal safety or property damage due to incorrect operation**
- Warning symbol note: Do not use in the bath or shower.
- **WARNING:** Do not use this product near water-filled bathtubs, shower basins, washbasin or other water-containing utensils. Please keep this product and power adapter dry; if the product is damaged, please do not touch any internal parts to avoid electric shock;
- Do not insert any objects into the product, and do not place the product of will;
- Do not use this product around your eyes. Wearing protective glasses is a normal operation that protects the user. If you operate around your eyes, protective glasses may not completely block the energy, and the light will be absorbed by the iris melanin in your eyes at close range;
- We recommend that you wear protective glasses every time you use this product.
- Although light will not directly affect your eyes, high-frequency flashing light may affect your vision;
- When you need back hair removal, you may need an assistant or partner. This is our recommendation, because we are concerned that you may be at risk if you use it alone;
- If your hair is longer than 1 mm, please shave it first. If the excessive hair absorbs heat on the skin surface, it will cause skin redness and swelling, and the burnt hair will produce peculiar smell, which will affect your experience of using the product;
- If the skin in the irradiated area is not dry or clean, it will affect the use effect and may cause the product to not flash or accidentally flash;
- We lack data on the use of the product in children under the age of 15, so we recommend that you take care to prevent children under the age of 15 from contacting the product, and please stay away from children under the age of 15 when operating the instrument.
- please use different light intensity gears for different body parts as their complexions are different;
- If the power cord is damaged, to avoid danger, it must be replaced by a specialist from the manufacturer's maintenance department or similar department.



Unsuited community and situation

- Not suitable for pregnant women, lactating women, women's menstrual periods, body scars, diabetes, fragile skin and light-sensitive skin;
- Severe heart disease, moderate to severe hypertension, malignant tumors, folliculitis, infectious skin disease, burned skin and anticoagulant disorders are prohibited;
- Do not use if you have abnormal collagen metabolism, including history of scars and poor wound healing;
- Patients with vascular diseases, such as varicose veins and vasodilation in the irradiation site;
- Do not use if your skin is sensitive to light and prone to rashes or allergies;
- Please do not use if your irradiated area is infected, patient with herpes simplex, skin lesions and hematomata is prohibited;
- Suppression of immune response diseases (including HIV infection, lupus erythematosus, porphyria, etc.) is prohibited;
- Any patient with bleeding disorder is prohibited, and those who use isotretinoin, fruit acid (AHA), salicylic acid (BHA) or azelaic acid are prohibited;
- Please use this product with caution for those who use chemical hair removal within 30 days, those who have used skin whitening products in a short time, or those who have used vitamin A ingredients;
- Do not use with essential oil products at the same time. Please ensure that the skin in the light area is clean and dry, and no lotions, lubricants and other cosmetics are used.
- If the skin surface is not clean and dry, it may affect the use experience of the product and may even damage the product;
- prohibited during taking aspirin drugs, using immunosuppressive drugs, isotretinoin and other drugs. The recent use of light-sensitive drugs is prohibited (such as tetracycline and riflax).

Do not use for the following body parts or the following conditions

- 1/Intra-oral, genitals, mucosa, etc.
- 2/Plastic surgery
- 3/throat knot
- 4/eyeball, eyelid, around eyes
- 5/head
- 6/Acne area
- 7/Tattoo, body painting
- 8/eyebrows, hair
- 9/ear
- 10/Nipple, Areola
- 11/inside the navel
- 12/large vein
- 13/Tattoo, body painting

(Note: This list is not exhaustive, if you are not sure whether you can use this product we recommend that you consult a doctor)

Notes

1. Do not use the instrument with naturally dark skin
 2. Don't expose yourself before hair removal
 3. Before shaving, please shave first (do not use chemical hair removal method)
 4. Be sure to start in first gear when using for the first time
 5. When preparing hair removal, make sure the hair removal window is close to the skin and then perform the hair removal operation.
 6. Prohibited to use in dark environment
 7. Wear protective glasses when using
 8. Do not use it in direct sunlight to avoid machine failure.
 9. Do not violently impact or shake this product
 10. Do not unplug directly without turning off the machine to avoid damage to the machine.
- Notice: When using the hair removal device for a break, please make sure that the inlet is facing up and the LCD is facing down, which is beneficial to thermal dissipation and cooling.

V. Product specifications

Principle of instrument		IPL 光子美容技术	
Net weight of instrument		260g	
Irradiation area		3cm ²	
Instrument spectrum		≥530nm	
Rated voltage		12V DC	
Input power		18W	
Transportation and storage environment		Operating environment	
Storage temperature	-20~50℃	working temperature	5~40℃
humidity	10%~90%RH	humidity	10%~90%RH
Relative pressure	700~1060hPa	Relative pressure	700~1060hPa

Thank you for purchasing the product. For convenience, please read and keep this manual carefully before use.

VI Introduction to use

6.1 Use preparation

1. Please remove the hair from the surface of the skin first. Be sure to wet the hair, use the method of wet scraping to do a good skin treatment in advance, so that the hair below the surface of the skin can be absorbed into the light source to ensure the effect of the course of treatment. If the hair does not exceed 2 mm, you can operate the instrument without shaving.
2. Clean the skin and make sure there is no hair on its surface, completely dry and free of any oily substances. Before using the photon depilator, it is strictly forbidden to apply conditioner, aloe gel and other products, otherwise it will lead to no light (only large equipment in beauty salons need to be applied, but not for household use).

6.2 Steps for usage

1. Take out the instrument, please clean the head filter of the machine with a soft cloth.
2. Before using this instrument to depilate, please wet the hair and scrape it off with a razor.
3. Connect the power adapter and the machine display displays the word "OFF".
4. Press "Power on key" for a long time until the display screen is turned on (the machine is turned on when the content of the display screen changes).
5. When the content of the display changes from "flashing" to "always bright", it means that the hair removal device is ready (the first flash is bright, the eyes may not adapt, please wear goggles to reduce the discomfort of the flash to the eyes).
6. The content of the display screen is "always on", and the instrument is in the state of "flash ready"; if the content of the display screen has not changed, try to unplug the power adapter and restart the machine.
7. Press the depilator vertically on the skin at an angle of 90° to ensure that the skin is in full contact, and the "number" in the display screen will "flash".
(The use of hair removal instrument for the first time, there will be a "micro-current" feeling because of uncomfortable of skin, which is the normal reaction of the human body, not caused by machine leakage.)
8. After the machine is glued to the skin, gently press the "flash button" of the instrument, and the instrument will send out a flash, accompanied by a clear "crackle sound", and you will feel the warm feeling produced by the flash at the same time (if the temperature is too high, it means that the gear you have chosen is too high, it is recommended to start from the first gear).

9. After using the instrument to send out a flash, change the depilation part, when the energy is ready, the "gear indicator" in the display screen will be "always on", and you can repeat the above "step 7".
10. Use the first gear to test the skin. If there is no adverse reaction after 24 hours, you can adjust the energy gear through the "power on key", and the energy gear can be divided into 5 gears.
11. Follow these steps, you can depilate the part of your choice. One hair removal course is usually lighted three times. First take off the whole body, and then do it for the second and third time.
12. Please flash as close to the last flash as possible to fully deal with the desired part.
13. After using, please press "Power on key" for a long time to shut down (after shutdown, the product will enter the countdown to the shutdown), the screen will enter the word "OFF" and cut off the power supply. Finally, the head and fuselage are cleaned.

Model introduction

- Continuous play mode (for large areas such as arms/legs): Place the light outlet on the area to be treated on the skin and fully fit the skin. At this time, the hair removal sensor lamp will light up, indicating that the lamp head is in full contact with the skin. Long press the flash button, and immediately start moving the instrument slowly and uniformly on the skin, keeping the fit with the skin and the flash will be emitted at a certain interval. At this time, please keep the hair removal moving, and slide the hair back and forth twice in the area where hair removal is needed. To ensure that the hair removal area is covered evenly and avoid missing light.
Tip: The skin's movement speed can be adjusted according to the flash speed, so that the coverage is more uniform, the flash frequency is faster in the first gear, the flash frequency is slower in the second gear, and the higher the gear, the slower the light output speed.

- Single play mode (suitable for small areas such as fingers and toes): Place the light outlet on the area to be treated on the skin and fully fit the skin. At this time, the hair removal sensor lamp will light up, indicating that the lamp head is in full contact with the skin. Short press the flash button to start hair removal. Move the lamp head to the next piece of skin, wait for the white light of the hair removal sensor lamp to light up, and then press the flash button to repeat the above operation in order.
Tip: At any time during the hair removal process, if your skin begins to feel overheated, just pause for a few seconds, allow the light outlet to cool down again, and then perform hair removal.

Therapeutic cycle reference table

1-4 weeks	5-6 weeks	7-8 weeks	More than 8 weeks
2-3 times a week	Twice a week	Once a week	1-2 times a week
More than 75% of the hair begins to fall off naturally.	New hair falls off and inhibits growth.	The hair becomes thinner and less, no longer repeated.	Use it selectively according to the growth status of the hair (when the hair is growing, it

VII Care after use

precautions after hair removal

- It is recommended to avoid contact with water within 6 hours after the hair removal operation, and then moisturizing water can be used to moisturize the irradiated area.
- Do not use perfumes, masks, or cosmetics containing frosted beads after skin depilation. Avoid UV exposure for 24 hours after depilation, and apply sunscreen if necessary.
- Sunbathing is likely to cause skin darkening or sunburn. Irradiation at this time will cause redness, swelling, and inflammation. Therefore, during the use of this product, be careful not to be exposed to the sun, and the time to repair the sunburn is not less than 1 week.
- Artificial tanning is not allowed during the hair removal operation.
- If you need plastic surgery after hair removal (only for the hair removal area), consult your doctor.
- If you feel unwell, consult your surgeon.

Cleaning and maintenance

- In order to ensure the best results and the longest service life of the product, please clean the product before and after each depilation to avoid long-term accumulation of hair and dust. If you do not clean it properly, the product will not get its proper effect.
- Before starting cleaning, the device must be switched off and unplugged to allow the machine to cool.
- Use a soft, clean and dry cloth to wipe and clean the mirror surface of the light outlet.
- Store the device in a dry, dust-free environment at 10-60° C.
- Never wash the product and accessories in faucets and dishwashers.
- Do not use steel balls, abrasive cleaners (toothpaste, etc.) or corrosive liquids (gasoline, motor oil, acetone, etc.) to clean the product.

Note: If the water does not clean well, you can use a few drops of high concentration alcohol (not more than 96%) to clean the above parts.

Environmental protection

When disposing of the product, please do not throw it with general household waste. You should send it to a designated recycling center. This is good for the environment.



VIII. Troubleshooting

This chapter outlines the most common problems encountered in using this instrument. If you cannot resolve the problem by using the following information, please contact the customer service center in your country.

Problem description	Possible reasons	Solutions
The device/power adapter becomes warm during use	This is a normal phenomenon	There is no need to take any action.
I put the depilator on the skin and the number flashes, but the cooling fan started running.	Overheat protection has been activated.	When the overheat protection starts, the fan will run continuously. Do not turn off the power of the device, please wait for it to cool for more than 10 minutes before continuing to use.
[Flash ready] flashes, but the device does not flash when the flash button is pressed.	The device needs to be reset.	Reset the device, turn the power off and back on.
I can't increase or decrease the flash intensity setting.	The device needs to be reset.	Reset the device, turn the power off and back on.
It smells bad when using the device.	Accessories or filters are too dirty	Please clean the accessories carefully, please remove the hair on the surface of the skin before using this instrument.
	You have not completely removed the hair from the treatment site. These hairs may burn and produce a bad smell.	Once the UV filter is damaged, do not continue to use this device. Please contact the customer service center in your country.
	The light intensity setting you are using is too high	Please check to make sure that the appropriate light intensity setting is selected. Please lower the setting if necessary.
	You have not completely removed the hair from the treatment site.	Please remove the hair on the surface of the skin before using this instrument.
	The ultraviolet filter at the outlet is damaged.	Once the UV filter is damaged, do not continue to use this device. Please contact the customer service center in your country.
	You use this device where it is not applicable.	Do not use this device in the following areas: nipple, areola, labia minora, vagina, anus, as well as the inside of nostrils and ears. Male users should not use the product on the face or scrotum.
	This device is not suitable for your skin.	The skin type is VI or V, brown-black or dark skin, please do not use this device.

The skin reaction was longer than usual at the end of the course of treatment. I am not satisfied with the result of the course of treatment.	The light intensity setting you are using is too high.	Next time, choose a lower light intensity. Please refer to the section "skin testing" in the "instructions for use" unit.
	The light intensity setting you are using is too low.	Please choose a higher setting next time.
	You have not done again around the area where you have been treated before.	Please emit a flash as close as possible to the previous part.
	This device has no effect on your hair color.	You need to successfully remove all hair. It is recommended that you follow the recommended treatment procedure.
	The frequency of using the instrument is insufficient.	It is recommended that you follow the recommended treatment procedure. You can shorten the treatment interval, but not more than once every two weeks.
	Your response to IPL treatment may be slow.	Continue to use the device for at least 6 months, because the hair will grow again during this period, but will continue to decrease.

IX Disclaimer

From the date of purchase, the warranty for the whole machine within 12 months any of the followings, it cannot be treated as a warranty even within the warranty period:

1. Failure to operate in accordance with the requirements of the Manual, maintenance and maintenance resulting in product damage
2. Product damage caused by during use misoperation (washing, immersion, etc.)
3. Product damage caused by exceeding the allowable use environment (such as high temperature, corrosion, etc.).
4. Product damage caused by abnormal external forces (fall, extrusion, impact, etc.).
5. Product damage caused by natural force majeure (fire, flood, lightning strike, etc.).
6. Damage caused by the use of non-factory parts and extensions.

If there are any changes listed above, the relevant regulations newly formulated by our company or the relevant laws and regulations updated by the state shall prevail.

X After-sales service card

Dear customer

First of all, thank you very much for choosing this product, so that we have the opportunity to provide you with the best service. In order to use the service of the instrument make you more satisfied, please read this description carefully and take good care of this warranty card after purchase

If you have any questions, please communicate with us in time and we will solve them for you as soon as possible!

1. If there is a quality problem causing the goods to be unusable, please contact the after-sales customer service of the seller of the products you buy in time to provide you with after-sales service.
2. We do not accept any surface mail and cash on delivery, and the special after-sales customer service will communicate with you to solve the postage.
3. Please fill in the after-sales service guarantee card patiently and return it with the goods so that we can deal with it as soon as possible.

Real name: (required) _____ Contact number: (required) _____

Purchase channels: _____ Purchase account: _____

Order number: _____

After-sales selection: (required) Return goods Exchange goods

Receiving address: (required) _____

Return instructions: _____