

the brush

interchangeable brush head


For optimal cleaning, replace your brush head every 8 weeks. When changing the brush head, remove it in a straight vertical motion. Avoid twisting or bending. The brush head is securely attached to the handle, with a small gap between the handle and the brush head.

battery

The battery has a 90-day battery life when used for 2x two-minute brushing sessions each day. Please note that the battery is not fully charged upon delivery. The indicator will turn red when the battery is running low, but the brush will continue to function for a few more days.

charging

The brush is IPX7 certified, making it waterproof. Ensure a dry connection while charging. The brush will only turn on when it has reached 67% charge after being completely empty. A full charge takes 10-12 hours, indicated by the switch from red to green before automatically turning off. The brush can't be used during charging. Use a 5V adapter, avoid fast chargers, and employ a USB to USB-C cable.

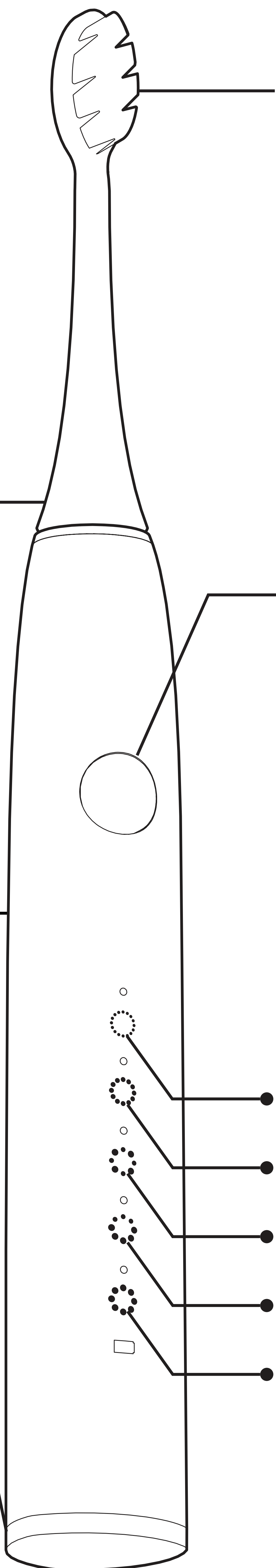

how to brush

Position the Boombrush on your teeth and start brushing using a regular brushing motion. Ensure that you hold the brush at a slight angle toward the gums ($\pm 45^\circ$) to ensure proper gum brushing. Avoid applying excessive pressure during brushing. It may take some time to get used to the brushing technique (typically 3-5 days), but afterward, you won't want anything else ;)

the button

Press the button briefly to activate the brush and switch between modes. The brush will remember your last setting once you're finished. After 4 cycles of 30 seconds, the brush will stop automatically. Vibrations briefly pause every 30 seconds to signal that you can switch quadrants.

the five modes

- General brushing mode
- Soft vibrations for sensitive teeth
- Massaging mode for tongue or gums
- Alternating vibrations
- Powerful mode to remove plaque



Scan the QR code if you want to read the instructions in other languages.