# auronic

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## Introduction

- Thank you for choosing a Auronic product! Please always follow the instructions for safe usage. If you have any questions about the product or experience any problems, please contact us by e-mail: [service@auronic.nl](mailto:service@auronic.nl).

### Intended use

- This product is only intended for the purposes described in the manual.
- The supplier is released from liability for any injury or damage caused by using the product outside of its intended use.
- Always observe the following precautions when using 'the massage pillow', hereinafter referred to as 'the product'.

## Warnings and safety instructions

- ⚠ Please ensure that the packaging is kept out of the reach of children! Risk of suffocation!
- ⚠ Make sure to keep this user manual. Read the manual carefully!

### General

- If this is your first time using the product, it is important to read the manual and instructions carefully.
- If the product is not working properly, has been dropped or damaged, do not use the product and contact our customer service.
- Before each use, inspect the parts and accessories for damage or wear (cracks/breaks) that could impair function and if so, replace them with original parts and accessories recommended by the manufacturer.
- Do not use accessories not offered by the supplier. They may pose a danger to the user and damage the product.
- Use the product only for the purposes indicated in this manual.
- Do not disassemble, repair or repair the product. This may damage the product. For your own safety, these parts should only be replaced by authorized service personnel. Always send the product to the service center for repairs to avoid damage and personal injury.
- Be aware of the risks and consequences associated with incorrect use of the product. Improper use of the product can lead to personal injury and damage to the product.
- The product is designed and intended for private use. Do not use the product for commercial purposes.
- Keep the product away from heat, areas where the temperature is relatively high and heat sources such as ovens.
- Turn off the product and remove the plug from the socket when you are not using or cleaning the product. Let the product cool down completely before cleaning the product.
- Make sure the product does not get wet.

### Persons

- Cleaning and maintenance of the product should not be carried out by children unless they are older than 8 years and supervised.
- Do not let children play with the product. Use and store the product out of the reach of children.
- The product should not be used by children. Use of the product is permitted only when a supervisor is present and responsible for the safety of the user.
- The product must not be used by people with physical or mental disabilities. Use of the product is allowed only when a supervisor is present, who is responsible for the safety of the user.
- Keep the product and cord out of the reach of children under 8 years old.
- Prolonged use on a high setting may result in burns.
- We recommend using the product only if you are 16 years or older.

### Use

- Do not cover the product when it is on. Never place pillows or blankets on or against the product. This could result in a fire hazard, electric shock, personal injury or damage to the product.
- Do not clean the product while it is still plugged in.
- Do not use the product near a bath, shower, swimming pool, or near a basin or sink filled with water. Do not touch the product and power cord with wet hands. Make sure that the product does not come into contact with water. Do not immerse the product in water. If water does get in or on the product, turn the product off immediately and contact our customer service.
- Use and store the product so that the cables cannot be damaged.
- Parts of the product become hot. Do not clean the product while it is still hot. Wait until the product has cooled down before cleaning it.
- Allow the product to cool completely before storing.
- Do not fall asleep while using the product.
- Do not use the product before bedtime. The massage has a stimulating effect and may make it difficult to fall asleep.
- After 15 minutes, the product turns off automatically. Let the product cool down before using it again.

### Medical notes

- If you experience pain or discomfort while using the product, stop immediately and consult your doctor.
- This product is not a substitute for treatment from a physician or physiotherapist. Chronic complaints can worsen.
- Do not use the product if you have circulation problems and do not use the product in areas that are sensitive to heat. This can lead to burns!
- Do not use the product or consult a doctor: if you are pregnant; if you have a pacemaker or other electronic implants, artificial joints, high blood pressure, heart problems, cerebral atherosclerosis, cerebral thrombosis, fractures, malignant tumors, circulation problems, varicose veins, open wounds, bruises, broken skin, inflammation of the veins, jaundice, diabetes, have nerve disorders, acute inflammation or other health problems; during or immediately after strenuous exercise; immediately after showering, taking a bath, or after the sauna; if you have drunk alcohol; during or immediately after eating; if you have a fever; if you have unexplained pain, are undergoing medical treatment or use medical equipment; if you have one or more of the following conditions: circulation problems, varicose veins, open wounds, bruises, broken skin, inflammation of the veins.
- Consult a physician before using the product if you have unexplained pain, are undergoing medical treatment or use medical equipment.

### Electronics

- If the power fails during use, turn off the product immediately and unplug the power cord.
- If the surface of the product is cracked, please turn off the product to avoid electric shock.
- Periodically inspect both power cord and plug for signs of damage or wear, and if any, do not operate the product and have the power cord and/or plug replaced or repaired by authorized service personnel.
- To avoid overheating, explosion, electric shock, fire, other damage and personal injury, use only the adapter and/or accessories supplied by the supplier.
- To avoid electric shock, do not use the product if you suspect that the power cord or plug is damaged. A damaged power cord or plug must be repaired or replaced by a professional (an authorized distributor).
- Unplug the product and turn off the product when not in use and before cleaning the product.
- Do not unplug the product by pulling on the power cord. This may cause fire or electric shock.
- Do not let the cable hang over the edge of a table or counter.
- Do not drop the product or subject it to large shocks.
- Do not place the cable on hot surfaces or near oil.
- Do not tuck the power cord under rugs or blankets. Make sure that the cable is never in a walkway.
- Do not move the product by pulling or twisting the cord.
- Before first use, remove all packaging material and any stickers.
- Before connecting the product, make sure that the mains voltage indicated on the type rating plate corresponds to that of your mains supply.
- Make sure that you do not trip over the cable and that the cable is not bent, twisted, or caught.

- If black smoke comes out of the product, unplug the power cord immediately. Wait for the smoke to clear before removing the contents of the product and contact our customer service.
- Do not insert objects into the openings of the product. Make sure that the moving parts are not blocked and that nothing is pinched.
- Do not leave the product unattended when in use. When not using the product, turn off the product and unplug the power cord.

### Receiving your shipment

- When receiving your shipment, make sure that the contents are consistent with the packing list. Notify your distributor of any missing items. If the equipment appears to be damaged, file a claim immediately with the carrier and notify your distributor at once, giving a detailed description of any damage. Save the damaged packing container to substantiate your claim.

### Contents of package

- Auronic Massage pillow - Model AU3030 - EAN 8720195254799


#### Included accessories

- 1x User Manual
- 1x Power adapter

### Product features

#### Product overview

1. Power button
2. Leather cover
3. Massage head
4. Mesh cover
5. DC connector
6. Euro plug
7. Power cable
8. DC plug
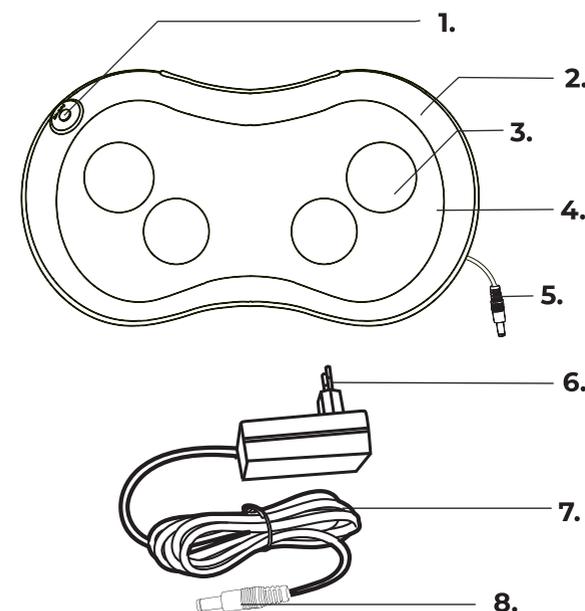
## Technical data

Dimension	36 x 17 x 12.5 cm	Output massage pillow	5A
Weight	1 kg	Frequency	50/60Hz
Color	Black	Current	24W
Material	PU leather and fabric	Automatic switch off	Yes, after 15 min
Cord length	1,8 m	Overheating protection	Yes
Input adapter	220 - 240V~	Operating temperature	-10°C to 45°C
Input massage pillow	12V	Storage Temperature	-20°C to 60°C

### Other features

- Suitable for legs, upper and lower back, neck, shoulders, feet, calves and arms;
- Shiatsu massage type;
- With 4 rotating massage heads;
- Includes heat function;
- With adjustable strap to attach to your headrest or chair.

## Operating instructions

Turn on the product by pressing the power button. You can switch between different functions with the power button.

- Press once: Turn the product on or off.
- Press twice: Change the rotation direction of the massage heads.
- Press three times: Turn the heat function on or off.

## Maintenance and cleaning

- Unplug the product from the electrical outlet before cleaning. Clean the product with a damp, clean cloth and dry it thoroughly. Do not immerse the product in water or other liquids. Do not use abrasive or harsh cleaning agents.

## Storage

- Store the product in a clean and dry place. Do not store the device in extreme high or extreme low temperatures.

## Disposal - recycle

- Think about the environment and contribute to a cleaner living environment! This product is subject to the European Directive (EU) 2012/19/EU. At the end of its service life, do not dispose of the product with normal household waste, but deliver it to a collection point designated by the authorities for recycling.

## Service and warranty

- If you need service or information regarding your product, please contact our Customer Support at [service@auronic.nl](mailto:service@auronic.nl).
- We provide a 2-year warranty on our products. To obtain service during the warranty period, the product needs to be returned with a proof of purchase. Product defects have to be reported within 2 years from the date of purchase. This warranty does not cover damages caused by misuse or incorrect repair; parts subject to wear; defects that the customer was aware of at the time of purchase; damages or defects caused by customer neglect; damages or defects caused by third parties.

## Disclaimer

- Changes reserved; specifications are subject to change without stating grounds for doing so.

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© LifeGoods B.V.