

TECHNICAL SPECIFICATIONS

Name: New concept massager
 Rated voltage: DC-12.6V
 Rated frequency: 24W(including heating 5W)
 Time setting: 10 minutes
 Size:45*25*9cm
 NW: 1.5 KGS
 Others: one charger, one user manual



WARNING

If the charger is damaged and needs to be changed over, you must replace it with the same specification charger to ensure safety in use (power DC-12.6V).

Manufacturer reserves the right to make changes to information published in this document, including without limitation specifications and product descriptions at any time and without notice.

NEW CONCEPT MASSAGER




Before using, please read this User manual carefully

PRODUCT DESCRIPTION

This product is designed for convenience and comfort. Our unique, heated quad-roller is designed to massage the desired area of the user. The Cordless Massage Pillow gives a comfortable cushion as a throw pillow on one side and provides a relaxing therapeutic massage on the other side.

Made with durable and tear resistant material combined with easily operated remote button on the side of the massage pillow.

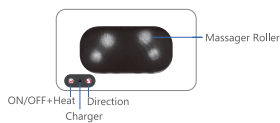

Massager Roller

ON/OFF+Heat
 Direction
 Charger

SAFETY PRECAUTIONS

Recommended Use

- It is suggested to keep the massage to 15-20 minute sessions, and gradually allow your body to get comfortable with longer sessions as needed and to extend the life of the massage pillow.
- If the massage ball intensity feels too strong, you can put a towel or blanket between the body and massage ball to reduce intensity.

Charger

- After fully charged, you should pull the plug out the jack by pulling the plug, do not pull the wire.

Locating place

- Do not operate the massage under where dusty, greasy, damp or smoggy conditions.
- Do not place near steam oven, open fire or heat.
- Do not place directly in the sun.
- Do not use under damp place(such as swimming pool, bathroom).
- Do not use the massager outdoor.
- Do not use in smog or oxygen lacking place.

Users

- If you have any injuries of any kind, consult a doctor before using.
- This massager not suit for children to use.
- People with long hair, please not put the hair between body and massager.
- When the massager is working, people shouldn't leave it alone.

HOW TO USE

- When charging, the charger light displays "red".
- Plug in Charger and allow the massage pillow to fully charger.
- When fully charged, the charger light displays"green".
- Hold on the "ON/OFF" switch for three seconds, the pillow will start to work. Press it one more time again for heat function ON/OFF.
- Hold on the "Direction" switch for mode one kneading. Press it one more time again for mode two kneading.
- Under normal use, the massager will automatically stop working in about 10 minutes, if you want to continue to use, hold the switch for 3 seconds again.

TROUBLESHOOTING GUIDE

Problems	Possible Cause	Trouble shooting method
Sounds like tapping or hum sounds	Normal phenomenon	Normal phenomenon
Massage sounds noise	Whether continues 30minutes	Turn power switch off and stop 1 hour
Hot or Wont cool down	Over heating	Turn off and let pillow rest