
Instructions EN

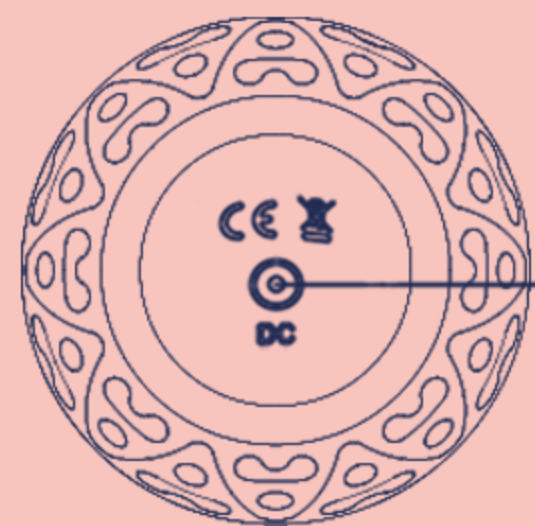
DC charging port

(A) Press and hold 1,5 sec to turn the Massageball ON/OFF + Vibration modes

(B) Press and hold 1,5 sec to turn on heating function (up to 42 °C)

USB charging plug in PC, or in 5V power supply

DC charging plug, in Massageball

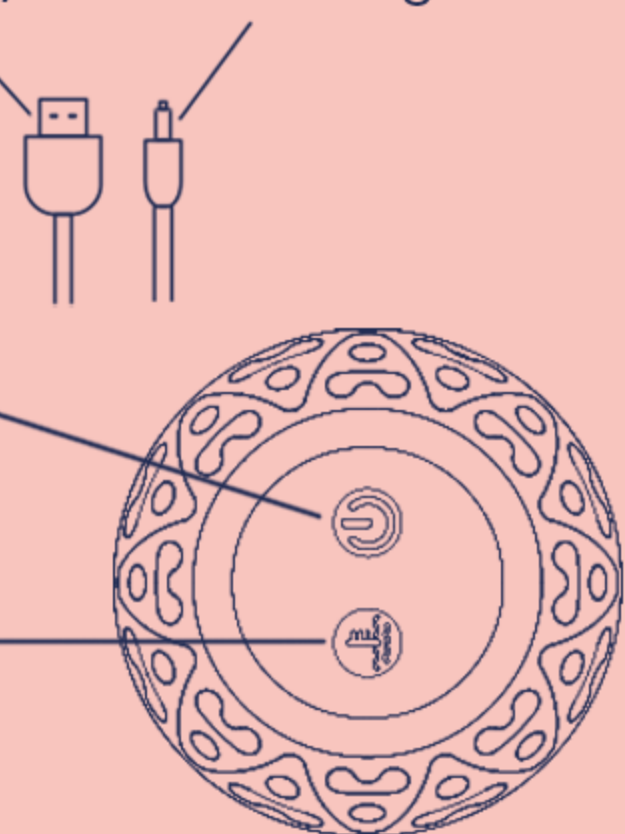Usage

- Fully charge the product prior to first use.
- Press and hold button A for 1,5 seconds to turn the product on/off.
- Press button A to navigate through the 10 different vibration modes.
- Press and hold button B for 1,5 seconds to switch the heating function on/off. When turned on, the middle area of the Massageball will start heating over a couple of minutes, up to 42 °C.

Recharge

- Put the DC charging plug in the port of the product to charge. Subsequently, put the USB charging plug in the USB-port of your PC or another 5V power supply. (Output 5V=1A or 5V=1.5A)
- When the product is charging, a small LED light will blink on the control panel.
- When the LED light remains on continuously, the product is fully charged.
- Do not use the product while it is charging. While charging, the product should remain off.

▲ Safety warnings

- Use this product consciously. In case you suffer from a medical (skin)disease or other physical circumstances which are not natural, consult a doctor prior to using this product.
- Stop using this product if any discomfort occurs or when the appearance of the product is no longer in original state.
- Charge this product on a safe location. Do not charge, for instance, near water or heat.
- Use this product with care. Do not throw or drop from high altitude.
- Not suitable for use underwater. IP65 verified: this means splash waterproof so the product can be cleaned, for instance, by holding it under water jets from the water tap.
- To clean the product, use warm water and antibacterial soap. Never clean it with organic cleaner that contains alcohol, gasoline or acetone. Avoid boiling water.
- When not using, keep this product in a cool, dry and clean place. Avoid direct sunshine.
- This product contains a lithium battery. Only dispose of it at designated waste collection points.

Disclaimer: This product is for personal use. Only use this product as intended. No medical claims are warranted or implied by the use of this product. Sold as an adult massage device only. Keep out of reach of children. All rights reserved.

Our mission

Life is too beautiful not to make the most of it. An important part of that is to feel good about yourself and this is what we want to contribute to for everyone.

We strive for a world with more love and intimacy.

Love starts with yourself, and to contribute to this we provide you with products that you can be proud of; user-friendly, innovative and overwhelming.

We are here to let you enjoy the best things in life.

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
Specifications EN

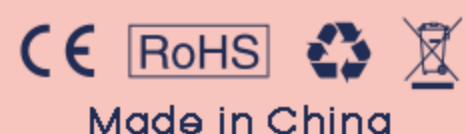Model: MFM001 / YAI88W 003
Material: bodysafe silicone, ABS
Weight: 415g
Size: 177 x 81 mm
Volume: <50db
Battery: 1400 mAh rechargeable built-in lithium battery
Power: USB 5V 1A/1.5A DC rechargeable
Charging time: approximately 3 hours
Working time: approximately 1-2 hours
IP65 splashwaterproof

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