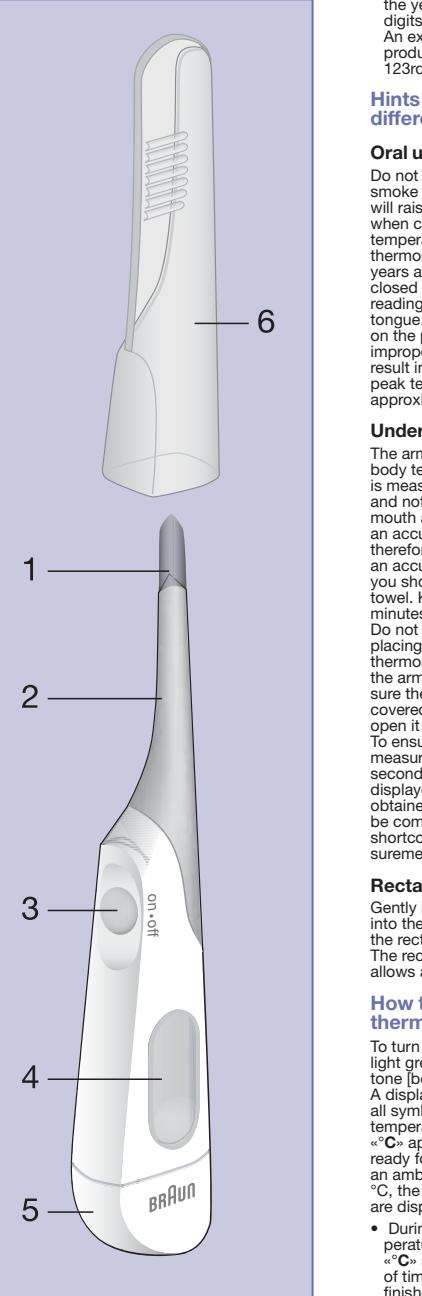
BRAUNHigh speed
thermometer

GB English
DE Deutsch
DK Dansk
FI Suomi
FR Français
GR Ελληνικά
IT Italiano
NL Nederlands
NO Norsk
SE Svenska

Protected by the following patent nos.:
 US 6522,388
 EU 294426-0001

This product is manufactured by
 Kaz Europe SA under a license
 to the 'Braun' trademark. 'Braun' is a
 registered trademark of Braun GmbH,
 Kronberg, Germany.

Kaz Europe SA
 Place Châtelier 18
 CH-1003 Lausanne
 Switzerland

©2011, Kaz Europe SA
 www.kaz-europe.com
 Designed in Germany, Made in PRC.

PRT1000CE
 REV 27/MAY/11 P/N: 31M1000205

C 0297

the year of manufacture. The last three digits are the manufacturer designation. An example: LOT 12311ONB – This product has been manufactured on the 123rd day of the year 2011.

Hints when measuring at different sites

Oral use

Do not drink hot or cold fluids, exercise, smoke or perform other activities that may affect your temperature readings when compared to your normal body temperature. We recommend using the thermometer orally only with children four years and older. The mouth should remain open during measurement. Place the probe well under the tongue. Instruct children not to bite down on the probe. Opening the mouth or improper placement of the tip could result in a longer measurement time. The peak temperature should be reached in approximately 20 seconds.

Underarm/axillary use

The armpit is not an ideal site to measure body temperature, because temperature is measured near the surface of the skin and not in a closed cavity (such as mouth and rectum). The time required for an accurate underarm temperature may therefore take 90 seconds or longer. For an accurate temperature measurement, the thermometer must be placed in a dry towel. Keep the armpit closed for several minutes before placing the thermometer. Do not open the armpit too long before placing the thermometer. Place the thermometer in the armpit and gently press the arm tightly against the body. Make sure the measuring sensor is completely covered under the armpit and do not open it during the measurement process. To ensure an accurate underarm (axillary) reading, wait for 90 seconds or longer 20 seconds before reading the temperature displayed. Please note that the reading obtained immediately after the beep might be compromised by the above mentioned shortcomings related to the axillary measurement site.

Rectal use

Gently insert the tip of the thermometer into the rectum. Do not force the tip into the rectum if resistance is encountered. The rectum is a closed body cavity and allows a faster reading.

How to use the Braun digital thermometer

To turn on the thermometer, push the light [beep] signal. Thermometer on, short signaling tone sounds [beep]. Thermometer off, no signal on the display. In an ambient temperature below 32 °C ---- and a «C» appear and the thermometer is now ready for use on different body sites. (In the temperature higher than 32 °C, the temperature and a flashing «C» are displayed).

• During measurement, the current temperature is displayed continuously. The «C» symbol flashes during this period of time. The measuring process is finalized by a short signaling tone sounds [beep... beep... beep 10 times] and the «C» is no longer flashing.

• To prolong the life of the battery, turn the thermometer off by pressing the on/off button briefly. Otherwise the thermometer will automatically turn off after about 10 minutes.

Memory mode

• If you press the on/off button for more than 3 seconds when turning the thermometer on, the last measured temperature, which was stored automatically, will be displayed. At the same time, a small memory will appear at the right of the display. After 3 more seconds, this value will disappear and the thermometer will switch off. The stored value will be replaced by the new value, which is also stored automatically.

As with any other thermometer, proper handling and storage is necessary in order to obtain accurate temperature readings. When used properly, this thermometer will give you years of accurate temperature readings that you can depend on. Please read the use instructions thoroughly. The use of this thermometer is not intended as a substitute for consultation with your doctor.

Safety instructions

• This thermometer is intended for household use only.

• The thermometer is only suitable for the measurement of body temperature!

• Always adhere to the minimum measurement time (until the signaling tone [beep] is heard).

• Keep the thermometer away from children.

• Do not expose the thermometer to mechanical shock. Use the protective cap to prevent damage.

• Do not bend the tip more than 45°!

• Do not expose to temperatures above 60 °C. Never boil the thermometer! Keep the thermometer away from direct sunlight or high temperatures.

• Only use commercial disinfectants to clean the thermometer. The thermometer must be cleaned when immersed in liquid disinfectants.

• Attention: The thermometer contains a CR1615 3.6V battery. The battery life is associated with the lifetime of the thermometer. For your safety the battery compartment lid is securely attached to the thermometer and cannot be removed without a standard tool. For disposal, please comply with local regulations for special waste.

• The thermometer is initially calibrated at the time of manufacture. If the thermometer is used according to the use instructions, the accuracy of measurements will not be affected. If at any time you question the accuracy of the temperature measurement, please contact an authorized Kaz Service Center.

• Manufacturing date is given by the LOT number located in the battery compartment. The first three digits after LOT represent the day of the year of manufacture. The second two digits give

the year of manufacture. The last three digits are the manufacturer designation. An example: LOT 12311ONB – This product has been manufactured on the 123rd day of the year 2011.

Product specifications

Displayed temperature range

32,0 °C to 42,9 °C

Operating temperature

+10 °C to +40 °C

Operating humidity

up to 85% r.h.

Measurement accuracy
 ± 0,1 °C between 35,5 °C and 42 °C body temperature at an ambient temperature of 18 °C to 28 °C

Display

Liquid crystal display (LCD) with three digits. Smallest unit of display: 0,1 °C

Recommended storage

-25°C to +55°C

Battery

1,55 V

If the device is not used within specified temperature and humidity ranges the technical accuracy of the measurement cannot be guaranteed.

Subject to change without notice.

This appliance conforms to the following standards: EN 60061-1:3/96 – Medical electrical equipment

- Part 1: General requirements for safety EN 60061-2:3/96 – Medical electrical equipment

- Part 1-2: General compatibility – requirements and tests EN 12470-3:2000

- Clinical thermometers – Part 3: Performance of contact thermometers for oral, rectal and non-invasive use with maximum device MEDICAL ELECTRICAL EQUIPMENT needs special precautions regarding EMC. For detailed description of EMC requirements please contact an authorized local Service Centre (See Package Insert).

Portable and mobile RF communications equipment can affect MEDICAL ELECTRICAL EQUIPMENT

ASTM Standard 1112-00 – Standard Specification for Electronic Thermometer for intermittent determination of patient temperature

This product conforms to the provisions of the EC directive 93/42/EEC (Medical Device Directive).

Please do not dispose of the product in the household waste at the end of its life cycle. Disposal can take place at your local retailer or at appropriate collection points provided in your country.

ThermoScan ear thermometer

For newborns and children in particular, we also recommend the Braun ThermoScan ear thermometer for a notably quick, gentle and accurate temperature measurement. The ear is an excellent site to measure body temperature, because the eardrum sends blood supply with the temperature control center in the brain. The Braun ThermoScan ear thermometer uses infrared technology to detect the infrared energy emitted by the eardrum and the surrounding tissues within seconds. The brand most used by doctors*.

How to use the Braun digital thermometer

To turn on the thermometer, push the light [beep] signal. Thermometer on, short signaling tone sounds [beep]. Thermometer off, no signal on the display. In an ambient temperature below 32 °C ---- and a «C» appear and the thermometer is now ready for use on different body sites. (In the temperature higher than 32 °C, the temperature and a flashing «C» are displayed).

• During measurement, the current temperature is displayed continuously. The «C» symbol flashes during this period of time. The measuring process is finalized by a short signaling tone sounds [beep... beep... beep 10 times] and the «C» is no longer flashing.

• To prolong the life of the battery, turn the thermometer off by pressing the on/off button briefly. Otherwise the thermometer will automatically turn off after about 10 minutes.

Memory mode

• If you press the on/off button for more than 3 seconds when turning the thermometer on, the last measured temperature, which was stored automatically, will be displayed. At the same time, a small memory will appear at the right of the display. After 3 more seconds, this value will disappear and the thermometer will switch off. The stored value will be replaced by the new value, which is also stored automatically.

As with any other thermometer, proper handling and storage is necessary in order to obtain accurate temperature readings. When used properly, this thermometer will give you years of accurate temperature readings that you can depend on. Please read the use instructions thoroughly. The use of this thermometer is not intended as a substitute for consultation with your doctor.

Safety instructions

• This thermometer is intended for household use only.

• The thermometer is only suitable for the measurement of body temperature!

• Always adhere to the minimum measurement time (until the signaling tone [beep] is heard).

• Keep the thermometer away from children.

• Do not expose the thermometer to mechanical shock. Use the protective cap to prevent damage.

• Do not bend the tip more than 45°!

• Do not expose to temperatures above 60 °C. Never boil the thermometer! Keep the thermometer away from direct sunlight or high temperatures.

• Only use commercial disinfectants to clean the thermometer. The thermometer must be cleaned when immersed in liquid disinfectants.

• Attention: The thermometer contains a CR1615 3.6V battery. The battery life is associated with the lifetime of the thermometer. For your safety the battery compartment lid is securely attached to the thermometer and cannot be removed without a standard tool. For disposal, please comply with local regulations for special waste.

• The thermometer is initially calibrated at the time of manufacture. If the thermometer is used according to the use instructions, the accuracy of measurements will not be affected. If at any time you question the accuracy of the temperature measurement, please contact an authorized Kaz Service Center.

• Manufacturing date is given by the LOT number located in the battery compartment. The first three digits after LOT represent the day of the year of manufacture. The second two digits give

the year of manufacture. The last three digits are the manufacturer designation. An example: LOT 12311ONB – This product has been manufactured on the 123rd day of the year 2011.

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