Body Fat Scale BFS18W

Instruction Manual

English / Nederlands / Deutsch

moa
color your life

www.moacolors.com

BEFORE USE

Before use, please read this manual carefully.



Only connect the appliance to an earthed wall socket.

This appliance is for household use only. When using electrical appliances, basic safety precautions should always be followed to avoid fires, electric shocks, burns, or other injuries and damage. Read these operating and safety instructions carefully.

SPECIFICATIONS

- Equipped with high precision strain gauge sensor system Capacity: -
- Division: 100g/0.2lb
- 150kg/330lb
- Height unit: cm/fit-in
- Unit:kg-lb
- BMI/Body fat/ Muscle/Water division: 0,1%;
- Body bone division:0,1kg
- Mini effective weight:8kg
- Age range for body fat: 10-80 years
- 5 mm tempered safety glass platform
- Height range: 100-220 cm/3 ft 3in-7 ft 3 in
- Automatic power generation control
- Body bone range; 0-15kg
- Battery-free technology
- Muscle range; 13.0-100.0%
- Auto zero resetting / Auto power off
- Fat range: 5-50%
- Low power / Over-load indication
- Water range: 35-75%

PARTS & LOCATION

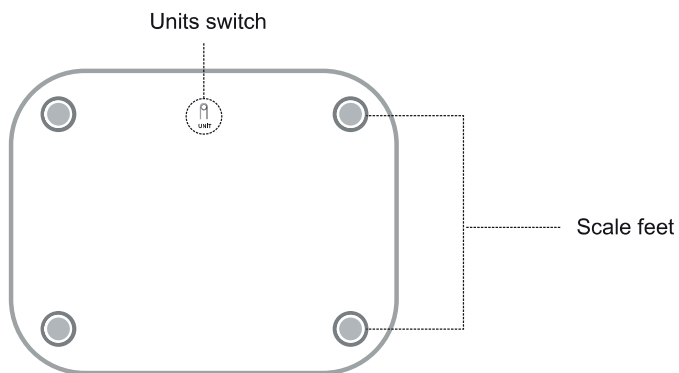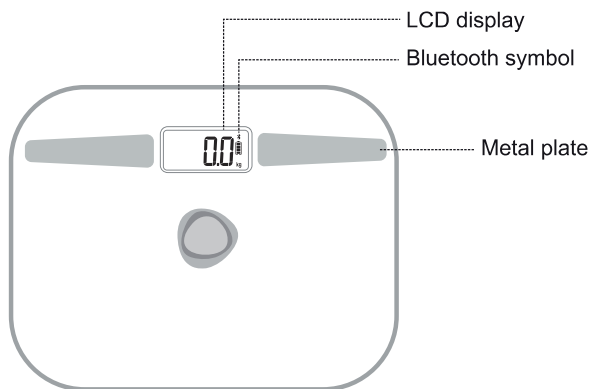APP

COMPATIBILITY

Battery free Bluetooth Smart Scale is made for the following models:

iOS: Above iPhone 5

Andriod: Above Android 6.0 verion

“AirScale U” DOWNLOAD

Download the app here:

Android

iOS

BLEUTOOTH SETTING

Open the Bluetooth of the iOS/Andriod device, go to Bluetooth setting.

USER SETTING

Press “+” and fill in Name, Gender, Birthday, Height, Goal, Weight and tick in the upper right cornor. Go to your registered mailbox for validation.

SETTING AND REVIEW


Press the hamburger-menu icon you can review your homepage, body's composition, health chart and set unit and slect language.

CHANGE LANGUAGE

“Settings” in Aircscale U app, select “Multilanguage”

HOW TO USE

"ECO-Power" One quick stamping on the power button by your feet and the scale charges for use.

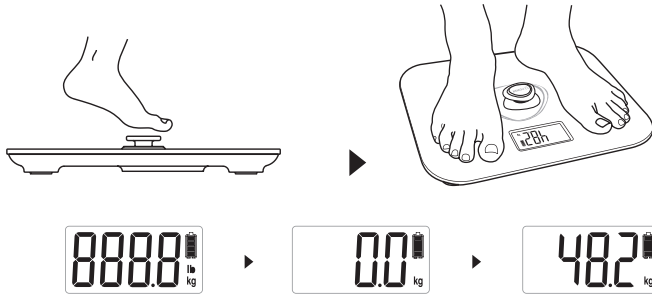1. Put the scale on a hard & flat surface.
2. One quick stamping on the ECO-power button will generate power for weighing, please wait until "0.0" appears.
3. The scale powers off automatically after one minute.

NOTE:

- Do not step on the scale until "0.0" appears.
- Please measure with bare feet and keep the scale clean
- For an accurate result, please keep your feet in a good contacting with the metal electrode.

FUNCTION

ERROR INDICATION

Over-load Indication

When the scale is over-loaded [more than 150kg], the screen will appear "Err".

Low battery indication

When the screen indicates "00" it will turn off automatically after few seconds. Simply reload it by stamping on the ECO-power button !

UNIT EXCHANGE (kg/lb)

Under weight mode, press the unit key to choose the unit, the key is on the back of the scale. Note: Please set the unit before weighing. It could not change the unit while weighing.

USE AND CARE

1. For accurate results, place scale on a flat surface.
2. Static pressure is effective, do not drop it or strike the scale, otherwise the scale is easy to damage.
3. In weighing, please keep balance. Do not allow water to get inside scale.
4. Please use a damp cloth to clean the scale. Do not allow water to get inside scale.
5. Do not use chemicals to clean the scale.
6. To ensure the life of the scale, do not put the scale under wet & hot condition.
7. Do not use for legal trade.
8. Do not overload, otherwise the sensor will be damaged.
9. Do not put heavy things on the scale when not use.
10. Do not store the scale in an upright position.
11. Due to the continuous improvements, this product is subject to change without notice.
12. Please do not jump or stamp on the scale and avoid dropping it, as the scale is highly precise electronic instrument.
13. Do not use it on the wet and slippery floor, do not stand on at the edge of scale to avoid any falling down.

ATTENTION FOR SAFETY

- The scale should not be used for anyone who has a pacemaker, or with other medical transplanting. The impedance passed through our body will have disturbance of the normal working of medical equipment.
- This scale is for family use only. We do not advise you to use it as medical equipment.
- Do not use it on a wet or slippery floor. Do not stand at the edge of the scale to avoid falling and injury.
- We advise pregnant women, elderly persons and disabled persons to use the scale with help.

DISPOSAL

Do not dispose of the device in normal domestic waste.
Dispose of the device over a registered waste disposal firm or through your communal waste disposal facility.
Observe the currently valid regulations. In case of doubt, consult your waste disposal facility.

WARRANTY

Dear Valued Customer,

Thank you very much for purchasing a MOA product.

We wish you to inform you that this product is covered by a warranty which complies with all legal provisions concerning existing warranty and consumer rights in the country where the product was purchased.

Should you find any defect or malfunction of your MOA product, please contact the appropriate Customer Care Center.

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