MEDISANA®

home of wellness

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
IMPORTANT INFORMATION! RETAIN FOR FUTURE USE!

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

Explanation of symbols

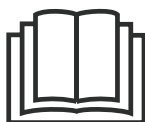This instruction manual belongs to this device. It contains important information about starting up and operation. Read the instruction manual thoroughly. Non-observance of these instructions can result in serious injury or damage to the device.

WARNING

These warning notes must be observed to prevent any injury to the user.

CAUTION

These notes must be observed to prevent any damage to the device.

NOTE

These notes give you useful additional information on the installation or operation.

Protection category II

Lot number

Manufacturer

about the power supply

Safety Information

- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- Only insert the mains plug into a socket when the appliance is switched off.
- Keep the cable and the appliance away from heat, hot surfaces, moisture and liquids. Do not touch the mains plug when you are standing in water and always take hold of the plug with dry hands.
- Appliance components carrying voltage may not come into contact with liquid.
- Do not take hold of a device if it falls into water. Pull out the power plug from the socket immediately.
- Connect the unit so that the power plug is freely accessible.
- Always switch off the appliance after use with the buttons on the control unit and unplug the power plug.
- To disconnect the unit from the power supply, always pull the plug out of the socket. Do not pull on the cable!
- Do not carry, pull or twist the device on the cable.
- You must not continue to use the unit if the lead or the plug are damaged.
- If the power supply cable is damaged, it may only be replaced by **MEDISANA**, an authorised dealer or by qualified personnel.
- Ensure that the cable will not cause tripping. The cable must not be kinked, trapped or twisted.

specific groups of people

- This device is not designed to be used by persons (including children) with limited physical, sensory or mental abilities, or by persons with insufficient experience and/or knowledge, unless under observation by a person responsible for their safety, or unless they have been instructed in the use of the device.
- Children must be supervised to ensure that they do not play with the device.
- Do not use this unit to supplement or replace medical treatment. Chronic complaints and symptoms may be made worse.
- Do not use the tapping massager **MNT**, or only after consulting your doctor, if:
 - you are pregnant,
 - you have a cardiac pacemaker, artificial joints or electronic implants,
 - you suffer from one or more of the following illnesses and complaints:
Circulation disorders, varicose veins, open wounds, bruises, broken skin, inflammation of the veins.
- Do not use the appliance near your eyes or near any other sensitive areas.

- If you experience any pain or discomfort during the massage, stop immediately and consult your doctor.
- Do not use the appliance on small children or vulnerable individuals.
- Before using the massager, consult your doctor regarding any unexplained pain, if you are undergoing medical treatment or using medical equipment.

before using the unit

- Check the plug, cable and massager carefully for damage before each use. Do not attempt to use a defective unit.
- Do not use the appliance if it, or the power cord show signs of damage, if it is not in perfect working order, or if it has been dropped or become damp. In order to prevent hazards, always send the unit in to the service centre.

operating the unit

The device is not to be used for commercial use or medical purposes. If you have health concerns, consult your doctor before using the device.

Only use the massager indoors.

Do not use the massager in damp rooms (e.g. when bathing or showering).

- The device must only be used for its intended purpose as described in the instruction manual.
- Using the unit for any other purpose invalidates the warranty.
- Only use the appliance in the proper place, that is around the neck and on the shoulders.
- Do not use the appliance while you are asleep or when you are lying in bed.
- Do not use the appliance before going to sleep. The tapping massage has a stimulating effect.
- Never use the appliance while driving a vehicle or operating a machine.
- The maximum operating time for one session is 10 minutes. Prolonged use reduces the service life of the appliance.
- Very long periods of use can lead to overheating. So after continuously operating for 10 minutes, user should stop it for **at least** 10 minutes before using it again.
- Never put the unit down next to an electric stove or other sources of heat.
- Never leave the unit unattended if it is connected to the power supply.

- Do not allow the appliance to come into contact with pointed or sharp objects.
- Do not cover up the appliance when it is switched on. Never use it under blankets or cushions. There is a risk of fire, electric shock or injury.

maintenance and cleaning

- The unit is maintenance free. If a fault does occur, just check whether the mains power plug is properly connected.
- You may only clean the device yourself. If the device malfunctions please do not repair it yourself as this leads to the automatic expiry of all guarantee rights. Ask your specialist dealer and have repairs carried out by the service only.
- Do not immerse the device in water or other fluids.
- If, nonetheless, liquid should penetrate the device, unplug the mains cable immediately.

Do not wash!

Do not dry clean!

2 Useful Information

Thank you very much

Thank you for your confidence in us and congratulations!

You have gained a high-quality product from **MEDISANA** with the purchase of this intensive tapping massager **MNT**. This appliance is intended for the tapping massage of the neck and shoulder area. To derive the most success and to enjoy your **MEDISANA** intensive tapping massager **MNT** for a very long time, we recommend that you read the following information on the use and care of the appliance carefully.

2.1 Items supplied and packaging

Please check first of all that the device is complete and is not damaged in any way. In case of doubt, do not use it and contact your dealer or your service centre.

The following parts are included:

- 1 **MEDISANA** Intensive tapping massager **MNT**
- 1 Instruction manual

The packaging can be reused or recycled. Please dispose properly of any packaging material no longer required.

If you notice any transport damage during unpacking, please contact your dealer without delay.

WARNING

Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!

The intensive tapping massager **MNT** is intended for targeted, deep-efficiently tapping massage treatment of the neck and shoulder area. The tapping massage is a relaxing massage applied to the nerve tracts to loose blockings and improve the blood circulation. The unit has 5 levels of intensity to choose from and 3 different program modes can be selected to match different purposes of treatment.

CAUTION

Make sure that the maximum operating time of 10 minutes is not exceeded.

3.1 Using

1. Insert the power plug **3** into the mains wall outlet.
2. Turn the knob of the position control **1** in clockwise direction to lower the distance of the two massage heads. Turning the knob in anti-clockwise direction increases the distance of the massage heads. Adjust the position of the massage heads according to your personal need.
3. Place the massager around the neck so that it has sufficient hold on the shoulders.
4. Push the ON-/OFF-button "⏻" **7** to switch on the massage function. The appliance starts working straight away at low intensity in mode 1.
5. To change the massage mode, press the selection button **3**, **4** or **5**. The 3 modes differ by the tapping frequency and the massage rhythm:

Mode 1:
Tapping massage mode to improve the deep-efficient relaxation; 1 to 16 tapping massage movements per second - alternatively left and right and then both sides. Intensity will increase after some time.

Mode 2:
Tapping massage mode to improve the blood circulation and to remove blockings; 1 to 5 tapping massage movements per second - the intensity and the massage speed vary constantly.

Mode 3:
Tapping massage mode for loosing tense muscles; 1 to 10 tapping massage movements per second with constant changing massage rhythms.
6. To increase the massage intensity press the button "+" on the intensity control **8**. Pressing the "-" button decreases the massage intensity. The intensity can be regulated in 5 levels. According to the selected level the LEDs will light up in the intensity control section **8**: 1 LED for level 1, 2 LEDs for level 2, 3 LEDs for level 3, 4 LEDs for level 4 and 5 LEDs for level 5 (highest level).
7. Press the ON-/OFF-button "⏻" **7** to switch off the massage. After 10 minutes of continuous use the appliance will switch off automatically.

4.1
Cleaning and maintenance

- Before cleaning the appliance, make sure that the appliance is switched off and the power supply is disconnected. Allow the unit to cool down.
- Clean the massager with a slightly damp sponge only. Never use brushes, aggressive cleaning agents, petrol, thinners or alcohol. Rub the unit dry with a clean dry cloth.
- Never immerse the appliance in water to clean it and ensure water cannot get into the appliance.
- Wait until the device is completely dry before using it again.
- Undo any kinks in the cable.
- Take care to roll up the mains lead to avoid damaging it.
- You should ideally store the device in its original packaging and keep it in a clean, dry place.

4.2
Disposal

This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner.

Please remove the batteries before disposing of the device / unit. Do not dispose of old batteries with your household waste, but at a battery collection station at a recycling site or in a shop. Consult your municipal authority or your dealer for information about disposal.

4.3
Technical specifications

Name and model	: MEDISANA Intensive tapping massager MNT
Voltage supply	: 230V~ / 50Hz
Power consumption	: 75 W
Automatic switch-off	: After approx. 10 minutes
Operating conditions	: Only use in dry rooms
Storage conditions	: In a cool, dry place
Dimensions	: Approx. 61 x 43 x 6 cm
Weight	: Approx. 2,6 kg
Power cable length	: Approx. 2,0 m
Article No.	: 88940
EAN Code	: 40 15588 88940 0

In accordance with our policy of continual product improvement, we reserve the right to make technical and optical changes without notice.

Warranty and repair terms Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.
The following warranty terms apply:

1. The warranty period for **MEDISANA** products is three years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

**MEDISANA AG**

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41468 NEUSS

Germany

eMail: info@medisana.deInternet: www.medisana.com

The service centre address is shown on the attached leaflet.

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