
MEDISANA®

home of wellness

- DE** Shiatsu-Massagekissen ^{MPD}
- GB** Shiatsu massage cushion ^{MPD}
- FR** Coussin de massage Shiatsu ^{MPD}
- IT** Cuscino per massaggio Shiatsu ^{MPD}
- ES** Almohadilla de masaje de Shiatsu ^{MPD}
- PT** Almofada de massagem de Shiatsu ^{MPD}
- NL** Shiatsu-massagekussen ^{MPD}
- FI** Shiatsu-hierontatyyny ^{MPD}
- SE** Shiatsu-massagekudde ^{MPD}
- GR** Μαξιλάρι Shiatsu μασάζ ^{MPD}

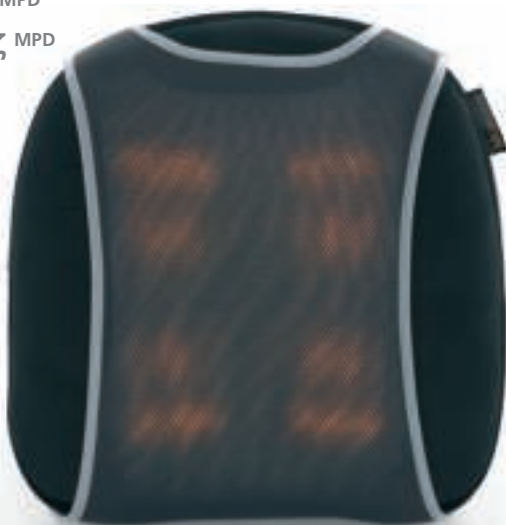Art. 88908

Gebrauchsanweisung
Instruction Manual
Mode d'emploi
Istruzioni per l'uso
Instrucciones de manejo

Bitte sorgfältig lesen!
Please read carefully!
A lire attentivement s.v.p.!
Da leggere con attenzione!
¡Por favor lea con cuidado!

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IMPORTANT INFORMATION! RETAIN FOR FUTURE USE!

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

Explanation of symbols

This instruction manual belongs to this device. It contains important information about starting up and operation. Read the instruction manual thoroughly. Non-observance of these instructions can result in serious injury or damage to the device.

WARNING

These warning notes must be observed to prevent any injury to the user.

CAUTION

These notes must be observed to prevent any damage to the device.

NOTE

These notes give you useful additional information on the installation or operation.

Protection category II

Lot number

Manufacturer

GB

1 Safety Instructions

MEDISANA[®]
home of wellness

Safety Instructions

about the power supply

- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- Only connect the mains plug to the wall socket when the device is switched off.
- Keep the power cord and the device itself away from heat, hot surfaces, damp and liquids. Avoid touching the plug or power switch with wet or even damp hands, or when you are standing in water.
- If the device falls into water, do not try to pick it up. Instead pull the plug out immediately.
- Disconnect the plug from the socket immediately after using the device.
- Never pull on the power cord to disconnect the appliance from the mains. Always pull out the mains plug!
- Do not carry, pull or rotate the device by its power cord.
- Ensure that the cable will not cause tripping. The cable must not be kinked, trapped or twisted.

specific groups of people

- This device is not designed to be used by persons (including children) with limited physical, sensory or mental abilities, or by persons with insufficient experience and/or knowledge, unless under observation by a person responsible for their safety, or unless they have been instructed in the use of the device.
- Children must be supervised to ensure that they do not play with the device.
- The surface of the appliance is hot. People who are insensitive to heat must be careful when using the appliance.
- Do not use this unit to supplement or replace medical treatment. Chronic complaints and symptoms may be made worse.
- You should refrain from using the Shiatsu massage cushion **MPD**, or should consult your doctor before using it, if
 - you are pregnant,
 - you have a cardiac pacemaker, artificial joints or electronic implants,
 - you suffer from one or more of the following illnesses and complaints: Circulation disorders, Diabetes, Thromboses, varicose veins, open wounds, bruises, broken skin, inflammation of the veins.
- If you experience any pain or discomfort during the massage, stop immediately and consult your doctor.
- Never leave the unit unattended if it is connected to the power supply.
- Take care when using the heater. Do not use if you have circulatory problems or on areas not sensitive to heat. **Risk of burning!**
- Consult your doctor if you have any questions concerning the therapeutic use of the Shiatsu massage cushion.
- If you have unexplained pain, are receiving treatment from a doctor and/or use medical equipment, consult your doctor before using the Shiatsu massage cushion.
- Use of the massage device is prohibited in the following cases:
 - not on children, - not on animals,
 - never while you sleep, - not in the car,
 - not while ingesting substances that reduce reaction times (e.g. painkiller medications, alcohol).

before using the unit

- Always carefully inspect the cable, the control unit and the Shiatsu massage cushion for damage before use. Do not attempt to use a defective unit.
- Do not use the unit if there are visible signs of damage on the unit or the cable components, the unit does not work properly or the seat or the control unit have been dropped or become damp. In order to prevent hazards, always send the unit in to the service centre.
- The zip fastener on the massage device must not be opened. It is only fitted for reasons relating to the manufacturing of the device.

operating the unit

The device is only intended for household use and must not be used for commercial or medical purposes.

If you have health concerns, consult your doctor before using the Shiatsu massage cushion.

The massage cushion must only be used in enclosed spaces!

Do not use the Shiatsu massage cushion in damp rooms (e.g. when bathing or showering).

- The Shiatsu massage cushion must only be used for its intended purpose as described in the instruction manual.
- Using the Shiatsu massage cushion for any other purpose invalidates the warranty.
- Never leave the unit unattended if it is connected to the power supply.
- The maximum operating time for one session is 15 minutes. Allow it to cool off for at least 15 minutes before using it again.
- Do not allow the Shiatsu massage cushion to come into contact with pointed or sharp objects.
- Never put the unit down next to an electric stove or other sources of heat.
- Never use the device under covers, such as blankets or pillows.
- For back massages operate the device in an upright position only. For legs massages operate the device in a flat position, to place the legs on it comfortably.
- Never sit or stand on the moving parts of the device with all your weight and do not place objects on the device.
- Do not jam or force any part of your body, in particular your fingers, between the rotating massage heads and their mountings in the device.
- In order to prevent yourself from being trapped between the heads, do NOT lean on the massage heads when changing your position or moving the massage device.
- Make sure that all moving parts can move freely at all times.

maintenance and cleaning

- The unit is maintenance free.
- You may only clean the device yourself. If the device malfunctions please do not repair it yourself as this leads to the automatic expiry of all guarantee rights. Ask your specialist dealer and have repairs carried out by the service.

GB

1 Safety instructions / 2 Useful Information

- Do not use the unit if there are visible signs of damage on the unit or the cable components, the unit does not work properly or the massage cushion have been dropped or become damp. In order to prevent hazards, always send the unit in to the service centre.
- Do not immerse the device in water or other fluids.
- If, nonetheless, liquid should penetrate the device, unplug the mains cable immediately.

Do not wash!

Do not dry clean!

2 Useful Information

Thank you very much

Thank you for your confidence in us and congratulations!
You have acquired a **MEDISANA** quality product with this Shiatsu massage cushion **MPD**.

This unit is designed for massaging back and legs.

We recommend that you carefully read the following information on use and care of the product so that you experience its benefits and enjoy using your **MEDISANA** Shiatsu massage cushion **MPD** for many years.

2.1 Items supplied and packaging

Please check first of all that the device is complete and is not damaged in any way. In case of doubt, do not use it and contact your dealer or your service centre.

The following parts are included:

- 1 **MEDISANA** Shiatsu massage cushion **MPD**
- 1 Cover
- 1 Instruction manual

The packaging can be reused or recycled. Please dispose properly of any packaging material no longer required.

If you notice any transport damage during unpacking, please contact your dealer without delay.

WARNING

Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!

2.2 What is Shiatsu massage?

Shiatsu is a form of finger pressure massage and one of the most important massaging techniques along with manual lymphatic drainage and reflexology. This form of physical therapy was developed in Japan but stems from traditional Chinese massage. Gentle stroking and beneficial pressure application techniques are used to awaken life forces and encourage them to flow through the body. The object of the massage is not just to ease tension, but also to optimise the flow of life force and harmonise body and soul.

2.3 How does the MEDISANA Shiatsu massage cushion work?

In the **MEDISANA** Shiatsu massage cushion **MPD** you have a modern device for genuine Shiatsu style back and legs massages. The **MEDISANA** Shiatsu massage cushion **MPD** has four three-dimensional rotating twin massage heads ❶ for Shiatsu massage of back and legs. The device is equipped with an additional heat function in the massage heads. The unique combination of intensive shiatsu massage with radiated heat makes the massage cushion an effective and relaxing treatment.

3 Operating

3.1 Usage

The use of the Shiatsu massage cushion will improve your wellbeing after a stressful day. The massage effect is particularly useful for treating tense muscles and tired tissue. The unit provides pleasant relaxation for overcoming everyday stress. The massaging can improve your well-being, and is also useful in the sport and fitness areas.

3.2 Operation

Start-up

- Fasten the protective cover ❷ to the surface of the massage cushion with the Velcro strips. Just as well the cover is detachable for easy cleaning. It is also possible to massage without the protective cover; this is generally found to be "harder".
- Place the massage cushion on the backrest of a high chair or armchair and attach it to this by means of the Velcro straps ❷ on the back. Ensure that the cushion is supported over its entire area and has a sufficiently high backrest. Alternatively place the cushion on an adequate base (e.g. couch, sofa, or similar) to massage your legs.
- Insert plug ❸ into the mains socket. Ensure that the power supply is easily accessible.
- Take a seat and gently press your back against the the Shiatsu massage cushion to familiarise yourself with the functions.
- Sit comfortably and upright. Ensure you sit centered and such that the massage heads ❶ can move to the left and right of your spinal column.

Operating the unit

The **MPD** shiatsu massage cushion has two functions for an effective massage: fingertip massage (shiatsu) and heat. The heat function can only be used together with the shiatsu massage. Do not use the unit for more than 15 minutes at a time. Always allow the unit to cool down fully before you use it again. Settings can be adjusted by means of switch **3**.

Switch **3**

Use the switch **3** on the side to operate the massage device.

The switch positions „I“ and „II“ determine the setting of function:

- Switch position **I**: Shiatsu massage ON
the massage heads swing
- Switch position **II**: Shiatsu massage with heat ON
the massage heads swing with heat radiated simultaneously by means of illuminated massage heads
- Switch position **0**: device OFF

Do not use the massage device for longer than 15 minutes. Using the device for longer can overstimulate the muscles and cause tension instead of relaxation.

To avoid overheating the device, the maximum operating time for one session is 15 minutes. Allow it to cool off for at least 15 minutes before using it again.

The massage can be finished at any time by moving the switch **3 to position 0.**

4.1 Cleaning and maintenance

- Before cleaning the appliance, make sure that the appliance is switched off and the power supply unit is disconnected. Allow the unit to cool down.
- Clean the massage cushion with a slightly damp sponge only. Never use brushes, aggressive cleaning agents, petrol, thinners or alcohol. Rub the seat dry with a clean dry cloth.
- Never immerse the appliance in water to clean it and ensure water cannot get into the appliance.
- Wait until the device is completely dry before using it again.
- Undo any kinks in the cable.
- Take care when rolling up the mains lead to avoid damaging it.
- You should ideally store the device in its original packaging and keep it in a clean, dry place.

4.2 Disposal

This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

4.3 Technical specifications

Name and model	: MEDISANA Shiatsu massage cushion MPD
Voltage supply	: 220 - 240 V~ 50 - 60 Hz
Power consumption	: approx. 30 W
Operating time	: max. 15 minutes
Operating conditions	: only use in dry rooms
Storage conditions	: in a cool, dry place
Dimensions (L x W x H)	: approx. 37.5 x 37 x 11.5 cm
Weight	: approx. 2.95 kg
Power cable length	: approx. 2.2 m
Article No.	: 88908
EAN Code	: 40 15588 88908 0

In accordance with our policy of continual product improvement, we reserve the right to make technical and optical changes without notice.

GB 5 Warranty**MEDISANA**
home of wellness**Warranty and
repair terms**

Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

1. The warranty period for **MEDISANA** products is three years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

 **MEDISANA AG**
Jagenbergstraße 19
41468 NEUSS
GERMANY
eMail: info@medisana.de
Internet: www.medisana.de

The service centre address is shown on the attached leaflet.

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