Waterflosser

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Inhoud

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Read the instructions before use

Contents

1. About the product	P3
2. In the box	P3
3. Important safety precautions	P3
4. Product overview	P4-5
5. Operating procedure	P6
6. User manual	P7
7. Recommended techniques	P8
8. Disassembly	P9
9. Demontage voorbeeld	P10
10. Maintenance and handling of the battery	P11
11. Troubleshooting	P12

2

1 About the product

Thank you for purchasing our TIGIOO Water Flosser! To ensure that this product is used correctly and safely, please read the user manual carefully before using it for the first time. Follow these instructions for safe use. Should you have any questions or problems with the product, please contact the seller/customer service.

2 In the box

1x	Waterflosser
2x	Nozzles
1x	Charging cable
1x	Charging plug
1x	User manual

3 Important safety precautions

When using electrical products, especially in the presence of children, the safety precautions listed below must always be followed. Read all instructions before using the water flosser.

Reducing the risk of electrocution

- Do not touch the charger with wet hands.
- Do not immerse the product or charger in water or other liquids.
- Do not use the product while using the bath or shower.
- Store the product in a place where it cannot be pulled into a bath, shower stall or wash basin.
- Do not touch the product if it has fallen into the water, unplug it immediately.
- Check that the charging cable is not damaged before using the product.

3

4 Warning

To reduce the risk of burns, electric shock, fire, serious injury or damage

- Do not plug the product into a power supply with a voltage other than that indicated on the product or on the charger.
- Do not plug the product into a power supply with a voltage other than that indicated on the product or on the charger.
- Do not use the product if the charging cable or plug is damaged. Please contact the seller or our Customer Service.
- Do not point water under the tongue, in the ears, nose or other sensitive areas. This product is capable of producing a high pressure that can cause serious injury in these areas. Refer to user instructions for the correct way to use the product.
- Use this product only as directed in this manual, as recommended by your dentist.
- Bleeding gums are common when using the TIGIOO Water Flosser for the first time due to the condition of your gums.
- Fill the reservoir with water or dentist recommended liquid only.
- Do not drop or insert foreign objects into the opening of the product.
- Keep the charging cable and plug away from heated surfaces.
- Do not use the product in areas where oxygen or any other aerosol is used.
- Do not use mouth rinses containing iodine, bleach or tea tree oil.
- Do not wear oral jewellery, remove it before use.

4


4 Warning

- Do not use this product if you have an open wound on your tongue or in your mouth.
- If your doctor has advised you to undergo pre-medical treatment with antibiotics prior to dental surgery, you should consult your dentist before using this product or any other hygiene aid.
- This product is for domestic use only and is not suitable for outdoor use.
- Instruct and supervise children and persons with special needs when using the product.
- Children should be supervised to ensure that they do not play with the product.
- The product should not be used by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they are supervised or have been instructed.
- The product is powered by a removable battery, the product can only be used with the charger <ACA1AAUEWH4>. Input: 100-240V ~, 50/ 60Hz, 0.8A, output: 5.0V, 1.0A
- Only the charger <ACA1AAUEWH4> can be used for charging the battery.
- This product contains batteries that may only be replaced by qualified personnel.
- This product contains batteries that cannot be replaced.
- This symbol means detachable charger.

5

5 Product overview

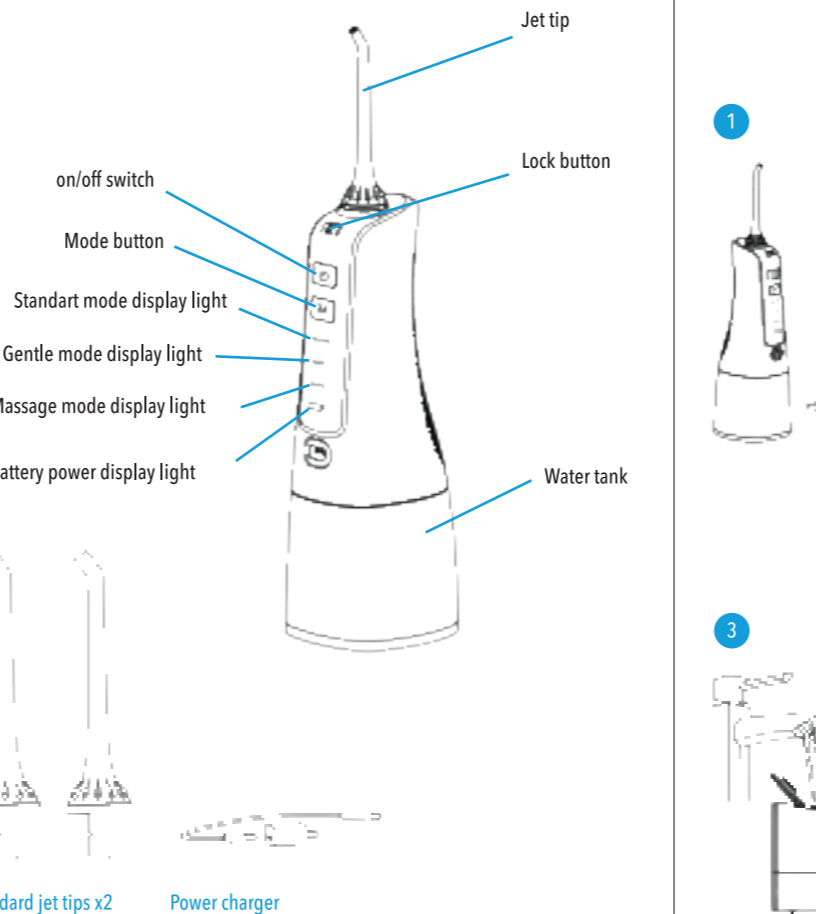
Standard jet tips x2 Power charger

6

6 Operating procedure

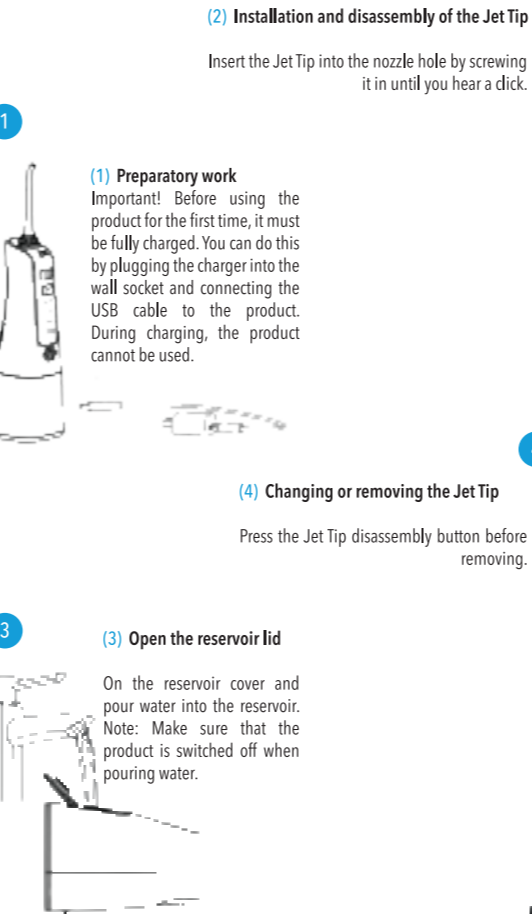7

7 User manual

- This product has a memory function and when the product is turned back on, it will be in the same mode as before the product was turned off.
- The Standard, Soft and Massage modes can be changed when the product is turned off.
- The mode the product is in can be changed before the product is turned on so that you can choose a mode you like.
- Regardless of the mode, the product turns itself off automatically after 120 seconds.
- When the massage mode is used for more than 3 seconds, the remaining energy will be displayed.
- During charging, the 3 display lights will flash. If all 3 remain on, this indicates that the battery is fully charged.
- During use, when all 3 display lights flash, it means that the battery is low and the product needs to be recharged.
- During charging, the product cannot be used.
- Turn the Jet Tip dial to change the direction of the water.

8

8 Recommended techniques

- Hold the product at a downward angle, put the Jet Tip in your mouth and point to your teeth before turning on the product.
- When the water flow and the gums are at a 90 degree angle to each other, you can close your mouth to avoid splashing water and achieve the best result. Please start by clenching your teeth and start at the incisors. Guide the water stream along the gums and pause between teeth until the area around the gums and the spaces between the teeth are completely cleaned.


9

9 Disassembly

Please follow the steps below.

- Turn off the product and empty the reservoir.
- Use both hands as shown in Figure 2 - use your left hand to hold the product and your right hand to hold the reservoir. At the same time, rotate the reservoir in a downward, counterclockwise motion to remove it.

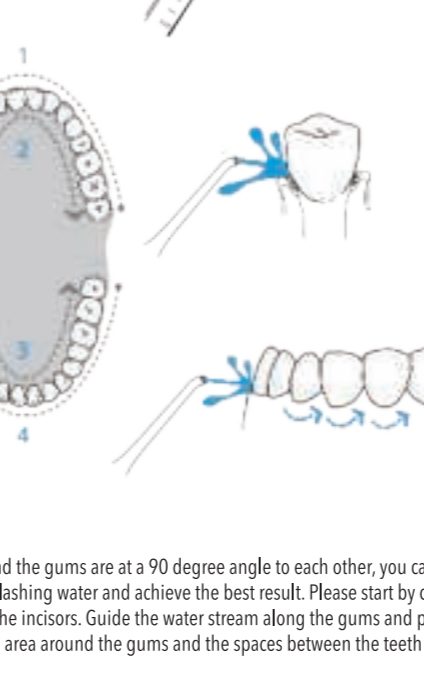10

10 Maintenance and handling of the battery

Battery maintenance

- Before using the product for the first time, it must be fully charged. This takes about 4 hours.
- When storing the product for a longer period of time, the product should be fully charged again before use.
- To keep the battery life as long as possible, the battery should be charged in time. A long time without use will negatively influence the battery life.

Handling the battery

The product uses a lithium battery and is not replaceable. At the end of the battery's service life, please return it to a collection point in your area.

11

11 Troubleshooting

Problems	Cause	Solution
Water leaking between the Jet Tip and the screw part	The sealing ring on the screw top is damaged and the Jet Tip can no longer be fixed properly	Please contact the retailer for repair, and install the Jet Tip properly
Reservoir is leaking	The sealing ring is loose	Reinstall the sealing ring correctly
Not enough water pressure	Low voltage	Charging the product
The product does not turn on	Low voltage	Charging the product
No water comes out of the Jet Tip	The filter of the reservoir is clogged	Clean the filter. Do this regularly

If the methods mentioned above do not solve the problems, you can always contact customer service.

Notice

- During the first use, only use the "Soft mode" to avoid bleeding gums, this is because the gums are not yet used to using this product. After using it daily for 2 weeks, progress can already be seen. If the gums still bleed a lot after two weeks, this may be a case of chronic gum inflammation. It is then recommended to contact a dentist.
- It is not recommended to use the product on children, people with dental implants or with open sores in the mouth.

12

Contact

Questions?

If you have questions about the product or are you experiencing problems with the product? Then please contact the seller or our customer service.

Let us know your opinion!

Are you satisfied with our product? Then we would like to hear from you!

By scanning the QR code on the right, you can easily rate our product. It would mean a lot to us!

Thank you and have fun with your Water Flosser!

Notice

- During the first use, only use the "Soft mode" to avoid bleeding gums, this is because the gums are not yet used to using this product. After using it daily for 2 weeks, progress can already be seen. If the gums still bleed a lot after two weeks, this may be a case of chronic gum inflammation. It is then recommended to contact a dentist.
- It is not recommended to use the product on children, people with dental implants or with open sores in the mouth.

12

Contact

Questions?

If you have questions about the product or are you experiencing problems with the product? Then please contact the seller or our customer service.

Let us know your opinion!

Are you satisfied with our product? Then we would like to hear from you!

By scanning the QR code on the right, you can easily rate our product. It would mean a lot to us!

Thank you and have fun with your Water Flosser!

Notice

- During the first use, only use the "Soft mode" to avoid bleeding gums, this is because the gums are not yet used to using this product. After using it daily for 2 weeks, progress can already be seen. If the gums still bleed a lot after two weeks, this may be a case of chronic gum inflammation. It is then recommended to contact a dentist.
- It is not recommended to use the product on children, people with dental implants or with open sores in the mouth.

12

Contact

Questions?

If you have questions about the product or are you experiencing problems with the product? Then please contact the seller or our customer service.

Let us know your opinion!

Are you satisfied with our product? Then we would like to hear from you!

By scanning the QR code on the right, you can easily rate our product. It would mean a lot to us!

Thank you and have fun with your Water Flosser!

Notice

- During the first use, only use the "Soft mode" to avoid bleeding gums, this is because the gums are not yet used to using this product. After using it daily for 2 weeks, progress can already be seen. If the gums still bleed a lot after two weeks, this may be a case of chronic gum inflammation. It is then recommended to contact a dentist.
- It is not recommended to use the product on children, people with dental implants or with open sores in the mouth.

12

TIGIOO
Edisonweg 10
1821BN Alkmaar

Model nr.: 0022.1
Versie: 04-2021