# MEDISANA®

**BS 430 CONNECT**  
**BS 440 CONNECT**

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### IMPORTANT INFORMATION! RETAIN FOR FUTURE USE!

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

- Use the device only according to its intended purpose as specified in the instruction manual.
- The warranty will be invalidated if the device is used for purposes other than those for which it is intended.
- This device is not suitable for pregnant women!
- Weighing only is no problem, but bodyfat analysis should be avoided for children because of the current impulse. This impulse is flowing through the body during the bodyfat measurement process.
- These scales are not suitable for people with pacemakers or other medical implants. The body fat information may be inaccurate in people with diabetes or other medical/physical limitations. The same also applies to very highly trained athletes.
- Do not use the device if it is not working correctly, if it has been dropped or has fallen into water or has been damaged.

Protect the unit against moisture. Keep the scale away from water. The surface of the scales is slippery when wet. Keep the surface dry.

**Risk of slipping!**  
Do not stand on the scale with wet feet. Do not stand on the scale wearing socks.

**Risk of tipping over!**  
Do not place the scale on an uneven surface. When weighing yourself, do not stand on one side or the edge of the scales.


We hereby declare that bathroom scales model BS 430 connect / BS 440 connect, art. 40422/40423 complies with the fundamental requirements of the European RTTE Directive 1999/5/EC. You can request the complete EC Conformity Declaration from Medisana AG, Jagenbergstrasse 19, 41468 Neuss, Germany, or you can also download it from the Medisana homepage ([www.medisana.com/bs430connect](http://www.medisana.com/bs430connect)) ([www.medisana.com/bs440connect](http://www.medisana.com/bs440connect)).

For further information und country-specific service addresses visit [www.medisana.com/bs430connect](http://www.medisana.com/bs430connect) / [www.medisana.com/bs440connect](http://www.medisana.com/bs440connect).

### Simple weight measurements with the „Step-On“ function


1 Simply stand on the scales.

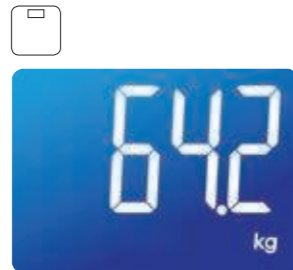2 Your weight is shown.


### Setting the personal data in the scales

The scales can save the data (gender, athletic form, height, age) of 8 people. The values preset in the scales are: Male, 165 cm, 30 years old, kg.

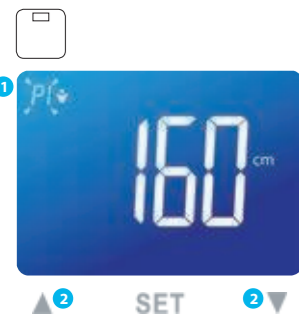
1 Press the SET button until the disk space blinks in the display.

2 Select the desired disk space using the ▲/▼ buttons.

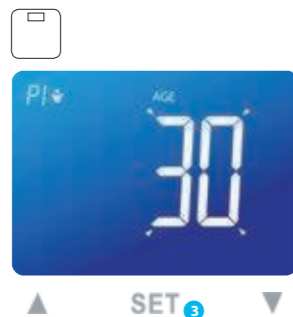
3 Press the SET button to save the selection. Continue accordingly to make the settings for gender, athletic condition, height and age.

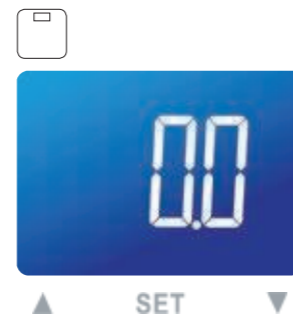4 After you have set your personal data, the display shows "0.0".

5 To finish the scales settings, simply stand on the scales.

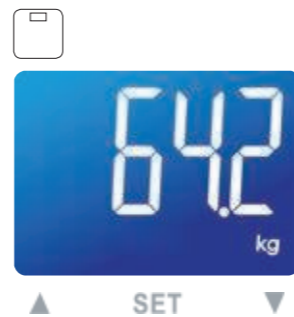
6 The display shows the following measurements:  
• Weight  
• Proportion body fat  
• Proportion body water  
• Proportion muscle mass  
• Bone mass  
• BMR (calory requirement)  
• BMI

### Connect your scales with your smartphone

The scales and the smartphone must be connected **once** using Bluetooth®. Please install the current VitaDock+ app and follow the instruction steps given in the following.


1 Free download of the VitaDock+ app in the app store or Google Play store.

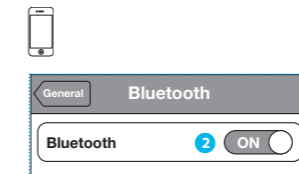2 Activate Bluetooth® on your iOS or Android appliance.

3 Start the app.

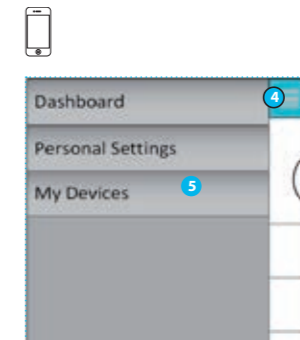4 Select the "Settings" area in the app.


5 Click there on "My appliances".


6 To set up the connection between the smartphone and scales, press the SET button, select your disk space and stand on the scales.

7 Select your scales to pair them with the app.

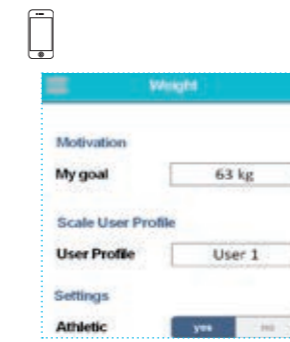
8 In order to be able to use all the functions of your scales, please update the personal settings in the app.

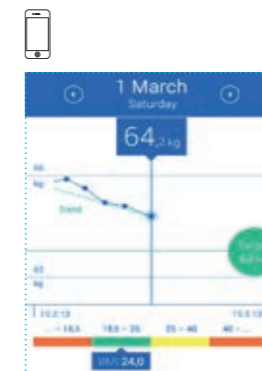9 The measurement data are automatically transferred to and stored by the app.

### Weight measurement using your smartphone and the VitaDock+ app

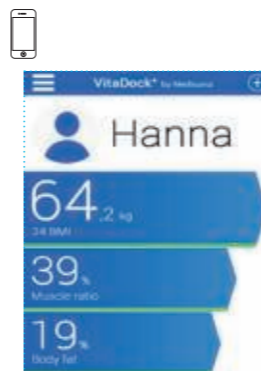
1 Start the app.

2 Stand on the scales.

6 The display shows the following measurements:  
• Weight  
• Proportion body fat  
• Proportion body water  
• Proportion muscle mass  
• Bone mass  
• BMR (calory requirement)  
• BMI

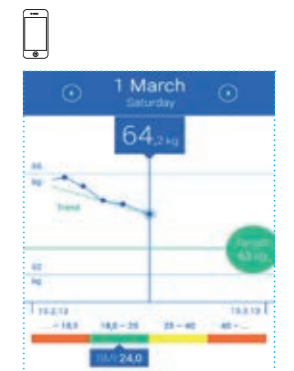4 The measurement data are automatically transferred to and stored by the app.

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