



Hydra.Pen
Skin Aesthetics Force



USER MANUAL

PLEASE READ THE MANUAL BEFORE USE



No.	Description
①	Power Interface
②	Speed Display
③	Power / Speed Button
④	Retaining Ring
⑤	Needle Cartridge



① Needle Cartridge

② Retaining Ring

③ Protective Cover

④ Soft Plug

Wireless Model


The instrument use both built-in battery and AC power


1.Plug in DC adaptor to the jack

2.Plug in DC adaptor to AC outlet

3.To charge the battery: Plug the adapter into an outlet and connect it to charger

Charging Status

 The power indicator stays red while charging

 The power indicator turns green after charging completed



- ① Put essence into the empty needle cartridge
- ② Put the soft plug at the bottom of needle cartridge
- ③ Remove the protective cover
- ④ Assemble the needle cartridge to the Hydra pen
- ⑤ Assemble the retaining ring to the hydra pen
- ⑥ Start the hydra pen. Press the power button for 2-3 seconds to start or shut down the instrument. It starts from slow to fast speed

Attention

In order to protect the instrument well, do not invert the instrument to prevent the liquid backflow



Needle cartridge introduction

4



Nano-HS



Nano-HR



H12



Packing list

Item	Quantity
Hand Piece	1 pc
Needle Cartridge	2 pcs
Manual	1 pc
Adapter	1 pc
Power Cable	1 pc
Retaining Ring	1 pc



Specification

5

Product`s name: Hydra pen H2
Adapter: 5V, 1000mA
Weight: 430g
Packaging size: 200*100*80mm



Maintenance

General maintenance

Keep it with its power off. Wipe its surface out with a dry cloth

Keep it away from direct sunlight or high temperature and humidity

Do not leave the measurement, repair or alter of the device to non-professional staff


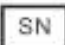





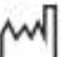








Make sure to inspect the device and adapter before use

When you use the device which has not been used for a long time, make sure to check if it works properly with safety

Warranty

The warranty period of Hydra pen H2 is 6 months. We offer repair to our products from the day they are purchased. If any of our products has damages or natural trouble. We offer repair or replacements



	Alternating current		Serial No
	Direct current		EU Representative
	Class I equipment		CE mark
	Attention, consult accompanying documents		Manufacturer
	Date of manufacture		Floating applied part
	Please be sure to read instructions before use and follow as described.		
	Do not apply it to the man as follow: 1. Who has acute disease or bloody disease over 38°C 2. Who has allergy and hemophilia 3. Who has melanosis of the face 4. Who has high internal heat 5. Who regularly visits hospital 6. Who takes medicines 7. Who is under 12 years old		
	Do not use it on the power other than AC 220V. It may cause fires or malfunction		
	Do not use it near the place with high temperature, dust or humidity (Near the humidifier). It may cause malfunction.		
	Do not repair, disassemble or alter the controller and the body. It may cause malfunction.		
	Do not touch the device or detach the power cords with wet hands. It may cause electric shock.		
	Do not use the stripped cord or propped cord. It may cause fires(Be sure to replace it with another).		
	Do not use the propped cord which has crack or scares. It may cause scares on the skin.		
	When you don't use it for a long time. unplug it from the outlet.		
	Be sure to have an inspection for test. If any troubles are found, get it repaired.		
	No one other than the operator can touch it. It may cause accidents or malfunction.		



Hydra.Pen
Skin Aesthetics Force